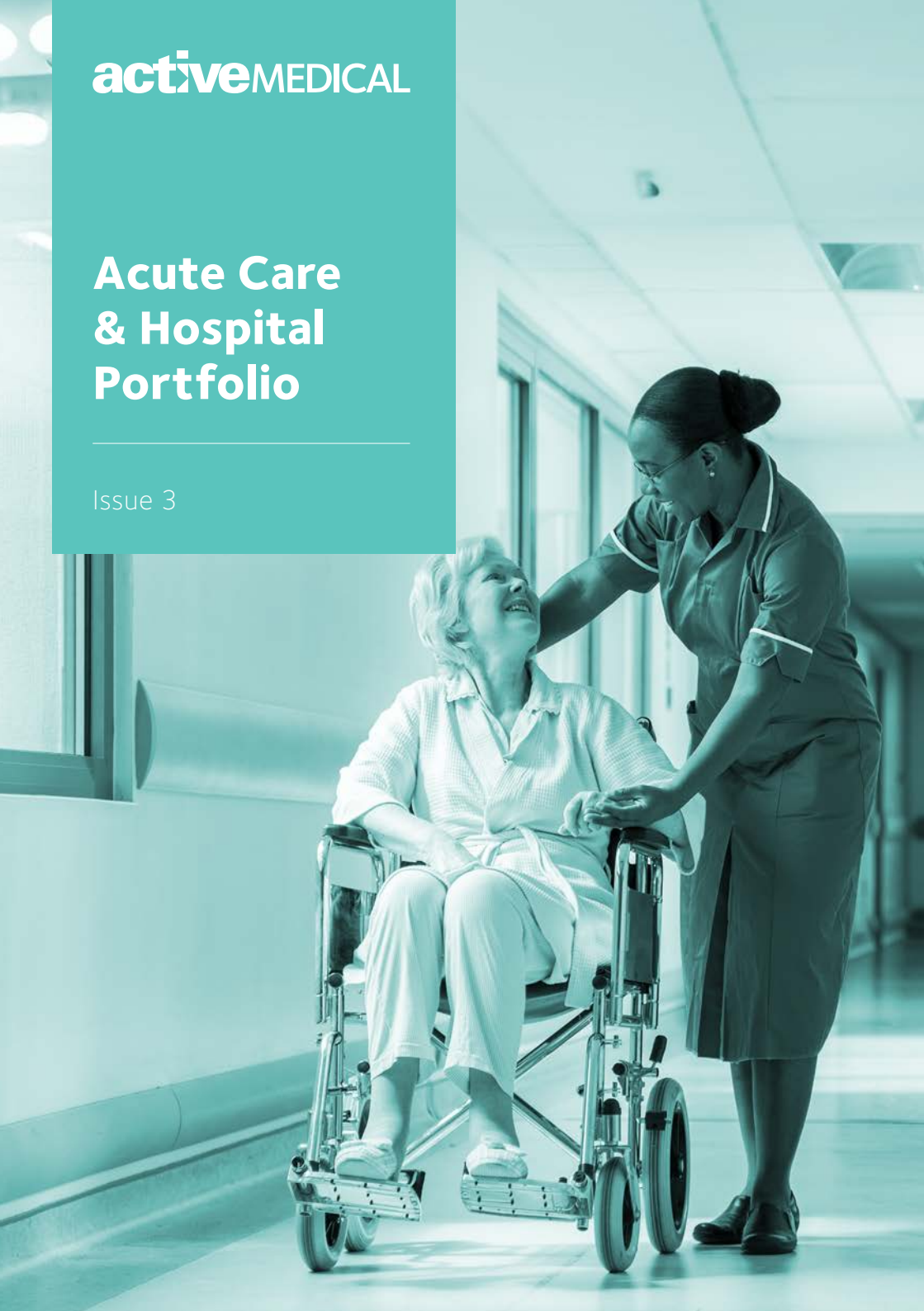


activeMEDICAL

Acute Care & Hospital Portfolio

Issue 3





Contents

01. Bath, Shower & Toilet	7
02. Patient Handling	37
03. Pressure Care	73
• Cushion Matrix	100
• Mattress Matrix	102
04. Wheelchairs	105
05. Seating	116
06. Mobility Solutions	138
07. Bedroom Solutions	146
08. Patient Safety	179
09. Treatment Tables	187
10. Mental Health Range	194

We are Here to Help

For over 30 years, Active Medical Supplies has delivered safe and innovative solutions for the care of clients in the Aged Care, Hospital, School and Home Care sectors. We strive to provide the highest quality equipment for Special Schools that provides a safe environment and will withstand the test of time.

We partner with experts in the field of design and manufacturing to provide the very best resources to ensure that your loved ones are provided with the correct equipment to receive optimum care.

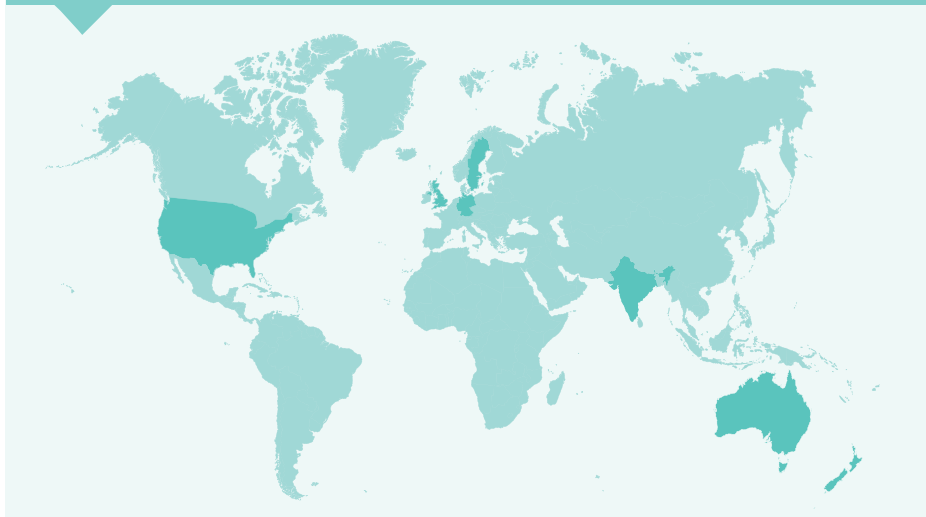
We have a commitment to lower our carbon footprint and protect our environment. We utilise our own delivery vehicles and use sustainable materials for packaging to reduce environmental damage. We are also proud to provide beds and equipment with world leading, energy-saving electronics.

For your peace of mind our product is manufactured to very high quality and we have a secure supply chain. 90% of our products are sourced from Australia, Europe, USA, Canada, United Kingdom, Spain, India, Sweden, Vietnam and Taiwan.



Procurement & Distribution

Active Medical has developed a network of manufacturers in Australia and around the world – from Germany and Sweden to the United States of America and Canada. We consult directly with industry professionals and ensure that all our products are ethically sourced and supplied and meet environmental standards.



*Disclaimer

Active Medical Supplies does not take responsibility for incorrect application of products. This publication is a guide only to assist in selecting products suitable to your clients. Please use in conjunction with professional advice. Products and information listed in this catalogue are subject to change. Website and digital information supercedes all print. For the most up to date information, please contact us directly.

Hospital & Aged Care Supplies

From sophisticated intensive care beds with alternating therapy systems, to a fundamental wheelchair, Active Medical Supplies a comprehensive range of hospital and aged care equipment that is sourced, designed and engineered to meet uncompromising standards for the care and safety

of patients and carers alike. The quality of product and operating simplicity of our equipment prioritises the wellbeing of patients, work place efficiency and a greater whole of life return on your investment.

Bariatric Services

An early appreciation and understanding of the complexities faced by bariatric patients and carers has seen Active Medical operating at the forefront of innovation in this specific care sector. Our comprehensive range of products and system solutions, backed by specific product knowledge has made us the preferred supplier to many of the country's leading providers of bariatric care. Our bariatric products are ideal for use in Hospitals, Aged Care and Home Care and meet the highest standards for efficient and safe use by patients and carers alike.



Dementia

With the right help and care, people who suffer with dementia can lead a happier, fuller life. Working hand in hand with international suppliers, such as Memoriana Furnishing Concepts, Active Medical is actively sourcing innovative products and holistic models of care that respond to the

demands of patients with dementia and their carers. All of our products are quality tested and are ideal for creating a space that is both relaxing and efficient, keeping the patient comfortable and giving the carer or staff member space to carry out their duties efficiently.



Home Based Care

Active Medical understands that home-based care sometimes demands a varied response to the hospital and aged care sectors. To support the needs of individual carers, Active Medical provides a range of specific needs equipment backed by the service and product knowledge of a team with a combined 145 years of experience. Our people can take you through the range of options available for your individual situation and demonstrate proper and practical use at our customised showroom.

Pressure Care

Pressure relieving and pressure distribution devices are high-accepted solutions for patients who experience a great deal of pain and discomfort. Our array of products focus on pressure injury prevention and treatments and our pressure care products are clinically proven to prevent pressure sores and ease pain.

Active Medical is dedicated to supporting the healthcare industry and is motivated by the desire to work for our community and increase the quality of life. Our service and support provide hospitals, aged care and home care facilities products to safeguard patients from pressure injuries.

Manual Handling

As part of our commitment to patient care and safety, we provide manual handling solutions to improve mobility for patients and carers alike. Our array of product encompasses solutions that will aid in movement, transferring and repositioning patients while minimising injuries. Solutions comprise innovative and user-friendly based products for hospitals, aged care and other medical facilities.





Our Commitment to Quality

All our products at Active Medical Supplies are TGA listed and comply with the latest international standard's where applicable.

Active Medical Supplies is helping the resellers provide better advice to clients and professionals. This is achieved through product innovation, cost reduction activities and compliance with ISO 9001 requirements.

This Quality Policy seeks to actively encourage stakeholders and members of the company to continually communicate about effective means of improvement for product quality and customer requirements. We will sustain this commitment by the involvement and further training and development of all our employees.

To assist us to do this, Active Medical Supplies undertakes to operate and maintain documented quality management systems which satisfy the appropriate controlling requirements linked to the inspection or certification service being provided.

This quality management system is designed to meet the requirements of AS ISO 9001:2015 Quality Management Systems.

01.

Bath, Shower & Toilet

Shower Chairs	8
Transfer Benches	9
Over Toilet Frames	11
Shower Trolleys	13
Shower Commodes	16
Accessories	28
Commode Seats	34

DVA

NDIS

MASS

NDIS

MASS

NDIS



Staebel Aluminium Shower Chair with Arms

This Shower Chair provides a comfortable seat for those unable to stand whilst showering. The armrests provide security when entering or exiting the shower, and assist the user on and off the seat.

- Strengthened frame for users needing extra room and weight capacity
- Seat height and angle adjusts to suit the individuals requirements
- One piece moulded supportive seat and backrest

Staebel Extracare Shower Chair with Arms

The Extracare Shower Chair provides a comfortable seat for those unable to stand whilst showering. The armrests provide security when entering or exiting the shower, and assist the user on and off the seat.

- Strengthened frame for users needing extra room and weight capacity
- Seat height and angle adjusts to suit the individuals requirements
- One piece moulded supportive seat and backrest
- Small drainage slots allow water to escape and remove risk of entrapment

Bariatric Shower Chair with Padded Seat

Padded seat and curved backrest provide comfort and postural support

- Powder coated steel frame
- Zinc treated to prevent rusting
- Made in Australia
- Suitable for Bariatric Users



Specifications:

Product Code	1315J
Seat Width	440 mm
Seat Depth	430 mm
Seat Height	400 - 550 mm
Bath Rim Clear	460 mm
Backrest Height	350 mm
Overall Width	570 mm
Overall Depth	540 mm
SWL	200 kg

Specifications:

Product Code	1320J
Seat Width	440 mm
Seat Depth	430 mm
Seat Height	400 - 550 mm
Width Between Arms	500 mm
Backrest Height	350 mm
Overall Width	620 mm
Overall Depth	570 mm
SWL	200 kg



Specifications:

Product Code	1579
Seat Depth	450 mm
Seat Width	600 mm
Seat Height	520 mm
SWL	200 kg



Staebel Shower Stool with Padded Seat

The aluminium shower stool with arms, provides security when entering or exiting the shower and assist the user on and off the seat.

- Aluminium construction for long life
- Plastic Seat
- Compact size for the smaller shower recess
- Adjustable seat height for the individual's requirements

Staebel Shower Stool with Plastic Seat

The aluminium shower stool with arms, provides security when entering or exiting the shower and assist the user on and off the seat.

- Aluminium construction for long life
- Plastic Seat
- Compact size for the smaller shower recess
- Adjustable seat height for the individual's requirements

Staebel Shower Stool with Arms & Back

This shower stool with arms and backrest, provides security when entering or exiting the shower and assist the user on and off the seat.

- Compact size for the smaller shower recess
- Padded seat and backrest for user comfort
- Adjustable seat height for the individual's requirements



Specifications:

Product Code	1567
Seat Width	370 mm
Seat Depth	290 mm
Seat Height	520 - 670 mm
Width Between Arms	460 mm
Armrest Height	640 - 840 mm
Overall Width	550 mm
Overall Depth	370 mm
SWL	125 kg



Specifications:

Product Code	1565J
Seat Width	370 mm
Seat Depth	290 mm
Seat Height	520 - 670 mm
Width Between Arms	460 mm
Armrest Height	640 - 840 mm
Overall Width	550 mm
Overall Depth	370 mm
SWL	175 kg

Also Available:

Shower Stool with Padded Seat	
Product Code	1567



Specifications:

Product Code	1569
Seat Width	370 mm
Seat Depth	290 mm
Seat Height	520 - 670 mm
Width Between Arms	460 mm
Armrest Height	640 - 840 mm
Overall Width	570 mm
Overall Depth	450 mm
SWL	125 kg

DVA
NDIS
MASS

NDIS
MASS



Standard Transfer Bench

Designed for individuals with restricted mobility, the Standard Transfer Bench provides safety and independence while transferring in and out of the bath tub. This durable aluminium transfer bench is adjustable in height.

- 3 piece plastic non-slip seat with drain holes
- Hand rail for added support while transferring
- Reversible backrest to suit every situation
- Welded frame

Padded Transfer Bench

Designed for individuals with restricted mobility, the padded Transfer Bench provides safety and independence while transferring in and out of the bath tub. This durable aluminium transfer bench is adjustable in height.

- 3 piece padded vinyl seat
- Hand rail for added support while transferring
- Reversible padded backrest to suit every situation
- Welded frame



Specifications:	
Product Code	1018J
Seat Width	670 mm
Seat Depth	410 mm
Seat Height	445 - 570 mm
Bath Rim Clear	370 - 495 mm
Seat-Back Height	410 mm
Overall Width	790 mm
Overall Depth	540 mm
SWL	140 kg

Options Available:	
Leg Extension With Suction Cups	
Product Code	1018J - 3
Leg Extension With Stoppers	
Product Code	1018J - 4



Specifications:	
Product Code	1018J - 1
Seat Width	670 mm
Seat Depth	410 mm
Seat Height	460 - 585 mm
SWL	140 kg

Options Available:	
Leg Extension With Suction Cups	
Product Code	1018J - 3
Leg Extension With Stoppers	
Product Code	1018J - 4

NDIS

MASS



Juvo Heavy Duty Transfer Bench

The Juvo Heavy Duty Transfer bench is both compact and sturdy providing safe and easy access for bathing or showering.

With a weight limit of 200kg this suits a vast range of clients whilst remaining lightweight for easy manoeuvrability. The Heavy Duty bench features height adjustable legs, non slip feet, the option of suction cups and reversible backrest to suit all bath types. With a textured non-slip seat with drainage holes and side rail the transfer bench provides stable seating for our clients.

- Leg Extensions available.
- Suitable for Bariatric Users



Specifications:

Product Code	1785
Seat Width	680 mm
Seat Depth	410 mm
Seat Height	470 - 570 mm
Bath Rim Clear	390 - 490 mm
Back Height	360 mm
Overall Width	780 mm
Overall Depth	550 mm
SWL	200 kg

NDIS



Staebelextracare Over Toilet Frame

This Extracare over toilet frame, one of four available sizes, features a zinc treated, epoxy coated frame with armrests and plastic seat and splash guard.

- Ideal for people with hip or knee injuries or those who need an increased seat height
- Bowl available as an accessory for commode use.
- Splash guard included



Specifications:

Product Code	1770J
Width Between Arms	500 mm
Seat Height	475 - 625 mm
Overall Width	620 mm
Armrest Height	700 - 850 mm
Footprint (W x D)	620 x 610 mm
SWL	150 kg

Also Available:

Product Code	1771	1772	1773
Width Between Arms	550 mm	600 mm	650 mm
Seat Height	500 mm	500 mm	500 mm
SWL	175 kg	195 kg	215 kg



DVA
NDIS
MASS



NDIS

Staebel Aluminium Over Toilet Frame

This aluminium over toilet frame is adjustable in height and fits over the toilet. This frame has been designed to assist the user who has difficulty sitting down or standing up by raising the toilet seat height.

- Aluminium construction for long life
- No risk of rust
- Removable splash guard included
- Lift up seat that removes for cleaning
- Ideal for those with hip or knee injuries or who require an increased seat height
- Unique stacking design saves space when storing more than one frame
- Multi D leg adjustment ensures stability on uneven surfaces
- Push button height adjustment
- Non slip guides on base of legs

Staebel Extra Care Over Toilet Frame

This Staebel Extra Care over toilet frame, one of four available sizes, features an Aluminium frame with armrests and plastic seat and splash guard.

- 500mm wide seat width
- Ideal for people with hip or knee injuries or those who need an increased seat height
- Bowl available as an accessory for commode use.
- Splash guard included
- 150kg Safe Weight Limit



Specifications:	
Product Code	1755J
Width Between Arms	445 mm
Seat Height	460 - 610 mm
Overall Width	550 mm
Armrest Height	700 - 840 mm
Footprint (W x D)	550 x 630 mm
SWL	150 kg



Specifications:	
Product Code	1750J
Seat Width	500 mm
Seat Height	500 mm
Width Between Arms	500 - 600 mm
Armrest Height	700 - 850 mm
SWL	150 kg



Staebel Junior Shower Trolley

- Shower trolley with one-section trolley top
- Trolley top length: 150 cm, trolley top width: 60 cm
- Maximum dimensions: 167 cm (length) and 76 cm (width)
- Adjustable 5° trolley top inclination
- Hydraulic lowering/raising of trolley top via a console from 54 cm to 97 cm
- Castors with diameter of 150 mm
- Central braking system
- Easy steering system (for ease of maneuvering trolley's front left castor can be locked in direction of steering)
- Side rails – foldable, height: 18,5 cm with safety release mechanism (material – powder coated stainless steel)
- Metal base cover
- Water-resistant shower mattress (height: 27 cm)
- Mattress integrated drain from 645 – 2000 mm long flexible hose
- Water-resistant head pillow – semi-cylinder



Specifications:

Product Code	1214ST
Length	1670 mm
Width	760 mm
Castors	150 mm
SWL	170 kg

Staebel Shower Trolley

- Shower trolley with one-section top
- Central braking system
- Easy steering system (for ease of maneuvering trolley's front left castor can be locked in direction of steering)
- Side rails – foldable, height: 185 mm with safety release mechanism (Material – powder coated stainless steel)
- Metal base cover
- Water-resistant shower mattress (height: 270 mm)
- Mattress integrated drain with 645 – 2000 mm long flexible hose
- Water-resistant head pillow – semi – cylinder



Specifications:

Product Code	1209ST
Shower Trolley Top Width	650 mm
Shower Trolley Top Length	1900 mm
Shower Trolley Top Height	540 – 970 mm
Height with Castors	600 – 1030 mm
Trolley Top Inclination	0° – 5°
Length	2070 mm
Width	1070 mm
SWL	170 kg

Also Available:

Staebel Shower Trolley with Backrest	
Product Code	1211ST



Staebel Bariatric Shower Trolley

- Hydraulic lowering/raising of the stretcher top via a console lift
- Lowering/raising range: 540 mm to 970 mm
- Castors with a diameter of 150 mm
- Central braking easy steering system (for ease of maneuvering the trolley's front left castor can be locked in the direction of steering directional wheel)
- Side rails – foldable, height: 185mm with safety release mechanism (Material – powder coated stainless steel)
- Metal base cover
- Water resistant shower mattress (height: 270mm)
- Mattress integrated drain with a 645 up to 2000mm long flexible hose
- Water resistant head pillow - semi - cylinder

Staebel Bariatric Electric Shower Trolley

- Electric lowering/raising of the stretcher top via a console lift
- Lowering/raising range: 540 mm to 970 mm
- Castors with a diameter of 150 mm
- Central braking easy steering system (for ease of maneuvering the trolley's front left castor can be locked in the direction of steering directional wheel)
- Side rails – foldable, height: 185 mm with safety release mechanism (Material – powder coated stainless steel)
- Metal base cover
- Water resistant shower mattress (Height: 270 mm)
- Mattress integrated drain with a 645 up to 2000mm long flexible hose
- Water resistant head pillow semi cylinder



Specifications:	
Product Code	1216ST
Shower Trolley Top Width	900 mm
Shower Trolley Top Length	2000 mm
Shower Trolley Top Height	540 - 970 mm
Height with Castors	600 - 1030 mm
Trolley Top Inclination	0° - 5°
Length	2170 mm
Width	1070 mm
SWL	250 kg



Specifications:	
Product Code	1217ST
Shower trolley top Width	900 mm
Shower trolley top Length	2000 mm
Shower trolley top Height	540 - 970 mm
Height with Castors	600 - 1030 mm
Trolley top Inclination	0° - 5°
Length	2170 mm
Width	1060 mm
SWL	250 kg



Staebel Shower Trolley with Backrest

- Shower trolley with two-section top
- Maximum dimensions: 207 cm (length) and 81 cm (width)
- Trolley top length: 190 cm, trolley top width: 65 cm
- Adjustable 5° trolley top inclination
- Electrically lowering/raising of trolley top via console from 54 cm to 97cm
- Electrically adjustable head section inclination 45° (head section length: 70 cm)
- Castors with diameter of 150 mm
- Central braking system
- Easy steering system (for ease of maneuvering trolley's front left castor can be locked in direction of steering)
- Side rails – foldable, height: 18.5 cm with safety release mechanism (material – powder coated stainless steel)
- Metal base cover
- Water-resistant shower mattress (height: 27 cm)
- Mattress integrated drain with 645 – 2000 mm long flexible hose
- Water-resistant head pillow



Specifications:

Product Code	1211ST
Length	2070 mm
Width	810 mm
Castors	150 mm
SWL	170 kg



Staebel Profiling Shower Trolley Electric

Profiling shower trolley has been designed to facilitate shower bathing and hygiene routines, for patients and service users with a wide range of postural and cognitive challenges.

It's three section profiling mattress platform can be adjusted to provide specific postural support in the four key regions of the anatomy, head, back, seat, and legs providing maximum safety, support, and comfort.

- Shower trolley with four-section top
- Electrically adjustable Trendelenburg 0° till 2° , antitrendelenburg 0° till 10°
- Foldable side rails
- Castors with diameter of 150 mm with central braking system
- Easy steering system (for ease of maneuvering trolley's front left castor can be locked in direction of steering)
- Metal base cover
- Water-resistant shower mattress
- Mattress integrated drain with 100cm long flexible hose
- Water-resistant head pillow – semi-cylinder




Specifications:


Product Code	1213ST
Shower trolley top Width	650 mm
Shower trolley top Length	1900 mm
Trolley Head Section Inclination	0° - 60°
Electric Hi-Lo	520 - 920 mm
Length	2070 mm
Width	810 mm
SWL	200 kg

Shower Commodes Innovation by Design

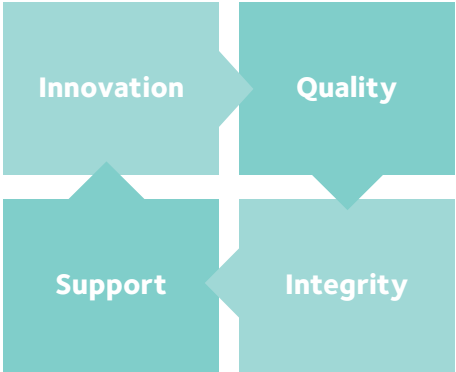Leading designer, manufacturer & supplier of Stainless Steel Shower Commodes



Build your Standard Shower Commode!

Choose your Commode	<input type="checkbox"/>	Mobile	<input type="checkbox"/>	Static		
Fixed or Tilted Commode	<input type="checkbox"/>	Fixed	<input type="checkbox"/>	Tilted		
Attendant or Self Propelled	<input type="checkbox"/>	Attendant Propelled	<input type="checkbox"/>	Self Propelled		
Footplate Options	<input type="checkbox"/>	Sliding Footplate	<input type="checkbox"/>	Swing-Away Footplate		
<div></div> Choose your seat width	<input type="checkbox"/>	360 mm	<input type="checkbox"/>	410 mm	<input type="checkbox"/>	460 mm
	<input type="checkbox"/>	510 mm	<input type="checkbox"/>	560 mm	<input type="checkbox"/>	610 mm
	<input type="checkbox"/>	660 mm	<input type="checkbox"/>	710 mm	<input type="checkbox"/>	910 mm (Trix)



Why not enhance your clients situation with high quality accessories! Options include customisable armrests, backrests, footplates & so much more!

See page 28 for all accessories!

NDIS

MASS



Fully adjustable height for improved client positioning



Fully adjustable depth to assist with client positioning

Trix 3 In 1 Shower Chair

- Assists patient with showering & toileting
- Clear rear open space enabling positioning over the toilet
- Easy access for hygiene with open armrest design
- Deeper access front trough for easy cleaning
- Integral pan carrier glides. Back saver design, assists with emptying the bowl, can be removed half way, self-supporting while the bowl/pan is removed
- Grade 304 stainless steel heavy wall thickness tubing forming robust frame to withstand extensive wet area use
- Reinforced seat/back frame and armrests – weight bearing
- Easy removable seat, using the Juvo Registered Design "Klipps" System
- Transportation wheels for effortless manoeuvring. Interchangeable so chair can be pulled by back or armrest
- Adjustable seat height for improved client positioning
- Fully adjustable seat depth to assist with client positioning
- Selection of seamless Polyurethane padded seats & pressure care seats (open or closed fronts)
- Armrests project beyond the seat to assist patient in sitting down and standing up
- Splayed armrest design to eliminate pinch point to patient when sitting down and standing up
- Flexible impervious non-woven back straps conforming to client shape for support and comfort. Antibacterial and nonsuction designed straps to avoid skin tear when client is wet
- Curved Backrest design to allow for natural prominence of the gluteal shelf



Specifications:

Product Code	1780
Seat Width	610 mm
Seat Height (25mm increments)	410 - 560 mm
Seat Depth (25mm increments)	470 - 570 mm
Width between Arms	660 mm
Overall Width	780 mm
Overall Depth	650 - 750 mm
Toilet Max Clearance H x W	540 x 500 mm
SWL	400 kg

Also Available:

1783	1779	1781	1782
460 mm	510 mm	710 mm	910 mm
410 - 560 mm	410 - 560 mm	410 - 560 mm	410 - 560 mm
470 - 570 mm	470 - 570 mm	470 - 570 mm	495 - 595 mm
510 mm	560 mm	760 mm	960 mm
630 mm	680 mm	880 mm	1080 mm
650 - 750 mm	650 - 750 mm	650 - 750 mm	710 - 810 mm
390 x 500 mm	440 x 500 mm	640 x 500 mm	840 x 500 mm
400 kg	400 kg	400 kg	400 kg



- DVA
- NDIS
- MASS



Juvo Attendant Propelled Shower Commode with Sliding Footplate

The ultimate in strength, quality and stability this Juvo mobile shower chair is ideal for the Petite client! The wheelbase design prevents forward tipping and has proven to be an effective solution in all areas from home to aged care to hospital.

- Attendant propelled
- Quality stainless steel frame for longevity in a demanding environment
- Flip up armrests for ease of side transfers
- To utilise as a commode you must source the slide-in pan rack

- Armrests can be flipped so that they are positioned in front of occupant to reduce risk of slipping forward when being wheeled
- Smooth, slatted backrest design to facilitate easy cleaning
- Generous size sliding stainless steel footplate, weight bearing

***Please note that the seat for this product is sold separately**
Optional accessories are available, refer to page 28



Specifications:	
Product Code	1904J
Seat Width	410 mm
Seat Depth	435 mm
Width Between Arms	440 mm
Overall Width	490 mm
SWL	200 kg

Seat colours can be varied:

Lemon TwistLimeLipstickRegalEmerald

Sizes Available:					
1910J	1914J	1912J	1930J	1936J	1971J
460 mm	510 mm	560 mm	610 mm	660 mm	710 mm
435 mm	435 mm	435 mm	435 mm	435 mm	435 mm
490 mm	540 mm	590 mm	640 mm	690 mm	740 mm
540 mm	590 mm	640 mm	690 mm	740 mm	790 mm
200 kg	300 kg	300 kg	300 kg	300 kg	400 kg



Juvo Mobile Shower Commode with Swing Away Legrests

The ultimate in strength, quality and stability the Juvo mobile shower commode is built for heavy use having a 200kg SWL! The long wheelbase stops forward tipping and has proven to be an effective solution in all areas from home to hospital.

- Attendant propelled
- Assists transportation of the patient for showering and toileting
- Can be positioned over the toilet
- Durable total stainless steel frame
- Selection of seamless PU padded seats & pressure care seats (open or closed front)
- Removable seat
- Reinforced seat/back frame

- Swing away, removable and height adjustable legrests with flip up footplates
- Patented design, no protrusions
- Stainless steel 125mm castors (4 braking)
- Swing back independent armrests to assist in side transfer and patient access
- Two position height adjustable armrests
- Safety arms provide user security and safety
- Stainless steel pan rack guides
- This product is tested and certified to Australian Standards

***Please note that the seat for this product is sold separately**

Optional accessories are available, refer to page 28

Seat colours can be varied:



Lemon Twist



Lime



Lipstick



Regal



Emerald



Specifications:

Product Code	1915J
Seat Width	460 mm
Anterior Seat Height	590 mm
Posterior Seat Height	565 mm
Width Between Arms	465 mm
Overall Width	540 mm
SWL	200 kg

Sizes Available:

1905J	1916J	1917J	1935J
410 mm	510 mm	560 mm	610 mm
595 mm	590 mm	590 mm	590 mm
505 mm	565 mm	565 mm	565 mm
440 mm	520 mm	570 mm	660 mm
490 mm	610 mm	640 mm	710 mm
200 kg	300 kg	300 kg	300 kg



NDIS
MASS



Minimum Height



Maximum Height

Juvo ‘Tilting in Space’ Shower Commode

- Assists transportation of the patient for showering and toileting
- Gas strut tilt mechanism for smooth and easy tilting
- Multi-position 35 deg. tilt
- Can be positioned over the toilet
- Durable total stainless steel frame
- Selection of seamless PU padded seats & pressure care seats (open or closed front)
- Removable seat
- Reinforced seat/back frame
- Swing away, removable and height adjustable legrests with flip up footplates
- Patented TEARSAFE design, no protrusions
- Stainless steel 125mm castors (4 braking)
- Swing back independent armrests to assist in side transfer and patient access
- Adjustable seat height
- Two position height adjustable armrests
- Safety arms provide user security and safety
- Stainless steel pan rack guides

***Please note that the seat for this product is sold separately**
Optional accessories are available, refer to page 28

Electric upgrade available upon request. Call the team to discuss



Specifications:	
Product Code	1951J
Seat Width	460 mm
Overall Width	570 mm
Armrest Height - High	240 mm
Armrest Height - Low	200 mm
Backrest Height from Seat	900 mm
Width Between Arms	470 mm
SWL	200 kg

Seat colours can be varied:



Also Available:				
1950J	1952J	1956J	1961J	
410 mm	510 mm	560 mm	610 mm	
560 mm	590 mm	660 mm	590 mm	
240 mm	240 mm	240 mm	240 mm	
200 mm	200 mm	200 mm	200 mm	
900 mm	900 mm	800 mm	800 mm	
440 mm	520 mm	590 mm	610 mm	
200 kg	300 kg	300 kg	300 kg	

NDIS

MASS



Juvo Tilt & Recline Shower Chair

- Manual Reclining Shower Chair which can be positioned over the toilet
- Durable total stainless steel frame
- Selection of seamless PU padded seats & pressure care seats (open or closed front)
- Removable seat
- Reinforced seat/back frame
- Patented TEARSAFE design, no protrusions
- Stainless steel 125mm castors (4 braking)
- Swing away, removable and height adjustable legrests with flip up footplates
- Patented TEARSAFE design, no protrusions
- Swing back independent armrests to assist in side transfer and patient access
- Adjustable seat height
- Two position height adjustable armrests
- Safety arms provide user security and safety
- Stainless steel pan rack guides

***Please note that the seat for this product is sold separately**
Optional accessories are available, refer to page 28

Seat colours can be varied:



Lemon Twist

Lime

Lipstick

Regal

Emerald



Specifications:

Product Code	1976J
Seat Width	410 mm
Overall Width	560 mm
Overall Height	1465 – 1615 mm
Overall Depth	940 – 1600 mm
Width Between Arms	490 mm
Armrest Height	200 – 240 mm
Backrest Height from Seat	800 mm
Seat Tilt Angle	-5 ° - 38 °
SWL	200 kg

Sizes Available:

1977J	1978J	1980J
460 mm	510 mm	610 mm
560 mm	610 mm	730 mm
1465 – 1615 mm	1465 – 1615 mm	1465 – 1615 mm
940 – 1600 mm	940 – 1600 mm	940 – 1600 mm
490 mm	540 mm	610 mm
200 – 240 mm	200 – 240 mm	200 – 240 mm
800 mm	800 mm	800 mm
-5 ° - 38 °	-5 ° - 38 °	-5 ° - 38 °
200 kg	200 kg	300 kg



NDIS
MASS



Juvo Electric Hi-Lo Commode with Swingaway Legrests

- Independently operated electric control for both height and tilt
- Non-woven body contouring back straps with multiple backrest options including contour headrest
- Swing-away armrests with locking safety arms
- Seamless soft, easily removable polyurethane seat or enhanced pressure care impress seat option, including the Innovative Juvo Seat Klipps® also allows for custom fabricated seating options in addition to the standard Juvo seat options
- Standard side cut-outs allow for multiple seat access options without modification to the seat frame
- Universal pan carrier glides and carrier allow for multiple styles of bedpans and bowls to be utilised
- Robust and corrosion-resistant stainless steel frame with an open back for positioning over the toilet
- Tearsafe® swing away leg rests - no sharps remain to prevent skin tears
- New Linak slim-line LIFT50 battery and control system handheld remote
- Detachable light weight lithium battery IPX6 water-resistant, showersafe electric components
- Tilt range from -5 ° to 38 °
- 4 total locking castors for safe manual handling

***Please note that the seat for this product is sold separately**
Optional accessories are available, refer to page 28

Seat colours can be varied:



Lemon Twist Lime Lipstick Regal Emerald

Specifications:		Also Available:		
Product Code	1982J	1983J	1984J	1985J
Seat Width	460 mm	510 mm	560 mm	610 mm
Anterior Seat Height	680 - 875 mm	680 - 875 mm	680 - 875 mm	680 - 875 mm
Posterior Seat Height	430 - 640 mm	430 - 640 mm	430 - 640 mm	430 - 640 mm
Overall Width	635 mm	635 mm	635 mm	635 mm
Armrest Height	200 - 240 mm	200 - 240 mm	200 - 240 mm	200 - 240 mm
Overall Depth	1030 - 1600 mm	1030 - 1600 mm	1030 - 1600 mm	1030 - 1600 mm
Width Between Arms	485 mm	535 mm	585 mm	635 mm
SWL	200 kg	300 kg	300 kg	300 kg



Juvo Self Propelled Shower Commode with Sliding Footrest

- Assists transportation of the patient for showering and toileting
 - Can be positioned over the toilet
 - Durable total stainless steel frame
 - Selection of seamless PU padded seats & pressure care seats (open or closed fronts)
 - Removable seat
 - Reinforced seat/back frame
 - Sliding footrest
 - Patented TEARSAFE design, no protrusions
 - Now with braking stainless steel 125mm castors
 - Swing back independent armrests to assist in
 - Side transfer and patient access
 - Two position height adjustable armrests
 - Safety arms provide user security & safety
 - Stainless steel pan rack guides
 - This product is tested and certified to Australian Standards
- *Please note that the seat for this product is sold separately**
Optional accessories are available, refer to page 28

Seat colours can be varied:



Lemon Twist

Lime

Lipstick

Regal

Emerald



Specifications:

Product Code	1920J
Seat Width	460 mm
Anterior Seat Height	590 mm
Posterior Seat Height	565 mm
Width Between Arms	456 mm
Overall Width	520 mm
SWL	175 kg

Also Available:

1921J	1922J	1938J	1940J
510 mm	560 mm	660 mm	610 mm
590 mm	590 mm	590 mm	590 mm
565 mm	565 mm	565 mm	565 mm
530 mm	580 mm	680 mm	630 mm
580 mm	630 mm	730 mm	680 mm
300 kg	300 kg	300 kg	300 kg



Juvo Self-Propelled Shower Commode with Swing Away Legrests

- Assists transportation of the patient for showering and toileting
 - Can be positioned over the toilet
 - Durable total stainless steel frame
 - Selection of seamless PU padded seats & pressure care seats (open or closed fronts)
 - Removable seat
 - Reinforced seat/back frame
 - Swing away, removable & height adjustable legrests with flip-up footplates.
- Patented TEARSAFE design, no protrusions
 - Now with braking stainless steel 125mm castors
 - Swing back independent armrests to assist in
 - Side transfer and patient access
 - Two position height adjustable armrests
 - Safety arms provide user security & safety
 - Stainless steel pan rack guides
 - This product is tested and certified to Australian Standards

***Please note that the seat for this product is sold separately**
Optional accessories are available, refer to page 28



Seat colours can be varied:

Lemon TwistLimeLipstickRegalEmerald

Specifications:		Sizes Available:			
Product Code	1925J	1926J	1927J	1945J	1939J
Seat Width	460 mm	510 mm	560 mm	610 mm	660 mm
Anterior Seat Height	590 mm	595 mm	595 mm	595 mm	595 mm
Posterior Seat Height	565 mm	565 mm	565 mm	565 mm	565 mm
Width Between Arms	465 mm	540 mm	540 mm	540 mm	540 mm
SWL	200 kg	300 kg	300 kg	300 kg	300 kg



Juvo Height Adjustable Self Propelled with Swing Away Legrests

- Height adjustable up to 90mm front and rear axle
- Assists transportation of the patient for showering and toileting
- Can be positioned over the toilet
- Removable seat
- Reinforced seat/back frame
- Swing away, removable & height adjustable legrests with flip-up footplates.
- Patented TEARSAFE design, no protrusions
- Now with braking stainless steel 125mm castors
- Swing back independent armrests to assist in side transfer and patient access
- Two position height adjustable armrests
- Safety arms provide user security & safety
- Stainless steel pan rack guides
- This product is tested and certified to Australian Standards
- N.B. The legs elevated legrests shown do not come as standard

***Please note that the seat for this product is sold separately**
Optional accessories are available, refer to page 28

Seat colours can be varied:



Lemon Twist Lime Lipstick Regal Emerald



Specifications:

Product Code	1928J
Seat Width	460 mm
Anterior Seat Height	505 - 595 mm
Posterior Seat Height	475 - 565 mm
Width Between Arms	465 mm
Total Bowl Clearance (WxH)	360 x 500 mm
Overall Width	520 mm
SWL	200 kg

Also Available:

1929J	1946J	1974J	1997J
510 mm	560 mm	610 mm	660 mm
505 - 595 mm	505 - 595 mm	505 - 595 mm	505 - 595 mm
475 - 565 mm	475 - 565 mm	475 - 565 mm	475 - 565 mm
540 mm	590 mm	640 mm	660 mm
410 x 500 mm	410 x 500 mm	410 x 500 mm	410 x 500 mm
710 mm	760 mm	810 mm	860 mm
200 kg	300 kg	300 kg	300 kg

Height
Adjustable



Juvo Height Adjustable Attendant Propelled with Swing Away Legrests

- Attendant propelled
 - Height adjustable up to 90mm front and rear axle.
 - Assists transportation of the patient for showering and toileting
 - Can be positioned over the toilet
 - Selection of seamless PU padded seats & pressure care seats (open or closed front) (NB Pictured with hand padded closed front option)
 - Reinforced seat/back frame
 - Swing away, removable and height adjustable legrests with flip up footplates
 - Patented design, no protrusions
- Stainless steel 125mm castors (4 braking)
 - Swing back independent armrests to assist in side transfer and patient access
 - Safety arms provide user security & safety
 - Stainless steel pan rack guides
 - This product is tested and certified to Australian Standards

***Please note that the seat for this product is sold separately**
Optional accessories are available, refer to page 28



Seat colours can be varied:

Lemon TwistLimeLipstickRegalEmerald

Specifications:		Also Available:			
Product Code	1906J	1913J	1918J	1919J	1933J
Seat Width	410 mm	460 mm	510 mm	560 mm	660 mm
Anterior Seat Height	505 - 595 mm	505 - 595 mm	505 - 595 mm	505 - 595 mm	505 - 595 mm
Posterior Seat Height	475 - 565 mm	475 - 565 mm	475 - 565 mm	475 - 565 mm	475 - 565 mm
Width Between Arms	440 mm	490 mm	540 mm	590 mm	690 mm
Overall Width	490 mm	540 mm	590 mm	640 mm	740 mm
SWL	200 kg	200 kg	300 kg	300 kg	300 kg

DVA

NDIS

MASS



NDIS

MASS



Juvo Folding Self Propelled Commode with Swingaway Footplates

This Juvo Folding Self Propelled Commode with Swingaway Footplates Assists transportation of the patient for showering and toileting. It is made from Durable total stainless steel frame and can be positioned over the toilet

- Selection of seamless PU padded seats & pressure care seats (open or closed fronts)
- Reinforced seat/back frame
- Patented design, no protrusions
- Swing away, removable & height adjustable legrests
- Braking stainless steel 125mm castors
- Two position height adjustable armrests with safety arms
- Stainless steel pan rack guides
- Australian Standards Approved
- Patented folding design into standard suitcase size

***Please note that the seat for this product is sold separately.**



Specifications:

Product Code	1923J
Seat Width	460 mm
Anterior Seat Height	500 - 590 mm
Posterior Seat Height	475 - 565 mm
Width Between Arms	465 mm
Overall Width	520 mm
SWL	200 kg

Also Available:

Juvo Folding Self Propelled Commode, 510 mm

Product Code	1924J
--------------	-------

Juvo Paediatric Shower Commode

- Adjustable Backrest design - Adjust the backrest angle from 90 to 120 degrees to improve client comfort and access for the carer.
- Depth adjustable seat - Adjust the seat depth to suit smaller client by up to 75mm.
- Height Adjustable seat up to 90mm
- Self-propelled option available - frame ready to accept self-propelling kit
- Width adjustable
- This product is tested and certified to Australian Standards

***Please note that the seat for this product is sold separately**

Optional accessories are available, refer to page 28



Specifications:

Product Code	1903J
Seat Width	360 mm
Width Between Arms	390 mm
Anterior Seat Height	595 mm
Posterior Seat Height	565 mm
Overall Width	440 mm
SWL	136 kg

Bathroom Accessory Case Studies

Juvo Tilt & Recliner Commode with Contoured Headrest, Lateral Supports, Padded Armrests & elevated legrest



Juvo Attendant Propelled Commode with Sliding Footplate and one padded Armrests and one Arm Trough



Juvo Attendant Propelled Commode with Sliding Footplate and Padded Armrests, Pan/Bowl Carrier & Commode Bowl



Juvo Attendant Propelled Commode with Padded Armrests and a One Piece Footplate



Could not find the information you were after?

Contact our team for a consultation for more bath, shower & toilet accessories!



Juvo Angle Adjusted Backrest - 3 Strap

- Shower Commode Accessory
- Self Propelled with 3 Straps
- 340mm high to top strap
- 41cm, 46cm, 51cm, 56cm, 61cm and 66cm widths available

(Contact the team for more information)



Juvo Angle Adjusted Backrest - 4 Strap

- Shower Commode Accessory
- Self Propelled with 4 Straps
- 415mm high to top strap
- 41cm, 46cm, 51cm, 56cm, 61cm and 66cm widths available

(Contact the team for more information)



Backrest Extension - 4 Strap

- Suitable for Tilt in Space and Self Propelled with 4 Straps
- 415mm high to top strap
- 41cm, 46cm, 51cm, 56cm, 61cm and 66cm widths available

(Contact the team for more information)

Specifications:

Product Code 11991 - 46 - 3

Specifications:

Product Code 11991 - 46 - 4

Specifications:

Product Code 11945



Short Backrest Style

- Shower Commode Accessory
- Tilt in Space with 1 Strap
- 41cm, 46cm, 51cm and 61cm widths available

(Contact the team for more information)



Wheelchair Style Backrest- 3 Strap

- Shower Commode Accessory
- Self Propelled with 3 Straps
- 340mm high to top strap
- 41cm, 46cm, 51cm, 56cm, 61cm and 66cm widths available

(Contact the team for more information)



Wheelchair Style Backrest - 4 Strap

- Suitable for Tilt in Space and Self Propelled with 4 Straps
- 415mm high to top strap
- 41cm, 46cm, 51cm, 56cm, 61cm and 66cm widths available

(Contact the team for more information)

Specifications:

Product Code 11944 - 46 - 1

Specifications:

Product Code 11944-41-3

Specifications:

Product Code 11944 - 41 - 4



Backrest Extension - 5 Strap

- Suitable for Tilt in Space and Self Propelled with 5 Straps
- 490mm high to top strap
- 41cm, 46cm, 51cm, 56cm, 61cm and 66cm widths available

(Contact the team for more information)



Backrest Extension - 6 Strap

- Suitable for Tilt in Space and Self Propelled with 6 Straps
- 570mm high to top strap
- 41cm, 46cm, 51cm, 56cm, 61cm and 66cm widths available

(Contact the team for more information)



Waterfall Arm Rest

- Juvo waterfall armrest
- Comes in a pair

Specifications:

Product Code 11946

Specifications:

Product Code 11947

Specifications:

Product Code 19611 - 2



Juvo IV Pole

Designed for use with shower commodes, this heavy duty, adjustable stainless steel IV pole simply clamps onto the shower commode frame. Its height from floor adjusts 1200mm - 1800mm.



Stump Support Leg Extension

- Adjustable stump support for Juvo mobile shower commodes. Installed in replacement for normal design legrest on a needs basis.
- Adjustable angle and length to suit the individual user
- Left and Right side available (specify when ordering)



Contoured Headrest

- Contoured adjustable headrest to suit Juvo. Comes complete with mounting accessories.
- Fixing Clamp
- Small option available

Specifications:

Product Code 11942

Specifications:

Product Code 11943 - E

Specifications:

Product Code 19652



Chest / Pelvic Strap

- Adjustable Chest and Pelvic Strap suitable for Juvo Commodes.



Arm Trough

The arm trough is of benefit for the shower chair user who has upper limb weakness or paralysis. It attaches to the chair armrest via a sturdy clamp.



Lock Down Armrest

The lock down armrest option prevents the armrests from being flipped or pulled back and is useful in situations where this may pose a risk. The armrest is release by simply pulling the release button.

Specifications:

Product Code 11933

Specifications:

Product Code 19613

Specifications:

Product Code 19614



Transfer Bar / Safety Arm

- Supplied in a pair
- Handgrip
- Spring and button



Slip-On Lock Down Arm Kit

- Comes in a Pair
- No fixings required



One Piece Footplate

- One piece footplate for swing-away shower commode models
- Different widths for wider commodes

Specifications:

Product Code 11950 - T

Specifications:

Product Code 19620

Specifications:

Product Code 11963 - 46



NDIS

Elevated Legrest

Adjustable elevated legrest for Juvo mobile shower commodes. Installed in replacement for normal design legrest on a needs basis.

- Adjustable angle and length to suit the individual user
- Left and Right side available (specify when ordering)

Specifications:	
Product Code	Side
11940 - L	Left
11940 - R	Right

NDIS



Juvo Padded Armrest

- Shower chair accessory
- Padded urethane armrests including hardware to mount
- Sold singularly

Specifications:	
Product Code	19611



NDIS

Stump Support

Adjustable stump support for Juvo mobile shower commodes. Installed in replacement for normal design legrest on a needs basis.

- Adjustable angle and length to suit the individual user
- Left and Right side available (specify when ordering)

Specifications:	
Product Code	Side
11943 - L	Left
11943 - R	Right

NDIS



Commode Bowl with Handle

It is suitable for use with the current 7952PA bedside commode, and also with the current style of overtoilet frames and aids. Includes handle.

- Autoclavable
- Top diameter 290 mm
- Base diameter 210 mm
- Height 120 mm

Specifications:	
Product Code	1989



NDIS

Juvo Angle Adjusted Footrest

- Juvo Commode Part - Angle adjustable footrest
- Left and Right side available (specify when ordering)

Specifications:	
Product Code	Side
11931 - L	Left
11931 - R	Right

DVA

NDIS



Standard Bed Pan

These bedpans are made from polypropylene and are autoclavable. The standard bedpan is heavy duty and ideal for institutional use.

Specifications:	
Product Code	1990

NDIS



Lateral Support

- Quick release thoracic support
- Swing away
- Left and Right side available (specify when ordering)

NDIS



Calf Support Strap

Calf Support now has a colour matched, clip-on calf support for a complete, protective, hygienic solutions for adults and children using shower chairs.

NDIS



Juvo Oxygen Bottle Carrier

This heavy duty stainless steel oxygen bottle carrier clamps onto the rear of the shower commode frame, and can hold a 130 mm diameter oxygen bottle.

Specifications:

Product Code	Side
11999 - L	Left
11999 - R	Right

NDIS



Dropdown Pommel

- Height and depth adjustable



Specifications:

Product Code	11985
--------------	-------

Specifications:

Product Code	188216 - 22
--------------	-------------

NDIS



Juvo Commode Bowl & Lid

- Bowl with Lid and Handle
- Can be used with all Shower Commode types

Specifications:

Product Code	11936
--------------	-------

Specifications:

Product Code	11941
--------------	-------

NDIS

DVA



Pan/ Bowl Carrier

Pan and Bowl Carrier that can be used for different commodes.
Different sizing available.

Specifications:

Product Code	Size
19603	410 mm
19607	560 mm



Polyurethane Open Front Seat

This standard polyurethane seat is completely sealed for infection control and is ideal for the user who will not spend long on the shower chair.

- Seat is simply removed for ease of cleaning and/or replacement.
- Open Front

Polyurethane Closed Front Seat

This polyurethane seat is completely sealed for infection control and is ideal for the user who will not spend long on the shower chair.

- Seat is simply removed for ease of cleaning and/or replacement.
- Closed Front

Staebel Strato3 Commode Seat

This is a unique pressure management commode seat cushion combination specially developed for those at high risk of skin breakdown and pressure ulcers. It features a special hole positioning and seat depth, 60 mm cushioning to minimise pressure. The Strato3 is also available with a side bite.

- Completely sealed for infection control
- No sharp joins or seams
- Padded for increased pressure care
- Custom sizes are available in any size and shape.

* Contact us about customised seats.

Specifications:	
Product Code	Size
196430	410 mm
196440	460 mm
196560	510 mm
196560	560 mm
196610	610 mm
196660	660 mm

Specifications:	
Product Code	Size
19643C	410 mm
19644C	460 mm
19651C	510 mm
19655C	560 mm
19661C	610 mm
19666C	660 mm



Specifications:	
Product Code	19780



Staebel Deluxe Hand Padded Commode Seat - Open Front

The Hand Padded seat is completely sealed for infection control and has no sharp joins or seams. Deluxe model features generous padding for ultimate comfort. Easily removable for cleaning. Fits on to most commode types. Fits on to most commode types. Particularly ideal for users who have challenges with leg positioning.

- Completely sealed for infection control
- No sharp joins or seams
- Padded for increased pressure care
- Front opening
- Available in a range of different colours (Contact the team for more information)

Staebel Deluxe Hand Padded Commode Seat - Closed Front

The Hand Padded seat is completely sealed for infection control and has no sharp joins or seams. Deluxe model features generous padding for ultimate comfort. Easily removable for cleaning. Fits on to most commode types. Fits on to most commode types. Particularly ideal for users who have challenges with leg positioning.

- Completely sealed for infection control
- No sharp joins or seams
- Padded for increased pressure care
- Closed opening
- Available in a range of different colours (Contact the team for more information)



Specifications:

Product Code	19648D
Width	460 mm

Also Available:

Product Code	Width
196770	650 mm
19673D	660 mm
19691D	710 mm
19695D	780 mm



Specifications:

Product Code	19647D
Width	460 mm

Also Available:

Product Code	Width
19669D	610 mm
19672D	660 mm
19692D	710 mm



Stabel Hand Padded Commode Seat - Open Front

The Hand Padded Closed Front seat is completely sealed for infection control and has no sharp joins or seams. Easily removable for cleaning. Fits on to most commode types. Particularly ideal for petite users who have challenges with leg positioning.

- Completely sealed for infection control
- No sharp joins or seams
- Padded for increased pressure care
- Front opening
- Available in a range of different colours (Contact the team for more information)



Stabel Hand Padded Commode Seat - Closed Front

The Hand Padded Closed Front seat is completely sealed for infection control and has no sharp joins or seams. Easily removable for cleaning. Fits on to most commode types. Particularly ideal for petite users who have challenges with leg positioning.

- Completely sealed for infection control
- No sharp joins or seams
- Padded for increased pressure care
- Closed opening
- Available in a range of different colours (Contact the team for more information)



Specifications:

Product Code	196490
Width	410 mm

Also Available:

Product Code	Width
19648	460 mm
196630	510 mm
196640	530 mm
196550	560 mm
19667	600 mm
19668	610 mm
19728	710 mm



Specifications:

Product Code	19649C
Width	410 mm

Also Available:

Product Code	Width
19647	460 mm
19663C	510 mm
19664C	530 mm
19655C	560 mm
19667	600 mm

02.




Patient Handling

Staebel Hoists	38
Staebel Hoist Accessories	43
Ceiling Hoists	46
Staebel Slings	52
Transfers	56
Positioning	60
Support	66
Weighing	71






What is a Hoist?

Hoists are designed to assist with the manual handling of residents & patients. They're here to help reduce the risk of injury to carers, nurses and doctors. They are also portable and easy to move which allows for the safe movement of patients from one place to another. They are commonly used for shorter transfers such as bed to chair.

Pick your Hoist!

 Choose Model		Staebel Forte 320		Staebel Salsa 200
		Staebel Alto 200		Staebel Neos 150
		Staebel Sonata 150		
 Choose attachment style		Four Point Yoke	Standard Yoke	Tri Hook Yoke
		Powered Pivot Frame	Standard Yoke with Weigh Scale	Bariatric 4 Point Yoke
		Pivot Frame with Weigh Scale	Powered Pivot Frame with Weigh Scales	Bariatric 4 Point Yoke with Weigh Scales
 Accessories		Plastic Leg Protector		Sling Bag
		Staebel Salsa Handle Kit		Quick Release Adaptor
		Low Leg Height		Weighing Device
		Battery Replacement		Handset Replacement

Pick your Sling!

 Choose your material		Mesh
		Fabric
		Parasilk
 Choose Sling style	Depending on the model of hoist, choose the appropriate sling style. Contact the team for more information on choosing the right sling!	



Could not find the information you were after?

Contact our team for a consultation for more hoist options!



Staebel Sonata 150 Lifter

A lightweight, compact mobile patient lifting hoist with innovative and ergonomic design the Allegro Sonata 150 is an easy to use everyday lifter. With a safe working load of 150kg the Sonata offers an outstanding lift capacity in a small patient lifter. It utilises the SPACELINK™ concept, a breakthrough for optimised patient space, that enhances floor lift capability and increases maximum lifting height. With maintenance free bushes in all the major moving parts, this hoist is an Australian Made example of quality and strength that is hard to find in other patient lifters.

The Allegro range of patient lifting equipment, including lifters and accessories, has been designed and tested to comply with: AS ISO 10535-2002 Hoists for the transfer of disabled persons – Requirements and test methods.

- Manual leg adjustment is a standard feature
- Supplied complete with battery, in-built battery charger and instruction manual
- The yoke is the standard attachment for general purpose slings. A unique sling retainer has been designed into the hook of the yoke to prevent the sling loop from coming detached during lifting.
- Other lifter and attachment models are available



Specifications:

Product Code	2880	Castor Size	100 mm
Base Width Min	605 mm	Under Bed Clearance	126 mm
Base Width Max	1010 mm	Lifter Weight	39 kg
Lift Range - Yoke	414 - 1730 mm	SWL	150 kg

- DVA
- NDIS
- MASS



Staebel Alto MKII 200 Lifter

Alto Lift 200: A medium sized heavy duty mobile patient lifting hoist with a safe working load of 200kg and a compact ergonomic design. The Alto Lift 200, whilst being compact, provides superior patient lift space through the use of SPACE LINK™ technology.

The Alto Lift 200 is the lifter of choice for institutional use or home care where a high capacity lift is required. Designed and tested to comply with: AS ISO 10535-2002 Hoists for transfer of disabled persons - Requirements and test methods.

- Electric leg adjustment is a standard feature
- Battery life indicator on lifter as well as handset
- Additional up and down buttons on control box
- Emergency stop button
- Supplied complete with battery, battery charger and instruction manual



Specifications:			
Product Code	2891	Maximum Height Yoke	1920 mm
Base Length	1150 mm	Average Lift Range	1620 mm
External Base Width Closed	700 mm	Minimum Storage Height	1350 mm
External Base Width Open	1130 mm	Under Bed Height	55 mm
Internal Base Width Closed	600 mm	Hoist Reach	610 - 1060 mm
Internal Base Width Open	1030 mm	Hoist Reach at Maximum Height	515 mm
Base Width Adjustment	Electric	Mass of Heaviest Part (Base Assembly)	24 kg
Underbed Clearance	130 mm	Lifter Weight	46 kg
Minimum Height Yoke	300 mm	SWL	200 kg

DVA

NDIS

MASS



Staebel Alto MKII 200 Lifter with Powered Pivot Frame

A powered pivoting frame makes this hoist the perfect choice for lifting larger patients. With the push of a button repositioning and tilting the patient is safe and strain free.

A medium sized heavy duty mobile patient lifting hoist with a safe working load of 200 kg and a compact ergonomic design. Alto™ Lifter, whilst being compact provides superior patient lift space through the use of SPACELINK™ technology. The Alto™ Lifter is the lifter of choice for institutional use or home care where a high capacity lift is required. Electric leg adjustment is a standard feature.

It utilises the new SPACELINK™ concept, a breakthrough for optimised patient space, that enhances floor lift capability and increases maximum lifting height.

This hoist is an Australian Made example of Quality and Strength that is hard to find in other hoists.

Featuring the Linak Jumbo Care advanced diagnostic control system, with the latest developments in data collection.

This includes service indicator light, and with service data tool a service technician can view the history of the lifters usage including number of lifts, overloads and also program service intervals.

The Allegro range of patient lifting equipment including lifters and accessories has been designed and tested to comply with: AS ISO 10535-2002 Hoists for the transfer of disabled persons – Requirements and test methods.

- Electric leg adjustment is a standard feature
- Battery life indicator on lifter as well as handset
- Additional up and down buttons on control box
- Emergency stop button
- Supplied complete with battery, battery charger and instruction manual



Specifications:

Product Code	2892
Base Width Min	645 mm
Base Width Max	1130 mm
Lift Range - Pivot	325 - 1600mm
Caster Size	100 mm
Under Bed Clearance	126 mm

Lifter Weight	42 kg
SWL	200 kg

Also Available:

Staebel Alto MKII 200 Powered Pivot Frame with Weigh Scale

Product Code	2893
--------------	------

- DVA
- NDIS
- MASS



Staebel Forte 320 Hoist

A hoist designed the diverse resident portfolio and to enhance future proofing for your facility with mitigating manual handling risk and reducing the risk of skin integrity complications.

Featuring a class-leading lift height with extra low base for better access and maximised patient lift space. The large open cradle design allows the patient to be comfortably transferred without feeling pinched or squeezed. The concealed electric leg spread actuator promotes better hygiene practices and more space to walk behind hoist.

A long lasting lithium battery pack is a standard addition to the Forte Lifter, enabling more lifts from a single charge.

- Smooth mobility and enhanced manoeuvrability.
- Utilises the SPACELINK™ concept, a breakthrough for optimised patient space, that enhances floor lift capability and increases maximum lifting height.
- Emergency stop button & easy to use manual lowering on actuator
- Easy to service and spare parts readily available
- Additional hand control buttons on control box. Can be used when the handheld handset can't be used



Specifications:			
Product Code	2898	Overall Height (Boom Raised)	2243 mm
Base Length	1420 mm	Lifter Weight	63 kg
External Base Width Closed	715 mm	Minimum Storage Height	1445 mm
External Base Width Open	1200 mm	SWL	320 kg
Under Bed Clearance	111 mm		



Staebel Sling Bag

- Sling carry pouch, is effective for infection control where individual patient slings can be kept separate or for storage.
- Made from heavy duty green vinyl
- One size fits all sling sizes.



Plastic Leg Protector

For Alto MKII 200 Passive Hoist.

- Sold as each
- Secure press fit plastic protectors to minimise leg damage
- Easy clean option



Staebel Weighing Device

Weighing device designed to integrate with the Allegro Sonata and Allegro Alto.

Safe Weight Limit is 300kg.

Specifications:

Product Code 28804

Specifications:

Product Code 28803

Specifications:

Product Code 2889



Staebel Strap Kit for Lifting Frame & Yoke

The floor lift strap kit is an excellent accessory for lifting a person on a lifting frame from off the floor. The floor lift strap kit is an excellent accessory for lifting a person on a lifting frame from off the floor.



Staebel Lifting Frame with Gliders

The Lifting Frame is an essential piece of equipment for the safe lifting and handling of patients requiring spinal or orthopaedic immobility.

The specialised frame made from a durable metal comes complete with 12 removable nylon gliders and can be dismantled for ease of storage. Strap kit sold separately. *Pictured being utilised with 4-point yoke. SWL 200kg.



Staebel Salsa Handle Kit

- Salsa Handle Kit to aid in increased patient security whilst lift is in motion
- Suitable for the Salsa Standing Lifter
- Optional Extra

Specifications:

Product Code 27397

Specifications:

Product Code 227201

Specifications:

Product Code 27THK01



Staebel Four Point Yoke

The Allegro four point yoke can also be used with the general purpose slings and is particularly useful in the care and comfort of larger patients. A unique sling retainer has been designed into the hook of the yoke to prevent the sling loop from coming detached during lifting.



Staebel Standard Yoke

The yoke is the standard attachment for general purpose slings. A unique sling retainer has been designed into the hook of the yoke to prevent the sling loop from coming detached during lifting.



Staebel Tri Hook Yoke

These three heavy duty hoists are available with a choice of 6 different sling attachment options. Achieve quality patient transfers by selecting the fitting which best suits your client or healthcare facility.

- Can be used with general purpose (loop attachment) slings

Specifications:

Product Code	2886
--------------	------

Specifications:

Product Code	2882
--------------	------

Specifications:

Product Code	2883
--------------	------



Staebel Pivot Frame with Weigh Scale

Complete as one attachment, this allows for standard lifting without losing any critical lift height with separate attachment. Weighing scale Safe Weight Limit 300kg. Suitable for allegro Sonata & Alto Hoists.



Staebel Standard Yoke with Weigh Scale

Complete as one attachment, this allows for standard lifting without losing any critical lift height with separate attachment. Weighing scale Safe Weight Limit 300kg. Suitable for allegro Sonata & Alto Hoists.



Staebel Bariatric Four Point Yoke

Complete as one attachment, this allows for standard lifting without losing any critical lift height with separate attachment. Allegro Patient Lifter Option – Bariatric 4-point Yoke spreader bar to suit Forte 320

Specifications:

Product Code	2895
--------------	------

Specifications:

Product Code	2894
--------------	------

Specifications:

Product Code	2873
--------------	------



Stael Powered Pivot Frame with Weigh Scales

A powered pivoting frame makes this hoist the perfect choice for lifting larger patients. With the push of a button repositioning and tilting the patient is safe and strain free. A medium sized heavy duty mobile patient lifting hoist with a safe working load of 200 kg.



Stael Bariatric 4 Point Yoke with Weigh Scales

Allegro Patient Lifter Option - Bariatric 4-point Yoke spreader bar with integrated weighing device and cover to suit Forte 320.



Stael Powered Pivot Frame

A powered pivoting frame makes the hoist the perfect choice for lifting larger patients. With the push of a button repositioning and tilting the patient is safe and strain free. A medium sized heavy duty mobile patient lifting hoist with a safe working load of 200kg.

Specifications:

Product Code 2897

Specifications:

Product Code 2899

Specifications:

Product Code 2896



Stael Pivot Frame

The pivot frame is used with the pivot sling. The frame spreads the sling attachment points to ensure the patients hips are not flexed excessively, avoiding discomfort some patients experience with standard slings. The ability to position a person in a sitting or lying position, assisting in providing correct seating and positioning.



Stael Pivot Frame with Weigh Scale - Self Powered

The Self Power Pivot Frame is designed for self-powered patient positioning with no power required from the hoist. The integrated weigh scale is useful where weighing of patients is difficult due to level of support required



Stael Pivot Frame - Self Powered

The Self Power Pivot Frame is designed for self-powered patient positioning with no power required from the hoist. The integrated weigh scale is useful where weighing of patients is difficult due to level of support required. The integrated weigh scale is useful where weighing of patients is difficult due to level of support required. SWL 320kg

Specifications:

Product Code 2884

Specifications:

Product Code 2884 - 1

Specifications:

Product Code 2884 - 2

What is a ceiling hoist?

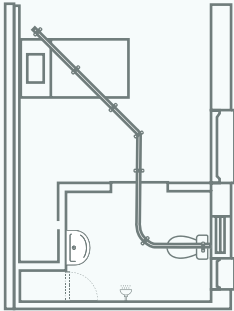
Ceiling Hoists are designed to assist with the manual handling of residents & patients. They're here to help reduce the risk of injury to carers, nurses and doctors. A ceiling hoist is an overhead hoisting system which can either be hung from the ceiling itself or hung from wall to wall. They are an effective alternative to a floorline hoist and are widely used in home care environments.

Build your Ceiling Hoist System

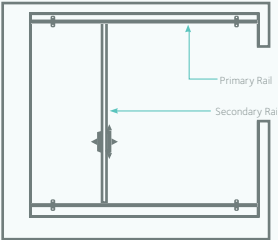
 Choose Transfer Type	Bed to Chair	Bed to Bathroom
	Bed to complete house	Full room access to bathroom
	Bariatric System	Access to other rooms

Choose Tracking System

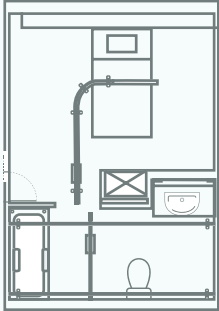
Basic Rail Systems






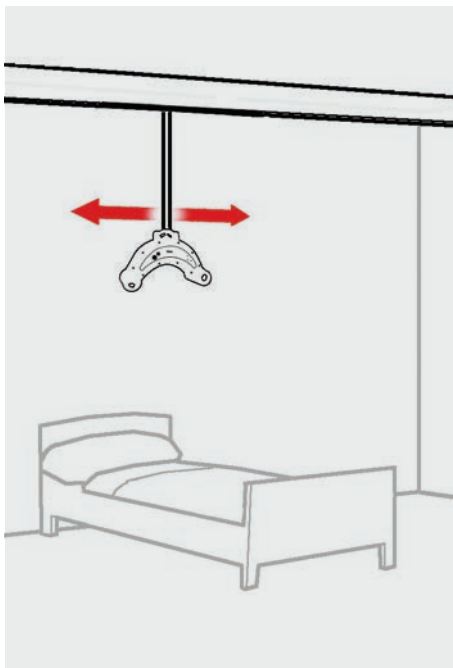
Full Room Coverage System - XY Traverse



Fully Integrated Rail System



 Weight Limit (SWL)	200 kg	270 kg
	350 kg	450 kg +
 Choose Motor Style	Portable Hoist	Fixed Powered Transverse
	Fixed Manual Traverse	Powered XY Motor
 Rail Options	Ceiling Mount	Wall Mount
		Freestanding



SugarGlider Nomad Portable Ceiling Hoist

The Nomad is a flexible and user friendly ceiling hoist. It is intuitive and easy to operate. The integrated handles make it truly portable. The curved shape is designed to afford user's head ample space while the integrated 4-point yoke gives ultimate safety, comfort and dignity during the lift and transfer. The user will feel safe and stable without being squeezed.

The combination with a compatible sling ensures that the user will achieve a correct and comfortable sitting position. If rail installations through door openings are not possible, the Nomad can climb from room to room, i.e. from one rail system to another. The Nomad has an emergency lowering and an emergency stop, as well as a soft stop-start function to prevent abrupt movements.

- Lift speed 35 mm/s with soft start/stop function
- High battery performance (70+ lifts per charge)
- Easy overnight charge system
- Built in electronic lift counter
- Splash proof casing
- Lifting Height 2m



Specifications:

Product Code	7324
Length	520 mm
Width	198 mm
Height	330 mm
Lift Range (Length of Lift Strap)	1980 mm
Lifter Weight	6.8 kg
SWL	230 kg

Also Available:

Product Code	7325
Length	520 mm
Width	198 mm
Height	330 mm
Lift Range (Length of Lift Strap)	1980 mm
Lifter Weight	6.8 kg
SWL	255 kg



SugarGlider Portable Motor

The Monarch Portable Ceiling Hoist is one of the lightest on the market weighing only 5kg, whilst lifting the industry standard of 200kg. Ideal for occasional or multiple room use, the product is robust and charges with an AC Charger.

- Weight: 5KG
- Lifts per Charge: 30 to 40 cycles
- Track Options: H system, Turntable, Gate, Curves, Room-to-Room, Freestanding Gantry: 2m, 3m, 4.5m, 6m, 8m
- Track Compatibility: Compatible with all Monarch tracks



Specifications:	
Product Code	272358000
Length	550 mm
Width	140 mm
Lifts per Charge	30 - 40 Cycles
Total Weight	4.8 kg
SWL	200 kg



SugarGlider In track Motor

The Monarch Fixed Ceiling Hoist is a complementary product to other products for individuals with mobility challenges. It is designed for use by caregivers and health professionals to lift and transfer a person in homecare, assisted-living or hospital environment.

- Unique Spreader Bar with Secure Sling Mechanism
- Fast - lifts up to 30% more quickly than other models
- Low lift clearance - slim profile means more room to manoeuvre patient
- Li-ion battery - recharges twice as fast as average competitor products
- Ultra lightweight - Weighing 8kg 35% lighter than competitors
- 200kg safe working load
- Ideal for Nursing Homes, Hospitals or Domestic Use



Specifications:	
Product Code	272358001
Length	640 mm
Width	300 mm
Lifts per Charge	30 - 40 Cycles
Total Weight	7.58 kg
SWL	200 kg

Also Available:	
Product Code	2074803
SWL	270 kg



SugarGlider Bariatric Ceiling Hoist

Designed from the ground up to better serve caregivers and patients, Monarch Bariatric ceiling lift uses a new-in-its-class powered unlimited rotation feature for superior manoeuvrability even in tight spaces. The uniquely designed powered 4-point tilting spreader bar provides enhanced patient positioning comfort during transfers.

The Monarch Bariatric ceiling lift is fully compatible with standard Savaria low profile track and the sleek motor design means the tight turning radius will be the same as other Monarch fixed lifts. The high efficiency Lithium-Ion battery delivers faster lifting speeds and built-in constant charging system compatibility means the lift is always ready to use.

When higher capacity lifting is needed, the Monarch Bariatric ceiling lift packages together a strong, fast lift that's easy to maneuver. Maintaining a patient's dignity and safety has never been simpler with this thoughtfully designed high-capacity ceiling lift.

- Powered track traversing
- Powered 360-degree unlimited rotation in both directions
- 4-point powered tilting spreader bar
- Compatible with constant charging system
- Fastest lift speed with true soft start and stop
- Superior battery recharge time
- Compatible with Monarch 4-Post Gantry System and track profiles and standard radius curves



Specifications:

Product Code	2074815
Lifting Clearance	560 mm
Lifting Speed	45mm / Second
Product Life Expectancy	10 Years
SWL	363 - 600 kg



Powered Bariatric 4 Point Tilting Spreader Bar

The uniquely designed motorised 4-point tilting spreader bar provides enhanced patient positioning comfort during transfers. For use with Monarch bariatric fixed hoists.

- SureClip System: Safety catches on spreader bar open outward to help stop the hoist loop from slipping away from the hoist during attachment
- Provides enhanced patient positioning: for comfort during transfers



Specifications:

Product Code	2074845
Spreader Bar Width	580 mm
Spreader Bar Length	607 mm
SWL	363 - 600 kg

NDIS

NDIS



SugarGlider Transport Trolley

The SugarGlider Transport Trolley makes it a simple task to transport the portable ceiling hoist unit, reach stick and slings from one location to another.

- Two large shelves
- Features a unique column to store reach stick
- Four high quality swivel castors for ease of manoeuvrability
- Robust, heavy duty design favoured by many aged care facilities


Flexi-Link Arm Reacher

Efficient, safe and simple to operate. Our patented Flexi-Link connection system is today's most advanced device for connecting the portable lift to the trolley. Its removable wand allows carers to move around freely as there are no handles in place and travel of the hoist is not impeded by lights, fans or doors.

Specifications:	
Product Code	7426
Width	520 mm
Depth	198 mm
Height	330 mm


Specifications:	
Product Code	7422FL
Length	600 mm

Also Available:	
Product Code	7423FL
Length	900 mm

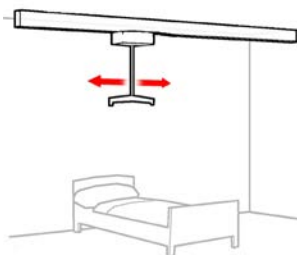

Single Track System

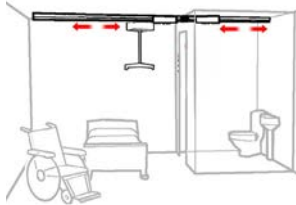
Single Track and Motor. For transfer from bed to mobile equipment e.g. wheelchair, commode.



Room To Room With Transition Gate

Room to room incorporating the Transition Gate.

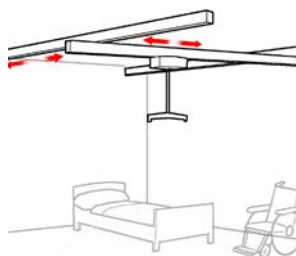Room to room transfers made easy with smooth entry to additional room via the Transition Gate. Ideal for easy access to amenities in ensuite.



X - Y System

X - Y System - Access all points in room.

The X-Y System allows access to all points within a room.



Specifications:

Product Code

2150400022



Specifications:

Product Code

2150400019

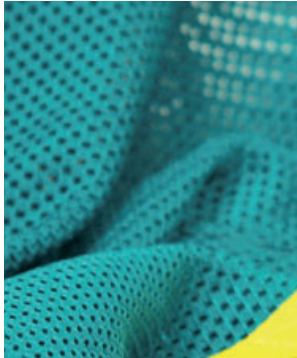
Specifications:

Product Code

2150400021

Unique Sling Features:


- Choice of mesh or fabric available (see options below).
- Complies with Australian Standards AS/ISO 10535:2006.
- Custom Slings & Standard Modifications: standard sling can be modified to include additional loops, straps or padding or a complete new sling can be designed. The trained sewing department review all custom slings to ensure they will achieve the best result possible. All that is required is a pattern or picture of the modification advising as many measurements and specifications as possible. The best results will be achieved by basing your modified sling on our existing templates.
- Label: Allegro Concepts slings have a quality woven identification label. This label will not fade over time or from the normal washing process. Important details such as sling description, serial number and date of manufacture are indelibly marked on the label. This is important for both lifetime traceability and auditing processes.
- E-Clip (optional): The patented E-clip is offered as an option on all Pivot slings and Standing slings. This revolutionary clip provides a mechanical means of ensuring the correct attachment and removes the possibility of inadvertent detachment of the sling. Depression of the side safety levers is the only way the E-clip can be removed. Visual control offers additional safety and wear over time is eliminated as the E-clip does not require a precisely engineered 'friction fit'. The E-clip is strongly recommended when prescribing our quality slings for use on lifters other than those manufactured by Allegro Concepts.



Mesh Material



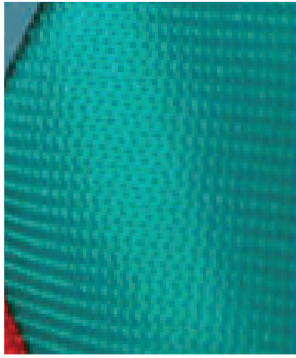Key Hole Buckle



Fabric Material



E-Clip (Optional)

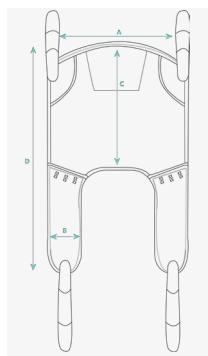


Woven Label





Choose
Sling style

Depending on the model of hoist, choose the appropriate sling style.
Contact the team for more information on choosing the right sling!



NDIS

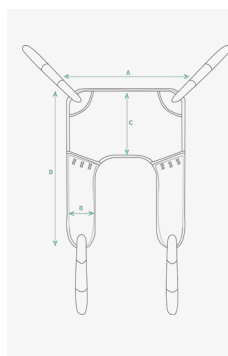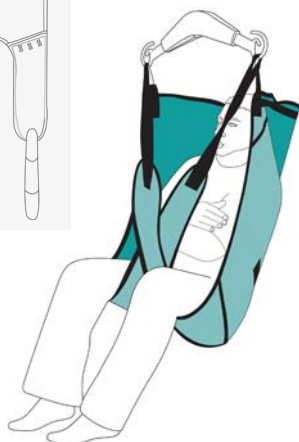MASS



NDIS

MASS



Staebel General Purpose with Head Support Sling

The General Purpose Head Support sling, as the name suggests, is a full body support style sling and is recommended for general lifting and in situations of lifting from the floor, and for the less secure patient. The loop style attachments have three lengths at each connection point and are used on a patient lifter with a yoke style attachment. Available in sizes from x-small through to XX-large, providing a popular choice for individuals and care facilities. Slings are manufactured from polyester fabric or mesh (suitable for bathing).

Also available:

Activecare General Purpose Hygiene Slings With Head Support

Staebel General Purpose Sling

The General Purpose sling is a standard sling providing good support for those with some upper body control. Suitable for general lifting for patient lifters with a yoke style attachment, and available in sizes from X-small to XX-large. Also with the option of fabric or mesh construction, this is a versatile sling that is simple and safe to use.



Mesh

XS	27GPHS02XS
S	27GPHS02S
M	27GPHS02M
L	27GPHS02L
XL	27GPHS02XL
XXL	27GPHS02XXL

Fabric

XS	27GPHS01XS
S	27GPHS01S
M	27GPHS01M
L	27GPHS01L
XL	27GPHS01XL
XXL	27GPHS01XXL

XS	S	M	L	XL	XXL
A 820	790	920	1110	1270	1360
B 210	230	240	260	270	290
C 750	780	920	1000	1020	1050
D 1310	1280	1410	1630	1680	1780
E 585	650	700	740	780	810




Mesh

XS	27GP02XS
S	27GP02S
M	27GP02M
L	27GP02L
XL	27GP02XL
XXL	27GP02XXL

Fabric


XS	27GP01XS
S	27GP01S
M	27GP01M
L	27GP01L
XL	27GP01XL
XXL	27GP01XXL

XS	S	M	L	XL	XXL
A 790	790	890	1110	1270	1400
B 210	230	240	260	270	290
C 470	460	530	650	730	750
D 1120	1130	1190	1390	1520	1590

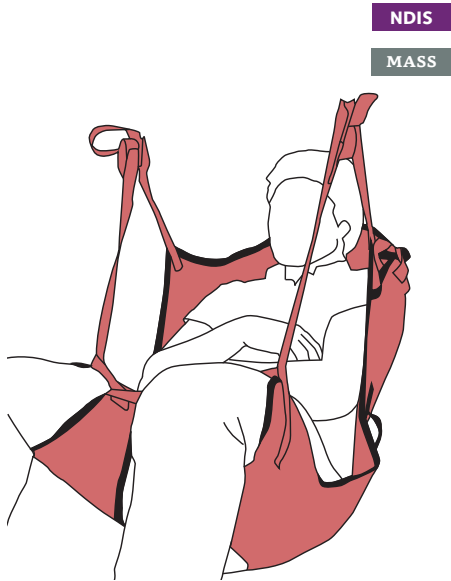

NDIS

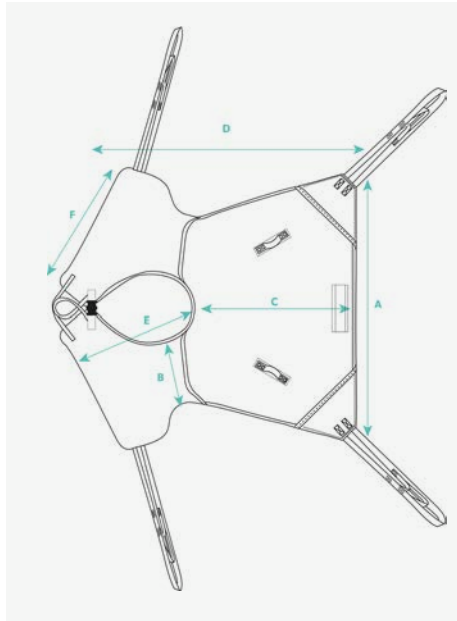MASS




Staebel Universal Deluxe Sling with Head Support

This sling is good for all standard transfers and is exceptionally comfortable for the patient. The leg sections are shaped to provide good support to the hip and thigh areas. It is designed to hold the leg more securely than standard sling leg sections, as the leg is held

from both sides instead of just the inside of the leg. The leg sections also spread the load over a larger area on the patient's legs so tend to reduce the pressure on the underside of the legs.




Mesh						
	XS	227430				
	S	227431				
	M	227432				
	L	227433				
	XL	227434				
	XXL	227450				
			Parasilk			
	XS	227440				
	S	227441				
	M	227442				
	L	227443				
	XL	227444				
	XXL	227212				
	XS	S	M	L	XL	XXL
A	715	800	900	1000	1105	1205
B	160	185	210	235	260	285
C	685	815	945	1075	1205	1270
D	1000	1175	1355	1535	1710	1820
E	295	340	385	430	475	520
F	375	440	500	560	625	690

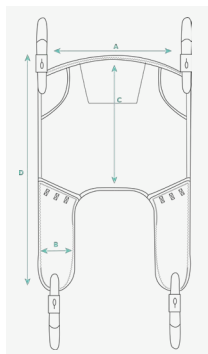

NDIS

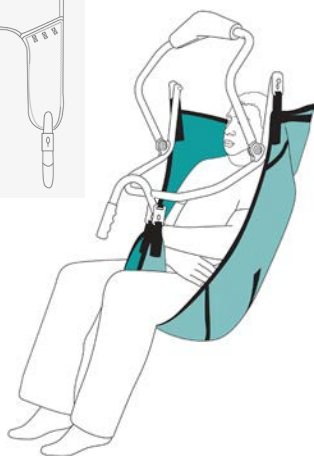MASS



Staebel Pivot Frame Sling

The Pivot sling is designed for use with the pivot/tilt frame. It attaches with the keyhole tags engineered to the highest standards in strength and size tolerances and feature auto clip on action (no safety straps required). It also features sling centre marking and positioning handles to assist correct positioning and woven label that will not wash out or fade.

Available in sizes from x-small through to XX-large and in fabric or mesh styles.



Mesh			Fabric			
XS	27P02XS		XS	27P01XS		
S	27P02S		S	27P01S		
M	27P02M		M	27P01M		
L	27P02L		L	27P01L		
XL	27P02XL		XL	27P01XL		
XXL	27P02XXL		XXL	27P01XXL		
XS	S	M	L	XL	XXL	
A	720	870	1000	1080	1170	1320
B	190	230	270	260	270	290
C	580	830	920	1010	1090	1180
D	870	1220	1310	1410	1540	1710

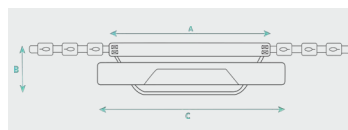
NDIS



Staebel Salsa Sling

The Salsa sling is for use with the Salsa stand up hoist. It features optimum back support and under arm padding for a safe and comfortable lift and also allows for the easy removal of clothing for toilet transfers.

The Allegro Salsa sling features a torso band for security and stability with Buckle (b) or Velcro (v) fastening.



Specifications:

Buckle Fastening

Product Code	227701	27702	27703	27704
Size	Small	Medium	Large	Extra Large
A	730 mm	880 mm	1060 mm	1220 mm
B	260 mm	290 mm	370 mm	410 mm
C (b)	900 mm	960 mm	1020 mm	1220 mm
C (v)	1140 mm	1180 mm	1320 mm	1440 mm

Velcro Fastening

Product Code	27720	27705	27706	27707
Size	Small	Medium	Large	Extra Large
A	730 mm	880 mm	1060 mm	1220 mm
B	260 mm	290 mm	370 mm	410 mm
C (b)	900 mm	960 mm	1020 mm	1220 mm
C (v)	1140 mm	1180 mm	1320 mm	1440 mm

DVA
NDIS
MASS



Staebel Salsa 200 Stand Up Lifter

The Salsa 200 was designed as an active standing aid which allows a single carer to raise a patient into a comfortable standing position. The Salsa is ideal for transferring between bed and chair, or chair and commode. Allied with the Salsa, our generously padded slings allow easy transfers whilst still providing superior access to client to remove clothing for toiletry purposes. This patient hoist has keyhole sling attachment points. The Salsa200 has a removable narrow foot plate (width 350 mm) which guarantees ease of use when utilised in conjunction with shower commodes etc. With a small footprint the Salsa is highly manoeuvrable and ideal for use in an institutional/ aged care environment.

Featuring the Linak Jumbo Care advanced diagnostic control system, with the latest developments in data collection. This includes service indicator light, and with service data tool a service technician can view the history of the lifters usage including number of lifts, overloads and also program service intervals.

The Allegro Salsa lifter has been designed and tested in house to comply with AS ISO 10535 – 2002 Hoists for the transfer of disabled persons – Requirements and test methods.

- The Salsa200 has the following features and benefits :
- Australian made quality product
 - Electric leg adjustment
 - Patient friendly lift arms enable multiple grip positions during lift
 - Removable footplate for ease of transfers
 - Height and angle adjustable contoured knee pad with removable calf strap
 - Comes complete with battery pack and charger
 - T.G.A. 136 065



Specifications:			
Product Code	27700	Castor Size	100 mm
Base Width Min	655 mm	Under Bed Clearance	130 mm
Base Width Max	1030 mm	Lifter Weight	51 kg
Storage Height	1200 mm	SWL	200 kg

NDIS

MASS



NDIS

MASS



Staebel Patient Transfer

The Staebel Patient Transfer Device is the perfect solution for safe transfer of a patient or resident within the home or facility. This is your all-in-one solution without the need for additional lifting straps and belts, and promoting independence. The width adjustable base allows good access to most seating, whilst also designed for good access beneath beds etc.



Specifications:

Product Code	212145
Seat Width	450 mm
Overall Width	640 mm
Overall Length	930 mm
SWL	150 kg

Also Available:

Product Code	212147
SWL	220 kg

Staebel Transfer Platform

Staebel Transfer Platform is a manually-operated device designed as a platform for use by healthcare professionals to assist in the physical moving (transfer) of a patient in an upright (standing) position over a short distance.

The device consists of a wheeled platform, a surrounding framework, kneepad and handlebar. It can include safety belt, to hold the person being transferred in the standing position.



Specifications:

Product Code	212150
Width	620 mm
Height	1180 mm
SWL	150 kg

NDIS

MASS



Staebel Transfer Belt

Staebel Transfer Platform Belt is used with the Staebel Transfer Platform when more support and stability is needed, for instance when the user's weight-bearing ability is questionable or varied. Staebel Transfer Platform Belt can also be used stand alone for sit-to-stand transfers.

- Anti-slip inner surface
- 170kg SWL
- Sizes available, XSmall, Small, Medium, Large, XLarge and XXLLarge
- Each size is colour coded for easy identification



Specifications:

Product Code	212152M				
Size	Medium				
Width	210 mm				
SWL	170 kg				

Also Available:

Product Code	212152XS	212152S	212152L	212152XL	212152XXL
Size	Extra Small	Small	Large	Extra Large	Extra Extra Large
SWL	170 kg	170 kg	170 kg	170 kg	170 kg

NDIS



Staebel Transfer Board Slide

The Transfer Board Slide is a simple low-cost device that allows for a sitting transfer across the board without risk of shearing/induced tissue damage.

This is a wonderful little product to use with any sitting transfer board and makes sure that the slide across is what it should be – effortless. Just pop it on top of your board, tuck it under your bottom a little and sliding from one surface to another will be plain sailing.


Simple and highly effective, this Transfer Board Slide is only as wide as the board itself and has a panel of antislip that grips the top of the board.

- Transferring from one sitting surface to another
- Helping people move, both independently or with assistance



Specifications:

Product Code	210360
Width	300 mm
Length	700 mm

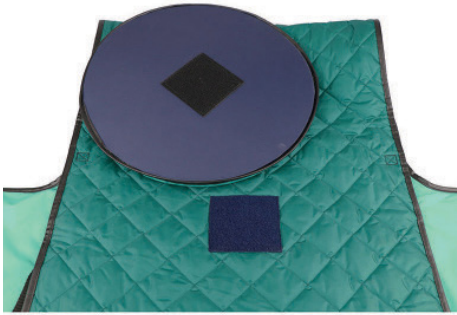

Staebel Bed Transfer Assist

The Bed Transfer Assistant has been designed to assist in supporting the lower limbs during functional transfers.

Simple and brilliant, our Bed Transfer Assistant allows one carer to safely raise a person's legs and swivel them around into bed and vice versa!

- Cotton upper for comfort.
- Raising legs and swivelling in a sitting position by a single carer



Specifications:

Product Code	210315
Width	600 mm
Length	640 mm
Thickness	30 mm
SWL	100 kg



Staebel Car Assist

Staebel Car Assist is simple and easy to use device to allow access into and out of the car seat. There are handles on each side of the Staebel Car Assist to allow the carer to provide assistance if required. Comes standard with a seat belt to secure to the Staebel Car Assist to the car seat, to prevent sliding while in motion.

- Easy to install.
- Robust handles enable carers to provide assistance.
- Secure fitting in the car seat prevents users from sliding forward.
- Soft, padded surface provides comfort



Specifications:

Product Code	210490
Width	450 mm
Length	450 mm



Staebel Slidesheet

The ideal solution when transferring or bed repositioning patients across, up or down a bed. Reduces risk injury by assisting sliding, rather than lifting. Made from strong, non flammable, nylon fabric our slide sheets will not use their slippery quality or colour even after numerous washes. Can be machine washed to 50 degrees Celcius, tumble dried to 60 degrees Celcius and even ironed on a cool setting.

Available Colours:

- Blue 2434B
- Green 2434G
- Orange 2434O
- Red 2434R



Staebel Disposable Rolling Slide

The Disposable Procedure Slide is a simple slide sheet that allows health professionals to easily move patients into the lithotomy position, ready for surgery. It is ideal for head-down procedures where the patient is positioned at an incline of up to 25 degrees. Those who work in Gynaecology will find the Disposable Procedure Slide hugely beneficial as it allows medical practitioners to:

- Slide a patient into a precise position (fine adjustments achievable), usually the Trendelenburg position
- Lock patient position with the sheet's unique and integral locking panels
- Position patients at inclines up to 25 degrees, facing head-down
- Adjust patient position mid-surgery
- Slide patients back up the operating table post-surgery
- Disposable slides help teams uphold excellent infection control
- Slippery material enables easy patient repositioning
- Integral locking pads ensure absolute safety and security
- Can remain under the patient
- 25 per pack



Specifications:

Product Code	2434
Width	1450 mm
Length	2000 mm



Specifications:

Product Code	210134
Width	1000 mm
Length	2000 mm



Staebel ComfortCare Bed Top Sheet

The ComfortCare is an in-bed sliding system for the care of the long term or 24 hour patient. The amazing ComfortCare offers patients 100% comfort and carers the ability to manoeuvre and reposition them in all directions with ease and then secure that position when needed.

- Moves a patient in all directions
- The locking panel will hold the patient steady
- Offers patients comfort and dignity
- Patients can remain on this indefinitely, either directly or under a cotton bed sheet
- Carers can re-position and manoeuvre patients without the need to disturb or handle directly
- Reduces risk of injury to carers
- Enhances patient care and clinical outcomes
- Suitable for use with hoists
- Easy to use
- Polycotton system, with low-friction inner linings
- Two panel system enables up/down manoeuvres, turning, swivelling, sitting up safely, transferring and turning into prone position
- Integral locking panel means safe position can be secured, even in sitting
- Suitable for single carers
- Washable at 40°C to 80°C and low temperature tumble dry



Specifications:

Product Code	210295
Width	1400 mm
Length	2000 mm
SWL	200 kg

Staebel ComfortCare Bed Base Sheet

The ComfortCare is an in-bed sliding system for the care of the long term or 24 hour patient. The amazing ComfortCare offers patients 100% comfort and carers the ability to manoeuvre and reposition them in all directions with ease and then secure that position when needed.

- Moves a patient in all directions
- The locking panel will hold the patient steady
- Offers patients comfort and dignity
- Patients can remain on this indefinitely, either directly or under a cotton bed sheet
- Carers can re-position and manoeuvre patients without the need to disturb or handle directly
- Reduces risk of injury to carers
- Enhances patient care and clinical outcomes
- Suitable for use with hoists
- Easy to use
- Polycotton system, with low-friction inner linings
- Two panel system enables up/down manoeuvres, turning, swivelling, sitting up safely, transferring and turning into prone position
- Integral locking panel means safe position can be secured, even in sitting
- Suitable for single carers
- Washable at 40°C to 80°C and low temperature tumble dry



Specifications:

Product Code	210296
Width	900 mm
Length	2000 mm



Staebel In Bed Turner 4-way with Handles

InBed Turner 4-Way with Handles allows safe and easy repositioning for people in bed – enabling them to turn from side to side and move up and down, with or without the aid of a hoist. They are designed for more dependent users who need assistance in bed, often from a single carer, and are designed to be used on top of one of our InBed Turner Base Sheets. They do not work on their own.

- Can remain under the patient at all times.
- Great for those who need multi positioning in bed –up / down and turning.
- Reduces risk of tissue and nerve damage.
- Helps to reduce incidence of pressure sores developing.
- The combination of both silky smooth satin panels working with each other mean easy, quiet and safe movements.
- Cotton and polyester mix for great comfort and repeated washing.
- Washable to 90C.
- Four-way movement in all directions for patient care.
- Allows more dependent people to be moved in comfort and safety.
- Cotton and satin fibres make for comfortable in-bed experience.
- Choice of torso or full length models.
- Main sizes are all 200cm wide, allowing the sheet to be tucked in after repositioning, increasing patient security.



Specifications:	
Product Code	210322
Width	2000 mm
Length	1400 mm
SWL	300 kg


Also Available:	
Product Code	210324
Width	2000 mm
Length	2000 mm
SWL	300 kg



Staebel In Bed Turner Base Sheet

The InBed Turner Base Sheet is part of our renowned satin range for in-bed repositioning and patient care. All InBed Turner Base Sheets allow safe and easy repositioning for people in bed enabling them to turn easily from side to side but not slip down. Extra safety is provided by the striped cotton edges on all of our Base Sheets. They can be used on their own or combined with one of our Satin Top Sheets for those who need more in bed repositioning care.

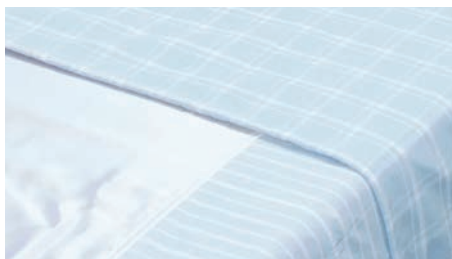- Reduces risk of tissue and nerve damage
- Helps reduce pressure sores developing
- No training required
- Easy turning and repositioning with little effort
- Great for those having difficulty turning in bed
- Washable
- Central satin panel is smooth for ease of movement
- Anti-slip cotton polyester on edges reduces risk of slipping
- Cotton polyester for enhanced comfort and repeat washing
- Single and king options are fitted sheets. The double is a flat sheet.



Specifications:

Product Code	210321
Size	Single
Width	900 mm
Length	2000 mm
SWL	300 kg

Also Available:

Product Code	210323	210325
Size	King Single	Double
Width	1050 mm	2000 mm
Length	2000 mm	2300 mm
SWL	300 kg	300 kg

NDIS

NDIS



Staebel One Way Slide & Lock

The One Way Slide & Lock are designed to allow carers to help a user slide to the back of a chair or wheelchair and not be able to slide down.

Handles are positioned equally along the length of the product and carers can use the handles directly or attach extension straps to the handles, to prevent them having to bend down.

- To help people to sit back properly in chairs or wheelchairs without sliding back down. People independently, or with assistance by one or two carers depending upon need
- Prevent people sliding down in wheelchairs, chairs and beds
- Improves people's seating posture and comfort
- Allows people to be positioned further back in a chair
- Increases people's independence
- Can be used unassisted or with carer assistance

Staebel One Way Slide & Lock Open Ended

The One Way Slide & Lock are designed to allow carers to help a user slide to the back of a chair or wheelchair and not be able to slide down.

Handles are positioned equally along the length of the product and carers can use the handles directly or attach extension straps to the handles, to prevent them having to bend down.

- To help people to sit back properly in chairs or wheelchairs without sliding back down. People independently, or with assistance by one or two carers depending upon need
- Prevent people sliding down in wheelchairs, chairs and beds
- Improves people's seating posture and comfort
- Allows people to be positioned further back in a chair
- Increases people's independence
- Can be used unassisted or with carer assistance



Specifications:	
Product Code	210190
Width	450 mm
Length	540 mm



Specifications:	
Product Code	210192
Width	450 mm
Length	1300 mm

Also Available:	
Product Code	210193
Width	550 mm
Length	1300 mm



Staebel Anti Shear Mattress Cover

A simple answer to the problem of constantly putting slide sheets in and out – lightweight and easy to use. You will love the Anti Shear Mattress Cover's simplicity as it makes it so much easier to adjust someone's position without disturbing them and really reduces the risk of injury for staff.

Thinner than a bed sheet, this product fits over the mattress and does not impinge on any pressure reducing qualities the mattress may have. You can sit the client/patient up in bed without the risk of them slipping down, and it only becomes fully slippery when you introduce a flat slide sheet on top to perform the manoeuvre that needs to be performed.

Just place the Anti Shear Mattress Cover over the bed mattress and, if you want to, put the normal bed sheet on top, tucking it in. When you need to reposition someone, untuck the bed sheet and unroll a flat slide sheet in between the Anti Shear Mattress Cover and bed sheet (without disturbing the patient) and use this to achieve a new position. Then slip out the slide sheet easily, retuck the bed sheet and your patient will be comfortable.

- Specialist technical fabric is fully breathable, lightweight and quiet
- Fabric is thinner than a bed sheet and will not impinge on pressure reducing qualities of mattresses
- Light-weight and compact
- Washable at 40 °C to 74 °C
- When using the base sheet to lie on, independent people can turn more easily in bed
- Breathable fabric means that people will not sweat or get too hot – as can happen with other systems



Specifications:

Product Code	210290
Width	900 mm
Length	2000 mm

Staebel Perfect Positioner

Do you hoist people into their wheelchairs and find it difficult to re-position them so that they are sitting properly?

Originally designed to help reposition a tetraplegic gentleman in a moulded wheelchair, this unique handling device has been redesigned and now enables carers to easily and comfortably reposition the user (once hoisted into the wheelchair) accurately at the back of the chair and ensure their pelvic tilt is correct.

- Especially useful in moulded and specialised wheelchairs
- Correct sitting position in wheelchairs
- Accurate pelvic tilt
- One carer or two



Specifications:

Product Code	210340
Width	900 mm
Length	2000 mm

NDIS

NDIS

NDIS



Staebel Slideezi with Handles

The Slideezi with Handles are manufactured using an ultra-slippy material developed by pioneering fabric specialists Slideezi®. While one side of the slide sheet is luxuriously soft to offer maximum comfort and quietness, technically advanced fibres give the other side a slippy finish, ideal for repositioning a patient in a bed or trolley.

The Slideezi with Handles allow health professionals to effortlessly manoeuvre patients/clients into a virtually endless list of positions so that they can remain comfortable at all times. Whether patients are using a bed or a trolley, health professionals can use the slide sheets to turn patients into comfortable positions and also transfer them from one surface to another.

- Superb slideability on both sides for effortless and easy repositioning
- Easy to insert and remove
- Superb comfort for during repositioning
- Quiet and soft
- Fabric is manufactured using advanced technical fibres
- Used in pairs for best performance



Specifications:

Product Code	210432
Width	900 mm
Length	2000 mm
SWL	300 kg

Also Available:

Product Code	210429
Width	1200 mm
Length	2000 mm

Staebel Safemove

The Safemove has been developed with practitioners and allows radiotherapists to position their patients far more precisely than has been possible until now.

Positioning patients for laser precision has always been difficult and achieving that specific position can be uncomfortable and awkward for the patient and extremely hard, tiring and heavy work for staff. It can result in potential skin damage to patients and injury risk to staff. However, by placing the Safemove under the patient (breast or pelvic areas), staff can use the straight and/or curved handles to much more easily reposition their patients to line up precisely with the lasers.

- No need to manually handle patients
- Reduced risk of tissue damage
- Less stress and discomfort for the patient
- Less strain and risk of injury to the staff



Specifications:

Product Code	210317
Width	370 mm
Length	620 mm

Staebel Sure Foot Pad

The Sure Foot Pad is a really helpful, inexpensive answer to the problem of feet slipping, on the floor or in bed.

Suitable for use on the floor or in bed, this Sure Foot Pad is designed as a safety aid that helps to prevent feet slipping. It can be used to stand on when sitting or rising from a chair, or it can be used with a slide sheet to aid movement in a bed or with a Slide & Lock to aid movement in a chair. A helpful aid for more independent living, it is hard wearing and easy to clean.

Simply place the foot pad on the surface you need to gain more grip on and your feet won't slip. You will love its simplicity because it gives you so much independence and is also a really useful aid for rehabilitation.

Ideal for use with small slide sheets to slide further up the bed or on the floor, on its own or with products such as Slide & Locks. In fact there are so many places that you can use this product, some people even have two; those for the floor and those for use in bed!

- Gaining more floor grip with your feet
- Helping people be more independent
- A rehabilitation aid



Specifications:

Product Code	210001
Width	1500 mm
Length	2000 mm





Staebel Knee Cushion

The Knee Cushion will provide a comfortable cushion for those who spend a long time kneeling as part of their job and/or find it painful to kneel on the floor.

Our Knee Cushions can be used on any floor surface, inside or out, and will provide a clean and comfortable surface to kneel on whilst protecting and insulating your knees.

They are completely sealed and have been welded to ensure that they are easy and safe to wipe down.

There is no need to send these Knee Cushions to the laundry or worry about cross-infection anymore.



Staebel Support Straps - Pair

These Support Straps are used to reduce the risk of over-stretching and reaching in frequent supine transfers. No need to launder reducing the chance of losses.

A new concept in support straps, these can be wiped down and/or washed in soapy water, rinsed and then left to dry. Used to reduce the risk of over-stretching and reaching in frequent supine transfers, they simply clip onto the handles of the flat Slideezi with handles and other patient handling equipment.

- For use in supine transfers
- Reducing the risk of overreaching and stretching



Staebel In Bed Sit Support

This easy-to-attach bed ladder allows patients to maintain control as they sit up and lie down. Patients who would otherwise struggle to manoeuvre themselves gain a whole new level of independence with the bed rope ladder, allowing them to reposition themselves safely and easily.

You will love how simple this multicoloured functional webbing ladder is to attach. No knots to tie or worry about. Just loop it through itself and then use the different coloured rungs to help you as you sit up or lie down. Ideal if you want to maintain your independence and just need a little more help or if you are going through a process of rehabilitation.

- To help people sit up/lie down independently and in a controlled manner



Specifications:

Product Code	210250
Width	450 mm
Thickness	40 mm

Specifications:

Product Code	210064
Width	30 mm
Length	750 mm
SWL	100 kg

Specifications:

Product Code	210440
Width	160 mm
Length	1880 mm
SWL	300 kg

NDIS



NDIS



NDIS



Staebel Easy Dress Arms

A pair of long slipper sleeves that can slide over a person's arm and ease the putting on of upper body garments.

For best practice: The sleeves work especially well with a flat slide sheet underneath the bed bound person to facilitate the slide.

- Help carers put a shirt, blouse or jumper on a bed bound person/child or someone who requires a lot of assistance

Staebel Easy Dress Legs

A pair of long slipper leg sleeves that can slide over a person's legs and ease the putting on of lower half garments.

For best practice: The sleeves work especially well with a flat slide sheet underneath the bed bound person to facilitate the slide.

- Help carers put trousers/pyjamas on a bed bound person/child or someone who requires a lot of assistance

Staebel Easy Hands

These silicon coated nylon slipper gloves are designed for carers and can be used to make patients/service users much more comfortable in many ways. For example: skin smoothing, positional adjustments and inserting slings etc.

- Technically designed fabric means the gloves slide easily on the outside but will not slide off the users hands
- Easy to slide under patients for improved patient care
- Can be used for a multitude of patient care tasks:
 - Skin smoothing
 - Tucking in sheets
 - Adjusting patient clothing
 - Moving pillows
 - Inserting slings etc
- Supplied in pairs
- Slippery outer and grippy inner material
- Fits to elbow



Specifications:

Product Code	210159
Size	Medium
Width	280 mm
Length	1070 mm

Also Available:

Product Code	210161
Size	Large
Width	270 mm
Length	1090 mm

Specifications:

Product Code	210160
Size	Medium
Width	340 mm
Length	1200 mm

Also Available:

Product Code	210162
Size	Large
Width	360 mm
Length	1300 mm



Specifications:

Product Code	210339
Width	200 mm
Length	520 mm



Staebel Easy Dress Kit

Getting dressed in bed has never been easier – the Easy Dress Kit enables dignified dressing for patients who need carer assistance.

This Easy Dress Kit is a complete set of in-bed dressing aids that enables one or two carers to dress a patient/client fully in bed. Where dressing is difficult or uncomfortable, they enable easy dressing and promote patient dignity.

- 2 x Easy Dress Legs
- 2 x Easy Dress Arms
- 1 x Flat Sheet (folded in a handy storage bag)
- Easy dressing aid for carers to clothe a person in bed
- Suitable as an aid for an individual who can clothe themselves but require a little help
- Super Slippery fabric for extra ease
- Available in Washable and disposable material
- Available in Standard (one size fits all)
- Washable 40° to 74° and low tumble dry
- Replacement parts available if lost



Specifications:

Product Code

210153 (Complete Kit)



Staebel Emergency Lifting Kit

An answer to the problem of handling a collapsed person in a number of situations – in a tight space, on the floor or in bed.

This specially designed pack contains all the equipment you will need, especially where there is no mechanical aid available. For example:

- Arrest in bed or in a chair
- Collapsed in an awkward corner or tight space.
- Collapsed on the floor and no mechanical hoist or lifting mattress available to lift off the floor.
- 2 x Slideezi + Handles
- 1 x Rolling Slide1
- x Lifting Sheet
- 2 x Support Straps



Specifications:

Product Code

210140 (Complete Kit)

NDIS

NDIS

NDIS



Staebel Lifting Sheet

The Lifting Sheet offers the ultimate falls management solution for transporting an injured person from the site of a car accident or other emergency.

No stretcher hoist? No way to lift a person safely? Look after your staff and patients by providing them with this inexpensive lifting sheet, designed to protect both patients and health professionals during lifting procedures.

The Lifting Sheet is a lightweight falls management solution, ideal for use when health professionals need to raise someone from the floor in a horizontal position. Designed in conjunction with the Guidelines for Safer Handling during Resuscitation in Healthcare Settings, the lifting sheet is perfect for post-resuscitation transfers where a stretcher hoist is not available to move the patient from the floor onto a bed or trolley.

- Lightweight solution when no stretcher hoist is available
- Robust design supports patients weighing up to 150kg
- Washable at 74°C
- Waterproof
- Fire-retardant
- Can remain under patient



Specifications:	
Product Code	210138
SWL	150 kg



Staebel Disposable Procedure Slide

The Disposable Procedure Slide is a simple slide sheet that allows health professionals to easily move patients into the lithotomy position, ready for surgery. It is ideal for head-down procedures where the patient is positioned at an incline of up to 25 degrees. Position patients at inclines up to 25 degrees, facing head-down

- Adjust patient position mid-surgery
- Slide patients back up the operating table post-surgery
- Disposable slides help teams uphold excellent infection control
- Slippery material enables easy patient repositioning
- Integral locking pads ensure absolute safety and security
- Can remain under the patient
- 25 per pack



Specifications:	
Product Code	210068
Width	410 mm
Length	410 mm



Staebel Octopannus – Pannus Support Belt

The Octopannus is a versatile pannus support belt that offers a wealth of uses in medical settings. Not only is the belt ideal for supporting bariatric patients during surgery, but it also provides the support patients often need while recovering from wound treatment. Plus, the belt helps caregivers access difficult-to-reach areas, which can be crucial for ensuring patient hygiene. During cleaning, carers can use the Octopannus to support a patient's pannus and limbs to make personal care easier.

- Assists with cleaning, personal care, catheterisation and probing
- Easy to insert and reposition
- Eliminates the need to reach over the patient or bed
- Soft and comfortable fabric
- Suitable for all bariatric patients
- Can support limbs during personal care
- Reduces musculoskeletal injuries for staff and carers
- 10 quantity per box



Specifications:	
Product Code	210271
Width	340 mm
Length	2900 mm
Box Qty	10



Staebel Support Walk Belt

A supportive Belt to help with transferring or walking. The ergonomically designed handles enable the care giver to hold the patient close without the need to grab clothes.

Ranging from Small to Large (adult sizes) this wonderful range of comfortable and supportive belts means that getting close to your patient has never been easier.

Simply place the belt around the person's waist, adjust to fit and then buckle. The ergonomically designed handles are perfectly located for user comfort and functionality. The belt is also washable.

- Helping a person from sitting to standing and vice versa
- Walking
- Supporting
- Helping whilst transferring on a sitting transfer board
- Comfortable to wear
- Give a feeling of added security
- Easily adjusted to fit



Specifications:

Product Code	210082
Width	250 mm
Length	730 - 910 mm

Also Available:

Product Code	210083	210084
Width	250 mm	250 mm
Length	880 - 1040 mm	1040 - 1170 mm



Staebel Chair Scale

The Staebel Chair Scale with its patented design maintains reliability and uncompromising precision. All of the user's care needs, including safety, comfort, caregiver's use and maintenance, are met. The compact design ensures the most cost efficient shipping method and reduces the risk of carer strain particularly when weighing an obese/bariatric patient.

- Twin wheel castors are standard and do not create excessive noise when wheeling
- One piece footplate flips up and has been designed to reduce skin tears
- Armrests flip back to allow scale to be utilised in bariatric care setting
- Base frame comes standard with a level indicator
- 5 digit LCD screen that can be tilted to accommodate different care staff heights
- Charger has updatable firmware and ports: RS232 and USB



Specifications:

Product Code	2957
Seat Width	420 mm
Seat Depth	420 mm
Overall Height	965 mm
SWL	300 kg

Sizes Available:

Staebel Chair Scale Charger	
Product Code	2957 - 1



Staebel Digital Chair Scale

The Staebel weigh chair is a stable well constructed Chair Scale which provides the operator and patient ease of use, especially small children, elderly, frail and handicapped patients. The arm rests move up on both sides of the chair for easy patient transfer to and from the chair scale and firmly fold down when the patient is seated, this reduces lifting strain on both the patient and staff.

- 300kg SWL
- 4x braking castors
- Swing up Armrests for easy patient transfer
- Hinged Footrest for easy patient transfer.
- Auto turn off Function



Specifications:

Product Code	2958
Seat Width	530 mm
Seat Depth	350 mm
Seat Height	550 mm
Overall Width	690 mm
Overall Depth	1000 mm
Overall Height	990 mm
SWL	300 kg





NDIS

NDIS

NDIS

Physicians Scale With Height

- Unique design
- Heavy duty construction
- Large platform to suit bariatric users
- Integrated height measure
- Suitable for outpatient and bariatric clinics

Digital Wheelchair Scale

With its high capacity of 360 kg, this scale makes it easy to weigh bariatric patients. Thanks to its extremely low platform and the self lowering ramps, this scale is easy to access. The integrated pre-TARE function determines in seconds the new weight of the patient utilizing internal weight memory. The swivelling LCD display with large letters is easy to read both for the patient and for the carer. The scale can be folded up in no time and can be transported anywhere thanks to its unique roller castors.

- Power pack is included
- Functions : Automatic switch off, pre-TARE with 3 memory cells, TARE, HOLD, BMI, damping, kg/lbs/ sts switch-over
- Optional : RS232 adapter set
- 24 month warranty

Portable Wheelchair Scale With BMI

- One-piece and multifunctional
- Large weighing area
- Fitted castors for easy transport
- Easy to store
- 2 ramps
- Column for display available as an accessory
- Handrails are an optional accessory
- Tare to deduct weight of wheelchair



Specifications:	
Product Code	2770
Overall Width	440 mm
Overall Height	1240 mm
Overall Depth	460 mm
SWL	300 kg

Specifications:	
Product Code	2664
Platform Width	760 mm
Platform Depth	930 mm
Platform Height	55 mm
Graduation	50 g
Device Weight	24.5 kg
SWL	360 kg



Specifications:	
Product Code	2637
Platform Width	1150 mm
Platform Depth	8800 mm
Platform Height	66 mm
SWL	300 kg

03.

Pressure Care

LOW RISK CUSHIONS

Action Polymer Pad	74
Staebel Deluxe Foam Cushion	74
Low Resilient Foam Cushion	75
Relax Trio Cushion	75
Staebel Bariatric Cushion	76

MEDIUM RISK CUSHIONS

Static Air Cushion	76
Action Gel Professional Cushion	77
Action Gel Pilot Cushion	77
Action Centurian Cushion	78
Action Commuter Cushion	78
Relax Easy Cushion	79
Relax Easy High Profile Cushion	79

HIGH RISK CUSHIONS

Relax Duo Profile Cushion	80
Relax DuoGel Cushion	80
Relax Gelcell Cushion	81
Relax Gelcell G2 Cushion	81
Starlock Cushion	82
Standard Contour Cushion	82
Staebel Action Gel - Segment Pad (3 Per Bed)	83

MEDIUM RISK MATTRESS

Staebel Zonecare Mattress	83
---------------------------	----

HIGH RISK MATTRESS

Metzeler Ultratube FirmFeel Mattress	84
Staebel Soft EuroFlex Hygiene Mattress	84
CuroCell® A4 CX16 Air Mattress	85
CuroCell® A4 CX20 Air Mattress	85
CuroCell IQ CX15 Air Mattress	86
CuroCell iA Automatic C17 Mattress	86
CuroCell iA Automatic C20 Mattress	87
Staebel OptiCell SMART R6 Mattress with Pump	87
Staebel OptiCell AIR1S R7 Mattress with Pump	88
Staebel OptiCell SMART R8 Mattress with Pump	88
CuroCell A4 Cirrus Air Mattress	89
Optimal 5zon Pressure Care Mattress	90
Optimal 5zon Mattress with Evac Handles	90
CuroCell Area Zone Mattress	91
CuroCell Sam Pro Air Mattress with Pump	91
PrimaCare 10	92
ViskoMatt 30	92
Precioso	93
Air2Care	93
CliniCare 100 HF Mattress	94
OptiCare Intergrated Mattress	94
Positioning & Support Products	95
Heel Lift Suspension Boot	99

NDIS

NDIS

MASS



Action Polymer Pad

Made from polymer gel, this product is effective in preventing the build up of pressure, promoting even pressure distribution. The polymer gel does not displace from areas requiring pressure relief and protects areas where shear forces and pressure points are prevalent. This gel pad can be used in many situations, and for spot protection.

Easily cleaned, ideal for multi client use.

Staebel Deluxe Foam Cushion

A comfortable cushion for the low risk user, the deluxe pressure relief foam cushions under pressure to provide a fluid like support. It helps to eliminate pressure points by transferring support uniformly. A top layer of memory foam is supported by 50 mm of standard foam to ensure no bottoming out. The vinyl cover is zipped with a soft stretch vinyl top.

- Suitable for low risk clients
- Zipped, removable cover included
- SWL is 100kg



Specifications:	
Product Code	340402
Width	230 mm
Length	280 mm
Thickness	95 mm

Also Available:		
Product Code	340403	340404
Width	280 mm	550 mm
Length	430 mm	500 mm
Thickness	60 mm	103 mm

Specifications:	
Product Code	3218
Width	410 mm
Length	410 mm
Thickness	85 mm
SWL	100 kg

Also Available:				
Product Code	3220	3221	3223	3225
Width	450 mm	510 mm	550 mm	450 mm
Length	400 mm	450 mm	450 mm	450 mm
Thickness	85 mm	85 mm	85 mm	85 mm



Low Resilient Foam Cushion

A very popular choice for the low risk user, the low resilient foam cushion flows under pressure to provide a fluid like support. It helps to eliminate and reduce pressure points by transferring support evenly.

- 2-way stretch top cover
- Zipped cover included

Relax Trio Cushion

The Relax Trio is an ultralight cushion that benefits from four sections of FloGel being positioned into a polyurethane base, with a viscoelastic foam upper layer. It is specially designed to conform to the shape of the user, increasing surface area and relieving pressure. Suitable for medium risk users, the cushion can be rotated, it has no designated front or back.

- 1.5" (38mm) thick
- Clear protective PU film coating is easy to wipe clean and is suitable for incontinent users
- Lightweight and easy to transport
- Includes Outer cover- Can be machine washed up to 60° C



Specifications:

Product Code	3210
Width	450 mm
Length	400 mm
Thickness	90 mm
SWL	100 kg

Specifications:

Product Code	35330
Width	410 mm
Length	410 mm
Thickness	38 mm
SWL	80 kg

Also Available:

Product Code	35332	35333
Width	450 mm	450 mm
Length	400 mm	450 mm
Thickness	38 mm	38 mm
SWL	80 kg	80 kg

NDIS

NDIS

MASS



Staebel Bariatric Cushion

A comfortable bariatric cushion for the low risk user, Helps to eliminate pressure points by transferring support uniformly. 2 foam layers, top layer of 50mm memory foam for comfort and is supported by 100mm of standard foam to ensure no bottoming out. Stretch cover with waterproof zip for easy cleaning.

- Suitable for low risk clients
- Zipped, removable cover included
- 260kg SWL
- Other sizes available upon request

Static Air Cushion

CuroCell® S.A.C. is a static air seat cushion for comfort, prevention and as an aid in the treatment of pressure injury up to category 2.

The cushion pre-filled and very easy to handle. The durable soft material for long term use provides a superior combination of comfort and support.

- Wipe with cleaning agent and/or disinfectants
- Air holes for venting heat and moisture
- Pre-inflated low profile design
- Fits all chairs
- Durable soft material for long term use
- Patient ready

***Pressure Mapping Reports and successful Case Studies Available on request**



Specifications:

Product Code	3270
Width	660 mm
Length	520 mm
Thickness	150 mm
SWL	260 kg

Also Available:

Product Code	3271	3272
Width	710 mm	610 mm
Length	550 mm	460 mm
SWL	260 kg	260 kg



Specifications:

Product Code	3417
Width	450 mm
Length	450 mm
SWL	300 kg

Also Available:

Product Code	3418
Width	500 mm
Length	500 mm
SWL	300 kg



Action Gel Professional Cushion

The Professional provides exceptional pressure relief for users at risk of skin breakdown or with a history of pressure ulcers. A full 4.1cm thickness of Akton polymer provides effective weight distribution on a contoured or flat surface. The Professional provides an additional 50% cushioning against pressure, shear and vibration. This is significant for those active clients using a power wheelchair who have concerns about bottoming out and potential effects of increased temperatures.

- Suitable for very high risk clients, including medium risk as a preventative measure
- Not restricted by any weight limit
- Easy to care for
- Basic cover included



Specifications:

Product Code	3540
Width	400 mm
Length	400 mm
Thickness	41 mm
SWL	Unlimited

Also Available:

Product Code	3541	3542	3543
Width	400 mm	450 mm	450 mm
Length	450 mm	450 mm	500 mm



Action Gel Pilot Cushion

The Pilot is a low-profile cushion that provides exceptional pressure relief for most users with a pressure injury history. 2.54cm Akton polymer provides immersion for bony prominences and protects clients with fragile skin integrity. This multi-use cushion is recommended for use on manual and power wheelchairs, truck or taxi seats, or for other uses requiring significantly reduced high-pressure areas, vibration, and shear. The low profile allows users to effectively reach the floor for more efficient foot-assisted propulsion. It also conforms to pre-contoured foam or rigid cushions to maximize stability and positioning when required. Custom sizes available upon request.

- 1" (25mm) thick
- Suitable for medium-high risk clients
- Not restricted by any weight limit
- Easy to care for
- Basic cover included



Specifications:

Product Code	3530
Width	400 mm
Length	400 mm
Thickness	25 mm
SWL	Unlimited

Also Available:

Product Code	3531	3532	3533
Width	450 mm	450 mm	500 mm
Length	400 mm	450 mm	450 mm

NDIS
MASS

NDIS
MASS



Action Centurian Cushion

The Centurian offers pressure and shear reduction for low-to-moderate risk users. Combined layers of 1.6 cm Akton[®] polymer and 1.6 cm foam are sealed in waterproof film for incontinence protection. Offering significant comfort for active clients who need a light to medium weight cushion with low maintenance, the Centurian is also the first choice of long term care facilities for decubitus prevention and treatment.

- 32 mm thick
- Suitable for medium and low risk clients
- Easy to care for
- Not restricted by any weight limit
- Basic cover included

Action Commuter Cushion

The Action Commuter offers basic pressure and shear reduction for low risk users. Combined 1.27cm Akton polymer and 4.45cm foam are sealed in a water proof film for incontinence protection.

- 57mm thick
- Suitable for medium risk clients
- Not restricted by any weight limit
- Easy to care for
- Basic cover included



Specifications:	
Product Code	3570
Width	400 mm
Length	400 mm
Thickness	32 mm
SWL	Unlimited

Also Available:			
Product Code	3571	3572	3573
Width	450 mm	450 mm	500 mm
Length	400 mm	450 mm	450 mm

Specifications:	
Product Code	3560
Width	400 mm
Length	400 mm
Thickness	57 mm
SWL	Unlimited

Also Available:			
Product Code	3561	3562	3563
Width	450 mm	450 mm	500 mm
Length	400 mm	450 mm	450 mm



Relax Easy Cushion

As the name suggests, this cushion from the Relax range is easy to use and is an ideal pressure relieving cushion for use in the nursing home environment. It has no designated up or down, or front or back, so can be simply placed into position. It features a layer of FloGel sandwiched between an upper and lower layer of viscoelastic foam. This combination of pressure relieving FloGel and viscoelastic foam facilitates effective pressure re-distribution while allowing the patient to immerse themselves comfortably into the cushion. It conforms to the shape of the user, increasing surface area and relieving pressure.

- 2" (50mm) thick
- Suitable for medium - high risk clients
- Clear protective PU film coating is easy to wipe clean and is suitable for incontinent users
- Lightweight and easy to transport
- Includes waterproof wipeable PU stretch outer cover
- Outer cover can be machine washed (60°C)

* High profile cushion available

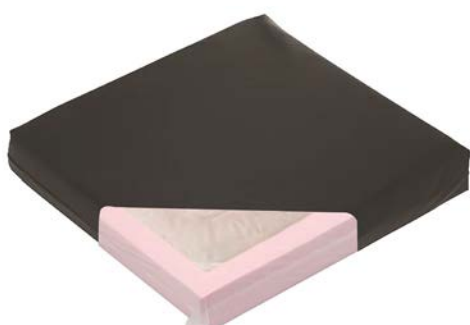


Specifications:

Product Code	35100
Width	400 mm
Length	400 mm
Thickness	50 mm
SWL	130 kg

Also Available:

Product Code	35312	35313	35315
Width	450 mm	450 mm	500 mm
Length	400 mm	450 mm	450 mm



Relax Easy High Profile Cushion

As the name suggests, this cushion from the Relax range is easy to use and is an ideal pressure relieving cushion for use who are over 100kgs and are in a long term care environment. It has no designated up or down, or front or back, so can be simply placed into position. It features a layer of FloGel sandwiched between an upper and lower layer of viscoelastic foam. This combination of pressure relieving FloGel and viscoelastic foam facilitates effective pressure re-distribution while allowing the patient to immerse themselves comfortably into the cushion. It conforms to the shape of the user, increasing surface area and relieving pressure.

- 85mm high and suited for those up to 190kgs
- Composition reduces the likelihood of hammock effect
- Suitable for medium - high risk clients
- Clear protective PU film coating is easy to wipe clean and is suitable for incontinent users
- Lightweight and easy to transport
- Includes waterproof wipe down outer cover
- Outer cover can be machine washed (60°C)



Specifications:



Product Code	35392
Width	400 mm
Length	400 mm
Thickness	85 mm
SWL	190 kg

Also Available:

Product Code	35394	35395	35396	35397
Width	450 mm	450 mm	500 mm	500 mm
Length	400 mm	450 mm	450 mm	500 mm

NDIS

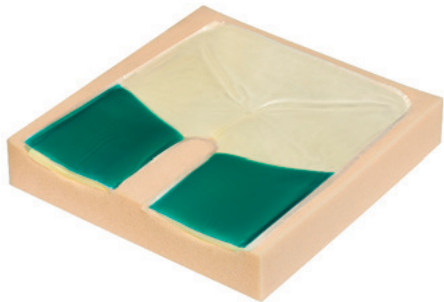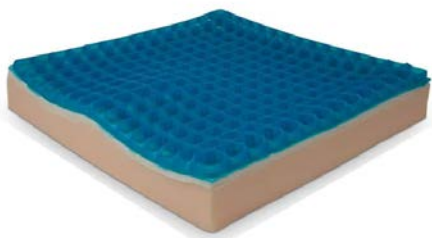NDIS



Relax Duo Profile Cushion

The Duogel is a unique combination of TruGel and FloGel that work together effectively to provide excellent pressure relief. This gel and foam combination allows for greater circulation and weight distribution. The Duogel provides for medial and lateral thigh positioning as well as support for the greater trochanter of the thigh bone. The thigh troughs are lined with TruGel and the ischial area is lined with FloGel which is set on a viscoelastic foam giving greater comfort and seating stability.

- 3" (75mm) thick
- Suitable for high risk clients
- Clear protective PU film coating is easy to wipe clean and is suitable for incontinent users
- Lightweight and easy to transport
- Includes machine washable cover

Relax DuoGel Cushion

The Duogel is a unique combination of TruGel and FloGel work together effectively to provide excellent pressure relief and weight distribution. This gel and foam combination allows for greater circulation and weight distribution. The Duogel provides for medial and lateral thigh positioning as well as support for the greater trochanter of the thigh bone. The thigh troughs are lined with TruGel and the ischial area is lined with FloGel which is set on a viscoelastic foam giving greater comfort and seating stability.

- 3" (75mm) thick
- Suitable for high risk clients
- Clear protective PU film coating is easy to wipe clean and is suitable for incontinent users
- Lightweight and easy to transport
- Includes machine washable cover




Specifications:	
Product Code	35380
Width	400 mm
Length	400 mm
Thickness	75 mm
SWL	130 kg



Specifications:	
Product Code	35300
Width	400 mm
Length	400 mm
Thickness	75 mm
SWL	130 kg

Also Available:			
Product Code	35306	35303	35304
Width	450 mm	450 mm	508 mm
Length	450 mm	450 mm	450 mm

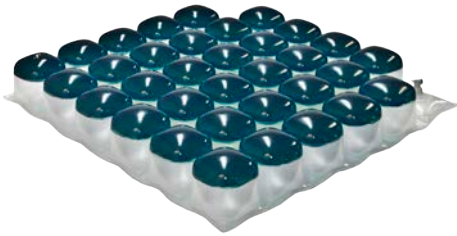

Relax GelCell Cushion

A unique combination of AIR and GEL technology effectively provides pressure relief and comfort to the high risk user. The Gelcell consists of 34 interconnected air-filled cells, each covered with a layer of pressure relieving silicone-based TruGel. The action of the air in the cushion enhances capillary function in the skin, maintaining blood flow and oxygen uptake. When pressure is applied, air circulates to each of the cells, distributing weight evenly, thus reducing high pressure concentrations around bony prominences. It ensures enhanced pressure re-distribution and also has a positive effect on skin temperature management whilst minimizing the effects of shear forces.

- 2" (50mm) thick
- Suitable for high risk clients
- Easy to wipe clean and is suitable for incontinent users
- Lightweight and easy to transport
- Easy to inflate and adjust
- Includes outer cover and pump



Specifications:

Product Code	35320
Width	400 mm
Length	400 mm
Thickness	50 mm
SWL	130 kg

Also Available:

Product Code	35323	3532320
Width	450 mm	500 mm
Length	450 mm	500 mm



Relax GelCell G2 Cushion

The Gelcell G2 consists of 48 interconnected cells, each covered with a layer of pressure relieving silicone-based TruGel

- Enhances pressure re-distribution
- Minimises the effects of shear forces
- Has a positive effect on skin temperature
- Silicone gel is impact dampening

The cells are contoured towards the rear, promoting good postural tilt. When pressure is applied, air circulates to each cell in both chambers distributing weight evenly. This process reduces high pressure concentrations around bony prominences, preventing skin breakdown.

The Gelcell is made up of 3 layers of polyurethane (PU) film making it highly durable and difficult to puncture. It has also easy to wipe clean and suitable for incontinent users.

- 1.5" (38mm) thick
- Suitable for high risk clients
- Easy to wipe clean and is suitable for incontinent users
- Lightweight and easy to transport
- Easy to inflate and adjust
- Includes outer cover and pump



Specifications:

Product Code	35324
Width	355 mm
Length	355 mm
Thickness	38 mm
SWL	130 kg

Also Available:

Product Code	3532315	35325	35328	35346
Width	380 mm	400 mm	400 mm	450 mm
Length	380 mm	400 mm	450 mm	450 mm



Starlock Cushion

First in the development of complete CELL LOCK technology, Star offers you the StarLock Cushion. This is our premier wheelchair cushion offering you the ultimate in aggressive positioning, unparalleled stability and the highest level of therapeutic prevention and treatment of pressure ulcers you've come to expect from a Star Cushion Product.

- No Weight Limit
- Approximately 4 lbs/1.36 kg
- Made with high-quality, flame resistant, latex-free neoprene rubber
- 24 Month Manufacturers Limited Warranty
- Available in 2", 3" & 4" cell heights.

Standard Contour Cushion

The 4" high cell cushion, or Standard Contour, offers deep immersion to provide a high level of pressure reduction for tissue protection and treatment. The Standard Contour's cells move and shape to the contact area providing a low shear, low friction surface ideal for the prevention for and healing of up to and including stage IV ulcers.

- No Weight Limit
- Light weight, Approximately 1.36 kg
- Made with high-quality, flame resistant, latex-free neoprene rubber
- Available in Single or Dual Valves
- Fitted cover, hand pump and repair kit is included



Specifications:				
Product Code	3324 - 1515			
Width	380 mm			
Length	380 mm			
Thickness	100 mm			
SWL	Unlimited			

Also Available:				
Product Code	3324 - 1715	3324 - 1717	3324 - 1917	3324 - 1919
Width	430 mm	430 mm	480 mm	490 mm
Length	380 mm	430 mm	430 mm	490 mm



Specifications:		
Product Code	3322 - 1715	
Width	380 mm	
Length	430 mm	
Thickness	100 mm	
SWL	Unlimited	

Also Available:		
Product Code	3322 - 1717	3322 - 1917
Width	430 mm	430 mm
Length	430 mm	480 mm

NDIS

MASS



Staebel Action Gel - Segment Pad (3 Per Bed)

The Action bed overlays (pictured) are made of Akton polymer which mimics human fatty tissue in consistency, does not migrate like gel and provides equally distributed pressure management, with varying sizes available, these overlays can be used to cover the entire surface or isolated areas of an existing standard mattress. Patients can use one, two or three segments in combination on their areas of need e.g. shoulders and hips, feet and hips or full body etc.

- Easily cleaned, no maintenance
- Not restricted by any weight limit
- 3 Segments cover length of single bed. Each sold separately




Specifications:

Product Code	3853
Width	630 mm
Length	690 mm
Thickness	18 mm
SWL	Unlimited

NDIS

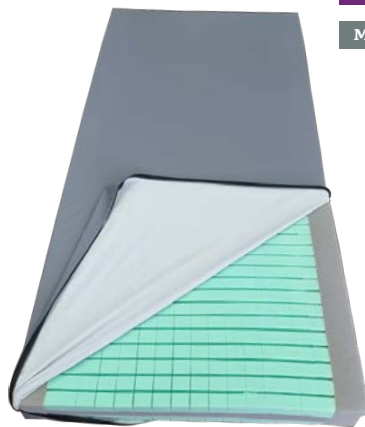MASS



Staebel ZONEcare Mattress

The ZONEcare foam mattress is a high quality foam mattress featuring a profile-cut pressure relieving foam insert. It provides pressure reduction and assists in the treatment of pressure ulcers, whilst enhancing patient comfort levels in homecare and residential care environments. The Eurocare's specialised inner cover aids in shear reduction on a profiling bed.

- Unique reinforced edges designed with longevity in mind
- Mattress turning is not required
- Waterproof zip for infection control
- Mattress Cover is vapour permeable, water resistant and fire retardant
- Mattress cover has welded seams to meet strict infection control standards
- Inner anti-shear cover



Specifications:

Product Code	3728Z
Width	850 mm
Length	1980 mm
Thickness	150 mm
SWL	150 kg

Also Available:

Product Code	3729Z
Width	950 mm
Length	1980 mm
SWL	150 kg

NDIS



Metzeler Ultratube Firm Feel Mattress

The heart of the Metzeler Ultratube Mattress is the newly developed Rubex Tubes System. The cold formed foam tubes in the mattress form a vertical tube system for perfect air circulation. The tubes, which differ in colour according to their degree of hardness, are integrated into the complete surface. This guarantees the user an ergonomically perfect lying position. Up to 200 litres of perspiration a year leave the body and go into the mattress. The Rubex Tubes System dispels this moisture through the vertical air channels.

- The tubes make up a system of s for perfect air circulation during sleep. The tubes of varying degrees of firmness are arranged to provide an ergonomically ideal sleeping posture reducing the risk of injury caused by shear and friction.
- Reinforced edge makes it easy to get in and out of bed and helps prevent roll-out.
- Standard with a bi-elastic, two-way stretch, waterproof, cover, designed to reduce shear and friction forces.
- Stain resistant, fire-resistant, waterproof, antiallergenic, antimicrobial, hygienically clean. Can be machine washed up to 95°C. Antibacterial treated.
- The Rubex Tubes System is symmetrically structured so that you always lie straight.
- Fully washable in temperatures up to 60°C

Staebel Soft EuroFlex Hygiene Mattress

Soft Euroflex is a hygiene mattress developed as an aid in the treatment of pressure ulcers up to category 3. The mattress is designed so that the entire mattress has relief for protruding body parts such as shoulders, heels and elbows. Soft Euroflex is made up of three different foam cores for good comfort and pressure relief. The mattress, , has a laminated surface layer of Lentex® and is completely devoid of seams and zippers, which provides increased conditions for good hygiene.

- Top White surface layer also makes it easier for staff to detect contamination on the surface of the mattress, which can be quickly cleaned by wiping and disinfecting.
- The surface layer is liquid-repellent and tight. Active Membrane technology allows the mattress to breathe while reducing the risk of liquid ingress.
- Approved for bed and mattress washing
- Lentex mattresses have been tested and approved completely without remark for use with bed and mattress cleaners
- The foam core consists of three layers with carefully selected foam types to provide good comfort and good pressure-relieving ability.
- Lentex® Soft Euroflex has been tested by Swerea IVF according to measurement method SS 876 00 13



Specifications:	
Product Code	3892
Width	1050 mm
Length	2000 mm
Thickness	180 mm
SWL	250 kg

Also Available:	
Metzeler Ultratube Soft Feel Mattress	
Product Code	3893

Specifications:	
Product Code	3271361
Width	850 mm
Length	2000 mm
Thickness	130 mm
SWL	150 kg

Sizes Available:	
Product Code	3271363
Width	1050 mm
Length	2000 mm



CuroCell® A4 CX16 Air Mattress

The CuroCell® A4 CX16 provides simplicity together with user comfort while effectively aiding in prevention and treatment of pressure injuries up to and including category 4, and unclassified pressure injuries. It is designed for excellent user comfort, reduction of shear forces and improvement of the microclimate.

The air mattress is 16cm thick with a foam base for added comfort. The alternation of the cells on the CX16 is every second cell down ensuring 50% of the body is being supported at one time. The pump is intelligent and fully automatic without need for manual settings. It identifies the user and automatically manages all settings to optimise the pressure according to the user's needs.

All mattresses are supplied with a four way stretch hygienic cover which meets standards for hygiene, quality and safety. All covers are manufactured from materials which reduce the risk of skin maceration.

- Unique CPR-function with transport function
- Heel Function
- Silent running and vibration free
- Three program modes: Pulsating, Alternating and CLP (constant low pressure)
- Individually replaceable cells with safety mattress in HR-foam base
- Audible and visual notifications
- Smart cable management



Specifications:

Product Code	3627
Width	1000 mm
Length	2000 mm
Thickness	160 mm
SWL	200 kg

Also Available:

Product Code	3626	3628
Width	850 mm	1050 mm
Length	2000 mm	2000 mm
SWL	200 kg	200 kg

CuroCell® A4 CX20 Air Mattress

The CuroCell® A4 CX20 provides simplicity together with user comfort while effectively aiding in prevention and treatment of pressure injuries up to and including category 4, and unclassified pressure ulcers. It provides features for excellent user comfort, reduction of shear forces and improvement of the microclimate.

The CX20 mattress is 20 cm thick and the alternation of the cells is every third cell down ensuring 75% of the body is being supported at one time. The pump is intelligent and fully automatic without need for manual settings. It identifies the user and automatically manages all settings to optimise the pressure according to the user's needs.

All mattresses are supplied with a four way stretch hygienic cover which meets standards for hygiene, quality and safety. All covers are manufactured from materials which reduce the risk of skin maceration.

- Unique CPR-function with transport function
- Heel Function
- Audible and visual notifications
- Silent running and vibration free
- Three program modes: Pulsating, Alternating and CLP (constant low pressure)
- Individually replaceable air cells with safety mattress
- Smart cable management



Specifications:

Product Code	3642
Width	1000 mm
Length	2000 mm
Thickness	150 mm
SWL	250 kg

Also Available:

Product Code	3641	3643
Width	850 mm	1050 mm
Length	2000 mm	2000 mm
SWL	250 kg	250 kg



CuroCell IQ CX15 Air Mattress

The CuroCell® A4 CX15 provides simplicity together with user comfort while aiding in the prevention and treatment of pressure ulcers/pressure injuries up to and including category 4, and unclassified pressure ulcers(1). It is designed for simplicity and provides features for good user comfort, reduction of shear forces and improvement of the microclimate.

CuroCell® A4 CX15 has an automatic control unit with a user-friendly interface. The mattress system regulates independently and without manual adjustment the inner pressure of the mattress to different values based on the user's weight and position. No manual action needs to be done to adjust the internal pressure of the mattress.

All mattresses are supplied with a four way stretch hygienic cover which meets standards for hygiene, quality and safety. All covers are manufactured from materials which reduce the risk of skin maceration.

- Individually replaceable cells with safety mattress
- Static head cells
- Heel function
- Multifunctional inner cover
- CPR-function with transport function
- Smart cable management
- Easy handling and usage
- Well protected zipper



Specifications:	
Product Code	3847
Width	850 mm
Length	2000 mm
Thickness	150 mm
SWL	220 kg

Sizes Available:	
Product Code	3848
Width	1000 mm
Length	2000 mm

CuroCell iA Automatic Ci17 Mattress

With Gentle Alternating Low Pressure (GALP), the CuroCell iA® system makes the step from mechanical to digital, from manual to autonomus, and from on-off comfort valves to a fine tuning of total support that we've never been able to deliver before in this segment. GALP has been developed to enable a system to sense, learn and take action whenever the high pressure that can cause ulcers/injuries arises.

CuroCell iA® uses artificial intelligence to constantly monitor pressure across the patient's contact area and calculates its response. The system identifies the patient's weight and height and can detect changes, such as when the patient adopts a new position.

All mattresses are supplied with a four way stretch hygienic cover which meets standards for hygiene, quality and safety. All covers are manufactured from materials which reduce the risk of skin maceration.

- Gentle Alternating Low Pressure (GALP) and Constant Low Pressure (CLP) and Maximum pressure (Caring mode)
- Panel Lock
- Information signals
- CPR connection with transport function
- Air Flow Control - saves energy and creates an even temperature in the support surface
- Low noise-and vibration level
- Individual and replacement cell construction
- Heel Function



Specifications:	
Product Code	3928
Width	850 mm
Length	2000 mm
Thickness	170 mm
SWL	210 kg

Sizes Available:	
Product Code	3929 3930 3931
Width	900 mm 1000 mm 1050 mm
Length	2000 mm 2000 mm 2000 mm

NEW

AUTO



NDIS

NEW

staebel



CuroCell iA Automatic Ci20 Mattress

With Gentle Alternating Low Pressure (GALP), the CuroCell iA® system makes the step from mechanical to digital, from manual to autonomous, and from on-off comfort valves to a fine tuning of total support that we've never been able to deliver before in this segment. GALP has been developed to enable a system to sense, learn and take action whenever the high pressure that can cause ulcers/injuries arises.

CuroCell iA® uses artificial intelligence to constantly monitor pressure across the patient's contact area & calculates its response. The system identifies the patient's weight and height and can detect changes, such as when the patient adopts a new position.

All mattresses are supplied with a four way stretch hygienic cover which meets standards for hygiene, quality and safety. All covers are manufactured from materials which reduce the risk of skin maceration.

- Gentle Alternating Low Pressure (GALP). Constant Low Pressure (CLP), Maximum pressure (Caring mode)
- Panel Lock
- Information signals
- CPR connection with transport function
- Air Flow Control - saves energy and creates an even temperature in the support surface
- Low noise-and vibration level
- Individual and replacement cell construction
- Heel Function



Specifications:

Product Code	3935
Width	850 mm
Length	2000 mm
Thickness	200 mm
SWL	230 kg

Sizes Available:

Product Code	3936	3937	3938
Width	900 mm	1000 mm	1050 mm
Length	2000 mm	2000 mm	2000 mm

Staebel OptiCell SMART R6 Mattress with Pump

In order to improve the efficiency of the care staff, time has been spent on developing an easy-to-use and intuitive panel. In addition, learning time is drastically reduced as the same pump is used regardless of the mattress. The pump has built-in features such as:

- Automatic adjustment to the user's weight and position in bed.
- Automatic return to normal mode after nursing mode.
- Battery saves pump settings in case of power failure

You can easily and quickly decouple individual cells without affecting the mattress's function. Each individual cell has a built-in check valve that can be released with a simple hand movement. The cell can then remain in the mattress case or be removed completely. It allows all switching to take place at the top - True Alternating - to reduce the risk of nausea.

- A SMART pump – three compatible series with a total of eleven mattresses.
- The user-friendly display has a layout that is easy to understand.
- The pump has a smooth surface without cavities to facilitate daily cleaning and disinfection
- The mattress system provides good comfort, good support and pressure relief through immersion and envelopment.
- Micro Air Loss air cells provide good ventilation and help reduce heat and moisture.
- The vapor permeable OptiTex™ mattress cover is both water resistant and fire rated



Specifications:

Product Code	3271091
Width	850 mm
Length	2000 mm
Thickness	160 mm
SWL	230 kg



Staebel OptiCell AIR1S R7 Mattress with Pump

In order to improve the efficiency of the care staff, time has been spent on developing an easy-to-use and intuitive panel. In addition, learning time is drastically reduced as the same pump is used regardless of the mattress. The pump has built-in features such as:

- Automatic adjustment to the user's weight and position in bed.
- Automatic return to normal mode after nursing mode.
- Battery saves pump settings in case of power failure

You can easily and quickly decouple individual cells without affecting the mattress's function. Each individual cell has a built-in check valve that can be released with a simple hand movement. The cell can then remain in the mattress case or be removed completely. It allows all switching to take place at the top - True Alternating - to reduce the risk of nausea.

- The user-friendly display has a layout that is easy to understand.
- The pump has a smooth surface without cavities to facilitate daily cleaning and disinfection
- The mattress system provides good comfort, good support and pressure relief through immersion and envelopment.
- Micro Air Loss air cells provide good ventilation and help reduce heat and moisture.
- The vapor permeable OptiTex™ mattress cover is both water resistant and fire rated

Staebel OptiCell SMART R8 Mattress with Pump

In order to improve the efficiency of the care staff, time has been spent on developing an easy-to-use and intuitive panel. In addition, learning time is drastically reduced as the same pump is used regardless of the mattress. The pump has built-in features such as:

- Automatic adjustment to the user's weight and position in bed.
- Automatic return to normal mode after nursing mode.
- Battery saves pump settings in case of power failure

You can easily and quickly decouple individual cells without affecting the mattress's function. Each individual cell has a built-in check valve that can be released with a simple hand movement. The cell can then remain in the mattress case or be removed completely. It allows all switching to take place at the top - True Alternating - to reduce the risk of nausea.

- The user-friendly display has a layout that is easy to understand.
- The pump has a smooth surface without cavities to facilitate daily cleaning and disinfection
- The mattress system provides good comfort, good support and pressure relief through immersion and envelopment.
- Micro Air Loss air cells provide good ventilation and help reduce heat and moisture.
- The vapor permeable OptiTex™ mattress cover is both water resistant and fire rated



Specifications:	
Product Code	3271171
Width	850 mm
Length	2000 mm
Thickness	180 mm
SWL	240 kg



Specifications:	
Product Code	3271111
Width	850 mm
Length	2000 mm
Thickness	200 mm
SWL	250 kg

Sizes Available:		
Product Code	3231112	3231113
Width	900 mm	1050 mm
Length	2000 mm	2000 mm



CuroCell A4 Cirrus Air Mattress

CuroCell® A4 Cirrus is used as an aid in the prevention and treatment of pressure ulcers up to and including category 4 and unclassified pressure ulcers/pressure injuries⁽¹⁾. It is a replacement system suitable for all healthcare environments.

CuroCell® A4 Cirrus combines air with a soft foam surface and provides good user comfort and pressure redistribution⁽³⁾. The added heel function offers pressure redistribution on the heels.

- Automatic control unit which adapts individual to the user
- User friendly interface
- Audible and visual notifications
- Three program modes: pulsating mode, alternating mode and constant low pressure
- Silent running and vibration free
- Air flow control; Constant temperature, good user comfort, energy reduction
- Soft foam surface
- Smart cable management
- Easy handling and cleaning
- Heel slope, for pressure redistribution of heels
- CPR-function with transport function



Specifications:

Product Code	3841
Width	850 mm
Length	2000 mm
Thickness	210 mm
SWL	230 kg

Also Available:

Product Code	3843	3844
Width	1000 mm	1050 mm
Length	2000 mm	2000 mm
SWL	230 kg	230 kg

Case Study

History

- Female, born 1932
- Operated breast cancer 2003, vertebral compression fractures, hip osteoarthritis with chronic pain
- Arterial fibrillation, hypertension and bilateral leg ulcers.
- Pressure ulcer category 4 in sacrum 7x2 cm with fistulas, R/T ongoing treatment leg ulcer.
- Living at home with help from home health care.
- Mentally clear, no incontinence, normal diet, sitting up during meals.
- Taking regular pain relief, Oxykontin and Oxynorm.
- Points 23, modified Norton scale.
- BMI 21



Action

- 2010–12–31 prescribed CuroCell Cirrus, alternating pressure mattress with foam core, already have a pressure relieving wheelchair cushion.
- Can manage repositioning by herself, aided move from bed to wheelchair. Wounds redressed daily. Daily supplements of nutrient drinks.
- 2011-01-21 admitted to hospital due to the leg ulcers, mattress goes with the patient.
- 2011-02-18 left leg amputated thigh bone After hospital stay the patient is accommodated short-term at a home and is later placed in a permanent nursing home, mattress goes with the patient. Wounds is decreasing in size.
- 2012-04-04 Wounds completely healed

NDIS
MASS

NDIS
MASS



Optimal 5zon Pressure Care Mattress

The Optimal 5zon® is used in all types of health care environments for pain relief and as an aid in the prevention and treatment of pressure ulcers up to and including category 2. A soft and flexible mattress replacement system. The Optimal 5zon® offers very good pressure relief properties and is now also equipped with the unique sit-function DuoCore®. This is an insert within the mattress core, that provides support when seated, while offering good comfort in the reclining position. It comes complete with a durable cover that is easy to clean and is equipped with divisible, drip shield zip. The soft and elastic material is vapour permeable, which reduces the risk of skin shear.

- Chlorine can be utilised for cleaning purposes
- Frequent rotation is not required
- Suitable for clients up to 180kg
- Patented sit function DuoCore® is designed to give support when sitting and also to provide exceptional comfort when lying
- Tested by Berlin Cert to prove pressure relief capability and hygienic standard.
- Mattress took on rigorous Clinical studies - ask for a copy of these, to peruse.



Specifications:	
Product Code	3771
Width	850 mm
Length	2000 mm
Thickness	140 mm
SWL	180 kg

Also Available:		
Product Code	3778	3775
Width	1000 mm	1050 mm
Length	2000 mm	2000 mm
SWL	180 kg	180 kg

Optimal 5zon Mattress with Evac Handles

The Optimal 5zon® is ideal for all healthcare environments and is used as an aid in the prevention and treatment of pressure ulcers up to and including category 2, and for pain therapy. It is a soft and flexible mattress replacement that is placed directly on the bed and replaces both the standard mattress and mattress pad.


- Highly elastic, high-quality material and a unique zone and module system gives effective pressure equalisation.
- Specially designed to minimise shear forces
- Simple and effective cleaning, cover has a very well documented hygiene standard
- Evacuation handles on the side of the mattress
- The hygienic cover is vapour permeable, which reduces the risk of skin maceration
- The patented sit function DuoCore® is designed to give support when sitting and also to provide good comfort when lying
- Heel Protection provides comfort for relief of heel. The user can move the leg and make position changes. Comes with a removable liquid resistant hygiene cover.



Specifications:	
Product Code	3770
Width	800 mm
Length	2000 mm
Thickness	140 mm
SWL	180 kg

NDIS

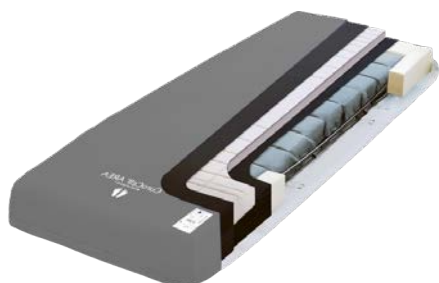MASS



CuroCell Area Zone Mattress

The AREA Zone mattress replacement system is designed with special valve technology to provide an enhanced mattress surface for the prevention and management of skin breakdown. The AREA provides clinically proven non-powered pressure redistribution and therapeutic support through the use of interconnected pressure redistribution air sectors. The valve system allows for communication between the air chambers which in turn provides a unique, custom low pressure profile for any patient. Upon ingress, a natural chain of events begins in each of the connected air chambers. The slow release valves allow for pressure changes to take place very gently. The mattress top is a layer of visco type memory foam which further enhances the low pressure profile.

- Suitable for high risk patients, and also medium risk as a preventative measure
- Static air replacement mattress system (not an overlay)
- Reinforced sides for improved lateral stability
- Unique integrated heel protection area
- Anti-shear inner cover
- No pump or power cord for hassle free use



Specifications:

Product Code	3678Z
Width	850 mm
Length	2000 mm
Thickness	180 mm
SWL	230 kg

Also Available:

Product Code	3679Z	3686Z	3680Z
Width	900 mm	1000 mm	1050 mm
Length	2000 mm	2000 mm	2000 mm
SWL	230 kg	230 kg	230 kg

NDIS

MASS



CuroCell Sam Pro Air Mattress with Pump

The CuroCell SAM PRO is a static air replacement mattress which has been proven to provide exceptional comfort, pressure and pain relief. It is used preventively and in the treatment of pressure ulcers up to and including sore degree 3. The mattress is divided into three zones, which amongst other advantages gives better support when sitting on the mattress. The support pump delivered together with the mattress is easy to use and very quiet. The CuroCell SAM PRO is delivered with a anti-shear inner cover and a soft, elastic, breathable outer cover which has a very well documented hygienic standard.

- Suitable for high risk patients, and also medium risk as a preventative measure
- Dynamic air mattress replacement system (lies directly on bed surface)
- Air cells produced from soft, durable polyurethane Not PVC. Significant noise reduction when cells move compared to other brands
- Comes complete with an anti-shear inner cover which has been proven to reduce pressure by up to 8%
- Easy to use dial pump with table to calculate most appropriate setting
- 2 year manufacturers warranty
- Quality product from Europe



Specifications:

Product Code	3683
Width	850 mm
Length	2000 mm
Thickness	160 mm
SWL	250 kg

Also Available:

Product Code	3670
--------------	------



PrimaCare 10

PrimaCare 10 is made of a monolithic foam core. The upper surface of the mattress has transverse and longitudinal profiling evenly distributed over its entire surface. Due to this profiling, there is better air circulation between the patient's skin and the base. They also help to optimally distribute pressure.

Increased density of profiling at the head and foot of the mattress: Due to the increased density of profiling in the head and foot area, the mattress ensures better comfort and distribution of pressure between the patient and the mattress.

Mattress cover made of vapour-permeable material: Helps reduce patient perspiration and also reduces skin moisture.

ViskoMatt 30

The body of the mattress consists of a combination of two layers: a lying surface with a viskoelastic foam and a support layer from cold foam.

Combination of viskoelastic and cold foam: Because of two different types of foam the mattress helps to provide pressure relief and also ensures patient comfort.

Wave Contour Cut: 2 foam layers not glued. This enables uniform distribution of the contact pressure without hard seams and helps to provide higher pressure relief.

Edge Zone Reinforcement: Stable reinforcement on both long sides helps to the relief of self-mobilization of the patient and prophylaxis of falls.

Tricot sock - a special textile layer between foam and cover: this feature reduces shear and friction forces by cover and tricot, which follow patient and not get stuck on the mattress. The patient doesn't feel turning movement during repositioning on the bed

- The mattress cover is made of PERLASTIC® Silber cover dark grey with concealed L-zipper and two handles on one side.
- Reinforced Edges
- Upper Layer is 5 cm Viskoelastic foam and the lower layer 9 cm Cold foam
- With the correct positioning of the patient, we would recommend to using it for patients with Low, Medium and High levels of risk of Pressure Injuries.



Specifications:	
Product Code	3760
Width	860 mm
Length	2000 mm
Thickness	140 mm
SWL	200 kg

Specifications:	
Product Code	30406VM30
Width	860 mm
Length	2000 mm
Thickness	140 mm
Weight	14 kg
SWL	165 kg



Precioso

The Precioso is a 24-hour system for wound care management. The system comes as an overlay or mattress replacement and can be used with a specific cushion. The Precioso will significantly contribute towards the healing and prevention of all grades of pressure ulcers. Adequate for use in general wards in hospitals and nursing homes.

The Precioso is fitted with 4 different pressure sensors and a microprocessor controlling the precise cell pressure. For patient comfort it is possible to increase or decrease the pressure by adjustment the levels on the membrane panel.

Alternating Pressure Therapy (3 cell technology, zero pressure, 7.5 min.cycle) is the fundamental therapeutic mode that simulates natural body movement and therefore helps to prevent and treat pressure ulcers.

Air2Care

The Air2Care alternating pressure system has been designed for the prevention and management of pressure ulcers in those at low to moderate risk. The Air2Care is suitable for use in various healthcare settings including long term care. For convenience the system is available as either a mattress replacement or as an overlay giving the flexibility to choose a solution suitable for the patient's and caregiver's individual requirements. Two options of system control unit are available. The addition of an active cushion to the Air2Care system gives 24 hour prevention and management of pressure ulcers which is simple to use and reliable.



Specifications:

Product Code	36021
Width	860 mm
Length	2000 mm
Thickness	170 mm
SWL	160 kg

Specifications:

Product Code	36010
Width	860 mm
Length	2000 mm
Thickness	140 - 170 mm
SWL	200 kg



CliniCare 100 HF Mattress

The CliniCare 100 HF hybrid mattress is comfortable and highly efficient for the prevention of pressure injuries. It uses a combination of an active system and passive mattress to provide many benefits across the spectrum of hospital care.

When a patient's condition requires active therapy, the mattress can be connected to a compressor for alternating pressure therapy in a 2 cell cycle to enhance and optimize pressure redistribution and pressure injury prevention.

The mattress provides pressure redistribution in non-powered mode, air moves from cell to cell as reaction to patient's body movement and weight.

A viscoelastic foam layers adapts to the body to provide comfort and proper support without exerting pressure on exposed area. The mattress has a 7° slope that helps to reduce pressure in vulnerable heel area.

The system inside the mattress keep blowing air and the mattress cover is automatically dried. In the context of pressure injuries, microclimate usually refers to skin temperature and moisture conditions at the skin support surface interface.

- With microclimate



Specifications:

Product Code	3817	3818
Width	860 mm	860 mm
Length	2000 mm	2080 mm
Thickness	200 mm	200 mm
SWL	230 kg	230 kg

OptiCare Intergrated Mattress

New smart technology for the prevention of pressure ulcers. The system detects when the patient has entered the bed. This automatically initiates the optimisation process, thereby saving the time normally required for nurses to set the mattress parameters. The mattress settings are automatically optimised for each patient based on the patient's weight, height, and even the patient's position on the mattress.

- The mattress is equipped with a whole range of sophisticated functions which enhance the therapeutic and preventive effect of the optimised pressure.
- OptiCare mattress is compatible with Multicare and Eleganza 5 beds. The controls for both the mattress and bed are integrated into a single control panel built into the bed. The customisable user interface provides the clinical staff with the option to modify the controls on the interface and allows control of only the required features.
- In the context of pressure ulcers, microclimate usually refers to skin temperature and moisture conditions at the skin-support surface interface. OptiCare helps maintain the natural thermoregulation by circulating air through the mattress surface.
- Helps control the heat and humidity of the patient's skin, which is an important part of maintaining tissue viability.
- The heel section has a 7-degree slope and is lined with a softer foam that gives the patient above-standard comfort. This significantly reduces the risk of heel pressure ulcers.
- Integrated systems to be used with Multicare and Eleganza 5 beds only.



Specifications:


Product Code	74200 - 1
Width	860 mm
Length	2140 mm
Thickness	220 mm
SWL	250 kg



Staebel Positioning & Support Wedge

This foam positioning support wedge is used to place a user in a side position in order to relieve the sacrum, which is an extra sensitive area. For the side position, two support wedges are used together. The support wedge is easy to use and easy to clean. It is supplied with a zipped, liquid resistant and vapour permeable cover. The cover can be wiped or machine washed.

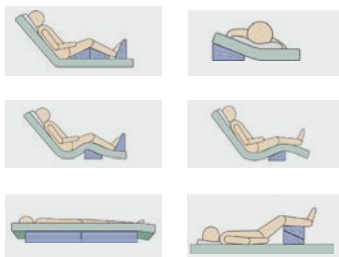- Soft comfortable foam, stretch vinyl covered



Specifications:

Product Code	3689
Width	260 mm
Length	700 mm
Height	50 - 150 mm



Back Pillow

The Back Pillow is a winged-pillow for back relief. The pillow is placed behind the back, and the wing rests under the user to set the position.



Specifications:

Product Code	32301
Width	500 mm
Length	800 mm

NDIS

NDIS



Positioning Roll

The positioning roll can give support to the resident in sitting and lying positions. It can be placed in various ways to provide support and relief.

- Supports shoulders, spine, hips, buttocks and heels
- All-round support in seated and reclining positions
- Relieves the knees and supports feet
- Supports in lateral position



Positioning Pillow

The Positioning Pillow range is an all-round pillow that can be used for different types of patient positioning.

- Support after a hip replacement
- Can be positioned to help breathing problems and general expansion of the lungs
- Positioning after a stroke
- Supports in lateral position

*** Available in different sizing**



Specifications:	
Product Code	32307
Width	150 mm
Length	600 mm

Also Available:		
Product Code	32308	32309
Width	200 mm	200 mm
Length	1250 mm	2250 mm

Specifications:	
Product Code	32304
Width	250 mm
Length	800 mm



Multi Pillow

The Multi Pillow is a versatile pillow that can be placed between the knees or under legs and heels for support and relief.

- Supports the feet, counters heel pressure against surfaces
- Locked position in connection with spasticity
- Supports in lateral position
- Relief for heels
- Supports under the knees

Molly Heel Protection

- Easy to use
- Equalizes pressure
- Faster healing
- Improves circulation
- Makes it easier to freely move the legs
- Prevents crossing the legs after hip surgery
- Soft and comfortable
- Wipe with cleaning agent or disinfect
- Machine wash max. 95°C
- Easy to clean
- Liquid resistant
- Long durability



Specifications:

Product Code	32302
Width	450 mm
Length	450 mm

Specifications:

Product Code	32311
Leg Length	780 - 850 mm
Thickness	75 mm

Also Available:

Product Code	32312	32313
Leg Length	860 - 940 mm	1020 - 1030 mm
Thickness	90 mm	110 mm

NDIS

NDIS

NDIS



Positioning Curved Pillow

If a person makes too few spontaneous movements, there is an increased risk of pressure ulcers. Regular changes of position reduce the risk of pressure ulcers. Position as required in both reclining and seated position.

- Curera Positioning Neck Pillow offers support and pressure relief of head/neck.
- Machine washable

Positioning Triangle Pillow

A Triangle Pillow under the knees can stabilize the pelvis to prevent sliding in combination with an elevated backrest.

The Triangle Pillow is also suitable for user's with shortened hip flexors and/or hamstrings. The increased flexion in hips and knees can relieve lower back strain. In some cases, the Triangle Pillow can aid in alleviating pressure from under the heels.

Positioning S Pillow

The 30° side lying is commonly used to provide pressure relief to immobilized users. For users unable to lie in prone positions, a supported semi-prone position can be used as an alternative.

The S-Shaped Pillow can be used to simulate 30° side lying position (when placed behind the user) and a semi-prone position (when placed in front of the user), instead of a Positioning Roll or Support Wedges. The S-shaped Pillow is available in two sizes, depending on the length of the user.



Specifications:

Product Code	32321
Width	750 mm
Length	1200 mm
Thickness	300 mm

Sizes Available:

Product Code	32325
Width	750 mm
Length	400 mm
Thickness	250 mm



Specifications:

Product Code	32318
Width	280 mm
Length	470 mm
Thickness	280 mm

Sizes Available:

Product Code	32319
Width	350 mm
Length	620 mm
Thickness	350 mm



Specifications:

Product Code	32322
Width	280 mm
Length	1700 mm



Positioning Torso Support Pillow

The Torso Support Pillow can be placed either behind the user for support when seated on the bedside or in front of the user, preferably in combination with a Curved Pillow behind.

The Trunk Support Pillow can also be used as a support in a half sitting position when in bed or while seated in a chair or wheelchair.

Positioning Neck Horseshoe Pillow

If a person makes too few spontaneous movements, there is an increased risk of pressure ulcers. Regular changes of position reduce the risk of pressure ulcers. Position as required in both reclining and seated position.

- Curera Positioning Neck Pillow offers support and pressure relief of head / neck.
- Machine washable

Heel Lift Suspension Boot

The Heel Lift Suspension Boot is designed to provide the most effective solution for prevention and assistive healing of pressure ulcers at the heel, while delivering positioning and alignment of the lower limb.

Heel Lift Suspension Boot Convoluted:

By using an inbuilt wedge of soft, Convoluted medical grade foam the heel is ensured a no contact area while giving good support over the calf.

- The base of the Heel Lift is made up of a tricot low shear/friction surface which allows the boot to slide easily over the mattress surface.
- The Heel Lift comes with an extra foam piece which can be fitted to prevent footdrop or prevent hip rotation.
- The Heel Lift fits comfortably onto the foot and is secured with Velcro® straps.



Smooth



Convoluted



Specifications:

Product Code	3140
Style	Smooth Standard

Sizes Available:

Product Code	Style
3141	Convoluted Standard
3142	Smooth Petite
3143	Convoluted Petite
3146	Smooth Bariatric

Specifications:

Product Code	32324
Width	650 mm
Length	850 mm
Thickness	300 mm

Specifications:

Product Code	32303
Width	380 mm
Length	500 mm

Mattress Matrix						
Mattress Name	Size	Comments	Waterlow scale / Risk Score			
			At Risk (10+)	High (15+)	Very High (20+)	Skin Intact
<div>Zonecare Pressure Care Foam Mattress</div> <div></div>	Single	Standard mattress for all replacements	✓			✓
	King Single					
<div>Optimal 5zon Pressure Care Mattress</div> <div></div>	Single	Highly elastic, high quality material and a unique zone and module system gives effective pressure. Specifically designed to minimise the effects of shear forces. Duocore is designed to give support when sitting and provide good comfort when lying. The cover is vapour permeable, which reduces the risk of skin maceration. The covers are easy to clean and handle and has a very well documented hygienic standard.	✓			
	King Single					
<div>Curocell S.A.M F16 Pro Static Air Mattress Replacement & Pump</div> <div></div>	Single	A protective / preventative strategy for residents with intact skin but identified as at risk of skin integrity issues of pressure injuries. Beautifully comfortable due to the static nature plus it comes with the comfort cover so you do not create Cat 1-3 pressure injuries.		✓	✓	
	King Single					
<div>Curocell A4 CX16</div> <div></div>	Single	Cat 1 - 4 Pressure Injuries. CX16 is a 3 cell cycle, opposed to the regular 2 cell. 2 cells up, 1 down means the body has 66% support not 50%, however the pulsating is the best program for the greatest comfort.		✓	✓	
	King Single	Great for persons who need the pressure care, but don't like the sensation. Excellent for where high incontinence issues are present.				
<div>CuroCell Cirrus</div> <div></div>	Single	Cat 1 - 4 Pressure Injuries. Uniquely designed to help in pain management, for end of life care. Thick top layer of foam to help create sensory recognition of the user with pulsation mode cells below that create a wave like relief effect opposed to the traditional alternation.		✓	✓	
	King Single	Heel slope, silent pump with sit position function, and option to adjust the timing of the cell change.				

Current Pressure			May suit people with particular needs:				
Skin Broken	Underweight	Bariatric	Client In Pain	Bed Mobility		Bed Transfers	Heels at risk
				Independent	Assisted / Static		
		 0-125KG		✓		Firm Transfers	
	(~30KG)	 0-180KG	✓	✓	✓	YES - the Duo Core adds stability	Zero pressure ability
	(~30KG)	 0-250KG	✓	May be difficult	✓	Staebel as it is a static mattress	Inbuilt heel slope
	(~1KG)	 0-200KG	✓	May be difficult	✓	Max function	Inbuilt heel slope
	(~1KG)	 0-230KG	✓	May be difficult	✓	Foam surround	Inbuilt heel slope

Cushion Matrix					
		Weight Capacity	Cushion Profile Height	Size Range	
				Width - inches	Width - mm
3210 Budget		100 kg	50 mm		
SAC		160 kg	25 mm		
Action Pilot		Unlimited	30 mm		
Action Professional		Unlimited	45 mm		
Relax Easy Low Profile		130 kg	50 mm	14 x 16" to 22 x 18"	360 x 400mm to 550 x 450mm
Relax Easy High Profile		190 kg	85 mm	16 x 16" to 20 x 20"	400 x 400mm to 500 x 500mm
Duogel		130 kg	85 mm	14 x 14" to 24 x 18"	360 x 360mm to 600 x 450mm
Duo Profile		130 kg	85 mm	16 x 16" to 18 x 18"	400 x 400mm to 450 x 450mm
Gel Cell		130 kg	40 mm	15 x 14" to 20 x 18"	380 x 360mm to 500 x 450mm

PI Rating	Thigh Support	Contoured seat	Skin Integrity			Suitable for Lounge chairs
			Skin Intact	Skin Broken	Boney Prominence	
Prevention	✗	✗	✓	✗	✗	✗
Stage 2	✗	✗	✓	✗	✗	✓
Stage 3	✓	✗	✓	✗	✗	✓
Stage 4	✓	✗	✓	✓	✓	✗
Stage 2	✗	✗	✓	✗	✓	✓
Stage 3	✗	✗	✓	✓	✓	✗
Stage 3	✓	✓	✓	✓	✓	✗
Stage 2	✓	✓	✓	✗	✓	✗
Stage 4	✗	✓	✓	✓	✓	✓

WATERLOW PRESSURE ULCER PREVENTION/TREATMENT POLICY
RING SCORES IN TABLE, ADD TOTAL. MORE THAN 1 SCORE/CATEGORY CAN BE USED

BUILD/WEIGHT FOR HEIGHT	SKIN TYPE VISUAL RISK AREAS	SEX AGE	MALNUTRITION SCREENING TOOL (MST) (Nutrition Vol.15, No.6 1999 - Australia)	
AVERAGE BMI = 20-24.9 ABOVE AVERAGE BMI = 25-29.9 OBESE BMI > 30 BELOW AVERAGE BMI < 20 BMI=Wt(Kg)/Ht (m) ²	HEALTHY TISSUE PAPER DRY 1 OEDEMATOUS CLAMMY, PYREXIA 2 DISCOLOURED GRADE 1 3 BROKEN/SPOTS GRADE 2-4	0 MALE 1 FEMALE 1 14 - 49 1 50 - 64 2 65 - 74 2 75 - 80 3 81 +	A - HAS PATIENT LOST WEIGHT RECENTLY YES - GO TO B NO - GO TO C UNSURE - GO TO C AND SCORE 2	B - WEIGHT LOSS SCORE 0.5 - 5kg = 1 5 - 10kg = 2 10 - 15kg = 3 > 15kg = 4 unsure = 2
C - PATIENT EATING POORLY OR LACK OF APPETITE 'NO' = 0; 'YES' SCORE = 1	NUTRITION SCORE If > 2 refer for nutrition assessment / intervention	SPECIAL RISKS		
CONTINENCE	MOBILITY	TISSUE MALNUTRITION	NEUROLOGICAL DEFICIT	
COMPLETE/ CATHETERISED URINE INCONT. FAECAL INCONT. URINARY + FAECAL INCONTINENCE	FULLY 0 RESTLESS/FIDGETY 1 APATHETIC 2 RESTRICTED 3 BEDBOUND e.g. TRACTION CHAIRBOUND e.g. WHEELCHAIR	0 1 2 3 4 5	8 8 5 5 2 1	4-6 4-6 4-6 5 5 8
SCORE		MAJOR SURGERY or TRAUMA		
10+ AT RISK		ORTHOPAEDIC/SPINAL		
15+ HIGH RISK		ON TABLE > 2 HR#		
20+ VERY HIGH RISK		ON TABLE > 6 HR#		
		MEDICATION - CYTOTOXICS, LONG TERM/HIGH DOSE STEROIDS, ANTI-INFLAMMATORY MAX OF 4		

© J Waterlow 1985 Revised 2005*
Obtainable from the Nook, Stoke Road, Henlade TAUNTON TA3 5LX
* The 2005 revision incorporates the research undertaken by Queensland Health.
Scores can be discounted after 48 hours provided patient is recovering normally
www.judy-waterlow.co.uk

REMEMBER TISSUE DAMAGE MAY START PRIOR TO ADMISSION, IN CASUALTY, A SEATED PATIENT IS AT RISK ASSESSMENT (See Over) IF THE PATIENT FALLS INTO ANY OF THE RISK CATEGORIES, THEN PREVENTATIVE NURSING IS REQUIRED A COMBINATION OF GOOD NURSING TECHNIQUES AND PREVENTATIVE AIDS WILL BE NECESSARY ALL ACTIONS MUST BE DOCUMENTED

PREVENTION
PRESSURE
REDUCING AIDS
Special
Mattress/beds:

10+ Overlays or specialist foam mattresses.
15+ Alternating pressure overlays, mattresses and bed systems
20+ Bed systems: Fluidised bead, low air loss and alternating pressure matting mattresses

Note: Preventative aids cover a wide spectrum of specialist features. Efficacy should be judged, if possible, on the basis of independent evidence.
No person should sit in a wheelchair without some form of cushioning. If nothing else is available - use the person's own pillow. (Consider infection risk)

Cushions:

10+ 100mm foam cushion
15+ Specialist Gel and/or foam cushion
20+ Specialised cushion, adjustable to individual person.
Avoid plastic draw sheets, inco pads and tightly tucked in sheet/sheet covers, especially when using specialist bed and mattress overlay systems
Use duvet - plus vapour permeable membrane.

NURSING CARE

General

HAND WASHING, frequent changes of position, lying, sitting. Use of pillows

Pain

Appropriate pain control

Nutrition

High protein, vitamins and minerals

Patient Handling

Correct lifting technique - hoists - monkey poles
Transfer devices
Real Sheepskin - bed cradle

Patient Comfort Aids

Operating Table

100mm(4ins) cover plus adequate protection

Theatre/A&E Trolley

Skin Care

General hygiene, NO rubbing, cover with an appropriate dressing

WOUND GUIDELINES

Assessment

odour, exudate, measure/photograph position

GRADE 1

Discolouration of intact skin not affected by light finger pressure (non-blanching erythema)

GRADE 2

This may be difficult to identify in darkly pigmented skin
Partial thickness skin loss or damage involving epidermis and/or dermis
The pressure ulcer is superficial and presents clinically as an abrasion, blister or shallow crater

GRADE 3

Full thickness skin loss involving damage of subcutaneous tissue but not extending to the underlying fascia
The pressure ulcer presents clinically as a deep crater with or without undermining of adjacent tissue

GRADE 4

Full thickness skin loss with extensive destruction and necrosis extending to underlying tissue.

Dressing Guide

Use Local dressings formulary and/or www.worldwidewounds

IF TREATMENT IS REQUIRED, FIRST REMOVE PRESSURE

04.

Wheelchairs

StaebeL Emperor Wheelchair Transit	106
StaebeL Emperor Wheelchair	106
StaebeL Liberation Lightweight Wheelchair	107
StaebeL Thrifty Wheelchair	107
StaebeL Sovereign Wheelchair	108
StaebeL Deluxe Alloy Wheelchair	109
StaebeL Deluxe Alloy WardWheelchair	110
Eclipse Bariatric Manual Wheelchair	110
StaebeL Gusto Heavy Duty Wheelchair	111
StaebeL Patient Transport Chair	112
Wheelchair Accessories	112



Staebel Emperor Wheelchair Transit

The Staebel Emperor is a comfortable, lightweight aluminium wheelchair, with a really helpful range of high quality features. With its sleek modern frame and lightweight foldable design, the wheelchair is perfect for travel. All features of this chair have been designed to make life easier for the user.

- Foldable backrest
- Lap belt
- Heel strap
- Height adjustable armrests and footplates
- Removable swing-away footrest
- Available in orange or blue.



Specifications:	
Product Code	40413
Seat Width	450 mm
Overall Width	670 mm
Overall Height	950 mm
Product Weight	12 kg
SWL	115 kg

Also Available:	
Product Code	40414
Seat Width	500 mm
Overall Width	720 mm
Overall Height	950 mm
Product Weight	15 kg
SWL	135 kg

Staebel Emperor Wheelchair

The Staebel Emperor is a comfortable, lightweight aluminium wheelchair, with a really helpful range of high quality features. With its sleek modern frame and lightweight foldable design, the wheelchair is perfect for travel. All features of this chair have been designed to make life easier for the user.

- Foldable backrest
- Lap belt
- Heel strap
- Height adjustable armrests and footplates
- Removable swing-away footrest
- Available in orange or blue



Specifications:	
Product Code	40411
Seat Width	450 mm
Overall Width	670 mm
Overall Height	950 mm
Product Weight	12 kg
SWL	115 kg

Also Available:	
Product Code	40412
Seat Width	500 mm
Overall Width	720 mm
Overall Height	950 mm
Product Weight	15 kg
SWL	135 kg



Staebel Liberation Lightweight Wheelchair

The Staebel Liberation Lightweight Wheelchair is a great travel companion. It is a sleek modern design that has been kept compact for easy transport and storage.

- Fully foldable backrest
- Dual brake
- Parking brake push handle
- Fold-up armrests
- Lightweight foldable swing-way footrests



Specifications:

Product Code	40400
Seat Width	450 mm
Overall Width	540 mm
Overall Height	920 mm
Product Weight	7 kg
SWL	100 kg



Staebel Thrifty Wheelchair

The Thrifty is a sturdy and reliable standard weight wheelchair, with removable armrests and leg rests, built for individuals, nursing homes and non-prescription requirements. This wheelchair meets Australian Standards for Manual Wheelchairs in Product Performance. All of these features combine to make the Thrifty of unbeatable value in standard wheelchairs.

- Swing away, nonremovable armrests
- Swing away, detachable legrests
- 24" rear wheels with urethane tyres (maintenance free)
- Australian Standard AS/NZS 3695.1:2011 approved



Specifications:

Product Code	4246
Seat Width	450 mm
Seat Depth	400 mm
Overall Width	650 mm
Overall Depth	1080 mm
Overall Height	880 mm
Fold Width	300 mm
Fold Length	790 mm
SWL	115 kg



NDIS
MASS



Staebel Sovereign Wheelchair

The Sovereign wheelchair by can be adjusted to suit almost any user. With much consideration in the design from the perspective of OTs, this chair is suited to the vast range of needs required by our clients. Rather than numerous script forms for various users, this chair meets the needs of almost any patient. All wheelchair parts suit any adjustments made - no additional parts required.

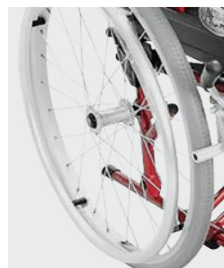
- Adjustable & removable legrests
- Flip back & removable armrests
- Adjustable armrest height & length
- Seat depth, height & angle adjustment
- Backrest height, tension & angle adjustment

See all Sovereign Accessories on page 112



Specifications:						
Product Code	4340C - 43	4340C - 46	4340C - 49	4352C - 52	4352C - 55	4352C - 58
Seat Width	430 mm	460 mm	490 mm	520 mm	550 mm	580 mm
Anterior Seat Height	420 mm	445 mm	450 mm	470 mm	475 mm	500 mm
Posterior Seat Height	410 mm	430 mm	450 mm	470 mm	490 mm	510 mm
Seat Depth	400 mm	425 mm	450 mm	475 mm	475 mm	500 mm
Backrest Height (2cm increments)	400 - 460 mm	400 - 460 mm	400 - 460 mm	400 - 460 mm	400 - 460 mm	400 - 460 mm
Backrest Angle (5° increments)	69° - 94°	69° - 94°	69° - 94°	69° - 94°	69° - 94°	69° - 94°
Lower Leg Length (2cm increments)	400 - 480 mm	400 - 480 mm	400 - 480 mm	400 - 480 mm	400 - 480 mm	400 - 480 mm
Armrest Height (2cm increments)	180 - 260 mm	180 - 260 mm	180 - 260 mm	180 - 260 mm	180 - 260 mm	180 - 260 mm
Armrest Length (4cm increments)	330 - 450 mm	330 - 450 mm	330 - 450 mm	330 - 450 mm	330 - 450 mm	330 - 450 mm
Footplate Width	105 - 155 mm	105 - 155 mm	105 - 155 mm	105 - 155 mm	105 - 155 mm	105 - 155 mm
Gross Weight	17 kg	17 kg	17 kg	21 kg	21 kg	
SWL	125 kg	125 kg	125 kg	170 kg	170 kg	170 kg

NDIS



Staebel Deluxe Alloy Wheelchair

The Deluxe Ultralight Alloy Wheelchairs have been designed to accurately meet the needs of our clients. They are an ultralight weight wheelchair with the benefit of many additional features. The improvement compared to previous designs (ranging from 25 - 30 kg gross weight) is considerable. Much input from clients and nurses has assisted in the development of these superior quality wheelchairs.

- Flip back armrests
- Swingaway, detachable legrests
- Fold down backrest
- Quick release rear wheels
- 24" rear wheels with urethane tyres
- Aluminium frame
- Padded nylon upholstery
- Complies with AS3695:2008



Specifications:


Product Code	4629
Seat Width	460 mm
Seat Depth	400 mm
Seat Height	480 mm
Backrest Height	430 mm
Armrest Height	230 mm
Overall Width	650 mm
Overall Depth	1060 mm
Overall Height	960 mm
Fold Width	320 mm
Net/Gross Weight	11 / 14 kg
SWL	125 kg

Also Available:

Product Code	45129	45629
Seat Width	510 mm	560 mm
Seat Depth	400 mm	400 mm
Seat Height	490 mm	490 mm
Backrest Height	430 mm	430 mm
Armrest Height	210 - 280 mm	210 - 280 mm
Overall Width	700 mm	750 mm
Overall Depth	1000 mm	1000 mm
Overall Height	960 mm	960 mm
Fold Width	300 mm	300 mm
Net/Gross Weight	14.5 / 17 kg	15 / 17.5 kg
SWL	150 kg	150 kg

NDIS

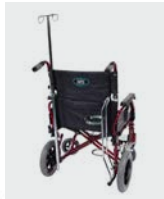NDIS



Deluxe Alloy Ward Wheelchair

The Deluxe Alloy Wheelchairs have been designed to accurately meet the needs of our clients. Much input from clients and nurses has assisted in the development of these superior quality wheelchairs.

- Flip back, height adjustable armrests
- Swingway, detachable legrests
- Height adjustable, flip back armrests this assists in side Transfers of the client
- Transit wheels with solid tyres
- Aluminium frame
- Padded nylon upholstery
- Complies with AS3695:2008



Specifications:

Product Code	45128
Seat Width	510 mm
Seat Depth	400 mm
Seat Height	490 mm
Backrest Height	430 mm
Armrest Height	210 - 280 mm
Overall Width	700 mm
Overall Depth	1000 mm
Overall Height	960 mm
Fold Width	350 mm
Net/Gross Weight	15 / 17 kg
SWL	150 kg

Eclipse Bariatric Manual Wheelchair

The Eclipse is a rigid chair designed by 7 Canadian engineers for real bariatric use, and ensures complete safety for the care staff attending the occupant.

Ultra-light in design and extended castor positioning for perfect weight balance, results in manual handling excellence for care staff when transporting the patient or resident through the facility.

- Removable armrest and legrest to assist in preferred transfer.
- Rigid frame for strength and reducing frame movement and fatigue.
- Overall width is only 19cm extra on the seat width, keeping it streamlined for going through doorways
- Solid seat pan to overcome slung seat sagging issues.
- Anti-tip wheels for complete user safety.
- Colour options for easy identification of bariatric equipment.



Specifications:

Product Code	4865
Seat Width	500 - 1200 mm
Seat Depth	400 - 550 mm
Seat Height	350 - 510 mm
Backrest angle	85 - 100°
Product Weight	20 kg
SWL	450 kg



Staelbel Gusto Heavy Duty Wheelchair

The Staelbel Gusto heavy duty wheelchair caters for the bariatric client. Featuring padded vinyl upholstery, heavy duty folding frame, 24 inch rear wheels and inset foot plates for client comfort, this wheelchair is ideal in the aged care and hospital environment. It comes standard with urethane tyres , eliminating the need to maintain tyre pressure. The swing away leg rests are removable and also adjust in length. The armrests will remove to enable side transFers if necessary. Footplates are alloy for longevity and ease of cleaning.

- Removable, swingaway, adjustable length legrests
- Removable armrests
- Padded vinyl upholstery
- 24" rear mag type wheels, urethane tyres



Specifications:	
Product Code	46108
Seat Width	610 mm
Seat Depth	460 mm
Seat Height	520 mm
Backrest Height	430 mm
Armrest Height	220 mm
Overall Width	850 mm
Overall Depth	1100 mm
Overall Height	930 mm
Fold Width	350 mm
Net/Gross Weight	26 / 30 kg
SWL	240 kg

Also Available:		
Product Code	46608	47108
Seat Width	660 mm	710 mm
Seat Depth	520 mm	550 mm
Seat Height	520 mm	520 mm
Backrest Height	480 mm	520 mm
Armrest Height	235 mm	220 mm
Overall Width	880 mm	930 mm
Overall Depth	1150 mm	1170 mm
Overall Height	1000 mm	990 mm
Fold Width	350 mm	350 mm
Net/Gross Weight	36 / 43 kg	36 / 46 kg
SWL	270 kg	290 kg



Staebel Patient Transport Chair

Staebel Patient Transport Chair is a handy, easy to use patient transport chair with flip-up armrests, generous 52cm side open seat and a retractable foot rest.

- 200kg SWL
- PVC document holder
- Adjustable braking system
- Silicone sealed seat for infection control
- HPL footplate
- UV protected Powdercoating
- Dual front castors for easy maneuvering.



Specifications:	
Product Code	4750
Width	660 mm
Length	870 mm
Height	1090 mm
Seat Height	620 mm
Seat Width	520 mm
Seat Depth	430 mm
Weight	28.5 kg
SWL	200 kg




Staebel Sovereign One Arm Drive Kit


The Staebel Sovereign One Arm Drive Kit assists one-handed wheelchair users in controlling the wheelchair. One arm drive systems for manual wheelchairs are designed for people who prefer the advantages offered by lightweight, folding manual wheelchairs.

Double hand rim one arm drives have both two hand rims on the same side of the wheelchair that control the rear wheels on both sides.


One arm drive system with double push rim (Must specify left or right at time of order)




Easy to Use



Intelligent



Extra Mobility



Additional Grip

Outer Hand Rim control

The outer (larger) hand rim controls the wheel on that side of the wheelchair where the kit is located.

Inner Hand Rim Control

The inner (smaller) hand rim controls the wheel on the opposite side of the wheel chair. The inner hand rim is connected to the wheel by way of a bar between the hand rim and the wheel.

Specifications:	
Product Code	43425



Stael Sovereign Arm Trough One Piece

The arm trough is of benefit for the wheelchair/ shower chair user who has upper limb weakness or paralysis. It attaches to the chair armrest via a clamp.

- One piece moulded polyurethane trough
- Supplied with hardware
- Can be mounted left or right side



Stael Sovereign Seat Belt

This seat belt is a simple device for safety and positioning. It has a push button buckle and is adjustable in length. It has metal eyelets to enable attachment via existing seat upholstery screws on the wheelchair.

- Adjustable in length
- Simply and easily fitted to wheelchairs



Stael Sovereign Drum Brake Kit

- Drum Brake Conversion Kit for the Sovereign Wheelchair
- Includes Brakes, Cables & Wheels

Specifications:

Product Code 4985

Specifications:

Product Code 4974

Specifications:

Product Code 43410



Stael Sovereign Stump Support

The stump support is used on a wheelchair in place of the regular legrest, and is useful for amputees with a below knee amputation.

- Adjustable pad height
- Left or right side configured – specify when ordering



Stael Sovereign Elevating Legrest

This seat belt is a simple device for safety and positioning of the wheel chair user. It features a push button buckle and has metal eyelets to enable attachment via existing seat upholstery screws on the wheelchair.

- Adjustable in length



Stael Sovereign Pelvic Strap

- Plastic side release buckle belt, 2-piece design
- Loop tie attachment around wheelchair/ shower chair frame
- Length adjustable 80–155 cm

Specifications:

Product Code 43417 - L

Product Code 43417 - R

Specifications:

Product Code 43416 - L

Product Code 43416 - R

Specifications:

Product Code 19636



Staebel Sovereign Power Assist

- Only suitable for Sovereign Wheelchair
- Aluminium frame with PU grey wheels.
- Includes controller, battery and charger.
- Two sizes available: 180W (43428) and 320W (43429).

Specifications:

Product Code	43428
--------------	-------



Staebel Anti - Tip Bars


The anti-tip bars are fitted to the wheelchair to prevent the wheelchair tipping backwards. It is particularly beneficial for amputees or active users.

- Attaches to wheelchair tilt bar
- Includes small plastic rollers

Product codes may vary depending on wheelchair model. Please call 1800 267 267 for the correct code. Please note all chairs are supplied with Anti-tip Bars. These are replacements.

Specifications:

Product Code	4982
--------------	------

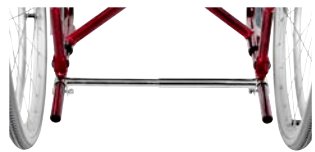

Staebel Anti Fold / Anti Theft Device

The anti-fold/ anti-theft device is used on a wheelchair to prevent the frame being folded, and to minimise the risk of theft. It is easily fitted to the frame between the rear wheels.

- Adjustable in width
- Larger sizes available

Product codes may vary depending on wheelchair model. Please call 1800 267 267 for the correct code.

Specifications:

Product Code	4978	4978L	4978XL
Sizing	400 - 490 mm	510 - 560 mm	610 - 710 mm



Staebel Elevating Legrest

The elevating legrest gives the wheelchair user the ability to support their leg in the elevated position. It attaches to the wheelchair in place of the existing legrest.

- Adjustable angle and length to suit the individual user
- Configured left or right - specify when ordering

Product codes may vary depending on wheelchair model. Please call 1800 267 267 for the correct code.

Specifications:

Left Leg	4980 - L
Right Leg	4980 - R



Staebel Handy Bar

Simple to use and install/remove Transfer Bar wheelchair/car transfer bar.

- Firm secure handle
- Can be used to break glass in emergency
- Inbuilt seat belt cutter
- Fits most vehicles - no modification needed
- Recommended by Doctors and Health Professionals
- Body point belts available. Contact the team for information

Specifications:

Product Code	4998
--------------	------



Staebel Calf Strap

This padded calf strap provides comfortable support and keeps the feet safely on the footplates. The Velcro straps simply and easily secure to the legrests of the wheelchair.

- Generous size (suitable for 430-490mm chair width)
- Calf section measures 460mm x 89mm

Specifications:

Product Code	4979
--------------	------



Staebel Oxygen Bottle Carrier

The oxygen bottle carrier is designed to be attached to the wheelchair to enable the oxygen dependent user to be transported in the wheelchair.

- Suitable for type C oxygen cylinders
- Removes easily for storage and transportation
- Slides onto existing wheelchair tilt bar

Product codes may vary depending on wheelchair model. Please call 1800 267 267 for the correct code.

Specifications:

Product Code 4983



Staebel IV Pole

The IV pole is attached to the wheelchair and allows the intravenous drip dependent client to be transported in a wheelchair.

- Easy height adjustment with triwheel
- Slides onto existing wheelchair tilt bar

Product codes may vary depending on wheelchair model. Please call 1800 267 267 for the correct code.

Specifications:

Product Code 45129 / 45629



Staebel Brake Lever Extension

The brake lever extension is used on the wheelchair and is ideal for those who find it difficult to operate standard brakes.

- Lengthens brake handle by approx 10cm
- Increases leverage for easier brake application
- Supplied as each

Specifications:

Product Code 4987



Staebel Seat Belt

This seat belt is a simple device for safety and positioning of the wheel chair user. It features a push button buckle and is adjustable in length. It has metal eyelets to enable attachment via existing seat upholstery screws on the user's wheelchair.

- Adjustable in length
- Simple and easy fit to most wheelchairs

Specifications:

Product Code 4974



Staebel Stump Support

The stump support is used on a wheelchair in place of the regular legrest, and is useful for amputees with a below knee amputation.

- Padded stump support
- Adjustable pad height
- Left or right side configured - specify when ordering

Product codes may vary depending on wheelchair model. Please call 1800 267 267 for the correct code.

Specifications:

Left Leg 4981 - L

Right Leg 4981 - R



Staebel Arm Trough

The arm trough is of benefit for the wheelchair/ shower chair user who has upper limb weakness or paralysis. It attaches to the chair armrest via a sturdy clamp.

- One piece moulded trough
- Supplied with hardware
- Can be mounted left or right side

Specifications:

Product Code 4985

05.

Seating

Regency Care Chair	117	Oasis Mobile Recliner	129
Regency Care Chair- Electric	118	Oasis Mobile Recliner with Poly Arm Caps	129
Sertain Care Chair Electric Hi-Lo	119	Comflex Windsor Chair	130
Sertain Bariatric Care Chair Electric Hi-Lo	120	Carson Recliner	130
Sertain Care Chair Electric Hi-Lo	121	Verity Mobile Recliner	131
Regency Result3 Recliner with Vienna Armrests - Electric	122	Shoalhaven Aircell Recliner	131
Regency Result3 Electric Recliner	122	Armchair Blues Comfort Pull Out Couch	132
Regency Bahama Comfort Chair	123	Staebel Bariatric Chair	132
Sertain Hi-Lo Day Chair with Euro Back	123	Embark Deluxe Day Chair	134
Regency & Sertain Accessories	124	Embark Superior Day Chair	134
Ashley Mobile Recliner Luxor Optima	126	Embark Utilite Day Chair	135
Ashley Mobile Recliner- Standard	126	Embark Deluxe Extracare Day Chair	135
Ashley Luxor Optima Lift & Recline Chair -One Motor	127	Embark Day Chair Legrest	136
Ashley Luxor Optima Lift and Recline Chair - Two Motor	127	Padded Footrest	136
Ashley Lift & Recline Chair - One Motor	128	Staebel Blood Donor Chair	137
Ashley Lift & Recline Chair - Two Motor	128	Blood Collection Chair	137



Regency Care Chair

The Regency Care Chair is designed to give the patient comfort by combining the original tilt system with an independent multi-positional back rest. Its many features take into account the needs of both carer and user of the chair. The Regency R2900 derives its status from a proven and successful track record across a broad spectrum of the industry.

- Reclining backrest – full recline to supine position
- Pressure management seating – gel cushion in seat, self-air adjusting back, profiled head cushion
- 2-position armrests
- Forward seat tilt for ease to stand, legrest angles under seat

- Elevating contoured legrest complete with fold out, extendible footplate
- Low skin shear 2-way stretch fabric on seating surfaces
- Drop down armrests and swingaway headwings – side transfer friendly
- Covers replaceable
- Patient lifter friendly castor base
- 5 year limited warranty
- 175mm castors, braking on rear, directional lock on front
- Quality Australian made construction



Tangerine



Purple



Green



Pink



Lagoon



Mustard Pattern



Koala Pattern



Arcadia Marina




Bouquet Navy



Specifications:

Product Code	5050G
Seat Width	490 mm
Seat Depth	520 mm
Seat Height (Tilt/Normal)	570 - 620 mm
Backrest Height	700 mm
Armrest Height	175 mm
Legrest Length	480 - 650 mm
Overall Width	660 mm
Overall Depth	1040 mm
Overall Height	1335 mm
SWL	180 kg

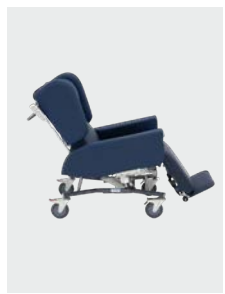NDIS



Regency Care Chair - Electric

The Wide Regency Care Chair follows the success of the R2900 Series model and is designed to give the large client unprecedented comfort by combining the tilt system with an independent multi-positional back rest and pressure area management seating.

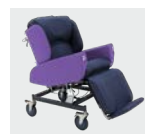
Its many features take into account the needs of both carer and user of the chair. The Regency derives its status from a proven and successful track record across a broad spectrum of the industry.

- Reclining backrest – full recline to supine position
- Gas strut assisted Tilting seat
- Pressure management seating – gel cushion in seat, self-air adjusting back, profiled adjustable head cushion

- 2-position armrests which drop down for sling insertion
- Forward seat tilt for ease to stand, legrest angles under seat
- Elevating contoured legrest complete with fold out, extendible footplate
- Low skin shear 2-way stretch fabric on seating surfaces
- Drop down armrests and swingaway headwings – side transfer friendly
- Covers are replaceable
- Patient lifter friendly castor base
- 150mm Twin wheel castors




Tangerine



Purple



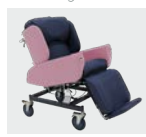Green



Pink



Lagoon



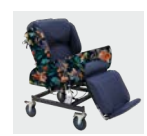
Mustard Pattern



Koala Pattern



Arcadia Marina



Bouquet Navy



Specifications:

Product Code	5050GE - SR
Seat Width	620 mm
Seat Depth	570 mm
Seat Height	550 mm
SWL	300 kg

Sizes Available:

Regency Care Chair - Electric - 520mm Wide	
Product Code	5050GE - S - SR
Regency Care Chair - Electric - 720mm Wide	
Product Code	5050GE - XL - SR



Sertain Care Chair Electric Hi-Lo

The Sertain® S4540 Series Pressure Care Seating System is a Team Player! As part of the Care Team the S4540 Series is there when you need it, to reposition the occupant, ease the transfer from bed to chair, increased access for wound dressing and procedures, maneuverable in small spaces – a relatively quiet contributor and achiever but a consistent performer as it adjusts to many of the challenges the Care Team faces – it doesn't take over, it joins in; a real Team needs a real Team Player!

- Electric hi-lo 540-940mm – low enough to stand from with ease, but high enough to achieve a comfortable working height for nursing staff
- Electric reclining backrest – full recline for supine position

- Electric tilting seat and legrest
- Pressure management seating – gel cushion in seat, air in backrest
- Fold out, extendable footplate
- Low skin shear 2-way stretch fabric on seating surfaces
- Removable two-position armrests and head wings
- Battery back up system
- Quality Australian made product
- Grey Frame
- Warranty – 5 year Structural Frame, 2 years Castors, 3 years Electric System



Tangerine



Purple



Green



Pink



Lagoon



Mustard Pattern



Koala Pattern



Arcadia Marina



Bouquet Navy



Specifications:

Product Code	5141
Seat Width	490 mm
Seat Depth	520 mm
Seat Height (Tilt/Normal)	540 mm / 940 mm
Backrest Height	700 mm
Armrest Height	150 mm
Legrest Length	440 mm
Overall Width	740 mm
Overall Depth	1080 mm
Overall Height	1200 mm – 1520 mm
SWL	160 kg

NDIS



Sertain Bariatric Care Chair Electric Hi-Lo

For the bariatric client, this model will support up to 300kg, with functionality, comfort and the needs of both user and carer in mind. Versatile in design and function, our Sertain HILO Care Chair is a popular choice for many wards including intensive care units. The HILO function of the chair is a nurse's dream, allowing you to adjust the height of the chair to a position comfortable to assist the patient.

- Electric 4 function (backrest, seat-tilt & leg-rest)
- Auto lay flat and auto chair function
- 2 position swing-away head wings
- 2 position standard removable armrests

- Deluxe leg rest including a foldaway footplate
- Multicush pressure management seating system
- Hi-Low range: 550 – 850mmH
- Quality Australian made product
- Available in seat widths 520mm, 620mm and 720mm



Tangerine



Purple



Green



Pink



Lagoon



Mustard Pattern



Koala Pattern



Arcadia Marina



Bouquet Navy



Specifications:

Product Code	5145
Seat Width	580 mm
Seat Depth	520 mm
Seat Height (Tilt/Normal)	550 / 950 mm
Backrest Height	700 mm
Armrest Height	160 mm
Legrest Length	440 mm
Overall Width	790 mm
Overall Depth	1080 mm
Overall Height	1335 mm
SWL	300 kg

NDIS



Sertain Care Chair Electric Hi-Lo

The Sertain® Care Chair Series is there when you need it, to reposition the occupant, ease the transfer from bed to chair, increased access for manual dressing and procedures, manoeuvrable in small spaces.

- Electric hi-lo ; low enough to stand from with ease, but high enough for a comfortable working height for carers
- Electric reclining backrest – full recline for supine position
- Electric tilting seat and elevating legrest
- Handset operated with preset auto chair and auto layflat functions
- Pressure management seating – gel cushion in seat, air in backrest in seat, air in backrest
- Fold out, extendible footplate

- Low skin shear 2-way stretch fabric on seating surfaces
- Removable two-position armrests
- Swingaway three position headwings
- Adjustable head rest cushion
- 150mm castors, 3 x braking, 1 x directional lock
- Battery back up system
- Quality Australian made product
- Available in seat widths 520mm, 620mm and 720mm



Tangerine



Purple



Green



Pink



Lagoon



Mustard Pattern



Koala Pattern



Arcadia Marina



Bouquet Navy

Specifications:

Product Code	5142
Seat Width	490 mm
Seat Depth	520 mm
Seat Height (Tilt/Normal)	540 / 940 mm
Backrest Height	700 mm
Armrest Height	150 mm
Legrest Length	440 mm
Overall Width	740 mm
Overall Depth	1080 mm
Overall Height	1200 - 1520 mm
SWL	160 kg

Sizes Available:

Sertain Bariatric Hilo Care Chair 450kg SWL

Product Code	5150
--------------	------

NDIS

NDIS



Regency Result3 Recliner with Vienna Armrests - Electric

Our Regency Result3 recliner chair is the perfect mix of style, comfort and function. With an electric Euro style backrest, elevating legrest and an auto layflat function, it makes delivering care that bit easier. Created with infection control and easy-clean in mind, this chair is a great addition to a range of healthcare settings.

- Electric 'tilt in space' seat and legrest action
- Electric Euro style reclining backrest and elevating legrest
- Incorporating pressure management Multicush gel seating and air cushion backrest system
- Auto lay flat function, auto chair function, M3 customer set memory position
- Including push handle
- 4x 100mm swivel braking castors
- 200kg SWL
- Manufactured in Australia
- 5 year limited warranty

Regency Result3 Electric Recliner

Packed full of features and with a 200kg weight capacity, this sleek recliner boasts a modern design appealing to wards, aged care facilities and clinics. The fixed lounge-style armrest has been uniquely designed to fit right up to most dining tables to increase social interaction.

- Electric 'tilt-in-space' seat and leg-rest action
- Electric Euro-style reclining backrest
- Elevating standard leg-rest
- Incorporating pressure management Multicush gel seating and air cushion backrest system
- Fixed lounge-style armrest
- Auto lay-flat function, auto chair function, M3 customer set memory position
- 4 x 100 swivel braking castors
- 5 year warranty
- Australian made



Specifications:	
Product Code	5027
Seat Width	550 mm
Seat Depth	520 mm
Seat Height	550 mm
Backrest Height	730 mm
Overall Width	810 mm
Overall Height	1280 mm
SWL	200 kg



Specifications:	
Product Code	5025GH
Seat Width	550 mm
Seat Depth	520 mm
Seat Height	550 mm
SWL	200 kg

NDIS



Regency Bahama Comfort Chair

Our Complex Care Chair is a popular choice for many healthcare environments. It's available in gas or electric formats, easy to maneuver and simple to clean for infection control. It's been carefully designed to keep your patients or residents social at meal times, fitting in closely to many dining tables.

- Gas Assisted design
- 'Tilt in space' seat and leg-rest action
- Euro style reclining backrest and Fixed lounge style armrest
- Elevating standard leg-rest
- Incorporates pressure management Multicush gel seating with air cushion backrest system
- Includes push handle
- Multi stretch material
- Available in a huge range of upholstery colours
- Manufactured in Australia



Specifications:

Product Code	5022GH
Seat Width	550 mm
Seat Depth	570 mm
Seat Height	550 mm
SWL	180 kg

Sizes Available:

Regency Bahama Comfort Chair with Lateral Supports	
Product Code	5023GH

NDIS



Sertain Hi-Lo Day Chair with Euro Back

Boasting a sleek Euro style backrest to add an extra touch of class to your ward, this Sertain HILO Care Chair is a popular choice for many hospitals wards. Available in a huge range of colours to compliment your ward, it promises an aesthetic touch without compromising on functionality.

- Electric 4 function (seat-tilt, backrest, Hi-Lo seat & legrest)
- Auto lay flat and auto chair function
- No head-wings – Euro-style backrest
- 160kg SWL (300kg SWL option also available)
- Standard removable armrest with no meal tray provision
- Standard push handle
- Hi-Low range: 550-850mmH
- Australian made

*Available in other colours and styles



Specifications:

Product Code	5026GH
Seat Width	520 mm
Seat Depth	520 mm
Seat Height	550 - 850 mm
SWL	160 kg



Regency Widening Kits

The widening kit is used when the user would benefit from a little more space inside the armrests of the Regency chair. It is a simple kit that is fitted to position the armrests further apart by 60mm.

- Enables the larger client to feel comfortable
- Often used in conjunction with the 5004W (wings only) widening kit for side wings/floor for standing transfers

Meal Tray

A popular option for the Regency & Sertain Chairs chair is the standard tray. Ideal for meals and activities. Simple locking device secures tray in various positions

- Standard flat table top
- Cut-out means tray can be positioned close to user
- T-Mould Edge

IV Pole

The IV Pole option allows the intravenous drip dependent client to stay seated in the comfort of their chair.

- Easy height adjustment with triwheel
- The relevant bracket to attach pole needs to be purchased as well

Required with purchase of 5009 IV Pole, to attach pole to Regency Chair.

Specifications:	
Product Code	5004 / 5004W

Specifications:	
Meal Tray to Suit 5050G	5007
Also Available:	
Meal Tray To Suit 5090	5007RR
Meal Tray To Suit 5030	5007N
Meal Tray To Suit 5120	5007S
Meal Tray To Suit 5131	5007S4
Meal Tray to Suit 5135	5007LG

Specifications:	
Regency IV Pole	
Product Code	5009
Sertain IV Pole	
Product Code	5109
Also Available:	
Product Code	5009B3 IV Pole Bracket



Sertain Oxygen Bottle Carrier

Oxygen Bottle Carrier, including bracket for mounting.



Sertain Layflat Bar Option

Supports chair at lieflat position, during CPR emergency procedures.

Specifications:

Product Code	5108
--------------	------



Sertain Procedure Arm With Tray

Procedure arm accessory allows for extra support & can be folded down when not in use.

Specifications:

Product Code	5103
--------------	------



Sertain Widening Kit Upgrade

Widening Kit Upgrade, at time of manufacture, increases available seat width by 50mm

Specifications:

Product Code	Right Hand	5105 - R
Product Code	Left Hand	5105 - L

Specifications:

Product Code	5106
--------------	------



NDIS



NDIS

Ashley Mobile Recliner
Luxor Optima

The Ashley Luxor range of chairs have been specially developed for those with pressure care needs or those needing increased comfort when seated. This is achieved with layered pressure reducing foam and stretch PU upholstered on the seating surfaces. All the favourite features of the standard Ashley are still included – a waterfall style backrest, sculptured armrests, side pockets and padded leg chaise. Cascading roll back with isolated head rest and lumbar roll to support both these areas

- Upholstered legrest with calf support
- Manual operated recliner action with lever
- Castors and push handle for mobility
- No sag springs in seat and lumbar for long lasting support
- Quality Australian made
- 12 month warranty on fabric, 24 month warranty on frame
- Adjustable head rest cushion
- 150 mm castors, 3 x braking, 1 x directional lock

Ashley Mobile Recliner - Standard

Economically priced, the Ashley features a waterfall style backrest, sculptured armrests, side pockets and padded leg chaise. Australian made. Cascading roll back with isolated head rest and lumbar roll to support both these areas.



Specifications:			
Product Code	5519MI	5519	5519MX
Size	Mini	Medium	Maxi
Seat Width	480 mm	480 mm	530 mm
Seat Depth	470 mm	500 mm	520 mm
Seat Height	460 mm	480 mm	500 mm
Backrest Height	600 mm	640 mm	680 mm
Overall Width	780 mm	780 mm	830 mm
Overall Height	1020 mm	1050 mm	1150 mm
SWL	120 kg	120 kg	120 kg

Specifications:	
Product Code	5521
Size	Medium
Seat Width	480 mm
Seat Depth	500 mm
Seat Height	480 mm
Backrest Height	640 mm
Overall Width	780 mm
Overall Height	1050 mm
SWL	120 kg



NDIS



NDIS

Ashley Luxor Optima Lift & Recline Chair -One Motor

The Ashley Luxor Optima range of chairs have been specially developed for those with pressure care needs or those needing increased comfort when seated. This is achieved with layered pressure reducing foam and a specialised 2-way stretch upholstery on the seating surfaces. All the favourite features of the standard Ashley are still included – a waterfall style backrest, sculptured armrests, side pockets and padded leg chaise. This model which features a single motor reclines to enable the occupant to rest and will also raise to the standing position.

- Upholstered legrest with calf support
- One motor wallsaver action – chair can be positioned close to within 100mm from wall and still recline
- Backrest reclines to half lying down position
- Chair raises to assist standing transfer
- No sag springs in seat and lumbar for long lasting support
- Tasmanian Oak hardwood frame for durability
- High quality Dunlop Enduro foams used
- Battery back up capability
- Small castors enable chair to be wheeled (not patient transport)
- 2 year electronic warranty, 5 year frame warranty

* One or two motor available



Specifications:

Product Code	5522MI	5522	5522MX
Size	Mini	Medium	Maxi
Seat Width	480 mm	480 mm	530 mm
Seat Depth	490 mm	500 mm	520 mm
Seat Height	470 mm	500 mm	530 mm
Backrest Height	600 mm	640 mm	680 mm
Overall Width	780 mm	780 mm	830 mm
Overall Height	1020 mm	1070 mm	1170 mm
SWL	120 kg	120 kg	120 kg

Ashley Luxor Optima Lift & Recline Chair - Two Motor

The Ashley Luxor Optima range of chairs have been specially developed for those with pressure care needs or those needing increased comfort when seated. This is achieved with layered pressure reducing foam and stretch PU upholstered on the seating surfaces. All the favourite features of the standard Ashley are still included – a waterfall style backrest, sculptured armrests, side pockets and padded leg chaise. The Optima chairs can be ordered in a variety of sizes to suit the individual.

The two motor action enables independent control of backrest and legrest operation for ultimate flexibility of chair position.

- Upholstered legrest with calf support
- Backrest reclines to flat position
- Chair raises to assist standing transfer
- No sag springs in seat and lumbar for long lasting support
- High quality Dunlop Enduro foams used
- Battery back up capability
- Small castors enable chair to be wheeled (not patient transport)
- Quality Australian made product
- Battery back up system

Specifications:

Product Code	5532MI	5532	5532MX
Size	Mini	Medium	Maxi
Seat Width	480 mm	480 mm	530 mm
Seat Depth	490 mm	500 mm	520 mm
Seat Height	470 mm	500 mm	530 mm
Backrest Height	600 mm	640 mm	680 mm
Overall Width	780 mm	780 mm	830 mm
Overall Height	1020 mm	1070 mm	1170 mm
SWL	120 kg	120 kg	120 kg

NDIS

NDIS



Ashley Lift & Recline Chair-
One Motor

The Ashley range of chairs have long been the favourite choice with their comfortable seating and range of options. Economically priced and upholstered in standard fabric, the Ashley features a waterfall style backrest, sculptured armrests, side pockets and padded leg chaise. The Ashley chairs can also be ordered in a variety of sizes to suit the individual. This quality powered rise and recline chair is both comfortable and durable, and provides the user assistance with hassle free push button operation.

This model features a single motor allowing the occupant to recline and be raised to the standing position.

Available in either a fabric or vinyl upholstery in a pattern or plain. Contact our team to discuss! Vinyl available in Austex Studio Encore colour range.



Ashley Lift & Recline Chair -
Two Motor

The Ashley range of chairs have long been the favourite choice with their comfortable seating and range of options. Economically priced, the Ashley features a waterfall style backrest, sculptured armrests, side pockets and padded leg chaise. Australian made, The Ashley chairs can also be ordered in a variety of sizes to suit the individual. This quality powered rise and recline chair is both comfortable and durable, and provides the user assistance with hassle free push button operation. The two motor action enables independent control of backrest and legrest operation for ultimate flexibility of chair position.

Specifications:			
Product Code	5524MI	5524	5524MX
Size	Mini	Medium	Maxi
Seat Width	480 mm	480 mm	530 mm
Seat Depth	490 mm	500 mm	520 mm
Seat Height	470 mm	500 mm	530 mm
Backrest Height	600 mm	640 mm	680 mm
Overall Width	780 mm	780 mm	830 mm
Overall Height	1020 mm	1070 mm	1170 mm
SWL	120 kg	120 kg	120 kg

Specifications:			
Product Code	5525MI	5525	5525MX
Size	Mini	Medium	Maxi
Seat Width	480 mm	480 mm	530 mm
Seat Depth	490 mm	500 mm	520 mm
Seat Height	470 mm	500 mm	530 mm
Backrest Height	600 mm	640 mm	680 mm
Overall Width	780 mm	780 mm	830 mm
Overall Height	1020 mm	1070 mm	1170 mm
SWL	120 kg	120 kg	120 kg

NDIS



Oasis Mobile Recliner

Thoughtfully designed and engineered for use in rigorous 24/7 healthcare environments. Generally padded all over with a gas cylinder mechanism allowing infinite recliner positioning.

- Ergonomically placed handles enable easy patient transport
- Polyurethane arm caps (black) and durable vinyl upholstery
- Heavy duty design for extra strength and product longevity
- Arm cap designed to assist client ingress and egress
- Side lever permits both patient and caregiver to operate leg extension
- Trendelenburg positioning function included
- Central castor locking system available
- Trendelenburg is easily activated from either side of the recliner



NDIS



Oasis Mobile Recliner with Poly Arm Caps

Thoughtfully designed and engineered for use in rigorous 24/7 healthcare environments. Generally padded all over with a gas cylinder mechanism allowing infinite recliner positioning.

- Ergonomically placed handles enable easy patient transport
- Polyurethane arm caps (black) and durable vinyl upholstery (Available in a huge range of finishes)
- Heavy duty design throughout for extra strength and product longevity
- Arm cap designed to assist client ingress and egress
- Infinite recline position
- Side lever permits both patient and caregiver to operate leg extension



Specifications:

Product Code	5824
Width	902 mm
Depth	914 mm
Reclined Depth	1880 mm
Height	1219 mm
Seat Height	521 mm
Seat Width	584 mm
Arm Height	692 mm
SWL	225 kg

Specifications:

Product Code	5814
Width	902 mm
Depth	914 mm
Reclined Depth	1880 mm
Height	1219 mm
Seat Height	520 mm
Seat Width	584 mm
Arm Height	692 mm
SWL	225 kg



NDIS



NDIS

Comflex Windsor Chair

This comfortable recliner is packed with features for both the patient and carer. The framework has been specifically engineered around stability for anti-tipping and it has a one piece backrest to keep cleaning easy. Choose to add a mobile base and slide-out footplate to increase functionality.

- Solid Tasmanian Oak frame
- Available as electric chair or manual (flush trigger pull operation)
- Unique design to limit tipping
- One piece backrest for easy clean & infection control
- Padded leg-rest chaise
- Available in a huge range of Healthcare grade upholstery colours
- Australian made



Carson Recliner

Carson Recliner provides a clean, modern aesthetic and is thoughtfully designed for healthcare environments, with a high level of comfort and ease-of-use for patients and caregivers.

- Available in 23" and 26" seat widths
- Field-replaceable upholstery covers
- Back construction features 14-gauge inner steel structure
- Lever positioning permits both caregiver and patient to operate footrest
- Ergonomically-positioned handles for easy transport
- Trendelenburg and Central Locking included
- Available in a huge range of finishes



Specifications:

Product Code	5850RB
Seat Width	580 mm
Seat Depth	660 mm
Seat Height	520 mm
Backrest Height	700 mm
Overall Width	740 mm
Overall Depth	910 mm
Overall Height	1220 mm
SWL	120 kg

Sizes Available:

Carson Recliner, 660 mm Seat Width	
Product Code	5851RB



NDIS



NDIS

Verity Mobile Recliner - Three Position Recliner

The Verity Recliner offers healthcare-grade durability and exudes an inviting, homelike aesthetic that brings residential warmth to clinical environments.

- 4 locking casters
- 3 position recline with easy to operate mechanism
- Ergonomically-positioned push handles for patient transport
- Stitching method reduces seams for increased infection resistance
- CAL TB133 compliance
- Moisture Barrier
- Optional arm caps available in polyurethane
- Contrasting fabrics are obtainable

Shoalhaven Aircell Recliner

Following on the success of the Shoalhaven comes the Shoalhaven Aircell. This has been engineered to ensure maximum comfort for those confined to a chair, whilst allowing the carer to operate the chair with maximum safety and care. Soft, stretch fabric, quality pressure sensitive foam and a shaped backrest giving lumbar and lateral support provide ultimate comfort.

- Non maintenance, air cell cushioning combined with slow release memory foam, for excellent pressure care.
- Adjustable / removable wings
- Compact postural design with tilt in space recline action, and high legrest elevation for reduced fluid retention
- Resilient, multiple stretch, water proof fabric
- Low seat height for fallout safety, and stand assist transfers
- Hospital quality 100mm locking castors
- Manufactured in Australia



Specifications:

Product Code	5266
Width	740 mm
Depth	740 mm
Height	1120 mm
Seat Height	480 mm
Seat Width	584.2 mm
Arm Height	660 mm
SWL	225 kg

Specifications:

Product Code	5230
Seat Width	520 mm
Seat Depth	500 mm
Seat Height	500 mm
Backrest Height	750 mm
Arm Height	220 mm
Legrest Length	400 mm
Overall Width	780 mm
Overall Depth	1000 mm
Overall Height	1170 mm
SWL	130 kg



Armchair Blues Comfort Pull Out Couch

Armchair Blues Comfort pull-out couch/sofa bed can be used as an extra bed in patients rooms in hospitals, hospice, etc. The simple unfolding mechanism made of metal construction ensures stability when lying on the armchair. Upholstered in vinyl leatherette welded, to prevent the water getting inside the foam.

- Mobile with four castors, front two are lockable
- Choose from 10 wooden decors and 14 vinyl colours! Contact the team to find out more information



Specifications:	
Product Code	78030
Chair Width	865 mm
Chair Depth	780 mm
Chair Height	990 mm
Bed Length	1950 mm
Mattress Platform Height	450 mm
SWL	150 kg



Staebel Bariatric Chair

Specially developed for Acute Care/Aged Care residence use, with a strong frame and high weight capacity, vinyl upholstery and tilt and move wheels to cater for the bariatric patient.

- Adjustable height
- High Contoured backrest incorporating lumbar support
- Contoured seat
- Padded drop down, removable arms for side transfer
- Adjustable seat depth
- Standard with tilt and move wheels for easy manoeuvrability



Specifications:	
Product Code	5966
Seat Width	780 mm
Seat Depth	600 mm
Seat Height	460 - 540 mm
Backrest Height	675 mm
Arm Height	170 mm
Overall Width	900 mm
Overall Height	1000 - 1085 mm
Width Between Arms	790 mm
SWL	425 kg

The Embark Range

Specially developed for use in various healthcare environments.



Choose your Colour!

Emerald

Periwinkle

Smoke

Black

Custom Colours available upon request. Contact the team to see custom colours!



Embank Deluxe Day Chair

An exceptional rehabilitation chair, the Embark Deluxe is available in sturdy quality vinyls for ease of cleaning. The adjustable seat height, seat depth and angle ensures Embark Deluxe can adapt to the widest range of people, making this chair ideal for hospital, home and nursing home. The padded seat and back rest have been designed to enhance the posture, whilst providing maximum comfort.

- High backrest incorporating lumbar and lateral support
- Padded armrests
- Comfortable vinyl upholstered seat and back
- Adjustable seat height and depth
- Powdercoated steel frame
- Standard with tilt and move wheels for easy manoeuvrability
- Available in a range of different colours. Contact the tem for more information



Specifications:	
Product Code	5930
Seat Width	500 mm
Seat Depth	450 - 530 mm
Seat Height	450 - 600 mm
Backrest Height	790 mm
Arm Height	230 mm
Overall Width	700 mm
Overall Depth	730 - 800 mm
Overall Height	1200 - 1350 mm
SWL	200 kg

Also Available:	
Embank Day Chair Legrest	
Product Code	5958



Embank Superior Day Chair

The Embark Superior Day Chair was designed and developed by a team of clinicians. It has a strong frame and a high weight capacity of 200 kg. The chair has low resilient foam and 2-way stretch upholstery on seating surfaces. For increased longevity in demanding environments it has 6cm tilt and move wheels and other small improvements.

There have been instances in the past, whereby patients, staff and visitors have sustained cuts and injuries from square tube frames. As a result, the Embark was not only designed for comfort, but to avoid these injuries.

- Low back option also available.
- Comfortable upholstered seat and back
- Padded armrests
- Adjustable seat height and depth
- Powder-coated steel frame
- High backrest incorporating lumbar and lateral support



Specifications:	
Product Code	5933
Seat Width	500 mm
Seat Depth	450 - 530 mm
Seat Height	450 - 600 mm
Backrest Height	790 mm
Arm Height	230 mm
Overall Width	700 mm
Overall Depth	730 - 800 mm
Overall Height	1200 - 1350 mm
SWL	200 kg

NDIS



NDIS



Embark Utilite Day Chair

The Utilite is a multi purpose adjustable seat height and angle chair. Ideal for use in an area where adjustable seating is required, for example, lounge, dining room, bedroom, activity room, outdoor area etc. This low back chair provides comfort, stability and practicality at an affordable price.

- Shallow curved low backrest
- Padded armrests
- Comfortable, durable upholstered seat, back and armrests
- Powdercoated steel frame
- Available in a range of different colours! Contact the tem for more information

Embark Deluxe Extracare Day Chair

The Extra Care Embark Deluxe has an impressive 250kg weight capacity and extra wide 60cm width, making it ideal for the heavier client. It still provides all the benefits of the Embark Deluxe, e.g. adjustable seat depth and height, and comes standard with tilt and move wheels.

- High backrest incorporating lumbar and lateral support
- Padded armrests
- Comfortable vinyl upholstered seat and back
- Adjustable seat height and depth
- Powder-coated steel frame with reinforcement for increased load



Specifications:

Product Code	5910
Seat Width	460 mm
Seat Depth	460 mm
Seat Height	485 - 640 mm
Backrest Height	455 mm
Arm Height	230 mm
Overall Width	610 mm
Overall Depth	550 mm
Overall Height	900 - 1075 mm
SWL	120 kg

Also Available:

Embark Day Chair Legrest	
Product Code	5958



Specifications:

Product Code	5931BN3
Seat Width	600 mm
Seat Depth	450 - 530 mm
Seat Height	470 - 620 mm
Backrest Height	765 mm
Arm Height	220 mm
Overall Width	800 mm
Overall Depth	740 - 820 mm
Overall Height	1180 - 1330 mm
SWL	250 kg

NDIS

NDIS



Embark Day Chair Legrest

This vinyl covered legrest provides comfort for resting legs, with its generously padded top set on a 10 degree angle. Height and angle adjustable using simple push button system. Fitted with non-slip rubber feet.

- Adjustable height 435-585mm
- Padded top measures 405mm x 510mm
- Sturdy powdercoated steel frame
- Complements the Embark chair range
- Vinyl upholstery



Padded Footrest

This is a padded footrest to provide support for the feet and assists in the correct posture when seated. The padded upholstered top provides user comfort.



Specifications:	
Product Code	5958
Height	435 - 585 mm
Width	550 mm
Lenght	430 mm
SWL	200 kg

Specifications:	
Product Code	5955
Width	240 mm
Length	400 mm
Height	125 mm

Also Available:	
Product Code	5956
Width	240 mm
Length	400 mm
Height	200 mm

NDIS



Stael Blood Donor Chair

An ergonomically designed chair, which enables a relaxed and comfortable sitting or lying during the procedure. The armrests are movable into all directions. Armrests suitable for blood donation procedures. All parts have quality upholstery with rounded corners with antibacterial artificial leather-skai, special for medical purposes.

- Patient maximum weight: 170 kg;
- Electrically adjustable seat section height
- Seamless padding 5 cm.
- Swivel double wheel casters with diameter of 75mm with brake
- Metal frame powder-coated with antibacterial paint
- The base is made of powder coated metal elements
- Backrest inclination adjustable via gas spring from -10° to +80°
- Seat section fixed 2,5° seat section inclination
- Inclination adjustable via gas spring from -90° to +10° footrest
- Upholstery colour by agreement. Contact our team for more colour options.



Specifications:

Product Code	5201
Seat Width	600 mm
Seat Length	550 mm
Overall Width	600 mm
Overall Length	1000 mm
Overall Length (Lying Position)	1850 mm
Seat Height Adjustment	550 - 750 mm
Backrest (Width x Length)	600 x 850 mm
Footrest (Width x Length)	450 x 400 mm
SWL	170 kg

NDIS



Blood Collection Chair

The Blood Collection chair has been designed in conjunction with leading pathology companies to ensure that sample collection is safe for both staff and patient, whilst affording the patient a reasonable degree of comfort. The sloping armrests are positioned to present the arm at the optimum angle and the step is fitted with a non slip surface.

- Elevated Seat height for ease of access
- Wide padded armrests
- Non-slip step
- Reinforced armrest support
- Range of hospital grade vinyl coverings
- A range of colours available, contact our team to discuss.



Specifications:

Product Code	5202G
Seat Width	550 mm
Overall Width	920 mm
Arm Rest Width	185 mm
Length	650 mm
Overall Length	840 mm
SWL	175 kg

06.

Mobility
Solutions

Walking Stick - Aluminium T-Handle	139	Galaxy Seat Walker	142
Walking Stick - Swan Neck	139	Staebel Server HD Seat Walker	142
Walking Stick - T-Bar Handle	139	Staebel Navigator Forearm Walker	142
Crutches - Aluminium Height Adjustable	140	byACRE Carbon Ultralight Walker	143
Staebel Anatomic Grip Elbow Crutches	140	byACRE Carbon Overland Fibre Walker	143
Lightweight Walking Frame	140	Rollz Motion 2 in 1 Wheelchair Walker	144
Staebel Cruiser Folding Walker	141	Stand Tall Bariatric Forearm Walker	145
Staebel Ligerio Seat Walker	141	Viking Forearm Walker	145
Staebel Petite Seat Walker	141	Staebel Stimulator Rehab Forearm Walker	145

NDIS

NDIS

DVA

NDIS



Walking Stick - Aluminium T-Handle

One of the best selling walking sticks favoured for its sheer versatility. The shaft is made of lightweight yet durable aluminium which folds into four sections and can be stored neatly in the plastic wallet supplied with it.



Specifications:

Product Code	6241
Height (Med)	790 - 840 mm
SWL	100 kg

Also Available:

Product Code	6240
Height (Lrg)	840 - 940 mm
SWL	100 kg



Walking Stick - Swan Neck

This adjustable height aluminium walking stick assists the user's mobility. The "swan neck" shaped handle helps to transmit the weight of the user vertically through the body of the stick. It provides exceptional support and comfort, eliminating slipping and encouraging correct posture. It also includes an anti rattle silencer and adjusts in 25mm increments.



Specifications:

Product Code	622J0
Height	730 - 990 mm
SWL	100 kg



Walking Stick - T-Bar Handle

This adjustable height aluminium walking stick assists the user's mobility. The T-shaped handle is strong and lightweight. It also includes an anti rattle silencer and adjusts in 25mm increments.

- Attractive timber look handle



Specifications:

Product Code	6230
Height	760 - 990 mm
SWL	100 kg



Crutches - Aluminium
Height Adjustable

The Aluminium Crutches offer push button height adjustment and are manufactured from a lightweight anodised aluminium. The solid one piece crutch frame gives strength and incorporates a contoured non rotating hand grip. A new anti rattle system makes these crutches silent.

- 3 sizes available (Codes 6011, 6012 and 6013)
- Supplied as a pair



Specifications:		
Product Code	6011J	
Crutch Height	930 - 1130 mm	
Patient Height	1372 - 1575 mm	
SWL	100 kg	

Also Available:		
Product Code	6012J	6013J
Crutch Height	1140 - 1345 mm	1335 - 1540 mm
Patient Height	1575 - 1778mm	1778 - 1981mm
SWL	100 kg	100 kg



Staebel Anatomic Grip
Elbow Crutches

The Anatomic Elbow Crutches feature contoured moulded handles that evenly distribute pressure across the palm of the hand, contributing to the users ease of operation and comfort. Added to this is the ability to adjust the crutch length from cuff to hand grip and also from hand grip to ground, enabling the user to adjust the crutch to their individual requirements. Manufactured from a larger tubing giving increased stability and weight capacity.

- Supplied as a pair
- 2 sizes available. (6151J, 6152J)



Specifications:		
Product Code	6151J	
Cuff to Hand Grip	930 - 1130 mm	
Hand Grip Height	1372 - 1575 mm	
SWL	180 kg	

Also Available:		
Product Code	6152J	
Cuff to Hand Grip	1335 - 1540 mm	
Hand Grip Height	1777 - 1981 mm	
SWL	180 kg	



Lightweight
Walking Frame

This folding lightweight walking frame provides assistance for those who have difficulty walking.

It folds in by pushing a simple release button. When folded this frame will fit into the smallest of cars and is ideal where storage is at a minimum. Other features are soft grip handles, simple push button height adjustment and anti rattle silencers.

- Height adjustable to suit the individuals requirements
- Lightweight aluminium construction
- Folds flat for easy storage and transportation
- Soft handles for comfort and grip



Specifications:		
Product Code	6380	
Height	690 - 870 mm	
SWL	120 kg	

Also Available:		
Product Code	6381	
Height	830 - 1010 mm	
SWL	120 kg	

DVA

NDIS

DVA

NDIS

NDIS

MASS



Staebel Cruiser Folding Walker

The Cruiser is a popular seat walker manufactured from durable lightweight aluminium. 6" castors make this an ideal outdoor walker. But is also compact enough for use indoors. Features include a padded seat, bag under seat, lockable brakes and padded, folding backrest for comfort while seated. The adjustable height handles are fitted with orthopaedic hand grips, therefore providing the user with excellent support. The walker also folds flat for transport

- Height adjustable ergonomic shaped handgrips
- Lightweight aluminium construction
- Lockable handbrakes
- 152mm castors



Staebel Ligero Seat Walker

The Ligero is a popular seat walker manufactured from durable lightweight aluminium. The larger castors make this an ideal outdoor walker, however is also compact enough for use indoors. Features include a padded seat and padded folding backrest for comfort while seated. The adjustable height handles are fitted with Orthopaedic hand grips therefore providing the user with excellent support. The walker also folds for storage and transport. Two seat heights available.

- Seat height of 55cm
- Height adjustable ergonomic shaped hand grips
- Lightweight aluminium construction
- Lockable handbrakes
- 200mm (8") castors



Staebel Petite Seat Walker

- Ideal for smaller people requiring mobility assistance
- Durable lightweight aluminium frame
- Folding for convenient transport and storage
- Height adjustable handles
- Ergonomically designed rubber
- Hand grips
- Fitted with rear braking and front swivel castors
- Locking brake mechanism with reflectors
- Padded seat
- Supplied with reinforced vinyl bag
- Removable backrest



Specifications:

Product Code	6730
Seat Width	440 mm
Seat Height	560 mm
Handle Height	820 - 930 mm
Walker Weight	6 kg
SWL	125 kg



Specifications:

Product Code	6865
Seat Width	440 mm
Seat Height	550 mm
Handle Height	790 - 910 mm
Walker Weight	7.5 kg
SWL	120 kg

Also Available:

Ligero Seat Walker - 510 mm	
Product Code	6869



Specifications:

Product Code	6771
Seat Width	380 mm
Seat Depth	300 mm
Seat Height	470 mm
Overall Width	560 mm
Overall Depth	670 mm
Overall Folded (WxDxH)	570 x 350 x 815 mm
Handle Height	740 - 850 mm
Walker Weight	5.5 kg
SWL	175 kg

NDIS

NDIS

NDIS

MASS

MASS



Galaxy Seat Walker

Constructed from heavy duty aluminium, this bariatric seat walker can cater for those up to 220kg. It features a removable padded backrest and ergonomically shaped handgrips for user comfort. The handbrakes are lockable for stability. A wire basket is included also.

- Extra wide seat width and reinforced frame to suit the bariatric user
- Folds for easy storage and transportation
- Height adjustable handles
- Handbrakes can be locked for added stability
- Supplied with removable backrest, basket and padded seat



Specifications:

Product Code	6736
Overall Width	710 mm
Width Between Handles	610 mm
Handle Height	770 - 920 mm
Walker Weight	11.1 kg
SWL	220 kg



Staebel Server HD Seat Walker

Server HD folds easily and stays that way with a simple locking system that doubles up as an ergonomically shaped carrying handle.

Frame and seat are tested and approved for maximum user weight of 200 kg. The brake mechanism is light yet positive. TPE covered, unbreakable wheels.

- The push handles are height adjustable with 10 positions, from 74 cm to 102 cm.
- User comfort is ensured with an optimal frame design and ergonomic push handles.



Specifications:

Product Code	6784
Width	690 mm
Depth	680 mm
Height	740 - 1020 mm
Seat Width	550 mm
Depth (When Folded)	680 mm
Height (When Folded)	860 mm
Width (When Folded)	225 mm
SWL	200 kg



Staebel Navigator Forearm Walker

The navigator is foldable, staying firmly closed making it is easy to transport. The frame and seat are tested and approved on users with a maximum weight of 150 kg. Light and reliable brakes act on the TPE covered unbreakable wheels. The optimal and rigid design of the frame and the fore arm support ensures the best possible user comfort (lightness, turning circle and steadiness).

- The Navigator is available in size large.
- The Navigator is easily height adjustable in the range from 900 - 1150 mm



Specifications:

Product Code	6782
Overall Width	550 mm
Overall Depth	610 mm
Height	740 - 1020 mm
Seat Height	620 mm
Handle Height	900 - 1150 mm
Walker Weight	10.8 kg
SWL	150 kg

NDIS

NDIS



byACRE Carbon Ultralight Walker

At just 4.8kg this carbon fibre walker / walking frame is the lightest in the world and uses space age technology and the best of Danish design.

- Soft inverted handles to support posture and comfort and reduce strain on fingers, hands and joints when steering
- Unbeatable in terms of the ease to fold, lift and grip when lifting it into the car or up stairs
- Award winning design providing a slick alternative to the traditional looking walking frame
- Brake cables hidden inside the frame for cleaner, less cluttered lines
- Shock absorbing frame and tyres made from a PU/rubber mix to reduce vibration
- Different colours available. Contact the team to see more!

byACRE Carbon Overland Fibre Walker

The byACRE Carbon Overland walker / walking frame is a sturdy, light weight and versatile multi-terrain walker

- Lightest off-road walker at 6.7 kg for easy lifting in y our vehicle
- Easy to fold, easy to grip and easy to lift
- Fully assembled, to make it hassle free to start exploring
- Extremely stable carbon fibre frame – 5 x stronger than steel
- Ultra-durable pneumatic rubber tires for better suspension
- Rigid front forks and sealed bearings
- Aircraft grade aluminium in gunmetal colour, for strength and lightness
- Different colours available. Contact the team to see more!



Specifications:	
Product Code	6884SR
Width	625 mm
Width when Folded	250 mm
Handle Height	790 - 910 mm
Width Between Handles	390 mm
Seat Width	420 mm
Seat Height	610 mm
Length	660 mm
Weight	4.8 kg
SWL	130 kg

Specifications:	
Product Code	6886DG
Handle Height	820 - 950 mm
Length	730 mm
Seat Width	470 mm
Seat Height	610 mm
Width	660 mm
Width when Folded	250 mm
Wheel Diameter	250 mm
Weight	6.7 kg
SWL	150 kg



Rollz Motion 2 in 1 Wheelchair Walker

The award winning Rollz Motion are one of a kind. Perfect for those who want to walk but are not sure how far they can go. With their large wheels for superior maneuverability and state of the art brakes, these are the ultimate in assisted walking. The Rollz Motion is an award-winning and beautifully designed rollator that converts simply into a wheelchair in seconds.

- Gives you confidence to go further
- Innovative and award-winning design
- Superior brakes suited to wet and dry conditions
- Folds into compact and transportable size
- Strong frame and wheels
- Allows you to stand and walk tall, instead of leaning forward
- Safe with superior manoeuvrability
- Adjustable footrests and handle heights
- Four colours available (Island Blue Matt Black, Dark Purple and Pebble White)

NEW

SAFE WORKING LOAD

725KG

2 in 1



Specifications:	
Product Code	688RMW
Seat Height	550 mm
Handle Height Range	870 - 980 mm
Seat Size Rollator	200 x 450 mm
Seat Size Wheelchair	420 x 450 mm
Overall Length	630 mm
Overall Width	310 mm
Width Folded	670 mm
Weight	11.6 kg
SWL	125 kg

NDIS



Stand Tall Bariatric Forearm Walker

Stand Tall is a multi-functional aid; it can be used for walking, mobilising from seated to standing, provides secure support for user and has an IV option.

Excellent rehabilitation support for bariatric patients. Provides a secure sit to stand process, then provides a stable support while enabling the patient to be mobile.

- Adjustable height to suit user requirements
- Easily manoeuvred and adjusted
- Simple construction allows easy cleaning and maintenance
- In-line braking system avoids skin tears



Specifications:

Product Code	6747
Overall Width	760 mm
Handle Height	1000 - 1300 mm
Walker Weight	710 mm
SWL	325 kg

NDIS



Viking Forearm Walker

The Viking Forearm Walker presents some unique features that particularly suit a rehab or post surgery environment. Its many adjustments also mean it can be adapted for a wide variety of patients.

- Adjustable angle forearm supports
- Adjustable angle hand grips
- Height adjustable
- Adjustable width base
- Large padded forearm pads for user comfort
- Strong aluminium frame
- Lockable rear castors



Specifications:

Product Code	6748
Overall Width	700 - 850 mm
Handle Height	880 - 1310 mm
Walker Weight	20 kg
SWL	200 kg

NDIS



Stael Rehab Forearm Support Walker

The best rehab program and the best rollator will guarantee you the best results – is this what you're looking for? Our Stimulator Rehab Forearm Walker has been specifically designed for the needs of the Aged and Acute Care/Hospital environment.

- Very lightweight, easy for the frail, weak patient to use
- High quality castors with super smooth rolling and minimal friction even on carpet, taking away any shuddering or wobbling
- Soft PU multi adjustable arm pads for any patient with any specific position requirement
- Multi adjustable hand grips, for the arthritic, ABI or stroke patients.
- Four swivelling castors for ultimate manoeuvrability
- Individual brakes on each castor, as well as hand brakes
- Modern streamlined design gives an up-to-date image to your rehab and physio area



Specifications:





Product Code	6755
Overall Width	700 mm
Height of Arm Support	905 - 1310 mm
Width Between Arms	310 - 470 mm
Overall Length	810 mm
SWL	200 kg

07.

Bedroom Solutions

Eleganza 1 Bed	148	Shelly Infant Incubator	168
Eleganza 2 Bed	151	Lilly Phototherapy Lamp	168
Eleganza 4 Bed	152	Tom 2 Cot	169
Eleganza 5 Bed	153	Nino Adolescent Bed	170
Multicare Bed	155	Staebel Maternity Bed with Adjustable Stirrups	171
Staebel Hospital Floorline Bed	157	Ave2 Birthing Bed	172
Essenza 300 Bed	158	Eleganza 2 With Junior Kit	174
Essenza 300 Lateral Tilt Bed	159	Staebel Deluxe Commode Bedside	175
Image 3 Bariatric XXL Bed	160	Staebel Standard Bedside Commode	175
Bariatric Tilt Hospital Bed	161	VT6 Mobile Bedside Cabinet - 2 Drawer	176
Sentida SC Bed	162	VT6 Mobile Bedside Cabinet - 3 Drawer	176
Staebel Floorline Bed	164	Staebel Deluxe Mobile Overbed Table	177
Luxa Homecare Bed	165	Staebel Basic Mobile Overbed Table	177
Bariatric Tilt Bed	166	Staebel Bed Rail Protector	178
Bassinets Trolley including Mattress & Tray	167	Blanket Support	178
Nataly Infant Radiant Warmer	167		

Build your Customisable Bed!

 Who is using the bed?	Infant		Child	
	Adult		Elderly	
 Choose your bed width	800 mm	850 mm	1000 mm	
	1050 mm	1200 mm	1350 mm	
	1500 mm	1800 mm	Companion Bed	
 Weight Limit (SWL)	80 kg	150 kg	200 kg	
	270 kg	320 kg	450 kg +	
 Lowest Bed Height Requirement	100 mm	230 mm	280 mm	
	380 mm	450 mm		

Let's stylise!

 Choose your bed style	Hospital Style		Homecare Style	
				
 Bed Frame Design	Colours	Stains	Laminate	
	Timber Head & Foot Boards	Upholstered Head & Foot Boards	Vinyl Head & Foot Boards	
 Accessories	Bed Functions	Castors	Controls / Hanset	
	Battery	Lighting	Rails / Support Rails	
	Self Help Pole	Clip on Extension	Handset Holders	



Could not find the information you were after?

Contact our team for a consultation for more bed options!



Eleganza 1 Bed

With full electrical adjustment, the Eleganza 1 bed is designed for the standard hospital ward. The bed has a simple construction, a wide range of features and functions and at a low cost, the Eleganza 1 provides customers with great value for money. Its sophisticated electric functions include the safety concept of the 'GO' button, protecting the patient against an injury caused by any unintentional movement of the bed; nurse control panel with pre-programmed cardiac chair, CPR and trendelenberg settings; ergonomic patient handset. The Eleganza 1 uses the Ergoframe, which enables the mattress platform to widen the space under the patients buttocks for optimum comfort and posture whilst eliminating the negative effects of compression, shear and friction forces which can lead to the development of pressure ulcers. Nosocomial infection prevention is achieved through joint-free construction and easy cleaning. Safety is paramount, with the Eleganza 1 bed meeting all of the requirements of the latest standard IEC 60601-2-52:2010.

- Electric auto regression of backrest and kneebend
- Electric forward and reverse trendelenberg
- CPR release
- Battery back up
- Nurse control panel as well as patient backlit handset
- Underbed lighting
- Cardiac protection
- Integrated bed extension 15cm
- Fold down 4-part side guards, soft drop
- 150mm anti static twin wheel central locking castors
- Corner buffers
- Standard moulded removable head and foot boards with safety locks
- Safe Working Load of 250kg



Specifications:	
Product Code	73014
External Dimensions	990 x 2180 mm
Bed Extension	150 mm
Mattress Platform	4 section electric adjustable
Mattress Surface	900 x 2000 x 140 mm
Electric Height Adjustment (125mm Castors)	370 - 735 mm
Electric Height Adjustment (150mm Castors)	395 - 768 mm
Electric Backrest Adjustment	0 - 70°
Electric Thighrest	0 - 39°
TR / ATR	15° - 15°
Ergoframe®	100 + 60 mm
Bed Weight	135 kg
Standard	IEC:60601-2-52
Cover Against Water and Dust	IP54
Clearance for Server (125mm Castors)	160 mm
Clearance for Server (150mm Castors)	1880 mm
SWL	250 kg

Control Options



Central Control Panel

The main central panel for the control and deactivation of the functions is located under the foot end and is integrated into the bed frame. This control unit contains the basic preset functions and the possibility of deactivating (locking) each function on the remote control unit. If the bed is equipped with this basic control panel, the nurse can use the remote control unit to operate the basic functions.



Supervisor Control Panel

Supervisor control panel has a coiled cable providing flexibility of movement and operation of the bed from an optimum position. For the nurse there is access to all of the standard functions + a useful pre-programmed function, the cardiac chair, available on the supervisor control panel only. The supervisor is an optional accessory and replaces the central control panel.



Hand Set

The ergonomic remote control unit allows the patient and nurse to control the basic positions of the bed (tilting of the backrest and thigh rest, adjustment of the mattress platform height) including autocontour, a combined function consisting of simultaneous positioning of the backrest and thigh rest. The control unit is equipped with an examination light. The backlit keyboard design would be of great assistance in the dark.



Patient Satellite

This control unit fixed on a flexible arm is always accessible and well within reach. It is designed particularly for elderly patients, so the symbols for controlling the functions are large and easily visible. The buttons are differentiated by contrasting colours and are of a raised design to suit patients with impaired vision or touch.



Side Rail Control Unit

Through extensive research it was found that controls located in the head side rails are more accessible, allowing the patient easier operation when adjusting the basic positions of the bed.



Plug & Play

In some cases, particularly when the patient is confused, it may be appropriate to remove the remote control unit. With the plug and play socket located on the accessory bar, the remote control unit can be easily removed.

Safe Siderails

The bed can be equipped with soft drop function of split plastic siderails or single $\frac{3}{4}$ siderails. In each case they comply with valid standards and all gaps on them are designed with regard to safety.



Infection Control



Nosocomial infection prevention

Assembly, dismantling and cleaning of the mattress platform is now much easier and more efficient due to the exceptional solution consisting of 8 selfcontained removable plastic components. Similar to the bed ends and side rails, these components are made using the joint-free technology of injection moulding. With a low weight (each part weighs only 100 grams) and dimensions, this solution is very popular. A figure pictogram can be used for quick identification of correct position of the bed.



Minimise Space for Infection

The latest technology is used during the manufacture of individual parts to minimise the space where unwanted infections may hide. In terms of easy cleaning, the bed is very safe and user friendly.



Easy Construction

From the beginning, in order to achieve easy cleaning, the Eleganza 1 bed was designed to maximise the simplification of the structure to include the individual parts and surfaces.



All Joints Covered

Cleaning of bolted joints can be very complicated, however, all joints on Eleganza 1 are covered, making cleaning much easier and more efficient.

Go Button



The safety concept of the GO button protects the patient against an injury caused by any unintentional movement of the bed. Use the GO button to activate the bed control units for a period of 3 minutes from the last depression. When the bed is not deliberately repositioned, the controls are automatically deactivated, which makes any unintentional activation impossible.

Ergoframe



More space for patients

When re-positioning the backrest and thigh rest, the Ergoframe function minimises the pressure and shear in the patient's pelvic region. This results in a significant reduction in the development of pressure ulcers. The effect is achieved by the automatic widening of the mattress platform by 16 cm.

Brake Alarm



Sound alarm for unbraked bed

Triggered when the bed is unplugged and unbraked. Available together with foot controls (pictured).



Eleganza 2 Bed

The LINET Eleganza 2 pushes the boundaries of bed design in terms of safety, quality of workmanship and innovation. The LINET Eleganza 2 is the first bed of its kind so has set the standard by defining the genuinely modern bed system. Containing sophisticated functions, the 'smart' features of the LINET Eleganza 2 are now accessible to the broad spectrum of healthcare environments.

The combination of practicality and beauty make the LINET Eleganza 2 a bed which will meet the needs and expectations of each hospital ward. The LINET Eleganza 2 makes demanding day-to-day patient care easier and thanks to its modern design, it shifts the perception of care provided to a far higher level.

- Integrated Mobilift innovation with support handles
- Highest possible level safety level to eliminate falls or injuries and to indicate safe bed positioning
- Attractive design with innovative column construction
- 'Smart' patient monitoring features: safety monitor, bed monitor and vital monitor
- Ease of transport with fifth castor Increased patient comfort and ease of cleaning
- Patient control of some functions
- Practical and ergonomic accessories available
- The Eleganza 2 can be ordered with the Junior
- Rail Kit, making it the perfect bed for safety with both adult and paediatric patients. For more information, please call 1800 267 267.



Specifications:

Product Code	73600
External Dimensions	995 x 2230 mm
Bed Extension	310 mm
Bed Height (Hi-Lo Adjustment)	395 – 775 mm
Internal Dimensions / (Compatible with Mattress Width)	2000 x 900 mm / (90 and 86 cm)
Max. Mattress Height	180 mm
Height of Split Plastic Side Rails	405 mm
Weight of Bed	140 kg
Clearance of Undercarriage	15 cm with 150 mm Tente castors
Max. Patient Weight	185 kg
ErgoFrame®	110 mm
Maximum Backrest angle	70°
Thigh/ Calf Rest angle	35° / 21°
TR/RTR Angle	16° / 16°
Lifting Pole Bearing Capacity	75 kg
SWL	250 kg

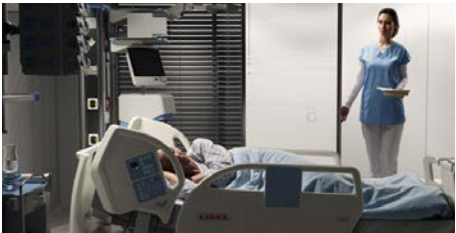


Eleganza 4 Bed

Eleganza 4 is an electrically positionable bed – designed especially for Critical Care settings. The bed platform gives you flexible and customised solutions to keep all care activities and all special patient needs together with an open architecture of different technical features. You can choose the right bed for each care and medical application, patient conditions and budget.

Most Critical Care treatments and tasks focus on the patient's head and chest. Eleganza 4 allows perfect access to the upper part of the patient's body. The unique siderails concept supports the staff in daily care activities as well as in critical situations – accurate, convenient and ergonomic at all times.

- The control panel is easily accessible, and displayed icons are simple to follow.
- The integrated control elements in both upper siderails, the caregiver has direct access to all functions and adjustment options of the bed
- Optimised work flow with Integrated Scales e X-ray translucent mattress and special bed construction and i-Drive Power
- Night Light and Bed Exit Alarm function
- The backrest stops automatically at a tilt of 30° and 45°.
- Safestop® The bed stops automatically if anything becomes trapped between the bed platform and undercarriage when moving to the lowest position.
- i-Brake® – The automatic brake can prevent falls caused by an unstable bed,
- The patient has immediate access to two separate controls for basic bed positioning
- Safe Working Load of 250kg



Specifications:	
Product Code	75101
External Dimensions	217.5 x 100 cm
Bed Extension	0 - 200 mm
Mattress Dimensions	2080 x 900 mm
Mattress Platform Height Adjustment	410 - 790 mm
Backrest Tilt	0 - 65°
Thighest Tilt	0 - 25°
TR / ATR	± 14°
SWL	250 kg



Eleganza 5 Bed

The Eleganza 5 is a symbiosis of top-level innovative technologies and the human approach used to create an optimum working and healing environment for ICU caregivers and patients.

- Eleganza 5 with integrated weighing system
- X ray cassette holder
- 4 castors Tente Integral double, 1 antistatic, 150 mm + 5th power drive wheel
- Centrally operated castors + i-brake with 4 pedals – head sided operation
- Foot controls (height adjustment, examination position) + foot controls (lateral tilt adjustment, GO button)
- Additional adaptor for lifting pole location (1 pc at the head part of undercarriage)
- Control panel iBoard in head section
- Illuminated patient controls in middle siderails
- Get up handles: Mobi Lift® handles; 1 pair with bed height adjustment buttons
- Mattress platform covers: Removable plastic covers (RAL 9006), backrest with X ray adaptation
- Integrated bed extension
- Head and foot ends with safety locks
- Plastic removable bedends design E5
- Backup battery with status signalisation
- Nurse controls: I Board, additional supervisor panel
- Handset with LED flashlight and with illuminated keyboard + plug – play adapter
- Night light and Bed exit alarm

- BedMonitor: EMR ready design (bed is equipped with sensors of the safe bed conditions; bed is ready for the field retrofitting of integration module and i-Track)
- USB charger on both bed sides
- “safe stop” – automatic stop of platform down movement if an object is detected



Specifications:

Product Code	76637
External Dimension (Shortest Position)	2190 × 1000 mm
External Dimension (Longest Position)	2410 × 1000 mm
Mattress Dimension	2080 × 900 mm
Max. Siderails Height	450 mm
Height Adjustment	435 – 815 mm
Lateral Tilt	±15°
TR / ATR Tilt	14° / 14°
Maximum Patient Weight	185 kg
SWL	250 kg

Art of Turning

Handle your patient with gentle care and save your back:

Innovative tilting technology with smooth columns helps to turn the patient gently with less physical effort.

- The nurse will carry out routine nursing procedures ergonomically and effortlessly.
- Improves pressure ulcer prevention – patient positioning is therefore easy.
- New technology such as smooth bed tilting for even more sensitive care

Hands-Free Tilt:

The hands-free operation using the foot controls increases the efficiency of the atraumatic patient turn method. Both of the nurse's hands are free to devote full attention to the patient without any effort whatsoever.



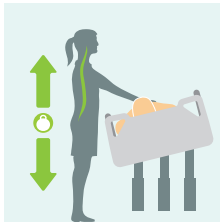
Art of Safety

Multiple Safety:


Sophisticated safety features controlled by the monitoring system.

- Safety under control
- Siderails concept
- I-brake
- Safestop®
- Multizone bed exit alarm
- 30 degree autostop
- Lowest height

Lateral Tilt 15 ° Left, 15 ° Right:



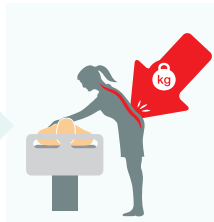The lateral tilt for ergonomically correct care



The normal way of turning a patient can lead to back injury

Art of Mobilisation



New progressive features work in symbiosis for early patient mobilization. Reduces bed exit effort by 50% Symbiosis of 3 elements: When exiting the bed the patient has 3 big aids to put him back on his feet. Mobilization can begin soon so the patient can gain in self-confidence.

1. Mobi-lift®: Unique! The patient controls the bed height using the mobi-lift® handle.
2. Mobirails: Ergonomically shaped handles in the siderails.
3. Lateral tilt: The height setting and side tilt will assist the patient's bed exit.

Art of Positioning



Multi-section frame design: allows unlimited patient comfort and improves nursing efficiency.

- Cardiac chair position
- Fowler position 30/30 with ergoframe®
- Tilttable position
- Vascular position
- Full vascular position

Art of Communication



Smart technology and bed connectivity for a new way of communication. iBoard – 1 control point for the staff

Bed Monitor System

Using sensors built into the bed, the system monitors the settings of the safety parameters.

- Allows long-term monitoring of the bed settings and adjustment to the needs of patients.
- The relevant information can be assessed and the effectiveness increased.
- An effective risk management tool in departments.



Multicare Bed

The Multicare is an accomplished member of the critical care team. It simplifies interventions in critical care situations improving patient care through treatment and rehabilitation whilst providing comfortable working conditions for nursing staff. A combination of advanced innovations with functional features ensures that the Multicare exceeds your expectations.

- Automatic Lateral tilting (both sides to 30 degrees)
- Multi Protect 4-part siderails
- Multiboard features an LCD touchscreen (in both siderails)
- Mattress platform tilt and extension adjustment
- Tilt indicators (backrest, lateral tilt, TR/ATR tilt)
- Programmed positions: CPR, orthopaedic chair, 30° HOB, mobilisation, Trendelenburg
- Activation GO button, safety STOP button
- Service data interface, user communication setup, quick reference guide on LCD touchscreen
- Battery back-up indicator (2 batteries)
- The Multicare also includes patient control panel in the side guards
- The Multicare can be customised to suit your specific needs: there are four different kits
- Mobi Lift handles – assists the patient to exit the bed safely.
- I-brake – an automatic brake that is activated after 10 seconds if the bed remains unbraked and is plugged in. It protects against uncontrolled movement, impacts and falls caused by leaning on an unbraked bed.
- I-Drive – fifth castor for safe transporting and manoeuvring of the bed, with automatic retraction of the castor to the chassis. If the bed is plugged in, it is automatically retracted. The fifth castor runs from its own battery.



Specifications:

Product Code	74200
External Dimensions (Shortest Mattress Platform Position)	1050 x 2150 mm
External Dimensions (Longest Mattress Platform Position)	1050 x 2370 mm
Mattress Dimensions	2080 x 860 mm
Maximum Mattress Height	230 mm
Height Adjustment	440 - 820 mm
Lateral Tilt	±30°
TR / ATR Angle	13° / 16°
Bed Weight	210 kg
SWL	250 kg



Auto Lateral Therapy



Automatic Lateral Therapy:

Maintains the bed and the patient in a permanent cycle of programmed lateral tilts.

Automatic Weighing (AWS)

Knowledge of the weight of a patient in critical condition: In the acute phase, this can involve, for example, determining an exact drug dosage. Historical information about the patient's weight will assist staff in visualizing the current state of nutrition and provide information on the fluid balance.

“FREEZE” Function:

The AWS also includes the “Freeze” function. This enables addition or removal of equipment from the bed (infusions, transport monitor, etc.) without any effect on measurement of the patient's weight.



Weighing:

The Advanced Weighing System integrated into the Multicare® bed is capable of calculating the exact weight of patients in all bed positions.

- Yields current value of the patient's weight
- Provides valuable information on changes over time.



Data Collection

The data-collection system automatically measures and regularly stores values of the patient's weight in its memory.

- Can be easily viewed in clear graphs showing daily, weekly or monthly reports.
- Simple export of the measured data to a PC.

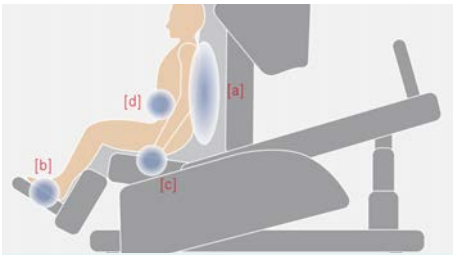
Integrated Mattress



Integrated mattress Symbioso:

- Dual modality
- Fully automatic
- Specially designed for heel protection
- Easy maintenance & operation
- Once touch CPR for bed & mattress
- Lying surface will suit most mattresses

Orthopneic Position



Orthopneic chair® position:

Is intended mainly for patients with resting dyspnea, facilitating respiration. It uses a high Fowler's position, pressure reduction in the abdominal area for deeper breathing and foot support, allowing the use of additional breathing muscles.

Orthopneic Chair®:

- sitting up
- supported feet
- supported arms
- low pressure in the abdominal area

Deep Breathing:

Regression of the back and thigh section reduces the pressure in the area of the pelvis and abdomen facilitating deep diaphragmatic respiration.



Staebel Hospital Floorline Bed

The Staebel Floorline Bed Single with Plastic Boards, CPR, 100mm castors, Battery Backup is one of the lowest beds available on the market and designed for intensive use. The bed lowers to floor level to prevent injuries and accidents caused by falls from beds. All actuators and mechanisms are concealed beneath the platform, allowing the bed to be constructed with no end pedestals! This allows visibility for users, and nursing staff access for procedures.

- The lifting range of the bed is not compromised by its low level – a platform height of 720 mm is achieved making this bed a superb companion in complex care.
- Electric 4-section profiling
- Auto-regression backrest for minimising shear forces
- Trendelenberg – 15° Reverse – 15°
- 75mm twin wheel castors. Optional central locking
- Various colours and styles available
- 10 year frame / 5 year electric warranty / 1 year Handset



Specifications:

Product Code	7779
External Dimensions	2080 x 905 mm
Sleeping Surface	2000 x 890 mm
Sleeping Surface Height	120 - 720 mm
Undercarriage Clearance	Clear View
Head Tilt / Foot Tilt	15° / 15°
Angle Back / Knee	62° / 35°
SWL	200 kg



Essenza 300 Bed

Universal low height bed with an exceptional safe working load of 300 kg improves safety for any patient. With a low height of 28 cm, Essenza 300 is a premium solution for any patient who may be at higher risk of falls. The many amazing bed features increase patients comfort during hospital stay lead to better satisfaction and health outcomes.

- Low height helps to reduce the severity of patients' falls.
- EasyDrive® decreases the necessary initial pushing force up to 50% in comparison with nonpowered bed.
- The safe mobilization concept supports the patient through
- specially designed features & can help to speed up the recovering process
- Single platform efficiency brought by one intuitive user interface. This means only one training for nurses and easier rotation across hospital wards.
- Incredible versatility that meets all the needs of each individual department. SafeSense 3 helps to save time by providing a remote overview of the patients' safe position in the bed or by smart alerting when the patient is at risk of pressure injuries or falls.
- SafeSense SmartTrack facilitates the work of caregivers by tracking hospital assets or patients at room level accuracy.
- Night Light
- Bed Exit Alarm
- USB Integration allows smart devices to be utilised safely from the bed and enables socialization.



Specifications:

Product Code	76771
External Dimensions (Width x Length)	995 x 2200 mm
Max Mattress Dimensions	900 x 2000 mm
Height Adjustment	280 - 780 mm
Split Siderails (Length x Height)	1065 x 420 mm
Backrest Angle	0 - 70°
Calf / Thigh Rest Angle	25° / 35°
TR / ATR Angle	14° / 14°
Bed Weight	235 kg
SWL	300 kg



Essenza 300 Lateral Tilt Bed

The Essenza 300 Lateral Tilt Bed turns into a new gold standard in acute care reducing physical strains on the caregiver and enhancing patient outcomes. With an universal low height bed with an exceptional safe working load of 300 kg improves safety for any patient. With a low height of 28 cm, Essenza 300 is a premium solution for any patient who may be at higher risk of falls. The many amazing bed features increase patients comfort during hospital stay lead to better satisfaction and health outcomes.

- Lateral tilt delivers premium ergonomics to caregivers and effective patients' mobilization.
- Low height helps to reduce the severity of patients' falls.
- EasyDrive® decreases the necessary initial pushing force up to 50% in comparison with nonpowered bed.
- The safe mobilization concept supports the patient through specially designed features & can help to speed up the recovering process
- Single platform efficiency brought by one intuitive user interface. This means only one training for nurses and easier rotation across hospital wards.
- SafeSense 3 helps to save time by providing a remote overview of the patients' safe position in the bed or by smart alerting when the patient is at risk of pressure injuries or falls.
- SafeSense SmartTrack facilitates the work of caregivers by tracking hospital assets or patients at room level accuracy.
- Night Light and Bed Exit Alarm
- USB Integration allows smart devices to be utilised safely from the bed and enables socialization.



Specifications:

Product Code	76772
External Dimensions (Width x Length)	995 x 2200 mm
Max Mattress Dimensions	900 x 2000 mm
Height Adjustment	280 - 780 mm
Split Siderails (Length x Height)	1065 x 420 mm
Backrest Angle	0 - 70°
Calf / Thigh Rest Angle	25° / 35°
Lateral Tilt	±15°
TR / ATR Angle	14° / 14°
Bed Weight	205 kg
SWL	270 kg



Image 3 Bariatric XXL Bed

Image 3, electrically adjustable bed respects the growing demands of the patient and nursing staff. Its attractive design will find universal use in a whole spectrum of hospital care. Image 3 deals with the patient's safety and mobilisation in a unique way: combining the SafeFree® siderails and reduced low height.

Image 3 was designed with focusing on 100% safe recovery of patient in the bed. Safety is the key characteristics of Image 3. Reduced height combined with telescopic SafeFree® siderails create absolutely safe environment for the patient, that can be adjusted according to their specific needs.

Also included is a complex system of supportive automatic features called 3-STOP STRATEGY. It contains functions for safe sleeping, independent mobilisation for quicker recovery, convenient care and cleaning procedures for hospital staff.

This bed contains extended features which prevent pressure injuries and contribute to their easier treatment. Mattress platform profiling called ErgoFrame® provides needed comfort for the patient.



Specifications:	
Product Code	76200
External Dimensions	1100 x 2200 mm
Mattress Dimensions	1000 x 2200 mm
Height of Siderails	490 mm including extender
Maximum Mattress Height	260 mm including extender
SWL	320 kg





Bariatric Tilt Hospital Bed

The Bariatric Hospital Bed offers impressive weight capacity with its maximum lift range. The increased weight capacity doesn't mean we compromise on comfort and quality – the Bariatric series provides ultimate comfort and is available in a range of widths to cater for all needs. The stress on the carers of heavier patients is reduced by relieving physical strain. The strong frame of this model has been enhanced to accommodate the high weight capacity. The heavy duty motors ensure reliable and efficient running of the bed for the safety and comfort of the patient.

- Electric Hi Lo, Electric Backrest and Electric Kneebreak

- Electric Forward & Reverse Trendelenberg
- 150mm Twin Wheel Central locking Castors
- Battery Backup
- Manual CPR release
- APO/Ironstone Powder Coated Frame
- Quality Australian made construction
- 10 year frame warranty, 5 year electrical warranty



Specifications:

Product Code	7622
External Dimensions	1090 x 2170 mm
Sleeping Surface Dimensions	1070 x 2000 mm
Height (Min / Max)	390 - 830 mm
Underbed Clearance	160 mm
Angle Back / Knee	77° / 40°
TR / ATR Angle	16° / 16°
SWL	400 kg





Sentida SC Bed


New concepts for living in old age are in demand. These include the desire to be provided with more homeliness, individuality, safety, comfort and flexibility in the equipment of corresponding living arrangements.

- With the new nursing bed platform sentida sc and customized concepts, we can respond to the different requirements even more flexibly, efficiently and demand-oriented.
- We want to support the maintenance and promotion of individual mobility as an overarching and essential goal in care in a solution-oriented manner.
- SafeFree® side rails – safe mobilization regardless of body size.
- MobiStick 2 – manual mobilization assistant
- SafeControl® handset – clear, ergonomic and intuitive
- Double retraction of the lying surface – safe and comfortable mobilisation in bed
- Automatic positions – quick, easy and ergonomic thanks to the one-button philosophy.
- SafeFree® in orientation position – fear reduction without unnecessary barriers.
- SafeSense® – autonomous, safe bed exit with individual call notification.
- 3-stop strategy – situation- and need-based care, provision and support
- Underbed clearance for mobile patient lifters at any position
- Quality German made construction
- Fully customisable with various models, decor and designs. Contact the team for more information.



Specifications:	
Product Code	77754C
External Dimensions	2080 x 1060 mm
Sleeping Surface	870/900 x 2000 mm
Backrest Adjustment	70° ± 5° (retraction by 120 mm)
Upper Backrest Adjustment	30° ± 5° (retraction by 50 mm)
Head/ Foot Swivel	17° / 17°
Height	230 - 800 mm
Bed Weight	115 kg (depending on model)
SWL	270 kg

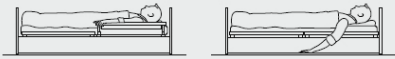Falls Prevention



Sleeping

The low-lying height of only 25 cm (without mattress) ensures safety during night rest. In this position, it is sufficient to set the head-end guard to the "orientation position", avoiding unnecessary barriers and reducing risks of falling.



Protection

Additional protection in the form of a fall mat is recommended for particularly agile residents, which cushions falls from a low position and reduces the risk of injury.



Mobilisation

The side guards, which can be used as mobilisation supports, ensure safety when getting in and out of bed. The bed level can be adjusted to the size of the person when doing this, even for very short residents.



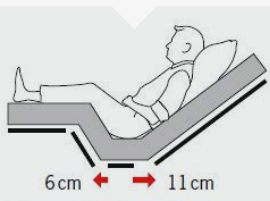
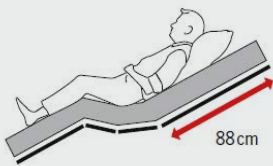
Nursing

So that all necessary nursing duties can be performed with a straight back, sentida can be adjusted to a working height that reduces back strain of up to 68 cm.

Patient Surfaces

Physiological Posture

A back rest with a length of at least 80 cm is an indispensable prerequisite for a physiological posture in bed. This reduces shearing forces, and breathing is aided by the natural bending in the area of the hips. In addition the head is optimally protected.



Free choice of patient surface

01 Aero – the wire mesh patient surface ensures good aeration.

02 Easyclean – the synthetic padding is particularly easy to clean.

03 Comfort – the Comfort patient surface, with its outstanding soft lying properties, provides an effective bedsores prevention. The patient surface partitions can be removed for ease of cleaning. The elements on the edges, with reinforced springs, are used as an orientation when sitting and while asleep.

01



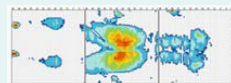
02




03



Dual retraction The dual retraction of backrest (11cm) and leg rest (6cm) helps to reduce pressure in the pelvic area and cut down friction and shearing forces. The risk of bedsores is reduced. In addition, 17cm more space is created in the bed and in the room – residents up to 180 cm in size do not require any patient surface extension that would take up space.



Standard patient surface

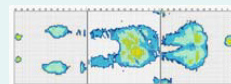

Dual retraction – 17 cm more space in bed

Integrated preparation for patient surface extension

With the patient surface extension prepared as standard, the nursing home bed can be adjusted quickly and without problem for taller patients. Two different types of patient surface extension are optionally available – a temporary extension of the bed (only for divided side guards) or a permanent extension of the patient surface.



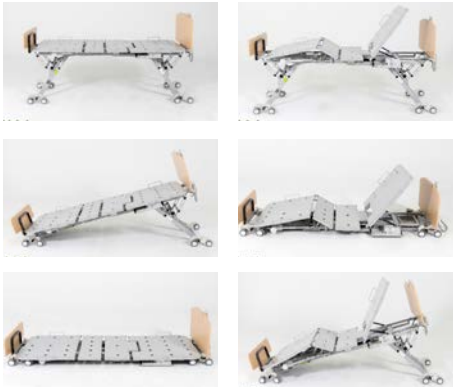
Scan the QR code
to watch the **video
demonstration!**



Staebel Floorline Bed

The Staebel Floorline Bed is one of the lowest beds available on the market and designed for intensive use. The bed lowers to floor level to prevent injuries and accidents caused by falls from beds. All actuators and mechanisms are concealed beneath the platform, allowing the bed to be constructed with no end pedestals! This allows visibility for users, and nursing staff access for procedures.

- The lifting range of the bed is not compromised by its low level – a platform height of 720 mm is achieved making this bed a superb companion in complex care.
- Electric 4-section profiling
- Auto-regression backrest for minimising shear forces
- Trendelenberg – 15° Reverse – 15°
- 75mm twin wheel castors. Optional central locking
- Various colours and styles available
- 10 year frame / 5 year electric warranty / 1 year Handset



Specifications:		Also Available:		
Product Code	7753	Product Code	7755	7758
Size	Single	Size	Wide	King Single
Overall Width	905 mm	Overall Width	1020 mm	1060 mm
Overall Length	2080 mm	Overall Length	2080 mm	2080 mm
Height	120 - 720 mm	Height	120 - 720 mm	120 - 720 mm
Sleeping Surface Dimensions	890 x 2000 mm	Sleeping Surface Dimensions	1000 x 2000 mm	1060 x 2000 mm
SWL	200 kg	SWL	200 kg	300 kg

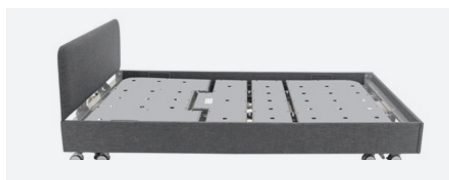


Luxa Homecare Bed

The Luxa is a premium homecare bed that is designed to give the user ultimate comfort and care, while maintaining dignity and style in the bedroom.

This bed is packed with premium features such as an auto-regression backrest, manual vascular support, a high height of 725mm, a low height of 150mm height and optional bluetooth handsets. It also has optional accessories including a sturdy horizontal rail, allowance for self-help poles, underbed lights and more.

- Extremely Low Minimum Position
- Vertical Lift
- Auto-Regression Backrest
- Ash or Ocean Padded Headboard with Full Padded Surround
- Proudly Australian made
- Optional features and styles can be added at time of manufacture (contact the team for more information)



Specifications:

Size	Single
External Dimensions	2160 x 1000 mm
Sleeping Surface	2040 x 900 mm
High Position	725 mm
Low Position	150 mm
Backrest / Knee Angle	62° / 35°
TR / ATR Angle	15° / 15°
SWL	250 kg

Also Available:

Size	King Single	Double	Queen
External Dimensions	2160 x 1170 mm	2160 x 1470 mm	2160 x 1600 mm
Sleeping Surface	2040 x 1070 mm	2040 x 1370 mm	2040 x 1500 mm
High Position	725 mm	725 mm	725 mm
Low Position	150 mm	170 mm	170 mm
Backrest / Knee Angle	62° / 35°	62° / 35°	62° / 35°
TR / ATR Angle	15° / 15°	15° / 15°	15° / 15°
SWL	250 kg	250 kg	250 kg



Bariatric Tilt Bed

The Alrick 2300 Tilt Series provides dignity and comfort for bariatric patients with a 300kg safe working load. This model (2300-T) includes the trendelenburg or tilt feature. With new features such as the auto regression backrest, back-lit handset and underbed light it gives the ultimate comfort and safety to both patients and carers. Available in a range of widths;

- Electric Hi Low and Electric Backrest and knee break
- Electric Forward & Reverse Trendelenburg
- Standard 125mm Braking Castors
- Curved Head & Foot Boards in Natural Oak Finish
- APO Grey & Ironstone Powdercoated Frame



Specifications:		Also Available:	
Product Code	7202004 - S	7202004 - KS	7202004 - Q
Size	Single	King Single	Queen
Overall Width	920 mm	1090 mm	1520 mm
Overall Length	2040 mm	2040 mm	2040 mm
Sleeping Surface Dimensions	900 x 2000 mm	1070 x 2000 mm	1500 x 2000 mm
Angle Back / Knee	77° / 40°	77° / 40°	77° / 40°
TR / ATR Angle	16° / 16°	16° / 16°	16° / 16°
Height	390 - 830 mm	390 - 830 mm	390 - 830 mm
SWL	300 kg	300 kg	300 kg



Bassinet Trolley including Mattress & Tray

Comes with a soft, covered mattress and tub. The frame is made from fully welded stainless steel with rounded corners for maximum infection control and ease of use.

- Rounded corners – minimise impact damage
- Angle adjustment to acrylic tub with flat bar support rails
- Supplied with mattress and storage trays



Specifications:

Product Code	70086
Height	1060 mm
Length	420 mm
Width	860 mm
Mattress Base Height	820 mm



Natally Infant Radiant Warmer

Safe and comfortable environment for newborn children in delivery rooms and during subsequent intensive care.

- Microprocessor-controlled infant warmer
- Measurement and control of manual mode, preheated mode
- Electric bed positioning (trendelenburg and anti-trendelenburg) 12°
- Rotating radiator for 90° on both sides
- 9' Lcd display – display of actual and set values and the time recording of temperature, humidity and oxygen concentration with displaying trends
- The possibility of inserting the x-ray cassettes
- Body temperature probe (multiple use) - 1 pc
- Built-in apgar score
- Undercarriage with fixed height and a shelf
- Guiding bar for accessories



Specifications:

Product Code	86600
Height	1730 - 1930 mm
Length	1024 mm
Width	560 mm
SWL	7 kg

NDIS

NDIS



Shelly Infant Incubator

The Shelly Incubator provides a protected and comfortable environment for children needing care.

- Acoustic Comfort with the noise level in the patient space not going over 38dB
- Thermal Comfort with an air curtain automatically activated to reduce heat loss when incubator door is open
- Plexiglass cover and door to enable high visibility of the child at all times.
- Foot pedals for comfortable height adjustment of the incubator.
- Integrated high precision scales to measure and record weight of the child.
- 9" LED Display can be positioned and rotated to any angle according to the needs of the operator.



Specifications:

Product Code	86500
Height	1450 mm
Length	1150 mm
Width	610 mm
SWL	7 kg



Lilly Phototherapy Lamp

Lilly Phototherapy Lamp - addition to standard blue light also white examining light, aiming points for precise light placement, auto or manual regulation of radiation intensity, equip programmable, patient body temp monitoring with alarm, low noise & easy maintenance



Specifications:

Product Code	86700
Height	119 mm
Length	588 mm
Width	2544 mm
LED Life	60 000 hrs



Tom 2 Cot

A state-of-the-art children's cot must be safe and be able to provide excellent access for nursing staff. The electric positioning of TOM 2 children's cot responds to these and further demands in its original concept of telescopic siderails, column construction mechanism and above standard safety parameters.

The Tom 2 bed fully complies with the new standard for hospital children's beds. DIN 32623:2009-11 and its construction guarantees perfect protection against children fall. Also complies with standards: EN 60601-1, EMC 60601-1-2, EMC 60601-1-6

- 3 height positions for side rails, 0cm, 40cm, 80cm
- Backrest retracts when it raises, meaning zero gap
- Integrated collapsible lockable side rails



Specifications:

Product Code	76765
Mattress Platform	1370 x 700 mm
Electric Hi- Lo	580 – 830 mm
Electric Tilt/Reverse Tilt	15°
Backrest Adjustment	40°
SWL	75 kg

Accessory Available:

Tom 2 Cot Mattress	
Product Code	76766

- 4 x Tente Motion 125mm castors
- Trendelenburg tilt buttons (Trendelenburg / Reverse Trendelenburg position)
- Central operated castors (one castor with steer system)
- 4 x holders for IV in corners of bed
- Undercarriage cover
- Quick removable head/foot ends
- Accessory holders: Pair of accessory rails (located on both sides under side rails)
- Battery Backup
- Horizontal Wall/corner bumpers
- Quality European product





Nino Adolescent Bed

The Nino Series allows children to perceive the world around them while at the same time providing protection and security. The transparent Perspex merges with the room allowing occupants to have an unrestricted view of their immediate surroundings yet still protects against falls and injuries. Also available in different powdercoat colours, or with higher door elements or doors on both sides.

- Electric Hi-low
- Electric Backrest and Kneebreak
- Standard Locking Castors
- Gates one side with Perspex infills
- Fixed Panel one side with Perspex Infills
- 750mm High Sides (from Mattress Base)
- Head & Foot Infills in Natural Oak Matt Finish
- Battery Backup



Specifications:	
Overall Dimensions	2200 x 1014 mm
Sleeping Surface Dimensions	2000 x 870 mm
NET Bed Weight	140 kg
Underbed Clearance	100 mm (at lowest point)
Low Position	390 mm
High Position	830 mm
SWL	200 kg



Configuration Possibilities:

Colour Options:

Poppy	Orange	Lemon	Smoke
Amethyst	Teal	Kiwi	Snow
Cloud	Frost	Pacific	Black



Nestie Maternity Bed with Adjustable Stirrups

The Nestie Double Maternity Bed with adjustable stirrups has been designed for quick and easy use in Birth Centre's and low risk birthing scenarios. The bed has a very high SWL to ensure both parents can be safely on the bed together and to guarantee stability in emergencies. Having the additional leg stirrups allows for easier observations or the ability to suture should it be required.

- Fully adjustable leg stirrups to cater for required position of the mother. These can be removed for easy storage if not in use
- CPR Emergency backrest release
- Removable head and foot boards
- Electric Functions – including forward and reverse trendelenburg, backrest and kneebreak.
- Battery backup
- Underbed and handset light



Specifications:

Product Code	7591
Overall Dimensiona	1390 x 2170 mm
Height	390 - 830 mm
Sleeping Surface Dimensions	1370 x 2000 mm
Underbed Clearance	160 mm
Angle Back / Knee	77° / 40°
TR / ATR Angle	16° / 16°
SWL	350 kg



Trendelenburg



CPR

Ave2 Birthing Bed

The AVE 2 birthing bed brings expectant mothers maximum safety and comfort. Thanks to the AVE 2, the birth takes place naturally according to their expectations. For healthcare staff the bed is a reliable and practical tool which improves the quality and facilitates the care of the woman in labour at all stages of birth. It creates a feeling of a pleasant and safe environment for the woman in labour and family members.

- Trendelenburg and CPR
- Extra battery backup
- ABS easily disinfect
- Lockout Control
- No scissor-lift column system
- Central brake easy to access
- Easy transport with directional castor and pushing handles
- Removable headboard for anesthesiologist (optional)
- Tilting holders for IV pole



Specifications:

Product Code	77120
Total Length of the Lying Area	2132 mm
Minimum Length in the default setting	1570 mm
Overall Width	980 mm
Castor Size	150 mm
Height	600 - 1050 mm
Back Section	-10° / +70°
Seat Section	0° / +18.5°
Trendelenburg Position	0° / -10°
Leg Rests Vertical	0° / +135°
Leg Rests Horizontal	0° / +60°
Foot Section Adjustment Angle	0° / -22°
Maximum Patient Load	240 kg
Overall Product Weight	236 kg
Lifting Pole Bearing Capacity	75 kg
SWL	250 kg

Siderails



Side rails with control panel easily integrated and hideable when transferring patients from one bed to another. Silent and easy turning with one hand. No additional space needed beside the bed. It serves as support to help patients sit up.

Simple Control

Integrated control panel with clear pictograms and symbols for woman in labor and nursing staff.

Nurses control locking of controlling movements.

Controls enabled:

1. High and Tilt Adjustment
2. Trendelenburg Position
3. Lockout Control
4. Under Bed Light (optional)
5. Nurse Call (optional)



Foot Section



Simple swivel mechanism to hide the foot section which is fast and easy for delivery. It remains attached to bed, no storage on the floor means infection prevention and prevention of dorsal injuries from being incurred by the personnel when manipulating with the product.



Hygiene & Safety



High hygienic standard and compact construction of particular sections of the rest area and the plastic cover of the bed base reduce the time needed for the cleaning of the bed.

1. Bacterial static leather
2. Seamless upholstery
3. Compact and separate space
4. Electrostatically painted with epoxy powder coating

The firm construction of the birthing bed contains no cutting points or sharp edges, which increases the safety of the attending personnel.



Legrest Concept

Simple swivel mechanism to hide the foot section which is fast and easy for delivery. It remains attached to bed, no storage on the floor means infection prevention and prevention of dorsal injuries from being incurred by the personnel when manipulating with the product.

Quick adjustment of the leg rests into the required position.

Independent Seat Tilting

- The unique oval shape creating space needed for perineal access fits closely to the birthing tub minimizing transfer gaps.
- There is independent motorized pelvic tilt for better coccyx elevation.



Hand Grips

Increase the comfort of a woman in the second stage of childbirth.

Column Concept

Innovative telescopic solution provides extra low height for accessing the bed safely and extra height for treatment.



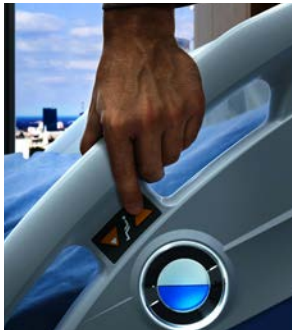


Eleganza 2 With Junior Kit

Maximum safety for small patients. Functionality and simple controls. This unique hospital bed, thanks to the side rail safety concept, is suitable for adult and paediatric patients. Patients from 90 to 220cm in height and can safely use the bed.

The combination of practicality and beauty make the Eleganza 2 with Junior Kit a bed which will meet the needs and expectations of each hospital ward.

- Integrated Mobilift innovation with support handles
- Highest possible level safety level to eliminate falls or injuries and to indicate safe bed positioning
- Attractive design with innovative column construction
- 'Smart' patient monitoring features: safety monitor, bed monitor and vital monitor
- Ease of transport with fifth castor
- Increased patient comfort and ease of cleaning
- Patient control of some functions
- Practical and ergonomic accessories available



Specifications:

Product Code	76618
External Dimensions	995 × 2230 mm
Bed Extension	310 mm
Bed Height (Hi-Lo Adjustment)	395 – 775 mm
Internal Dimensions / (Compatible with Mattress Width)	2000 × 900 mm / (900 and 860 mm)
Max. Mattress Height	180 mm
Height of Split Plastic Side Rails	405 mm
Clearance of Undercarriage	150 mm (150 mm tente castors)
Ergoframe®	110 mm
Maximum Backrest Angle	70°
Thigh / Calf Rest Angle	35° / 21°
TR/RTR Angle	16° / 16°
Lifting Pole Bearing Capacity	75 kg
Bed Weight	140 kg
Max. Patient Weight	185 kg
SWL	250 kg



DVA

NDIS



DVA

NDIS

MASS

Staebel Deluxe Commode Bedside

This bedside commode offers a full range of features including a slide back seat with durable toilet seat and pan rack. This comes standard with bowl and lid.

A skirt that surrounds the seat conceals the bowl etc. so the chair can be likened to a regular bedroom sitting chair

- Adjustable seat height
- Comfortable upholstered slide back seat
- Comfortable upholstered shallow curve backrest
- Plastic toilet seat, bowl and lid
- Padded armrests
- Powdercoated steel frame
- Vinyl upholstery



Specifications:

Product Code	7955
Overall Width	600 mm
Overall Height	870 - 1020 mm
Overall Depth	510 mm
Seat Width	460 mm
Seat Height	420 - 570 mm
Seat Depth	440 mm
Back Height	470 mm
SWL	120 kg

Also Available:

Product Code	7960	7961	7962	7963
Seat Width	500 mm	550 mm	600 mm	650 mm
SWL	135 kg	175 kg	195 kg	210 kg

Staebel Standard Bedside Commode

The standard bedside commode offers padded armrests for comfort and features adjustable height, a removable padded seat cover and a plastic toilet seat with bowl and lid.

- Value for money
- Stackable design
- Padded armrests
- Powdercoated steel frame
- Vinyl upholstery



Specifications:

Product Code	7952PA
Overall Width	575 mm
Overall Height	940 - 1090 mm
Overall Depth	600 mm
Seat Width	440 mm
Seat Height	500 - 650 mm
Seat Depth	440 mm
Back Height	475 mm
Toilet Seat	475 - 625 mm
SWL	125 kg

NDIS

NDIS



VT6 Mobile Bedside Cabinet - 2 Drawer

The VT6 mobile bedside cabinet stands out with its functional yet elegant design. This model features a top drawer, a shelf and a large dump drawer. With its curved surfaces and rounded handles this one is ideal in an aged care application. High quality swivel castors (2 with brakes) Integrated bottle holder in dump drawer. Available in a full range of decor to complement the beds.

- 510mm x 516mm x 770mm high
- 50mm double castors
- Option of top drawer lockable
- Option of middle drawer (instead of shelf)

Other design available upon request




VT6 Mobile Bedside Cabinet - 3 Drawer

Available in various woodstain laminates the Wissner-Bosserhoff VT6 mobile bedside cabinet is commonly chosen by health care facilities, including those where there are high falls risk residents. With curved surfaces this one is ideal in a high dependency residential care application.


- Styled to complement any bed any Movita, Carisma, Sentida bed
- This model has the attractive woodgrain decor in Beech B00 stain
- Generous size - 510mm x 516mm x 770mm high
- Features two standard drawers and one dump drawer
- Top drawer lockable - 2 keys supplied
- Quality swivel castors, two which are braking

Other design available upon request


Colour options:



Beech Natural



Lindberg Oak



Havana Cherry

Specifications:	
Product Code	77923C
Width	510 mm
Length	516 mm
Height	770 mm

Specifications:	
Product Code	77932C
Width	510 mm
Length	516 mm
Height	770 mm



Staebel Deluxe Mobile Overbed Table

The Staebel Deluxe Overbed Table offers the latest in smart and modern design, with an easy clean, seamless PE top to promote ease of use and satisfy infection control requirements. Featuring smooth lines, no joins or hard to reach contours, this innovative design makes previously hard to maintain items easier and quicker to clean.

- Easy gas lift height adjustment, non-tilt
- PE easy clean top, 34mm deep, max 20kg load.
- Base height 98mm
- Four castors, including two with brakes
- Supplied fully assembled



Specifications:

Product Code	4210
Top Width	900 mm
Top Depth	430 mm
Minimun Height	725 mm
Maximum Height	1050 mm
SWL	20 kg



Staebel Basic Mobile Overbed Table

A very popular mobile table with a number of useful features, making it an ideal choice in many situations. With its easy lever operated height adjustment and smooth curved U-design, it fulfils all needs.

- Spring assisted lever height adjustment – lever stays with top eliminating need to bend to adjust table height
- Powdercoated U-shape base, round tube minimises injury risk
- Swivel castors, two are lockable

*** Table top is sold separately. There are a range of colours on offer**

Other colours available upon request



Specifications:

Product Code	4215
Width	450 mm
Length	770 mm
Minimun Height	700 mm
Maximum Height	1080 mm
SWL	20 kg



Staebel Bed Rail Protector

This Bed Rail is made from a polyurethane material, coupled with a bioelastic incontinence cover, making it the ultimate answer for longevity and effectiveness.

No zips, clips, or flimsy foam means a zero risk of skin tears or wounds caused by the patient/resident hitting the side guards. Simply push it down onto the top rail and it is in place.

No time consuming, technical attachment points that foul up with the side guards.

- Skin friendly
- Anti Microbial safety
- Easy to clean and to disinfect
- Easy to use and install
- 2 x protectors can be used together to cover the complete length of a 4-part side guard or full length side rails.
- Sold singularly



Use with Wi-Bo Bed Mobilisation Aids



Specifications:	
Product Code	3983
Width	130 mm
Length	450 mm
Height	800 mm

Blanket Support

The bed cradle provides relief from the weight of bed clothes on delicate feet and legs

- Creates clearance between bed linen and the mattress
- Ideal for keeping bed linen off sensitive areas of the patient
- Slides between mattress and bed base
- Powder coated mild steel construction

Specifications:	
Product Code	7920
Overall Width	320 mm
Overall Height	470 mm
Overall Reach	450 mm

08.

Patient Safety

Safe Sense Kits	180	Cura1 Cordless Chair Pad with Transmitter	184
Staebel Socks	181	Cura1 Cordless Floor Mat with Transmitter	184
Posey Non Skid Slippers	181	Cura1 Pad Transmitter	185
Staebel Crash Mat	182	Under Mattress Cordless Bed Pad	185
Cordless Multiport Expander Bracket Mount	182	Cura1 Personal Emergency Transmitter	185
Ibeam Patient Monitor	183	Wireless Call Button with Pull Cord	185
EzyStand	183	Cura1 Wristband Alarm	185
Cura 1 Nine Channel Cordless Monitor	183	Cura1 Curly Cord Lead	185
Cura1 Cordless Bed Pad with Transmitter	184	Cura1 Replacement Break Cords	186



Safe Sense Nurse Call - Adaptor Wired

- SafeSense® = Safety + freedom - By day and by night
- Sustainable fall prevention for utmost safety. This is what the SafeSense® bed exit assistance system of wissner-bosserhoff stands for. It can considerably increase the level of safety, especially during night and weekend shifts with fewer staff. Especially at night, the bed exit assistant system can minimize the need for side guards. The freedom of the resident/patient is upheld.
- Faster response time in the event of falls
 - Greater safety during night and weekend shifts
 - Ideal in combination with sentida and image 3 low nursing beds and SafeFree® side guards
 - Sensor-switchable night light enables light e.g. on the way to the toilet at night
 - The nurse call timing is flexibly set from 0 to 30 min. according to the mobility of the resident/patient

Safe Sense Kit - Wireless

- SafeSense® = Safety + freedom - By day and by night
- Sustainable fall prevention for utmost safety. This is what the SafeSense® bed exit assistance system of wissner-bosserhoff stands for. It can considerably increase the level of safety, especially during night and weekend shifts with fewer staff. Especially at night, the bed exit assistant system can minimize the need for side guards. The freedom of the resident/patient is upheld.
- Faster response time in the event of falls
 - Greater safety during night and weekend shifts
 - Ideal in combination with sentida and image 3 low nursing beds and SafeFree® side guards
 - Sensor-switchable night light enables light e.g. on the way to the toilet at night
 - The nurse call timing is flexibly set from 0 to 30 min. according to the mobility of the resident/patient



Specifications:	
Product Code	77013

Specifications:	
Product Code	77015



Staebel Socks

Staebel socks provide safety and prevention of falls.

- Soft terry cloth upper. Latex-free.
- Colour helps care staff identify high fall risk patients
- Grip on top and bottom
- Mesh band on large socks comfortably fits swollen ankles and calves
- Standard size fits up to adult size 13

Available in Green, Orange, Purple, Red or Yellow in sizes Standard and Large

Posey Non Skid Slippers

- Posey fall management slippers provide safety and reduce the risk of falls.
- Soft terry cloth upper
- Latex-free rubber sole
- Colour helps care staff identify high fall risk patients
- Machine washable
- All colours available in five sizes
- One-pair per package

Available in Green, Orange, Purple, Red or Yellow

Please call 1800 267 267 for full range of colour and sizes



Specifications:

Small	9900
Medium	9902
Large	9904

Specifications:

Product Code	9PSFP	
Size	Ladies (Shoe Size)	Mens (Shoe Size)
Small	5 - 6	N/A
Medium	6 - 7	5 - 6
Medium / Large	7 - 8	6 - 7
Large	8 - 9	7 - 9
Extra Large	10 +	10 - 12

NDIS

NDIS



Staebel Crash Mat

- Foldable Falls Impact Reducing Mat.
- Easy to clean
- Made from heavy duty latex free material
- Chamfered edges to reduce the risk of tripping
- Can be easily folded for storage or transportation



Specifications:	
Product Code	9822
Width	910 mm
Length	1760 mm
Thickness	40 mm
Sizes Available:	
Staebel Crash Mat - 4buy	
Product Code	9822/4

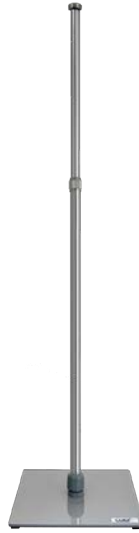


Cordless Multiport Expander Bracket Mount

Plugs into a nurse call system and can connect up to cura1 wireless devices (2 x pad & 3 x other transmitters eg PET), a 2638 Nurse Call Pendant and one other wired device (eg a Premium Floor Mat).

For nurse call panel configurations that are too small or crowded to accept the 2544 Cordless Multiport Expander, there is the 2549 Cordless Multiport Expander Bracket Mount, which has its nurse call plug on a Curly Cord, allowing the unit to be installed to the panel.

Specifications:	
Bracket Mount	
Product Code	92549
Cordless Multiport Expander	
Product Code	92549



iBeam Patient Monitor

The narrow beam creates a precise detection zone.

This is especially useful for preventing unwanted alarms if other people frequent the same room, or where only a narrow area requires monitoring.

Whether across a door, stairs or beside a bed, there are countless ways that our iBeam can be used to protect those in your care and to make the job of being a caregiver easier.

When used with our sturdy, adjustable-height EzyStand, precise positioning of the iBeam could not be simpler.

EzyStand

The EzyStand / iBeam combination allows a caregiver to quickly deploy targeted motion detection alarms in a variety of positions such as next to a bed, in front of a chair, near a door or stairs.

Our innovative 2723 EzyStand magnetically attaches to the 2721 iBeam and significantly adds to its ease of use.

The stand is sturdy, attractive and perfect for use in both professional and home environments.

Cura1 9 Channel Cordless Monitor

- Monitor up to 7 pads or mats simultaneously
- Comes with soft touch removable protective cover, fixing strap and batteries
- Adjustable alarm tone and volume
- Selectable 0, 2sec, 15 min delay setting
- Can be connected to a nurse call system
- Alarm memory and automatic reset
- Pause feature for caregiver to attend patient without setting off the alarm

Specifications:

Product Code 927413

Specifications:

Product Code 92723

Specifications:

Product Code 92700

NDIS

NDIS

NDIS



Cura1 Cordless Bed Pad with Transmitter

- Pad includes pouch that contains a transmitter
- Guaranteed 1 year life expectancy
- Easy to clean
- Non-allergenic soft PVC
- For use with cura1 cordless bed monitors
- 50 meter range (*Subject to environmental conditions)
- Vent holes in the neck reduce the risk of moisture -ingress
- Must be used with 92700 Monitor

Cura1 Cordless Chair Pad with Transmitter

- Simple operation and maintenance
- Includes transmitter
- For use with Cura1 cordless monitor
- Easy to clean
- Virtually no false alarms over the lifetime of the pad
- All bed/chair pads made from non-allergenic soft PVC
- 12 month life expectancy
- Must be used with 92700 Monitor

Cura1 Cordless Floor Mat with Transmitter

- This cordless floor mat is positioned beside a resident's bed to alert staff when the resident has put their feet on the floor, attempting to leave the bed.
- Simple operation
 - Easy to clean
 - Virtually no false alarms over the lifetime of the pad
 - Made from heavy duty non-slip material
 - Generous size compared to others on the market
 - Quality Australian made product
 - 12 month life expectancy
 - 50 meter range (*Subject to environmental conditions)
 - Includes transmitter
 - Must be used with 92700 Monitor



Specifications:	
Product Code	92662
Width	760 mm
Length	510 mm
SWL	30 - 200 kg



Specifications:	
Product Code	92661
Width	380 mm
Length	255 mm
SWL	30 - 200 kg



Specifications:	
Product Code	92663
Width	610 mm
Length	1220 mm
SWL	30 - 200 kg

NDIS



NDIS



NDIS



Cura1 Pad Transmitter

- Works with cordless bed pads, chair pads and floor mats
- Inserts into cordless pad and transmits alarms to 92700 Monitor
- Up to 50 wireless range
- 1 year battery life

Under Mattress Cordless Bed Pad

The ultimate in patient comfort and system reliability. This bed pad is discreetly installed underneath the mattress.

- Please be certain that the bed frame is compatible (has a flat area for the pad to sit on).
- Recommended weight range for reliable operation: 40kg to 200kg.

Cura1 Personal Emergency Transmitter

- Large rubber button for easy use by elderly patients
- Standard Lanyard for patient to wear around neck
- Works with single use hospital wristbands
- Changable battery with projected 1-year life
- Splash resistant

Specifications:

Product Code 92719

Specifications:

Product Code 93515
SWL 40 - 200 kg

Specifications:

Product Code 92588

NDIS



NDIS



NDIS



Wireless Call Button with Pull Cord

An IP 65 rated wireless call button for use with other cura[®] wireless equipment. Fully wireless – no cabling required.

Cura1 Wristband Alarm

- For use with cura[®] StaffAlert or LRCS Duress System
- Lightweight and small in size
- Low cost and low maintenance
- Replaceable battery
- Includes buckle strap. Optional lanyard available
- Medical Device Immunity: EN 60601-1-2

Cura1 Curly Cord Lead

RJ12 Plug to Right Angle 6.35mm Stereo Plug 300mm.

- Compatible with 92549 Cordless Multiport Expander Bracket

Specifications:

Product Code 92558

Specifications:

Product Code 92595

Specifications:

Product Code	Style
94160	T&S
94161	T&R
94162	R&S



Cura1 Nurse Call Cables

- Plugs into Cura1 Monitor to connect to other equipment e.g.. a Nurse Call system
- 90 degree 6.5 plug on one end
- 2 metre flexible cord

Contact the team to see more nurse call cables on 1800 267 267.

Specifications:

Product Code	
92643	



USB B

Cura1 Replacement Break Cord

USB B to right angle 6.5mm Mono Plug (T&S).

- For use with Bed / Chair Pads and Floor Mats with USB Breakcord connection.

Specifications:

Product Code	Length
92481	2 Metres
92456	3 Meters
92457	5 Metres



USB B

Cura1 Replacement Break Cord

USB B to right angle 6.5mm Stereo Plug T&R) 2m.

- For use with Bed / Chair Pads and Floor Mats with USB Breakcord connection.

Specifications:

Product Code	Type	Length
92494C	T&S	2 Metres
92494B	T&R	3 Meters



USB B

Cura1 Replacement Break Cord

USB B to right angle 3.5mm Mono Plug 2m.

- For use with Bed / Chair Pads and Floor Mats with USB Breakcord connection.

Specifications:

Product Code	
92482	



USB B

Cura1 Replacement Break Cord

USB B to right angle 6.5mm Stereo Plug 2m (T&S).

- For use with Bed / Chair Pads and Floor Mats with USB Breakcord connection.

Specifications:

Product Code	
92496	



USB B

Cura1 Replacement Break Cord

USB B to RJ10 Plug –(Wired Pin 1 & 4).

- For use with Bed / Chair Pads and Floor Mats with USB Breakcord connection.

Specifications:

Product Code	Length
92480	2 Metres
92458	3 Meters
92459	5 Metres

09.

Treatment & Physio Tables

Staebel 2 Section Physio Table - Electric Hi-Lo	188
Staebel 3 Section Physio Table - Electric Hi-Lo	188
Staebel 3 Section Physio Table - Fully Electric	189
Novak 3 Section Gynae Table - Electric Hi-Lo	189
Staebel Bariatric Treatment Table	190
Staebel Bobath / Neuro Treatment Table	190
Staebel Saddle Stool with Backrest	191
Staebel Saddle Stool	191
Staebel Therapy Stool	192
Staebel Gynecology Chair	192
Novak ZeroG Comfort Treatment Chair	193
Novak ZeroG Luxury Treatment Chair	193



Stael 2 Section Physio Table - Electric Hi-Lo

The preferred choice of physiotherapists this short head couch provides a strong platform for treatments. With a gas strut supported head and backrest coupled with a smooth hi-lo operation this table represents great value.

- Electric Hi Lo operation by 360 degree foot bar
- Sleek black frame with rounded frame legs for safety
- Head section has face hole and plug
- Adjustable feet for uneven surfaces
- 7.5cm retractable castors for safety and mobility
- Reliable TiMotion motor
- Supportive foam padding
- Certification: CE, FDA, ISO13485-2016
- Upholstery: Heavy duty PVC vinyl that is fire, stain, mildew, oil, water and wear resistant
- 5cm Thick Padding for user comfort



Specifications:	
Product Code	111004
Overall Width	660 mm
Overall Length	1950 mm
Head Tilt Angle	- 45 ° to 40 °
Back Section Angle	80 °
Height	475 - 930 mm
Head Section Length	400 mm
Mid Section Length	540 mm
Leg Section Length	930 mm
SWL	250 kg



Stael 3 Section Physio Table - Electric Hi-Lo

- Our most popular product this three equal section couch with dual gas strut supported sections & smooth hi-lo operation represents great value. Australia's number 1 choice for medical centres for many years.
- Electric Hi Lo operation by foot switch
- Sleek design with rounded frame legs for safety
- Three equal sections for universal use in any room
- Head section has face hole and plug
- Dual gas strut supported back & footrest
- Adjustable feet for uneven surfaces
- 7.5cm retractable castors for safety and mobility
- Reliable TiMotion motor
- Supportive 5cm foam padding
- Upholstery: Heavy duty PVC vinyl that is fire, stain, mildew, oil, water and wear resistant
- Certification: CE, FDA, ISO13485-2016



Specifications:	
Product Code	1133NB
Head / Foot Tilt Angle	- 40 ° to 65 °
Height	465 - 925 mm
Head Section Length	630 mm
Mid Section Length	630 mm
Leg Section Length	630 mm
Overall Width	710 mm
Overall Length	1950 mm
SWL	250 kg



Staebel 3 Section Physio Table - Fully Electric

- A high quality all electric treatment couch facilitating effective patient treatment all whilst looking after practitioner OH&S.
- Electric Hi Lo operation by foot switch
- Electric backrest & footrest for effortless operation
- Sleek new design with rounded frame legs for safety
- Head section has face hole and plug
- Adjustable feet for uneven surfaces
- 7.5cm retractable castors for safety and mobility
- Reliable TiMotion motors
- Supportive 5cm foam padding
- Upholstery: Heavy duty PVC vinyl that is fire, stain, mildew, oil, water and wear resistant
- Certification: CE, FDA, ISO13485-2016
- Warranty: 7 Years Construction, 3 Years Electrics



Specifications:

Product Code	11 - 156003E
Head / Foot Tilt Angle	- 40 ° to 75 °
Height	470 - 910 mm
Head Section Length	630 mm
Mid Section Length	630 mm
Leg Section Length	630 mm
Overall Width	710 mm
Overall Length	1950 mm
SWL	250 kg



Novak 3 Section Gynae Table - Electric Hi-Lo

Three-section gynaecology table S with its design and added elements is perfect for examinations in gynaecology, urology and general medicine.

Plastic base cover, smooth surfaces, seamless upholstery and PVC cover ensure easy cleaning and minimal use of disinfectants.

- Table width is 65 (75 cm option available)
- Thicker seamless padding (6 cm)
- Supplied with: retractable castors (Ø 75 mm) to move at short distances, gynaecological legrests on fixation rail, retractable inox bowl under the seat, support handles on fixation rail
- Electric High-Low Function



Specifications:

Product Code	111531C
Table Width	700 mm
Head Tilt Angle	-35 ° - 87 °
Min Height	430 mm
Max Height	910 mm
Head Section Length	685 mm
Mid Section Length	565 mm
Leg Section Length	680 mm
Overall Width	700 mm
Overall Length	1950 mm
SWL	255 kg

NDIS

NDIS



Staebel Bariatric Treatment Table

- 3 section support surface with electrical height adjustment specifically designed for Bariatric use
- Fitted with 4 retractable rubber castors
- Head section with ergonomically designed removable plug and facial opening
- Nicely curved cushions with premium padding
- MRSA Resistant Anti-bacterial & Fire Retardant upholstery
- 320kg safe working load
- 750mm wide x 2040mm long
- Head tilt angle -25 to +80 degrees



Specifications:

Product Code	111010
Table Width	760 mm
Head Tilt Angle	-25° - 80°
Min Height	460 mm
Max Height	910 mm
Head Section Length	620 mm
Mid Section Length	560 mm
Leg Section Length	800 mm
Overall Width	760 mm
Overall Length	2015 mm
SWL	320 kg



Staebel Bobath / Neuro Treatment Table

- Extra wide support surface ideal for neurological and bobath procedures
- Generously padded for maximum patient comfort
- Fitted with 4 retractable rubber castors
- Head section with ergonomically designed removable plug and facial opening
- Nicely curved cushions with premium padding
- MRSA Resistant Anti-bacterial & Fire Retardant upholstery
- Retardant upholstery
- Paper roll holder included
- Operated via convenient foot control (Can request with hand control as long as this is specified at time of order)



Specifications:

Product Code	111020
Table Width	1200 mm
Head Tilt Angle	-20° - 85°
Min Height	460 mm
Max Height	910 mm
Head Section Length	650 mm
Leg Section Length	1410 mm
Overall Length	2080 mm
SWL	320 kg



Staebel Saddle Stool with Backrest

This ergonomically designed saddle stool helps to reduce the discomfort of sitting for long periods and allows the user free range of upper body movement.

- High quality gas stool features a deluxe comfortable padded seat
- Anti-bacterial upholstery
- Five castor base offers greater stability and easy manoeuvrability
- Height adjustable for various positions

Staebel Saddle Stool

This ergonomically designed saddle stool helps to reduce the discomfort of sitting for long periods and allows the user free range of upper body movement.

- High quality gas stool features a deluxe comfortable padded seat
- Five castor base offers greater stability and easy manoeuvrability
- Height adjustable for various treatment positions
- Anti-bacterial upholstery
- Foam thickness: 6cm



Specifications:	
Product Code	111313
Min Height	550 mm
Max Height	760 mm
SWL	150 kg



Specifications:	
Product Code	111311
Min Height	560 mm
Max Height	630 mm
Seat Diameter	350 mm
SWL	160 kg

NDIS

NDIS



Staebel Therapy Stool

This high quality gas stool features a comfortable padded seat that is 9cm in thickness.

- Anti bacterial upholstery
- Five castor base offers greater stability and easy manoeuvrability
- Height adjustable for various treatment positions



Staebel Gynecology Chair

- Three section chair top
- Head inclination 12° (head length 30 cm)
- Back part inclination electrically adjustable from - 10° to +65° (length 60 cm)
- Seat inclination electrically adjustable from 0° to +30°
- Footrest manual adjustable
- Arm supports manual adjustable
- Height adjustment electrically adjustable - column lift from 56 to 94 cm
- Stable base on the leveling feets
- Seamless padding from artificial leather for medical purposes (5 cm) wear-resistant and disinfectant tolerant material
- Range of upholstery colours available



Specifications:		Specifications:	
Product Code	111305	Product Code	71160
Min Height	460 mm	Width	700 mm
Max Height	640 mm	Length	1380 mm
Seat Diameter	350 mm	Weight	74.5 kg
SWL	150 kg	SWL	150 kg





Novak ZeroG Comfort Treatment Chair

The ZeroG Comfort Treatment Chair.

- Electric chair adjustments with clear and understandable hand switch:
 - Independently adjustable back, seating and leg section.
 - Pre-saved positions
 - *Zero Gravity position, which offers pressure relief for the body & joints
 - *Seating position to sit and prepare for the therapy.
 - *Completely flat, lying position to rest or for reanimation.
 - *Trendelenburg position
- Soft padding made of high-quality foam, that provides full body support along the lying surface and prevents overheating and sweating
- Footrest, adjustable by pressing the button and supports patients in case of cramps
- Manually adjustable head section with 7 different levels
- Strong and simply adjustable armrests by pressing the lever
- Push handles and anti-static individually lockable castors
- Paper roll holder
- Superior design for hygiene and infection control



Specifications:

Product Code	111ZOG
Length	2200 mm
Width (Including Armrests)	1030 mm
Width of Rest Area	550 mm
Height	580 mm
SWL	200 kg



Novak ZeroG Luxury Treatment Chair

- Hand switch for adjustment of height and for back, sitting and leg section
- Extra soft padding (14cm)
- Adjustable head section
- Wide, comfortable and adjustable armrests by pressing the lever
- Adjustable footrest by pushing and releasing a built-in button
- Antistatic castors with central locking system (ø100mm)
- Paper roll holder and push handles
- Infusion pole
- Electric Hi-Lo



Specifications:

Product Code	111TCZ - LUX
Length	2200 mm
Width (Including Armrests)	1030 mm
Width of Rest Area	550 mm
Height	610 - 860 mm
SWL	200 kg

10.

Mental Health Range

Secura 1	195
Secura 2	196
Secura 3	197
Mental Health Fireproof Mattress	198
OASIS - 3-Position Sleeper	199
Empire Plastic Lounge Chair	200
Valiant Plastic Lounge Chair	201
Liberty Range	202
Fusion Range	203
Verity Patient Recliner	204



Secura 1

The Secura 1 is a robust bed for demanding environments, providing a sleeping surface which can reduce harm to both the patient and carers.

There are no electrics included, hence creating a safe solution for your facility, with a SWL of 200kg and durable framework. Suitable for your medium-high care patients.

- No electrics
- Fixed height solid frame
- Secure bolt on head and footboards
- Strong robust construction
- Enclosed mattress supports
- Anti-tamper screws
- Ironstone Powdercoat



Specifications:

Product Code	7691 - S
Size	Single
Overall Width	920 mm
Overall Length	2000 mm
Sleeping Surface Dimensions	900 x 2000 mm
Mattress Base Height	450 mm
SWL	200 kg

Also Available:

Product Code	7691 - KS
Size	King Single
Overall Width	1090 mm
Overall Length	2070 mm
Sleeping Surface Dimensions	1070 x 2000 mm
Mattress Base Height	450 mm
SWL	200 kg



Secura 2

The Secura 2 provides dignity and comfort for all mental health patients, with its impressive safe working load of 300kg, designed specifically to reduce stress on carers by relieving physical strain.

Secura 2 is fitted with heavy-duty actuators and concealed wiring to ensure reliable and efficient functioning. Suited for low to medium-care patients who may have clinical care needs with reduced anti-ligature features.

- Electric Hi-Low
- Enclosed anti-ligature area
- All bolts secured with Loctite

- Ironstone PowdercoatElectric Hi Low
- Braided electrical power supply
- Anti tamper screws
- Secure bolt on head & footboard
- Embedded footboard control
- Concealed wiring
- Curved Natural Oak head & footboard
- Enclosed mattress supports
- 125mm braking castors



Specifications:		Also Available:	
Product Code	7692S	Product Code	7692KS
Size	Single	Size	King Single
Overall Width	920 mm	Overall Width	1090 mm
Overall Length	2070 mm	Overall Length	2070 mm
Sleeping Surface Dimensions	900 x 2000 mm	Sleeping Surface Dimensions	1070 x 2000 mm
Height	390 - 830 mm	Mattress Base Height	390 - 830 mm
SWL	300 kg	SWL	300 kg



Secura 3

The Secura 3 is the ultimate safe and secure solution for your facility, eliminating access to the bed mechanisms, with the bed surround maintaining the concealment in both, Hi and Low positions, the lower bed section sitting stationary with the bed.

The sides of the bed are recessed to allow for nurses to kneel on the bed without the obstruction of a bed surround in the shin area. Secura 3 is the solution to your WHS concerns for the carers and cleaners.

(*Pictured with optional backrest)

- Concealed castors
- Embedded footboard control
- Brake lever access
- Braided electrical power supply
- Anti-tamper screws
- Lowered section for kneeling on bed
- Head & footboard lower than mattress height
- Natural Oak, 13mm Compact Laminate
- Ironstone powdercoat



Specifications:	
Product Code	7693 - S
Size	Single
Overall Width	930 mm
Overall Length	2030 mm
Sleeping Surface Dimensions	900 x 2000 mm
Height	400 - 700 mm
SWL	250 kg

Also Available:	
Product Code	7693 - KS
Size	King Single
Overall Width	1100 mm
Overall Length	2030 mm
Sleeping Surface Dimensions	1070 x 2000 mm
Mattress Base Height	400 - 700 mm
SWL	250 kg



Mental Health Fireproof Mattress

Fireproof Mattress is specially designed for demanding environments with an increased risk of fire. It has been tested and approved according to the world's most stringent fire safety standards.

The mattress cover resists open gas flames and glowing particles up to 900° C, which minimises the risk of fire and smoke generation. The cover is composed of an aramid-based composite material and sealed by high frequency (HF) welding, making the mattress waterproof and exceptionally durable. It is also easy to clean and disinfect.

Compliant fire safety standards

- SS 8760010 (test method NT-Fire 055)
- BS 6807:2006 Crib 7
- FAR 25.853
- FAP 915-1 (police custody)
- SWL: 250kg
- Easy to clean and disinfect
- Water resistant and extremely durable



Specifications:	
Product Code	3271050
Width	600 mm
Length	1900 mm
Thickness	100 mm
SWL	250 kg

Also Available:	
Product Code	327051
Width	900 mm
Length	1900 mm
Thickness	100 mm
SWL	250 kg



OASIS - 3-Position Sleeper

Sleep tight with the Oasis Sleeper, which opens from the side to conserve space in cramped patient rooms. Its central locking caster system – operated by a pedal – allows caregivers and cleaning staff to reposition the Sleeper with ease.

Features:

- Limited Lifetime Warranty
- Easily moved by caregivers or cleaning staff
- Central Locking Caster system activated by easy-access pedals
- Large, 28 x 80 inch sleep surface
- All upholstery is field-replaceable

Options:

- Choice of left or right hand side pull out.
- Custom finishes
- Choice of wood, polyurethane or solid surface arm caps
- Choice of upholstered or laminate front panel
- Contrasting fabrics
- CAL TB 133 Compliance
- Moisture Barrier

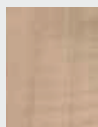
Healthcare Features:

- Heavy-duty all steel frame construction
- Seat and back foam meets or exceeds CAL 117-2013 flammability requirements


Cleaning-Vinyl:

Vinyl Clean with soft cloth dampened with warm water and mild soap solution. Stubborn stains and spills may be treated with a sponge or soft brush and mild non-abrasive cleaner.


Wood / Laminate finish selection : (Option)




Hardrock
Maple (HM)

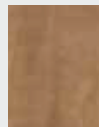


Candlelight
(CL)

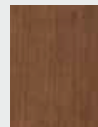

Rustic Cherry
(RC)

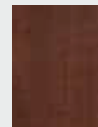


Summer Flame
(SF)

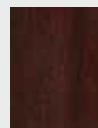

Mahogany
(MA)



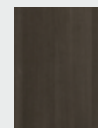Cherry
Blossom (CB)



Grenada
(GR)



Chocolate
(CH)

3-Positions



2 person sofa



Chaise position



Sleep position



Empire Plastic Lounge Chair

Empire is a one-piece, rotationally-moulded plastic lounge chair designed for intensive use environments. Durability and ergonomics come together in the Empire to provide maximum strength, comfort and safety.

Features:

- All fasteners are tamper-resistant
- 450kg SWL
- Cal TB 133 compliant
- Stain-resistant, seamless and easy to clean

Options:






- Can be used as stand-alone or be ganged together
- Can be bolted together, in place, as well as to the floor
- Ability to add weight



Specifications:

Seat Height	17 inch
Overall Width	24 inch
Overall Depth	32 inch
Overall Height	33 inch
SWL	450 kg

Colour Options:





Empire Armless



Empire with arms



Valiant Plastic Lounge Chair

Valiant is a one-piece, rotationally-moulded plastic lounge chair designed for intensive use environments. Durability and ergonomics come together in the Valiant to provide maximum strength, comfort and safety.

Features:

- All fasteners are tamperproof
- 680 kg static load weight rating
- Cal TB 133 compliant
- Stain-resistant, seamless and easy to clean

Options:

- Can be used as stand-alone or be ganged together
- Can be bolted together, in place, as well as to the floor
- Ability to add weight



Colour Options:



Specifications:

Seat Height 17 inch

Overall Width 31 inch

Overall Depth 32 inch

Overall Height 33 inch

SWL 680 kg



Liberty Range

The Liberty BH Collection exudes a warm, residential aesthetic without sacrificing durability and safety. Tops, doors and drawer fronts are surfaced with seamless thermofoil and feature softened edges to reduce risk of injury.

Options:

- Custom finishes and dimensions
- Can be bolted to wall if required
- Bed options include restraint loops and bolt-to-floor brackets
- Units can be supplied with various lock types

Features:

- Cabinets can be specified with drawers, door(s) or open fixed shelves
- Drawers are non-removable without a special tool
- Anti-ligature features include integral finger pulls, continuous hinges and J-Bar in wardrobe
- Cabinet tops are 1" thick and flush-mounted to inhibit handling of the unit
- Wardrobe doors hung on steel piano hinges
- Wardrobe top is sloped to eliminate non-visual surface area
- All joints reinforced with double-dowel construction

Construction - Case Goods:

- Thermoformed tops, drawer fronts, doors and toe kicks
Thermoformed components are constructed using 3/4" MDF core surfaced with wood grain thermofoil laminate. Reverse sides are finished in white melamine
- All drawer fronts and doors will be a minimum 3/4" thickness unless otherwise specified. Drawer sides will be epoxy coated steel with 5/8" white melamine bottoms and backs

- All cabinets feature 3/4" backs surfaced with matching laminate
- Wardrobe features steel piano hinge. Glides will be minimum 5/8" diameter and will be nonmarring, corrosionproof plastic

Construction of Bed:

- Platform Bed features double-wall construction consisting of a 3/4" 7-Ply poplar plywood ladder-style inner framework with (4) cross-members
- All joints are doweled and glued
- Exterior walls constructed of 3/4" laminated particle board, attached with concealed screws from the inside
- All edges are edged with 3mm PVC tape attached with high-strength adhesive
- All fasteners are tamper-resistant
- Glides are adjustable but non-removable



Cleaning Instructions

All of the surfaces listed below will withstand cleaning using a solution of 0.5 Hydrogen Peroxide or a 70% Isopropyl alcohol solution following the recommended cleaning frequency of a patient room environment. When cleaning with either of the above solutions, it is a requirement that the solution not be allowed to sit for an extended time on the surface but should be wiped thoroughly within a few moments of application.

Laminate Surfaces

Dust with a soft, clean cloth. To remove fingerprints etc. scrub with a damp soapy cloth followed by a rinse with cold water. Stubborn stains can be removed with hot water and a bit of vinegar. Avoid using abrasive cleaning tools as these may scratch the laminate surface. It may be necessary to use a mild non-abrasive cleaner.


Wood / Laminate finish selection : (Option)

							
Hardrock Maple (HM)	Candlelight (CL)	Rustic Cherry (RC)	Summer Flame (SF)	Mahogany (MA)	Cherry Blossom (CB)	Grenada (GR)	Chocolate (CH)



4 Shelf Wall Mount Unit

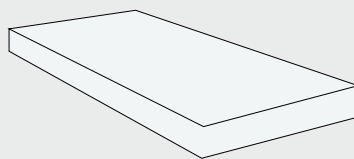The wall mount unit comes with 4 open fixed shelves and is attached to the wall to ensure ample storage. This cabinet is constructed using 3/4" MDF core surfaced with wood grain thermofoil laminate. The cabinet is made to ensure there are no anti-ligature features. This is done by ensuring the shelves are non-removable and minimum a 3/4" thick.



Behavioral Health Mattress

Behavioral Health mattress offers extreme durability, featuring a heavy-duty polyurethane coated nylon that is seamless, tear-resistant and treated to resist fungal growth and bacteria. Mattress size to suit Platform Bed 2.0.

With Draws



Open Fixed Shelves



Bedside Table

This Liberty bedside table has the option of being supplied with open fixed shelves or 3 draws. Like all items in this range, this cabinet is constructed using 3/4" MDF core surfaced with wood grain thermofoil laminate. The cabinet is made to ensure there are no anti-ligature features. This is done by ensuring the drawers are non-removable without a special tool and drawers will be minimum a 3/4" thick. This can be bolted to the floor or wall to inhibit handling of the unit.



With Draws



Open Fixed Shelves

Desk

The Liberty range desk can be supplied with open fixed shelves or 3 draws. Like all items in this range, the desk is constructed using 3/4" MDF core surfaced with wood grain thermofoil laminate. The desk is made to ensure there are no anti-ligature features. This is done by ensuring the drawers are non-removable without a special tool and drawers will be minimum a 3/4" thick. The desk can be bolted to the floor or wall to inhibit handling of the unit.

Interior shelf with J-Bar & door



Open Fixed Shelves



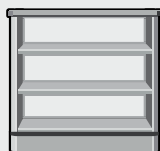
Single Wardrobe

The single wardrobe can have doors or be supplied with opened fixed shelves. The wardrobe is constructed using 3/4" MDF core surfaced with wood grain thermofoil laminate. The doors will be a minimum 3/4" thickness unless otherwise specified. The wardrobe will feature 3/4" backs surfaced with matching laminate. All wardrobe features steel piano hinge. The Wardrobe is made to ensure there are no anti-ligature features. This is done by including integral finger pulls, continuous hinges and J-Bar in wardrobe. The wardrobe tops are 1" thick and flush-mounted to inhibit handling of the unit. The top of the wardrobe is sloped to eliminate the non-visual surface area.



Platform Bed

The Platform Bed features a double-wall construction consisting of a 3/4" 7-Ply poplar plywood ladder-style inner framework with (4) cross-members. In ensure safety for every party, all joints are doweled and glued. The exterior walls constructed of 3/4" laminated particle board, attached with concealed screws from the inside. All edges are edged with 3mm PVC tape attached with highstrength adhesive. All fasteners are tamper-resistant. Its design ensures there are no ligature points for risk of harm to oneself or another. Additional weight can be added or it can be bolted to the wall or floor. It can come in your choice of a custom finish and dimensions.




Dresser Open Shelf

This dresser is supplied with 3 opened fixed shelves. It is constructed using 3/4" MDF core surfaced with wood grain thermofoil laminate. The dresser is made to ensure there are no anti-ligature features. This is done by ensuring the shelves are non-removable and minimum a 3/4" thick.



Fusion Range

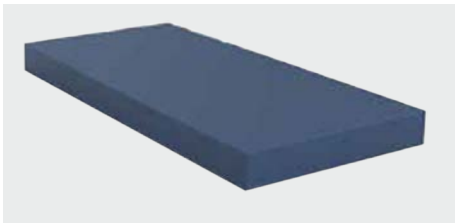Highly durable, safe and easy to maintain, the Fusion series includes beds, shelving and accessories for patient room spaces that are rotationally molded for maximum durability. That means all units are molded into one continuous, seamless piece of high impact polyethylene to achieve and sustain a consistency of strength.



Colour Options:



The Fusion collection is designed to fully furnish any behavioral health room. It is also ideal for use in detention, and correctional facilities.

- 10 Year Warranty
- Durable, one piece, rotationally molded design for increased Strength. Constructed of high impact "no break" polyethylene.
- No steel or wood components.
- Manufactured to ISO Quality Standards
- Resistant to stains.
- Easy to clean and maintain.
- Flame Retardant (passes Cal 133). Test reports available upon request.
- All pieces can be bolted to the floor and wall via concealed pilot holes.

Features:

- Corners and edges are rounded for increased safety
- Cal TB 133 compliant
- Stain-resistant, seamless and easy to clean

Options:

- Custom colors available for larger projects
- Can be bolted to floor or wall
- Bed options include storage compartment and bed riser

Cleaning Instructions:

Clean with a soft cloth dampened with a warm water and soap solution. Stubborn stains and spills may be treated with a sponge or soft brush and mild non-abrasive cleanser.

Mattress for Behavioral Health

Behavioral Health mattress offers extreme durability, featuring a heavy-duty polyurethane coated nylon that is seamless, tear-resistant and treated to resist fungal growth and bacteria.



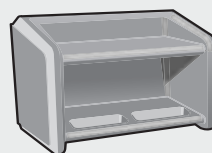
Wall Mount Desk

Bolt to the wall via concealed pilot holes with hardware of your choice.




4 Shelf Wall Mounted Unit

Shelving for open storage of personal items. Bolt to the wall via concealed pilot holes with hardware.




Wall Shelf

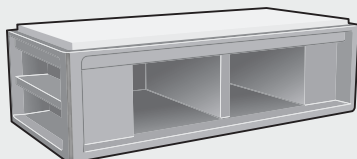Shelving for open storage of personal items. Bolt to the wall via concealed pilot holes with hardware of your choice.



Platform Bed

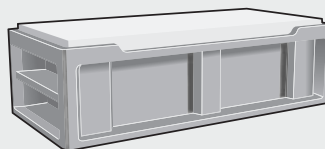End shelving for open storage of personal items. moulded side slots for provide increased storage and structural span strength. Bed can be reverse-mounted to eliminate storage feature. Can be bolted to the floor via concealed pilot holes with hardware of your choice.



Platform Bed 2.0

End shelving for open storage of personal items. Angled lowered side rail for easier entry and exit. Can be bolted to the floor via concealed pilot holes with hardware of your choice.

Options: Side slots for additional storage



Floor Mount Stool

Can be bolted to the floor



Waste Can

Virtually unbreakable and easy to clean and wipe down..



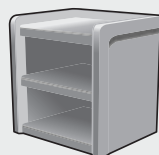
Mirror

Made from 100% virgin acrylic which is easy to install and hard to break.



Open Chest

Moulded in shelving for open storage of personal items. Can be bolted to the wall via concealed pilot holes with hardware of your choice.



Night Stand

Bedside storage of personal items. Can be bolted to the wall via concealed pilot holes with hardware of your choice.

Bed Riser (For Platform Bed 2.0)

The Fusion Bed Riser is designed to be combined with the Fusion 2.0 Bed for increased bed height for easier patient care, easier patient entry and exit, as well as allowing for mechanical lift device use. Bolts to bed through existing mounting holes.

Can be easily weighed down with sand on site. Can be bolted to the floor via concealed pilot holes with hardware of your choice.



Verity Patient Recliner

Verity Recliner offers healthcare-grade durability and exudes an inviting, homelike aesthetic that brings residential warmth to clinical environments.

Features:

- Limited lifetime warranty
- 4 locking casters
- 3 position recline with easy to operate mechanism
- Ergonomically-positioned push handles for patient transport
- Stitching method reduces seams for increased infection resistance

Options:

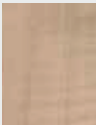





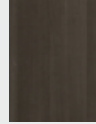

- CAL TB133 compliance
- Moisture Barrier
- Optional arm caps available in wood, polyurethane and solid surface
- Contrasting fabric

Cleaning-Vinyl:

Vinyl Clean with soft cloth dampened with warm water and mild soap solution. Stubborn stains and spills may be treated with a sponge or soft brush and mild non-abrasive cleaner.

Note that seat and back upholstery covers can be removed for more thorough cleaning.

Wood / Laminate finish selection : (Option)

			
Hardrock Maple (HM)	Candlelight (CL)	Rustic Cherry (RC)	Summer Flame (SF)
			
Mahogany (MA)	Cherry Blossom (CB)	Grenada (GR)	Chocolate (CH)

3-Positions



This image shows a full page of blank, lined paper. It features approximately 28 horizontal blue or grey lines spaced evenly apart, typical of notebook paper. The lines extend across the entire width of the page, leaving small margins at the top and bottom. There are no vertical lines, text, or other markings on the page.

We are
Here to Help



Let's have a chat!

If you see something you like, you can book in your consultancy with one of our friendly team members. See below information for ways you can get in contact...

1800 267 267 | www.activemedicalsupplies.com.au | hello@activemedicalsupplies.com.au



P. 1800 267 267 | **F.** 07 3899 8825

E. hello@activemedicalsupplies.com.au

www.activemedicalsupplies.com.au