

activeMEDICAL

Pressure Care Portfolio

Issue 2





Contents

01. Cushions	6
02. Cushion Matrix	17
03. Mattresses	20
04. Mattress Matrix	33
05. Ancillary	36
06. Seating	46

We are Here to Help

For over 30 years, Active Medical Supplies has delivered safe and innovative solutions for the care of clients in the Aged Care, Hospital, School and Home Care sectors. We strive to provide the highest quality equipment for Special Schools that provides a safe environment and will withstand the test of time.

We partner with experts in the field of design and manufacturing to provide the very best resources to ensure that your loved ones are provided with the correct equipment to receive optimum care.

We have a commitment to lower our carbon footprint and protect our environment. We utilise our own delivery vehicles and use sustainable materials for packaging to reduce environmental damage. We are also proud to provide beds and equipment with world leading, energy-saving electronics.

For your peace of mind our product is manufactured to very high quality and we have a secure supply chain. 90% of our products are sourced from Australia, Europe, USA, Canada, United Kingdom, Spain, India, Sweden, Vietnam and Taiwan.

Procurement & Distribution

Active Medical has developed a network of manufacturers in Australia and around the world – from Germany and Sweden to the United States of America and consults directly with industry professionals to ensure the products we provide responds to unique patient specifications.

The success of our business is supported by our ability to collaborate with manufacturers that not only provide the highest quality of products and services, but are aligned with our mission and commitment to

our social, corporate and environmental responsibilities.

Our efficient distribution systems are underpinned by our strong relationships with manufacturers and two centralised warehouses with a capacity to hold up to 50,000 items. This ensures the fast turn around and on time delivery demanded by the health care industry and the needs of people in care. We ensure that all our products are ethically sourced and supplied and meet environmental standards.



Did you know...



We sanitize **all products** that leave our warehouses



We offer full product **demonstrations and video trials**



We offer **contact free delivery**



Our staff are **double vaxxed** with boosters



Our showroom is a **controlled environment** that we can also offer trials


*Disclaimer

Active Medical Supplies does not take responsibility for incorrect application of products. This publication is a guide only to assist in selecting products suitable to your clients. Please use in conjunction with professional advice. Products and information listed in this catalogue are subject to change. Website and digital information supercedes all print. For the most up to date information, please contact us directly.

What is a Pressure Injury?

A pressure injury is an injury to a persons skin or underlying tissue; normally over a bony prominence, as a result of pressure. They take a long time to heal, negatively impacting the patients quality of life, causing severe pain, and lead to infection.

Fall Risks

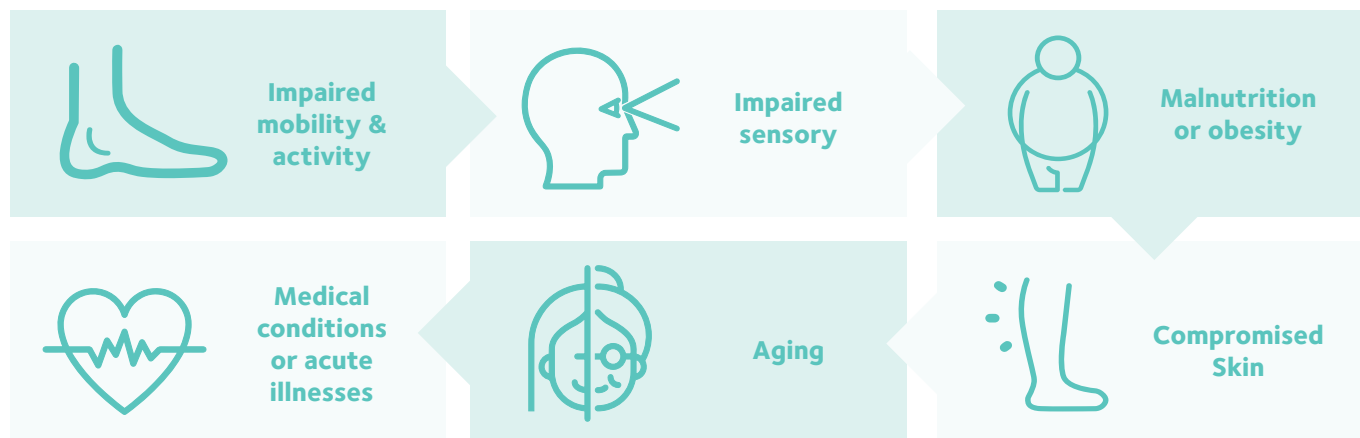

Did you Know...

The 3 most common pressure injury areas are:

34%

Sacrum/Coccyx

21%

Heels

5%

Nose



9.7 pressure care injuries occur every week in over **10k hospitals**

Stages of Pressure Injuries

01.

The area looks red and feels warm to the touch. It may burn, hurt, or itch.

02.

The area is damaged and may have an open sore. There is more pain and the skin around the wound may be discolored.

03.

The area now has a crater-like appearance due to damage under the skin's surface.

04.

The area is severely damaged with a large wound. The persons skin and underlying tissue is at a high risk of infection.

Selecting the correct product?

01.

Check residents weight

02.

Check resident cushion or mattress size

03.




Check residents pressure injury risk


04.

Use the icons in the catalogue as a guide

01.

Cushions

	LOW RISK			HIGH RISK	
Low Resilient Foam Cushion	7		Action Gel Professional Cushion	12	
Action Gel Commuter Cushion	7		Action Gel X-Act Classic Cushion	12	
Staebel Deluxe Foam Cushion	8		Relax DuoGel Cushion	13	
Relax Trio Cushion	8		Relax DuoProfile Cushion	13	
Staebel Bariatric Cushion	9		Relax GelCell Cushion	14	
	MEDIUM RISK		Relax GelCell G2 Cushion	14	
Static Air Cushion	9		Starlock Cushion	15	
Action Gel Pilot Cushion	10		Standard Contour Star Cushion	15	
Action Gel Shear Smart Pad	10		Stabil-Air Star Cushion	16	
Relax Easy High Profile Cushion	11		Galaxy Star Cushion	16	
Relax Easy Cushion	11				

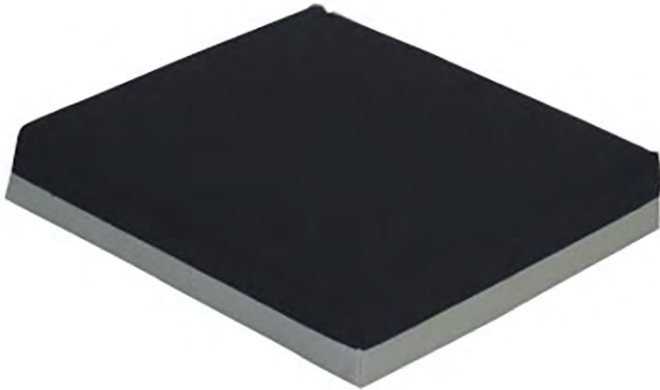

Low Resilient Foam Cushion

A very popular choice for the low risk user, the low resilient foam cushion flows under pressure to provide a fluid like support. It helps to eliminate pressure points by transferring support uniformly.

- 2-way stretch top cover
- Zipped cover included

Action Gel Commuter Cushion

The Action Commuter offers basic pressure and shear reduction for low risk users. Combined 13mm Akton polymer and 45mm foam are sealed in a water proof film for incontinence protection.

- 5.7cm thick
- Breathable double-layer cover included
- Easy to care for and wipe down



Specifications:

Product Code	3210
Width	450 mm
Length	400 mm
Thickness	90 mm
SWL	100 kg

Sizes Available:

Product Code	3219
Width	400 mm
Length	400 mm
Thickness	90 mm

Specifications:


Product Code	3561
Width	450 mm
Length	400 mm
Thickness	57 mm
SWL	Unlimited

Sizes Available:

Product Code	3560	3562	3563	3564
Width	400 mm	450 mm	500 mm	500 mm
Length	400 mm	450 mm	450 mm	500 mm
Thickness	57 mm	57 mm	57 mm	57 mm

NDIS

NDIS

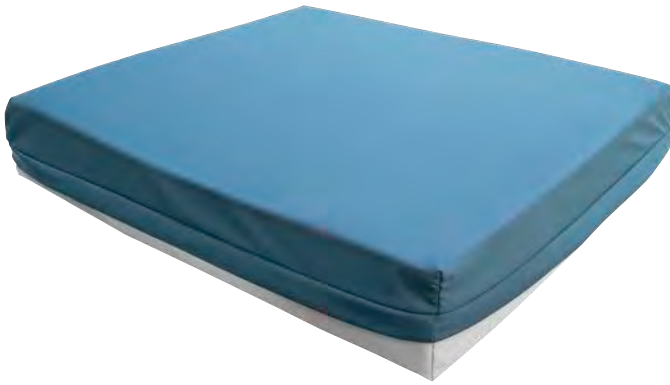

Staebel Deluxe Foam Cushion

A comfortable cushion for the low risk user, the deluxe pressure relief foam cushion flows under pressure to provide a fluid like support.

- Suitable for low risk clients
- 40mm top layer of memory foam. 50mm base of standard foam
- 2-way stretch top cover
- Zipped cover included with waterfall zip

Relax Trio Cushion

The Relax Trio is an ultralight cushion that benefits from four sections of FloGel being positioned into a polyurethane base, with a viscoelastic foam upper layer. It is designed to conform to the shape of the user, increasing surface area and relieving pressure.

- Viscoelastic slow
- Fast Recovery Foam
- Includes outer cover (not pictured)
- Pressure Relieving FloGel (4 sections)
- Internal clear protective PU film coating for easy clean and incontinent users

***Pressure Mapping Reports and successful Case Studies Available on request**



Specifications:

Product Code	3223
Width	550 mm
Length	450 mm
Thickness	90 mm
SWL	100 kg

Sizes Available:

Product Code	3225	3227
Width	450 mm	600 mm
Length	450 mm	450 mm
Thickness	90 mm	90 mm



Specifications:

Product Code	35332
Width	450 mm
Length	400 mm
Thickness	38 mm
SWL	80 kg

Sizes Available:

Product Code	35330	35333
Width	400 mm	450 mm
Length	400 mm	450 mm
Thickness	38 mm	38 mm
SWL	80 kg	80 kg

NDIS

NDIS

MASS



Staebel Bariatric Cushion

A comfortable bariatric cushion for the low risk user, Helps to eliminate pressure points by transferring support uniformly. 2 foam layers, top layer of 50mm memory foam for comfort and is supported by 100mm of standard foam to ensure no bottoming out. Stretch cover with waterproof zip for easy cleaning.

- Suitable for low risk clients
- Zipped, removable cover included
- 260kg SWL
- Other sizes available upon request

Static Air Cushion

CuroCell® S.A.C. is a static air seat cushion for comfort, prevention and as an aid in the treatment of pressure injury up to category 2.

The cushion pre-filled and very easy to handle. The durable soft material for long term use provides a superior combination of comfort and support.

- Wipe with cleaning agent and/or disinfectants
- Air holes for venting heat and moisture
- Pre-inflated low profile design
- Fits all chairs
- Durable soft material for long term use
- Patient ready

*Pressure Mapping Reports and successful Case Studies available on request



Specifications:

Product Code	3270
Width	660 mm
Length	520 mm
Thickness	150 mm
SWL	260 kg

Sizes Available:

Product Code	3271	3272
Width	710 mm	610 mm
Length	550 mm	460 mm
SWL	260 kg	260 kg



Specifications:

Product Code	3417	3418
Width	450 mm	500 mm
Length	450 mm	500 mm
SWL	160 kg	160 kg

NDIS

MASS

NDIS



Action Gel Pilot Cushion

The Pilot is a low-profile cushion that provides exceptional pressure relief for most users with a pressure ulcer history. 25mm of Akton polymer provides immersion for bony prominences and protects clients with fragile skin integrity. Incontinence covers are an optional extra

- Not restricted by a weight limit
- Easy to care for
- Quality, durable Pilot cover included - can be machine washed

*Custom sizes also available.

Action Gel Shear Smart Pad

The SHEAR SMART pad, combined with the breathable double-layer cover, significantly reduces the effects of shear while providing superb user comfort. This cushion will conform to contoured bases for additional support when required.

- Incontinence cover available on request.
- Breathable double layer cover included
- Cube construction allows air circulation and heat dissipations and comfort
- Action polymer gel
- Easy to care for and wipe down



Specifications:

Product Code	3531
Width	450 mm
Length	400 mm
Thickness	25 mm
SWL	Unlimited

Sizes Available:

Product Code	3530	3532	3533
Width	400 mm	450 mm	500 mm
Length	400 mm	450 mm	450 mm
Thickness	25 mm	25 mm	25 mm

Specifications:

Product Code	3524
Width	400 mm
Length	400 mm
Thickness	14 mm
SWL	Unlimited

Sizes Available:

Product Code	3526	3527	3528	35220
Width	450 mm	450 mm	500 mm	500 mm
Length	400 mm	450 mm	450 mm	500 mm
Thickness	14 mm	14 mm	14 mm	14 mm

NDIS
MASS



NDIS
MASS



Relax Easy High Profile Cushion

This cushion from the Relax range is easy to use and is an ideal pressure relieving cushion for use in a long term care environment. It features a layer of FloGel sandwiched between an upper and lower layer of viscoelastic foam. Relax High Profile has the same properties as the standard lax, therefore it conforms to the shape of the user, increasing surface area and relieving pressure.

- Wipe down outer cover, machine washable to 60° C (not pictured)
- 8.5cm thick for extra comfort and additional immersion
- Moving FloGel adds comfort and pressure relief
- Ideal for residents requiring stability or greater depth
- Viscoelastic foam on both sides reduces 'hammock effect'

*Pressure Mapping Reports and successful Case Studies Available on request

Relax Easy Cushion

This cushion from the Relax range is easy to use and is an ideal pressure relieving cushion for use in the nursing home environment. This combination of pressure relieving FloGel and viscoelastic foam facilitates effective pressure re-distribution while allowing the patient to immerse themselves comfortably into the cushion. The cushion has no specific front or back, top or bottom, making it simple to use.

- Conforms to shape of user, increasing surface area and relieving pressure
- Includes waterproof wipeable PU stretch outer cover
- Outer cover machine washable to 60° C
- Pressure Relieving FloGel - Segmented For Added Conformity

*Pressure Mapping Reports and successful Case Studies Available on request



Specifications:

Product Code	35392
Width	400 mm
Length	400 mm
Thickness	85 mm
SWL	190 kg

Sizes Available:

Product Code	35394	35395	35396	35397
Width	450 mm	450 mm	450 mm	500 mm
Length	400 mm	450 mm	500 mm	500 mm
Thickness	85 mm	85 mm	85 mm	85 mm

Specifications:

Product Code	35312
Width	450 mm
Length	400 mm
Thickness	50 mm
SWL	130 kg

Sizes Available:

Product Code	35100	35313	35315	353120
Width	400 mm	450 mm	500 mm	500 mm
Length	400 mm	450 mm	450 mm	500 mm
Thickness	50 mm	50 mm	50 mm	50 mm
SWL	130 kg	130 kg	130 kg	130 kg

NDIS

MASS

NDIS



Action Gel Professional Cushion

The Professional provides exceptional pressure relief for users at risk of skin breakdown or with a history of pressure ulcers. A full 4.1cm thickness of Akton polymer provides effective weight distribution on a contoured or flat surface. The Professional provides an additional 50% cushioning against pressure, shear and vibration. This is significant for those active clients using a power wheelchair who have increased concerns about bottoming out and the potential effects of increased temperatures (due to foam insulation properties). Custom sizes available upon request.

- 41mm thick
- Suitable for very high risk clients, including medium risk as a preventative measure
- Not restricted by any weight limit
- Easy to care for
- Basic cover included

Action Gel X-Act Classic Cushion

The Xact Classic offers outstanding pressure distribution and anti-shear properties. The 1.4cm Action polymer cube pad and laminated foam base help the user maintain a comfortable, therapeutic position. Rehab therapists most often recommend the Classic for users requiring forward and lateral positional support.

The recent re-design that incorporates viscoelastic pink foam offers improved support for gluteal areas and bony protuberances. Custom sizes available.

- 110mm thick for optimum comfort
- Suitable for clients up to very high risk, and medium risk as a preventative measure
- Easy to care for
- BodySmart cover included - is able to be machine washed
- Straps for attachment to wheelchairs.

SAFE WORKING LOAD
UNLIMITED

HIGH RISK

PREVENTION
P

Specifications:

Product Code	3541
Width	450 mm
Length	400 mm
Thickness	41 mm
SWL	Unlimited

Sizes Available:

Product Code	35422	3542	3543	3544
Width	450 mm	450 mm	450 mm	500 mm
Length	550 mm	450 mm	500 mm	500 mm
Thickness	41 mm	41 mm	41 mm	41 mm

SAFE WORKING LOAD
714KG

HIGH RISK


PREVENTION
P

Specifications:


Product Code	3552
Width	450 mm
Length	450 mm
Thickness	110 mm
SWL	114 kg



NDIS
MASS

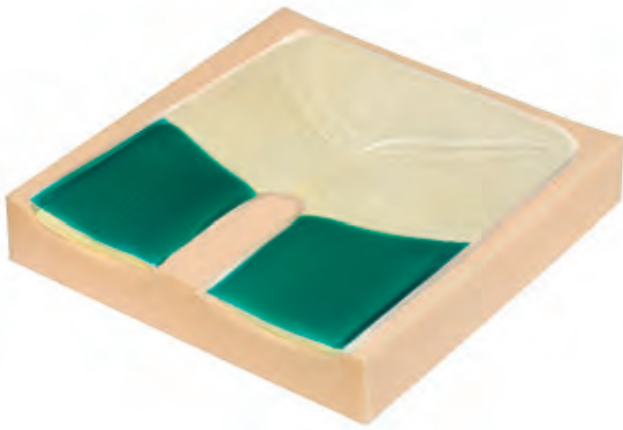

NDIS
MASS

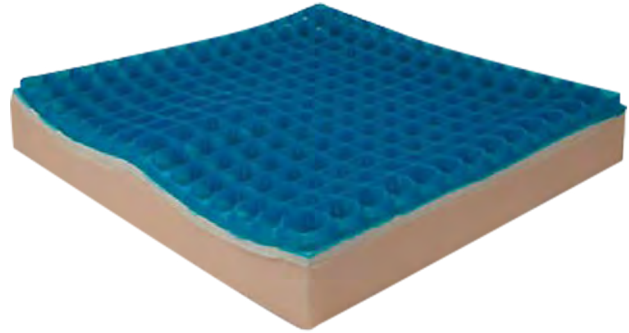

Relax DuoGel Cushion

The Relax DuoGel is a unique combination of TruGel and FloGel working together effectively to provide excellent pressure relief and weight distribution. The DuoGel's anatomical pre-contoured base provides for medial and lateral thigh positioning as well as support for the greater trochanter of the thigh bone. The thigh troughs are lined with TruGel and the ischial area is lined with FloGel which is set on a viscoelastic foam giving greater comfort and seating stability.

- Thigh troughs assist with seating stability
- Fast Recovery Foam (Anatomically Shaped)
- Includes machine washable outer cover (not pictured)
- Clear protective PU film coating - easy to clean and suitable for incontinent users
- Anatomical Pre-Contoured Foam Base

Relax Duo Profile Cushion

The Duo-Profile offers a pre-contoured surface giving excellent hip and leg alignment. The three layered design creates comfort for users with up to a medium pressure injury risk.

The Duo-Profile uses a collapsible gel topper. This unique design allows you to quickly immerse into the cushion, while it conforms to your shape.

- Top layer medical-grade silicone
- Includes machine washable outer cover (not pictured)
- Mid-section Slow Foam/ Memory Foam
- Base Fast Recovery Foam
- Clear protective PU film coating - easy to clean and suitable for incontinent users



Specifications:

Product Code	35302
Width	450 mm
Length	400 mm
Thickness	80 mm
SWL	130 kg

Sizes Available:

Product Code	35300	35303	35307	35308
Width	400 mm	450 mm	500 mm	550 mm
Length	400 mm	450 mm	500 mm	450 mm
Thickness	80 mm	80 mm	80 mm	80 mm
SWL	130 kg	130 kg	130 kg	130 kg



Specifications:

Product Code	35380
Width	400 mm
Length	400 mm
Thickness	80 mm
SWL	130 kg

Sizes Available:

Product Code	35382	35383
Width	450 mm	450 mm
Length	400 mm	450 mm
Thickness	80 mm	80 mm
SWL	130 kg	130 kg

NDIS

NDIS



Relax GelCell Cushion

A unique combination of AIR and GEL technology effectively provides pressure relief and comfort to the high risk user. The action of the air in the cushion enhances capillary function in the skin, maintaining blood flow and oxygen uptake. When pressure is applied, air circulates to each of the cells, distributing weight evenly, thus reducing high pressure concentrations around bony prominences.

- 1 air chamber
- 3 layer polyurethane film
- Pressure Relieving Silicone TruGel Surface
- 34 Air-Filled Interconnected Cells

Also available in the G2 with 48 Air-Filled Interconnected Cells and 2 separate air chambers. Contact Active Medical Supplies for more detail!

Relax GelCell G2 Cushion

The Gelcell G2 consists of 48 interconnected cells, each covered with a layer of pressure relieving silicone-based TruGel

- Enhances pressure re-distribution
- Minimises the effects of shear forces
- Has a positive effect on skin temperature
- Silicone gel is impact dampening

The cells are contoured towards the rear, promoting good postural tilt.

When pressure is applied, air circulates to each cell in both chambers distributing weight evenly. This process reduces high pressure concentrations around bony prominences, preventing skin breakdown.

The Gelcell is made up of 3 layers of polyurethane (PU) film making it highly durable and difficult to puncture. It has also easy to wipe clean and suitable for incontinent users.

- Suitable for high risk clients
- Easy to wipe clean and is suitable for incontinent users
- Lightweight and easy to transport
- Easy to inflate and adjust
- Includes outer cover and pump



Specifications:

Product Code	35322
Width	450 mm
Length	400 mm
Thickness	50 mm
SWL	130 kg

Sizes Available:

Product Code	35320	35323
Width	400 mm	450 mm
Length	400 mm	450 mm
Thickness	50 mm	50 mm
SWL	130 kg	130 kg

Specifications:

Product Code	35324
Width	360 mm
Length	360 mm
Thickness	40 mm
SWL	130 kg

Sizes Available:

Product Code	35327
Width	450 mm
Length	400 mm
Thickness	40 mm
SWL	130 kg



Starlock Cushion

The StarLock Cushion is a prime wheelchair cushion offering ultimate positioning, unparalleled stability and the highest level of therapeutic prevention and treatment of pressure ulcers.

Instead of locking air in sections, the StarLock allows the user to stop the flow of air between each cell in the cushion. The StarLock can be custom fit to the user's unique shape by locking the transfer of air between each cell in the cushion.

The StarLock also allows you to remove the air from individual cells and keep those cells from re-inflating once locked, creating a void, or no contact zone, where it is needed the most.

- Isolate individual air cells and customize 'no contact areas with cell lock technology
- 24 sizes available
- Complete with fitted cover, hand pump and repairs kit
- Made with high-quality, flame resistant, latex-free neoprene rubber

Standard Contour Star Cushion

The Standard Contour Cushion is Star Cushion Products premier cushion for the prevention and treatment of pressure ulcers. Flexible, interconnected 4' air cells allow for deep immersion into the cushion providing constant low pressure across the contact seating surface. In addition, the cell moves and shapes to the contact area providing a low shear, low friction surface.

- Comes with fitted cover, hand pump and repairs kit
- Air adjustable 4' cells offer deep immersion to provide a high level of pressure reduction for tissue protection and treatment
- Made with high quality, flame resistant, latex free neoprene rubber
- Also available in 50mm high

Additional sizes available



Specifications:

Product Code	3324 - 1515
Width	390 mm
Length	390 mm
Thickness	100 mm
SWL	Unlimited

Sizes Available:

Product Code	3324 - 1715	3324 - 1717	3324 - 1917	3324 - 1919
Width	440 mm	440 mm	490 mm	490 mm
Length	390 mm	440 mm	440 mm	490 mm
Thickness	10 mm	10 mm	10 mm	10 mm



Specifications:

Product Code	3322 - 1717
Width	440 mm
Length	440 mm
Thickness	100 mm
SWL	Unlimited

Sizes Available:

Product Code	3322 - 1715	3322 - 1917
Width	440 mm	490 mm
Length	390 mm	440 mm
Thickness	100 mm	100 mm

NDIS

NDIS



Stabil-Air Star Cushion

The Stabil-Air is a perfect combination of air and foam; air enables pressure equalization with foam bringing contoured stability. We do it through utilization of foam supporters internal to the air cells, further enhancing the dynamic weight distribution benefits of a cellular air cushion. It means reinforcing and shaping selected cells with internal support systems.

- Prevention and healing up to stage 4 pressure injury
- Made with high quality, flame resistant, latex free neoprene rubber
- Comes with fitted cover, hand pump and repairs kit
- Internal foam can be custom contoured to unique patient needs
- Internal foam is contoured depending on cushion height
- Also available in 75mm high

Galaxy Star Cushion

The Galaxy is a flame resistant neoprene cushion that combines multi-cellular air technology with a patented, contouring design, fitting the anatomical form. These detailed, contouring features, offer several benefits to aid in the patients' stability and posture, while greatly enhancing a true therapeutic advantage.

- Pelvic Supports:** A supporting area that gently adducts the trochanters creating a hydrostatic environment around the sitting tissues.
- Adductors:** Gently pushes femurs to the midline of the cushion.
- Coccyx/Perineum Relief:** Voided areas to keep excess pressures away from problematic areas.
- Pre-Ischial Ridge:** Posterior sloped air cells prevents forward migration.
- Abductors:** Built-in pommel for abduction of the knees to work in concert with the femur troughs.
- Femur Troughs:** A sculpted region to help position the thighs in a medial-lateral placement
- Prevention and healing up to stage 4 Pressure Injury
- Comes with fitted cover, hand pump and repairs kit

SAFE WORKING LOAD UNLIMITED HIGH RISK PREVENTION P TREATMENT T

Specifications:

Product Code	3320 - 1717
Width	440 mm
Length	440 mm
Thickness	100 mm
SWL	Unlimited

Sizes Available:

Product Code	3320 - 1715	3320 - 1917
Width	440 mm	490 mm
Length	390 mm	440 mm
Thickness	100 mm	100 mm

SAFE WORKING LOAD UNLIMITED HIGH RISK PREVENTION P TREATMENT T

Specifications:

Product Code	3334 - 1818
Width	460 mm
Length	460 mm
Thickness	127 mm
SWL	Unlimited

Sizes Available:

Product Code	3334 - 2020
Width	510 mm
Length	510 mm
Thickness	127 mm

02.

Cushion Matrix

The cushion matrix shows the relationship between our cushions and different variables (for example; weight limit, skin integrity or sizes) in a data set in grid format. Essentially, the cushion matrix is a quick guide which allows the reader to compare each cushion to help find the appropriate cushion for either their needs or patient needs in a time efficient matter.

Cushion Matrix

	Weight Capacity	Cushion Profile Height	Size Range	
			Width - inches	Width - mm
3210 Budget 	100 kg	50 mm		
SAC 	160 kg	25 mm		
Action Pilot 	Unlimited	30 mm		
Action Professional 	Unlimited	45 mm		
Relax Easy Low Profile 	130 kg	50 mm	14 x 16" to 22 x 18"	360 x 400mm to 550 x 450mm
Relax Easy High Profile 	190 kg	85 mm	16 x 16" to 20 x 20"	400 x 400mm to 500 x 500mm
Duogel 	130 kg	85 mm	14 x 14" to 24 x 18"	360 x 360mm to 600 x 450mm
Duo Profile 	130 kg	85 mm	16 x 16" to 18 x 18"	400 x 400mm to 450 x 450mm
Gel Cell 	130 kg	40 mm	15 x 14" to 20 x 18"	380 x 360mm to 500 x 450mm

PI Rating	Thigh Support	Contoured seat	Skin Integrity			Boney Prominence	Suitable for Lounge chairs
			Skin Intact	Skin Broken			
Prevention	✗	✗	✓	✗	✗	✗	
Stage 2	✗	✗	✓	✗	✗	✓	
Stage 3	✓	✗	✓	✗	✗	✓	
Stage 4	✓	✗	✓	✓	✓	✗	
Stage 2	✗	✗	✓	✗	✓	✓	
Stage 3	✗	✗	✓	✓	✓	✗	
Stage 3	✓	✓	✓	✓	✓	✗	
Stage 2	✓	✓	✓	✗	✓	✗	
Stage 4	✗	✓	✓	✓	✓	✓	

03.

Mattresses

LOW RISK

Visko Mattress	21
Silicone Fibre Mattress Overlay	21


MEDIUM RISK

Staebel Zonecare Mattress	22
Optimal Cura Mattress	22
Bariatric Foam Mattress	23

HIGH RISK

Soft EuroFlex Hygiene Mattress	23
Optimal 5zon Mattress	24
Optimal 5zon Mattress with Evac Handles	24
CuroCell® Area Zone Mattress	25
CuroCell® SAM BAS PRO Air Mattress	25

CuroCell® A4 CX16 Air Mattress	26
CuroCell® A4 CX20 Air Mattress	26
CuroCell iA Automatic Ci20 Mattress	27
CuroCell iA Automatic Ci17 Mattress	27
CuroCell® IQ CX15 Air Mattress	28
CuroCell® UNO Air Mattress	28
CuroCell® A4 Cirrus Air Mattress	29
Metzeler Ultratube Firm Feel Mattress	30
PrimaCare 10	30
Precioso	31
Air2Care	31
CliniCare 100 HF	32
OptiCare Intergrated Mattress	32

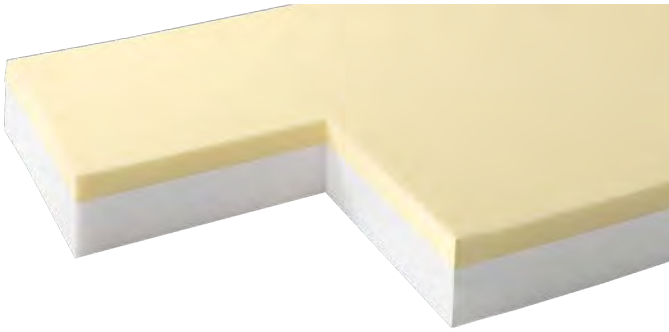

Visko Mattress

The Viskomatt includes a double layer foam mattress including a 4cm top layer of visco elastic foam, suitable for those who are still mobile but spending a long time in bed, providing optimum pressure distribution across all body area in contact with the mattress.

The mattress has moisture protection properties mean the cover has a high resistance to bacteria, viruses and fungi.

- Comes with two way stretch PU waterproof cover
- High Density foam base for stability

Silicone Fibre Mattress Overlay

This overlay is made of millions of silicone coated cylindrical hollow cored fibres that provide excellent pressure relief while the silicone coating acts as a lubricant to eliminate friction and reduce shear forces. These fibres do not absorb perspiration or urine, making this an excellent product for geriatric and palliative care with elastic corner retention straps.

- Suitable for low risk patients
- Soft breathable cotton cover on top, PVC on bottom
- Machine washable



Specifications:

Product Code	3755
Width	850 mm
Length	2000 mm
Thickness	140 mm
SWL	120 kg


Sizes Available:

Product Code	3755W
Width	1000 mm
Length	2000 mm
Thickness	140 mm
SWL	120 kg



Specifications:

Product Code	3740
Width	920 mm
Length	1900 mm

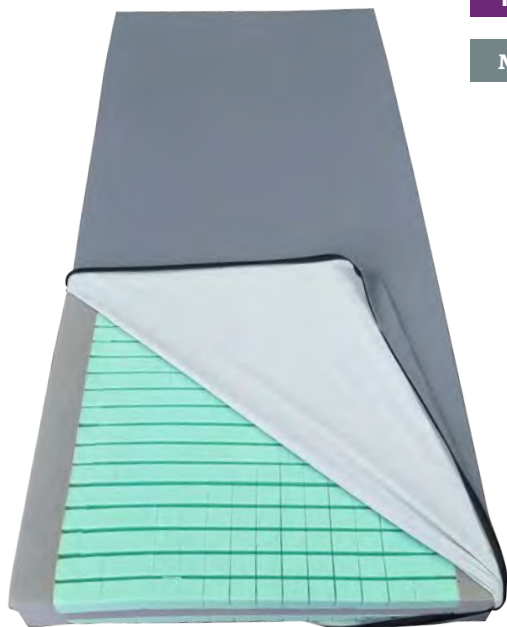

NDIS

MASS


NDIS



Staebel Zonecare Mattress

The Zonecare foam mattress is a high quality foam mattress featuring a profile-cut pressure relieving foam insert. It provides pressure reduction and assists in the treatment of pressure ulcers, whilst enhancing patient comfort levels in homecare and residential care environments. The Zonecare's specialised inner cover aids in shear reduction on a profiling bed.

- Unique re-reinforced edges designed with longevity in mind

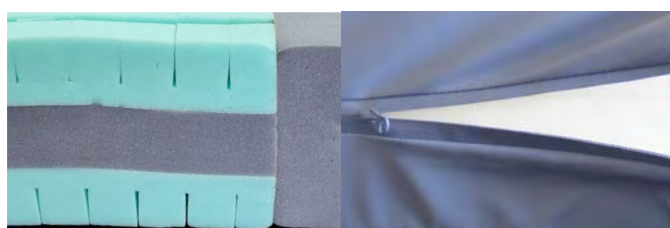

Specifications:

Product Code	3728Z
Width	850 mm
Length	1980 mm
Thickness	150 mm
SWL	150 kg

Sizes Available:

Product Code	3729Z
Width	950 mm
Length	1980 mm
SWL	150 kg

Optimal Cura Mattress

Optimal Cura is used as an aid in the prevention and treatment of pressure ulcers up to category 1.

Optimal Cura is a two layer mattress where the top layer is made from visco-elastic foam and the bottom layer is made from HR foam. The special design of the mattress and the materials gives very good pressure equalisation. The two layers allow the body to sink into the material, and thereby lower the pressure against the body.

- Easy to clean hygienic cover is vapour permeable, which reduces the risk of skin maceration
- Specially designed to minimise the effects of shear force
- Highly elastic, high quality foam with effective pressure equalisation that provides good comfort.

For more sizes please contact the Active Medical team on 1800 267 267




Specifications:

Product Code	3870
Width	850 mm
Length	2000 mm
Thickness	130 mm
SWL	130 kg

Sizes Available:

Product Code	3872
Width	1000 mm
Length	2000 mm
SWL	130 kg

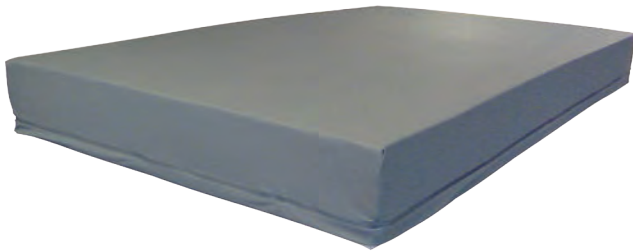

Bariatric Foam Mattress

The Bariatric foam pressure mattress offers the Bariatric patient comfort and pressure relief and has been designed to last. Manufactured from high quality foam of a higher density than standard mattresses and covered with two way stretch breathable vinyl, it is practical and affordable.

- Suitable for medium risk patients
- Full zip vinyl cover with flap over zip for infection control
- Cover made from heavy duty vinyl base, multi stretch vinyl top for ultimate patient comfort
- Does not require turning

Staebel Soft EuroFlex Hygiene Mattress

Soft Euroflex is a hygiene mattress developed as an aid in the treatment of pressure ulcers up to category 3. The mattress is designed so that the entire mattress has relief for protruding body parts such as shoulders, heels and elbows. Soft Euroflex is made up of three different foam cores for good comfort and pressure relief. The mattress, , has a laminated surface layer of Lentex® and is completely devoid of seams and zippers, which provides increased conditions for good hygiene.

- Top White surface layer also makes it easier for staff to detect contamination on the surface of the mattress, which can be quickly cleaned by wiping and disinfecting.
- The surface layer is liquid-repellent and tight. Active Membrane technology allows the mattress to breathe while reducing the risk of liquid ingress.
- Approved for bed and mattress washing
- Lentex mattresses have been tested and approved completely without remark for use with bed and mattress cleaners
- The foam core consists of three layers with carefully selected foam types to provide good comfort and good pressure-relieving ability.
- Lentex® Soft Euroflex has been tested by Swerea IVF according to measurement method SS 876 00 13
- Preventions of pressure sores up to and including category 3.
- Safe Working load of 150 kg



Specifications:

Product Code	3794
Width	1000 mm
Length	2000 mm
Thickness	200 mm
SWL	360 kg

Also Available:

Product Code	3791	3792	3793
Width	1050 mm	1180 mm	1350 mm
Length	1960 mm	1960 mm	1960 mm

Specifications:

Product Code	3271361
Width	850 mm
Length	2000 mm
Thickness	130 mm
SWL	150 kg

Sizes Available:

Product Code	3271363
Width	1050 mm
Length	2000 mm
SWL	150 kg



NDIS

MASS

NDIS



Optimal 5zon Mattress

The Optimal 5zon® is ideal for all healthcare environments and is used as an aid in the prevention and treatment of pressure ulcers up to and including category 2, and for pain therapy. It is a soft and flexible mattress replacement that is placed directly on the bed and replaces both the standard mattress and mattress pad.

- Highly elastic, high-quality material and a unique zone and module system gives effective pressure equalisation.
- Specially designed to minimise shear forces
- Simple and effective cleaning, cover has a very well documented hygiene standard
- The hygienic cover is vapour permeable, which reduces the risk of skin maceration
- The patented sit function DuoCore® is designed to give support when sitting and also to provide good comfort when lying
- Heel Protection provides comfort for relief of heel. The user can move the leg and make position changes. Comes with a removable liquid resistant hygiene cover.

Optimal 5zon Mattress with Evac Handles

The Optimal 5zon® is ideal for all healthcare environments and is used as an aid in the prevention and treatment of pressure ulcers up to and including category 2, and for pain therapy. It is a soft and flexible mattress replacement that is placed directly on the bed and replaces both the standard mattress and mattress pad.

- Highly elastic, high-quality material and a unique zone and module system gives effective pressure equalisation.
- Specially designed to minimise shear forces
- Simple and effective cleaning, cover has a very well documented hygiene standard
- Evacuation handles on the side of the mattress
- The hygienic cover is vapour permeable, which reduces the risk of skin maceration
- The patented sit function DuoCore® is designed to give support when sitting and also to provide good comfort when lying
- Heel Protection provides comfort for relief of heel. The user can move the leg and make position changes. Comes with a removable liquid resistant hygiene cover.



Specifications:

Product Code	3771
Width	850 mm
Length	2000 mm
Thickness	140 mm
SWL	180 kg

Sizes Available:

Product Code	3778	3775
Width	1000 mm	1050 mm
Length	2000 mm	2000 mm
SWL	180 kg	180 kg



Specifications:

Product Code	3770
Width	800 mm
Length	2000 mm
Thickness	140 mm
SWL	180 kg

NDIS

MASS



CuroCell® Area Zone Mattress

The AREA Zone is a mattress that combines air with foam for effective pressure redistribution and great comfort. The mattress is self-adjusting and non-powered. It is a flexible mattress replacement system that can be used in all kinds of healthcare settings, easy to use and handle. The mattress system is used as an aid in prevention and treatment of pressure injuries up to category 3.

- Zone castellations in top layer foam provides reduction of shearing forces
- Handle for transportation
- Flexible and self-adjusting internal air cells, no pump required
- Integrated heel function
- Reinforced cover with welded seams
- Non-powered mattress replacement system

*Clinical Test Report Available on request

NDIS

MASS



CuroCell® SAM BAS PRO Air Mattress

SAM BAS PRO is used as an aid in the prevention and treatment of pressure ulcers up to category 3, and for pain therapy. It is a mattress replacement system with an integrated foam core. The mattress is made of very soft and durable materials. The included support pump is easy to use and very quiet. The mattress comes with CPR (rapid deflation of air) and has a built-in heel function.

- This mattress is supplied with a four way stretch hygienic outer cover which conforms to relevant standards for hygiene, quality and safety.
- Silent, energy-saving pump
- Soft, elastic, vapour permeable hygienic cover
- Easy to handle and maintain
- Integrated heel function

*Clinical Test Report Available on request

For more sizes please contact the Active Medical team on 1800 267 267



Specifications:

Product Code	3678Z
Width	850 mm
Length	2000 mm
Thickness	180 mm
SWL	230 kg

Sizes Available:

Product Code	3686Z
Width	1000 mm
Length	2000 mm
SWL	230 kg



Specifications:

Product Code	3683
Width	850 mm
Length	2000 mm
Thickness	160 mm
SWL	250 kg

Sizes Available:

Product Code	3684	3685
Width	1050 mm	1000 mm
Length	2000 mm	2000 mm
SWL	250 kg	250 kg



NDIS

MASS



CuroCell® A4 CX16 Air Mattress

The CuroCell® A4 CX16 provides simplicity together with user comfort while effectively aiding in prevention and treatment of pressure injuries up to and including category 4, and unclassified pressure injuries. It is designed for excellent user comfort, reduction of shear forces and improvement of the microclimate.

The air mattress is 16cm thick with a foam base for added comfort. The alternation of the cells on the CX16 is every second cell down ensuring 50% of the body is being supported at one time. The pump is intelligent and fully automatic without need for manual settings. It identifies the user and automatically manages all settings to optimise the pressure according to the user's needs.

All mattresses are supplied with a four way stretch hygienic cover which meets standards for hygiene, quality and safety. All covers are manufactured from materials which reduce the risk of skin maceration.

- Unique CPR-function with transport function
- Heel Function
- Silent running and vibration free
- Three program modes: Pulsating, Alternating and CLP (constant low pressure)
- Individually replaceable cells with safety mattress in HR-foam base
- Audible and visual notifications
- Smart cable management

For more sizes please contact the Active Medical team on 1800 267 267



Specifications:

Product Code	3626
Width	850 mm
Length	2000 mm
Thickness	160 mm
SWL	200 kg

Sizes Available:

Product Code	3627	3628
Width	1000 mm	1050 mm
Length	2000 mm	2000 mm
SWL	200 kg	200 kg



NDIS

MASS



CuroCell® A4 CX20 Air Mattress

The CuroCell® A4 CX20 provides simplicity together with user comfort while effectively aiding in prevention and treatment of pressure injuries up to category 4, and unclassified pressure ulcers. It provides features for excellent user comfort, reduction of shear forces and improvement of the microclimate.

The CX20 mattress is 20 cm thick and the alternation of the cells is every third cell down ensuring 75% of the body is being supported at one time. The pump is intelligent and fully automatic without need for manual settings. It identifies the user and automatically manages all settings to optimise the pressure according to the user's needs.

All mattresses are supplied with a four way stretch hygienic cover which meets standards for hygiene, quality and safety. All covers are manufactured from materials which reduce the risk of skin maceration.

- Unique CPR-function with transport function
- Heel Function
- Audible and visual notifications
- Silent running and vibration free
- Three program modes: Pulsating, Alternating and CLP (constant low pressure)
- Individually replaceable air cells with safety mattress
- Smart cable management

For more sizes please contact the Active Medical team on 1800 267 267



Specifications:

Product Code	3641
Width	850 mm
Length	2000 mm
Thickness	200 mm
SWL	250 kg

Sizes Available:

Product Code	3642	3643
Width	1000 mm	1050 mm
Length	2000 mm	2000 mm
SWL	250 kg	250 kg

NDIS

NEW

AUTO



NEW

AUTO



NDIS

CuroCell iA Automatic Ci20 Mattress

With Gentle Alternating Low Pressure (GALP), the CuroCell iA® system makes the step from mechanical to digital, from manual to autonomus, and from on-off comfort valves to a fine tuning of total support that we've never been able to deliver before in this segment. GALP has been developed to enable a system to sense, learn and take action whenever the high pressure that can cause ulcers/injuries arises.

CuroCell iA® uses artificial intelligence to constantly monitor pressure across the patient's contact area and calculates its response. The system identifies the patient's weight and height and can detect changes, such as when the patient adopts a new position.

All mattresses are supplied with a four way stretch hygienic cover which meets standards for hygiene, quality and safety. All covers are manufactured from materials which reduce the risk of skin maceration.

- Gentle Alternating Low Pressure (GALP)
- Constant Low Pressure (CLP)
- Panel Lock
- Maximum pressure (Caring mode)
- Information signals
- CPR connection with transport function
- Air Flow Control - saves energy and creates an even temperature in the support surface
- Low noise-and vibration level
- Individual and replacement cell construction with safety mattress
- Heel Function

CuroCell iA Automatic Ci17 Mattress

With Gentle Alternating Low Pressure (GALP), the CuroCell iA® system makes the step from mechanical to digital, from manual to autonomus, and from on-off comfort valves to a fine tuning of total support that we've never been able to deliver before in this segment. GALP has been developed to enable a system to sense, learn and take action whenever the high pressure that can cause ulcers/injuries arises.

CuroCell iA® uses artificial intelligence to constantly monitor pressure across the patient's contact area and calculates its response. The system identifies the patient's weight and height and can detect changes, such as when the patient adopts a new position.

All mattresses are supplied with a four way stretch hygienic cover which meets standards for hygiene, quality and safety. All covers are manufactured from materials which reduce the risk of skin maceration.

- Gentle Alternating Low Pressure (GALP)
- Constant Low Pressure (CLP)
- Panel Lock
- Maximum pressure (Caring mode)
- Information signals
- CPR connection with transport function
- Air Flow Control - saves energy and creates an even temperature in the support surface
- Low noise-and vibration level
- Individual and replacement cell construction with safety mattress
- Heel Function

SAFE WORKING LOAD
240KG

HIGH RISK

PREVENTION
PCATEGORY
4

Specifications:

Product Code	3935
Width	850 mm
Length	2000 mm
Thickness	200 mm
SWL	240 kg

Also Available:

Product Code	3936	3937	3938
Width	900 mm	1000 mm	1050 mm
Length	2000 mm	2000 mm	2000 mm

SAFE WORKING LOAD
210KG

HIGH RISK

PREVENTION
PCATEGORY
4

Specifications:

Product Code	3928
Width	850 mm
Length	2000 mm
Thickness	170 mm
SWL	210 kg

Also Available:

Product Code	3929	3930	3931
Width	900 mm	1000 mm	1050 mm
Length	2000 mm	2000 mm	2000 mm



NDIS



NDIS

MASS




CuroCell® IQ CX15 Air Mattress

The CuroCell® IQ CX15 provides simplicity together with user comfort while effectively aiding in the prevention and treatment of pressure ulcers/pressure injuries up to and including category 4, and unclassified pressure ulcers. It is designed for excellent user comfort, reduction of shear forces and improvement of the microclimate.

The CuroCell® IQ CX15 identifies the user and automatically adapts itself to the individual. The pump is preset to our newly developed Pulsating mode, which offers a high level of comfort and effective pressure relief.

- Audible and visual notifications
- Silent running and vibration free
- Newly developed pulsating mode
- Unique CPR-function with transport function
- Heel Function
- Smart cable management
- Individually replaceable cells with safety mattress

For more sizes please contact the Active Medical team on 1800 267 267

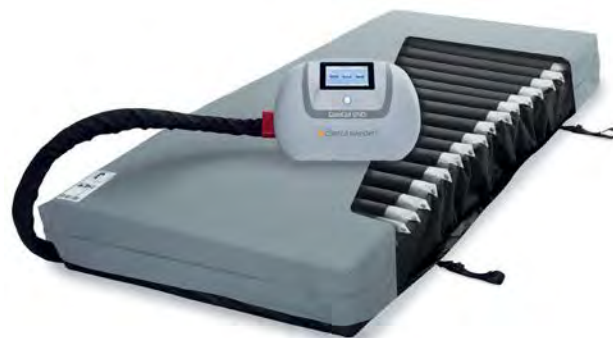

CuroCell® UNO Air Mattress

The CuroCell UNO® is used as an aid in the prevention and treatment of pressure injuries up to and including category 4. It is an alternating pressure mattress. It is equipped with a large touch screen panel with a user-friendly interface. The fully automated pump offers automatic adjustment, providing optimal pressure relief and comfort for each individual user while aiding simple use. The mattress has a four way stretch hygienic outer cover which conforms to relevant standards for hygiene, quality and safety. All covers are manufactured from modern, recyclable materials and are vapour permeable, which reduce the risk of skin maceration.

- Fully automated pump that adapts individually with large colour touch display
- Cell on Cell construction with safety mattress
- Unique design for easy handling and usage
- Silent running and vibration free
- Heel Function
- Smart cable management
- 3 modes: Alternating, Pulsating and Static



Specifications:

Product Code	3847
Width	850 mm
Length	2000 mm
Thickness	150 mm
SWL	200 kg

Sizes Available:

Product Code	3848
Width	1000 mm
Length	2000 mm
SWL	200 kg



Specifications:

Product Code	3651
Width	850 mm
Length	2000 mm
Thickness	200 mm
SWL	250 kg

Sizes Available:

Product Code	3652
Width	1050 mm
Length	2000 mm
SWL	250 kg



NDIS



CuroCell® A4 Cirrus Air Mattress

CuroCell® A4 Cirrus is used as an aid in the prevention and treatment of pressure ulcers up to and including category 4 and unclassified pressure ulcers/pressure injuries⁽¹⁾. It is a replacement system suitable for all healthcare environments.

CuroCell® A4 Cirrus combines air with a soft foam surface and provides good user comfort and pressure redistribution⁽²⁾. The added heel function offers pressure redistribution on the heels.

- Automatic control unit which adapts individual to the user
- User friendly interface
- Audible and visual notifications
- Three program modes: pulsating mode, alternating mode and constant low pressure
- Silent running and vibration free
- Air flow control; Constant temperature, good user comfort, energy reduction
- Soft foam surface
- Smart cable management
- Easy handling and cleaning
- Heel slope, for pressure redistribution of heels
- CPR-function with transport function



Specifications:

Product Code	3841
Width	850 mm
Length	2000 mm
Thickness	210 mm
SWL	230 kg

Sizes Available:

Product Code	3843	3844
Width	1000 mm	1050 mm
Length	2000 mm	2000 mm
SWL	230 kg	230 kg

Case Study

History

- Female, born 1932
- Operated breast cancer 2003, vertebral compression fractures, hip osteoarthritis with chronic pain
- Arterial fibrillation, hypertension and bilateral leg ulcers.
- Pressure ulcer category 4 in sacrum 7x2 cm with fistulas, R/T ongoing treatment leg ulcer.
- Living at home with help from home health care.
- Mentally clear, no incontinence, normal diet, sitting up during meals.
- Taking regular pain relief, Oxykontin & Oxynorm.
- Points 23, modified Norton scale.
- BMI 21




Action

- 2010-12-31 prescribed CuroCell Cirrus, alternating pressure mattress with foam core, already have a pressure relieving wheelchair cushion.
- Can manage repositioning by herself, aided move from bed to wheelchair. Wounds redressed daily. Daily supplements of nutrient drinks.
- 2011-01-21 admitted to hospital due to the leg ulcers, mattress goes with the patient.
- 2011-02-18 left leg amputated thigh bone After hospital stay the patient is accommodated short-term at a home and is later placed in a permanent nursing home, mattress goes with the patient. Wounds is decreasing in size.
- 2012-04-04 Wounds completely healed



NDIS

NDIS

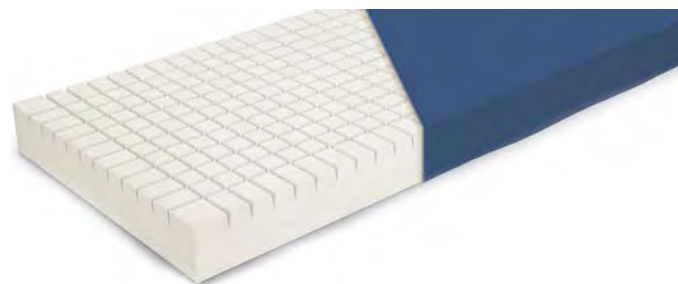

Metzeler Ultratube Firm Feel Mattress

The heart of the Metzeler Ultratube Mattress is the newly developed Rubex Tubes System. The cold formed foam tubes in the mattress form a vertical tube system for perfect air circulation. The tubes, which differ in colour according to their degree of hardness, are integrated into the complete surface. This guarantees the user an ergonomically perfect lying position. Up to 200 litres of perspiration a year leave the body and go into the mattress. The Rubex Tubes System dispels this moisture through the vertical air channels.

- The tubes make up a system of s for perfect air circulation during sleep. The tubes of varying degrees of firmness are arranged to provide an ergonomically ideal sleeping posture reducing the risk of injury caused by shear and friction.
- Reinforced edge makes it easy to get in and out of bed and helps prevent roll-out.
- Standard with a bi-elastic, two-way stretch, waterproof, cover, designed to reduce shear and friction forces.
- Stain resistant, fire-resistant, waterproof, antiallergenic, antimicrobial, hygienically clean. Can be machine washed up to 95°C. Antibacterial treated.
- The Rubex Tubes System is symmetrically structured so that you always lie straight.
- Fully washable in temperatures up to 60°C



Specifications:

Product Code	3892
Width	1050 mm
Length	2000 mm
Thickness	180 mm
SWL	250 kg

Also Available:

Metzeler Ultratube Soft Feel Mattress 1050x2000mm	
Product Code	3893

PrimaCare 10

PrimaCare 10 is made of a monolithic foam core. The upper surface of the mattress has transverse and longitudinal profiling evenly distributed over its entire surface. Due to this profiling, there is better air circulation between the patient's skin and the base. They also help to optimally distribute pressure.

Increased density of profiling at the head and foot of the mattress: Due to the increased density of profiling in the head and foot area, the mattress ensures better comfort and distribution of pressure between the patient and the mattress.

Mattress cover made of vapour-permeable material: Helps reduce patient perspiration and also reduces skin moisture.



Specifications:

Product Code	3760
Width	860 mm
Length	2000 mm
Thickness	140 mm
SWL	200 kg

Sizes Available:

Product Code	3761
Width	1000 mm
Length	2000 mm
Thickness	140 mm
SWL	200 kg



Precioso

The Precioso is a 24-hour system for wound care management. The system comes as an overlay or mattress replacement and can be used with a specific cushion. The Precioso will significantly contribute towards the healing and prevention of all grades of pressure ulcers. Adequate for use in general wards in hospitals and nursing homes.

The Precioso is fitted with 4 different pressure sensors and a microprocessor controlling the precise cell pressure. For patient comfort it is possible to increase or decrease the pressure by adjustment the levels on the membrane panel.

Alternating Pressure Therapy (3 cell technology, zero pressure, 7.5 min. cycle) is the fundamental therapeutic mode that simulates natural body movement and therefore helps to prevent and treat pressure ulcers.

Air2Care

The Air2Care alternating pressure system has been designed for the prevention and management of pressure ulcers in those at low to moderate risk. The Air2Care is suitable for use in various healthcare settings including long term care. For convenience the system is available as either a mattress replacement or as an overlay giving the flexibility to choose a solution suitable for the patient's and caregiver's individual requirements. Two options of system control unit are available. The addition of an active cushion to the Air2Care system gives 24 hour prevention and management of pressure ulcers which is simple to use and reliable.



Specifications:

Product Code	36021
Width	860 mm
Length	2000 mm
Thickness	170 mm
SWL	160 kg



Specifications:

Product Code	36010
Width	860 mm
Length	2000 mm
Thickness	170 mm
SWL	200 kg

NDIS



CliniCare 100 HF

The CliniCare 100 HF hybrid mattress is comfortable and highly efficient for the prevention of pressure injuries. It uses a combination of an active system and passive mattress to provide many benefits across the spectrum of hospital care.

When a patient's condition requires active therapy, the mattress can be connected to a compressor for alternating pressure therapy in a 2 cell cycle to enhance and optimize pressure redistribution and pressure injury prevention.

The mattress provides pressure redistribution in non-powered mode, air moves from cell to cell as reaction to patient's body movement and weight.

A viscoelastic foam layers adapts to the body to provide comfort and proper support without exerting pressure on exposed area. The mattress has a 7° slope that helps to reduce pressure in vulnerable heel area.

The system inside the mattress keep blowing air and the mattress cover is automatically dried. In the context of pressure injuries, microclimate usually refers to skin temperature and moisture conditions at the skin support surface interface.

- With microclimate



Specifications:

Product Code	3817	3818
Width	860 mm	860 mm
Length	2000 mm	2080 mm
Thickness	200 mm	200 mm
SWL	230 kg	230 kg

NDIS



OptiCare Intergrated Mattress

New smart technology for the prevention of pressure ulcers. The system detects when the patient has entered the bed. This automatically initiates the optimisation process, thereby saving the time normally required for nurses to set the mattress parameters. The mattress settings are automatically optimised for each patient based on the patient's weight, height, and even the patient's position on the mattress.

- The mattress is equipped with a whole range of sophisticated functions which enhance the therapeutic and preventive effect of the optimised pressure.
- OptiCare mattress is compatible with Multicare and Eleganza 5 beds. The controls for both the mattress and bed are integrated into a single control panel built into the bed. The customisable user interface provides the clinical staff with the option to modify the controls on the interface and allows control of only the required features.
- In the context of pressure ulcers, microclimate usually refers to skin temperature and moisture conditions at the skin-support surface interface. OptiCare helps maintain the natural thermoregulation by circulating air through the mattress surface.
- Helps control the heat and humidity of the patient's skin, which is an important part of maintaining tissue viability.
- The heel section has a 7-degree slope and is lined with a softer foam that gives the patient above-standard comfort. This significantly reduces the risk of heel pressure ulcers.
- Integrated systems to be used with Multicare and Eleganza 5 beds only.



Specifications:

Product Code	74200 - 1
Width	860 mm
Length	2140 mm
Thickness	220 mm
SWL	250 kg

04.

Mattress Matrix

The cushion matrix shows the relationship between our cushions and different variables (for example; weight limit, skin integrity or sizes) in a data set in grid format. Essentially, the cushion matrix is a quick guide which allows the reader to compare each cushion to help find the appropriate cushion for either their needs or patient needs in a time efficient matter.

Mattress Matrix

Waterlow scale / Risk Score

Mattress Name	Size	Comments	Waterlow scale / Risk Score			
			At Risk (10+)	High (15+)	Very High (20+)	Skin Intact
Zonecare Pressure Care Foam Mattress 	Single	Standard mattress for all replacements	✓			✓
	King Single					
Optimal 5zon Pressure Care Mattress 	Single	Highly elastic, high quality material and a unique zone and module system gives effective pressure. Specifically designed to minimise the effects of shear forces. Duocore is designed to give support when sitting and provide good comfort when lying. The cover is vapour permeable, which reduces the risk of skin maceration. The covers are easy to clean and handle and has a very well documented hygienic standard.	✓			
	King Single					
Curocell S.A.M F16 Pro Static Air Mattress Replacement & Pump 	Single	A protective / preventative strategy for residents with intact skin but identified as at risk of skin integrity issues of pressure injuries. Beautifully comfortable due to the static nature plus it comes with the comfort cover so you do not create Cat 1-3 pressure injuries.		✓	✓	
	King Single					
Curocell A4 CX16 	Single	Cat 1- 4 Pressure Injuries. CX16 is a 3 cell cycle, opposed to the regular 2 cell. 2 cells up, 1 down means the body has 66% support not 50%, however the pulsating is the best program for the greatest comfort. Great for persons who need the pressure care, but don't like the sensation. Excellent for where high incontinence issues are present.		✓	✓	
	King Single					
CuroCell Cirrus 	Single	Cat 1- 4 Pressure Injuries. Uniquely designed to help in pain management, for end of life care. Thick top layer of foam to help create sensory recognition of the user with pulsation mode cells below that create a wave like relief effect opposed to the traditional alternation. Heel slope, silent pump with sit position function, and option to adjust the timing of the cell change.		✓	✓	
	King Single					

Current Pressure			May suit people with particular needs:				
Skin Broken	Underweight	Bariatric	Client In Pain	Bed Mobility		Bed Transfers	Heels at risk
				Independent	Assisted / Static		
				✓		Firm Transfers	
	(~30KG)		✓	✓	✓	YES - the Duo Core adds stability	Zero pressure ability
	(~30KG)		✓	May be difficult	✓	Staebel as it is a static mattress	Inbuilt heel slope
	(~1KG)		✓	May be difficult	✓	Max function	Inbuilt heel slope
	(~1KG)		✓	May be difficult	✓	Foam surround	Inbuilt heel slope

05.

Ancillary

Staelbel Careturner	37	Positioning Roll	41
RollStop Coverlet	38	Molly Heel Protection	41
Rollstop Premium Foam Mattress	38	Staelbel Positioning & Support Wedge	42
Rollstop KS Premium Foam Mattress	38	Action Polymer Pad	42
Staelbel EvacSmart Evacuation Mat	39	Action Gel Segment Bed Pad	43
Staelbel Bed Rail Protector	39	Action Gel Elbow / Heel Protector	43
Back Pillow	40	Staelbel In Bed Turner 4-way with Handles	44
Multi Pillow	40	Staelbel In Bed Turner Base Sheet	45
Positioning Pillow	40		



Careturner works in conjunction with the beds existing movements.



Sides can be angled up to 60-80°. Sling can be used by 1 caregiver.



The programmable control box.

Staebel Careturner

Careturner is a unique and innovative medical device designed to reposition patients automatically and to support care staff with daily mobilization.

Careturner turns the patient gently from side to side throughout the night and provides pressure distribution to prevent development of pressure ulcers.

Careturner is optimised for single caregiver working environments. It optimizes the working environment and saves resources.

Patients

- Prevents development of pressure ulcers
- Better comfort and less pain
- Undisturbed sleep and more peace of mind

Caregiver

- Fewer heavy lifts
- Better working conditions
- Reduced attrition of care staff

Organisation

- Automated turning of patients saves resources
- More efficient mobilization (from 2 to 1 caregiver)
- Reduced health care spendings (pay back in 3-4 months)
-

Careturner turns the patient gently from side to side

Careturner is easy to install on a medical bed, and works in conjunction with the beds existing movements.

Manual mode: With a push on a button it turns the enduser from one side to the other to support the care staff in daily moving and handling procedures. In manual mode the wings can elevate up to 80°.

Automatic mode: The control box has a preprogrammed setting and option for two additional programs. Careturner operates very silent and without any vibration allowing the patient to sleep undisturbed. In automatic mode the wings can elevate up to 30°.



Specifications:

Product Code	78609
Width	660 mm
Length	1520 mm
Max Patient Weight	165 kg
Max Angle Auto	40°
Max Angle Manual	80°
SWL	185 kg

Want to see more?

Scan the QR code to watch our demonstration video on the Staebel Careturner online today!



NDIS



RollStop Coverlet

Another example of our commitment to innovation this cover can be placed on top of a mattress to safeguard the patients position in bed and prevent the patient being able to fall as a result of rolling to the edge of the bed. Ideal for those who get agitated about side rails being raised.

- Smooth material will prevent friction, shear and damage to patient skin
- Waterproof, breathable cover can be used over air mattresses
- Simple installation - with fastening straps to attach the cover to bed frame on each side
- Easy to clean with removable side foam inserts
- Easy patient transfer is achieved due to mid-section openings
- Removable heel protection that helps positioning patients feet



Specifications:

Product Code	3707
Width	920 mm

NDIS



Rollstop Premium Foam Mattress

The Rollstop mattress, as it sounds, is a high quality foam mattress with a concave lying surface to assist in preventing the user from rolling off the side of the bed.

This model consists of a five layer memory foam and has a hospital grade 2-way stretch vinyl cover with welded seams which contributes to optimum user comfort.


For more sizes please contact the Active Medical team on 1800 267 267



Specifications:

Product Code	30093
Width	860 mm
Length	1980 mm

NDIS

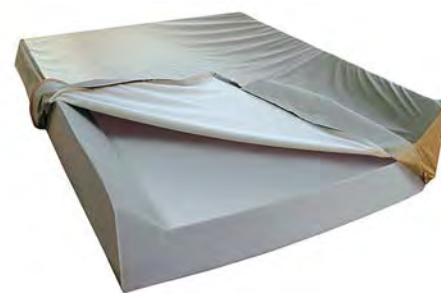

Rollstop Premium Foam Mattress King Single

The Rollstop mattress, as it sounds, is a high quality foam mattress with a concave lying surface to assist in preventing the user from rolling off the side of the bed.

This model consists of a five layer memory foam and has a hospital grade 2-way stretch vinyl cover with welded seams which contributes to optimum user comfort.

• King Single bed size is compatible with King Single size Wissner-bosserhoff Sentida Wide and Carisma Wide Beds



Specifications:

Product Code	30091
Width	1050 mm
Length	1960 mm
Thickness	230 mm



Staebel EvacSmart Evacuation Mat

The ActiveCare EvacSmart patient evacuation device is positioned under the mattress and allows a nurse/carer to remove the patient from their bed and escort them to safety in the event of a fire. The mat has four elastic straps to place over each corner of the mattress so that daily bed making can occur as normal.

The ultimate in strength and durability the EvacSmart patient evacuation device has three adjustable straps that, in the event of an evacuation, are secured around the patient. A lengthy strap at each end of the device allows the nurse/carer to gently lower the patient to the floor without causing them a manual handling risk trying to lift the patient.

The mattress and patient are then lowered onto the floor and can be slid to a safer area quickly and effectively. Once in a safe area the patient can remain on their mattress. Much thought has gone into the design of this device with durability and patient comfort our primary consideration.

- Strap retainers
- Long handles with multi positions
- Reduce the risk of manual handling injuries
- Can be remained under the patient



Specifications:

Product Code	99347
Width	900 mm
Length	2050 mm

Also Available:

Product Code	99347W
Width	1000 mm
Length	2050 mm



Staebel Bed Rail Protector

The newest bedside protector, made from a polyurethane material, coupled with a bioelastic incontinence cover, making it the ultimate answer for longevity and effectiveness.

No zips, clips, or flimsy foam means zero risk of skin tears. Simply Push it down onto the top rail and it is in place.

2x protectors can be used together to cover the complete length of a 4-part rail side guard or full length side rails, sold singularly.

- Easy to use and install
- Skin friendly
- Easy to clean and disinfect



Specifications:

Product Code	3983
Length	800 mm
Height	450 mm
Width	130 mm

NDIS

NDIS

NDIS



Back Pillow

The Back Pillow is a winged-pillow for back relief. The pillow is placed behind the back, and the wing rests under the user to set the position. It can be used for support in side lying, for beds without side guards. Use a back Support Pillow with a flap placed either under the user or the mattress overlay.

Multi Pillow

The Multi Pillow is a versatile pillow that can be placed between the knees or under legs and heels for support and relief.

- Supports the feet, counters heel pressure against surfaces
- Locked position in connection with spasticity
- Supports in lateral position
- Relief for heels
- Supports under the knees

Positioning Pillow

Pillows that provide support in multiple positions. If a person makes too few spontaneous movements, there is an increased risk of pressure ulcers. Regular changes of position reduce the risk of pressure ulcers. Position as required in both reclining and seated position.

Curera® Positioning Pillow can be used as an aid in the facilitation of lung expansion for users with breathing and/or ventilation difficulties. They can also be used to support side lying, Support for neck and arms, Support for users affected by stroke and Positioning for legs and feet.



Specifications:	
Product Code	32301
Width	500 mm
Length	850 mm



Specifications:	
Product Code	32310
Width	600 mm
Length	700 mm

Sizes Available:	
Product Code	32302
Width	450 mm
Length	450 mm



Specifications:	
Product Code	32304
Width	500 mm
Length	250 mm

Sizes Available:		
Product Code	32305	32306
Width	500 mm	500 mm
Length	600 mm	850 mm



Positioning Roll

The positioning roll can give support to the resident in sitting and lying positions. It can be placed in various ways to provide support and relief.

- Supports shoulders, spine, hips, buttocks and heels
- All-round support in seated and reclining positions
- Relieves the knees and supports feet
- Supports in lateral position

Molly Heel Protection

- Easy to use
- Equalizes pressure
- Faster healing
- Improves circulation
- Makes it easier to freely move the legs
- Prevents crossing the legs after hip surgery
- Soft and comfortable
- Wipe with cleaning agent and/or disinfect
- Machine wash max. 95°C
- Easy to clean
- Liquid resistant
- Long durability



Specifications:

Product Code	32307
Width	150 mm
Length	600 mm

Sizes Available:

Product Code	32308	32309
Width	200 mm	200 mm
Length	1250 mm	2250 mm

Specifications:

Product Code	32311
Leg Length	780 - 850 mm
Thickness	75 mm

Sizes Available:

Product Code	32312	32313
Leg Length	860 - 940 mm	1020 - 1030 mm
Thickness	90 mm	110 mm

NDIS

NDIS



Staebel Positioning & Support Wedge

The Support Wedge is used to place a user in side position in order to relieve pressure on the sacrum, which is an extra sensitive area. For the side position, two support wedges are used together.

The support wedge is easy to use and easy to clean. It is delivered with a liquid resistant and vapour permeable cover. The cover can be wiped off with an alcohol-based cleaning agent, as well as machine washed.

- Machine wash max. 95°C, tumble dry.
- Wipe off with cleaning agent, disinfect with alcohol-based cleaning agent.
- 100% HR-foam
- Polyurethane 2-way stretch polyester

Examples of different uses:

- When the user has problems to lie completely on the side
- To prevent pressure ulcers
- When the user has problems with position changes
- For other positioning



Specifications:

Product Code	3689
Width	260 mm
Length	700 mm
Thickness	50 - 150 mm

Action Polymer Pad

Made from polymer gel, this product is effective in preventing the build up of pressure, promoting even pressure distribution. The polymer gel does not displace from areas requiring pressure relief and protects areas where shear forces and pressure points are prevalent. This gel pad can be used in many situations, and for spot protection.

- Easily cleaned, ideal for multi client use.



Specifications:

Product Code	340402
Width	230 mm
Length	280 mm
Thickness	9.5 mm

Sizes Available:

Product Code	340403	340404
Width	280 mm	550 mm
Length	430 mm	500 mm
Thickness	6 mm	13 mm

NDIS

MASS

NDIS



Action Gel Segment Bed Pad

The Action bed overlays are made of Akton polymer which mimics human fatty tissue in consistency, does not migrate like liquid gel and provides equally distributed pressure management. With varying sizes available, these overlays can be used to cover the entire surface or isolated areas of an existing standard mattress. Clients can use one, two or three segments in combination on their areas of need e.g shoulders and hips, feet and hips or full body etc.

- Easily cleaned, no maintenance
- Not restricted by any weight limit
- 3 Segments cover length of single bed

Action Gel Elbow/ Heel Protector

These gel protectors are for use in long term care settings and for assistance in wound management. They allow clients taking medication that affects skin integrity or who are prone to skin tears or abrasions to continue daily activities or to sleep comfortably.

They need only a wipe over with antiseptic solution or mild detergent to be kept fresh and odour free.

- Unique polymer consistency (not liquid gel) so bottoming out cannot occur
- Helps prevent footdrop



Specifications:

Product Code	3853
Width	630 mm
Length	690 mm
Thickness	180 mm
SWL	250 kg

Specifications:

Product Code	3170
Gel Thickness	48 mm

Sizes Available:


Product Code	3171
Gel Thickness	95 mm



Staebel In Bed Turner Base Sheet

The InBed Turner Base Sheet is part of our renowned satin range for in-bed repositioning and patient care. All InBed Turner Base Sheets allow safe and easy repositioning for people in bed enabling them to turn easily from side to side but not slip down. Extra safety is provided by the striped cotton edges on all of our Base Sheets. They can be used on their own or combined with one of our Satin Top Sheets for those who need more in bed repositioning care.

- Reduces risk of tissue and nerve damage
- Helps reduce pressure sores developing
- No training required
- Easy turning and repositioning with little effort
- Great for those having difficulty turning in bed
- Washable
- Full length or midi length satin
- Central satin panel is smooth for ease of movement
- Anti-slip cotton polyester on edges reduces risk of slipping
- Cotton polyester for enhanced comfort and repeat washing

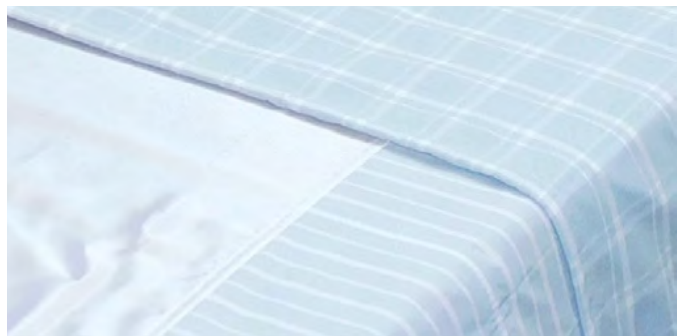

Specifications:

Product Code	210321
Size	Single
Width	900 mm
Length	2000 mm
SWL	300 kg

Also Available:


Product Code	210323	210325
Size	King Single	Double
Width	1050 mm	2000 mm
Length	2000 mm	2300 mm
SWL	300 kg	300 kg



Staebel In Bed Turner 4-Way with Handles

InBed Turner 4-Way with Handles allows safe and easy repositioning for people in bed – enabling them to turn from side to side and move up and down, with or without the aid of a hoist. They are designed for more dependent users who need assistance in bed, often from a single carer, and are designed to be used on top of one of our InBed Turner Base Sheets. They do not work on their own.

- Can remain under the patient at all times.
- Great for those who need multi positioning in bed – up / down and turning.
- Reduces risk of tissue and nerve damage.
- Helps to reduce incidence of pressure sores developing.
- The combination of both silky smooth satin panels working with each other mean easy, quiet and safe movements.
- Cotton and polyester mix for great comfort and repeated washing.
- Washable to 90C.
- Four-way movement in all directions for patient care.
- Allows more dependent people to be moved in comfort and safety.
- Cotton and satin fibres make for comfortable in-bed experience.
- Choice of torso or full length models.
- Available with or without handles for attachment to hoists.
- Main sizes are all 200cm wide, allowing the sheet to be tucked in after repositioning, increasing patient security.

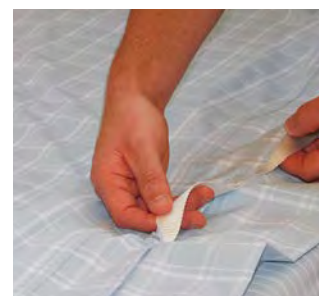

Specifications:	
Product Code	210322
Width	1400 mm
Length	2000 mm
SWL	300 kg

Also Available:	
Product Code	210324
Width	2000 mm
Length	2000 mm
SWL	300 kg

06.

Seating

Regency Care Chair	47
Regency Electric Care Chair	48
Sertain Care Chair Electric Hi-Lo	49
Regency Result3 Electric Recliner	50
Regency Bahama Comfort Chair	50
Ashley Mobile Recliner Luxor Optima	51
Ashley Mobile Recliner - Standard	51
Ashley Luxor Optima Lift & Recline Chair - One Motor	52
Ashley Luxor Optima Lift and Recline Chair - Two Motor	52
Ashley Lift & Recline Chair - One Motor	53
Ashley Lift & Recline Chair - Two Motor	53
Embark Superior Day Chair	54
Utilite Day Chair	55
Embark Deluxe Day Chair	55



Regency Care Chair

The 5050G Regency Care Chair is designed to give the client comfort by combining the original tilt system with an independent multi-positional back rest. Its many features take into account the needs of both carer and user of the chair. The Regency 5050G derives its status from a proven and successful track record across a broad spectrum of the industry.

- **Other sizes available.** Contact the team today to find out more.
- **Drop down armrests and side wings:** Regency armrests/wings fold right away to allow chair to be positioned close to bed for supine transfers. Also making sling insertion much easier.

- **Manoeuvrability:** Regency chairs have two front directional lock castors, two locking rear castors and a horizontal push handle reducing the effort to manoeuvre.
- **Gel cushion:** The guesswork has been removed ensuring the seat has the correct pressure to support the resident with the highest level of comfort.
- **Leg and footrest foldaway feature:** The uniquely designed non removable footplates prevents persons from sliding forward lowering the number of repositions required.
- **Manoeuvrable castor base:** 175mm castors reduces effort to move the chair.

Colours available:



Specifications:

Product Code	5050G
Seat Width	460 mm
Seat Depth	520 mm
Seat Height	600 mm
Backrest Length	700 mm
Armrest Height	175 mm
Legrest Length	480 - 650 mm
Overall Width	660 mm
Overall Depth	1040 mm
Overall Height	1310 mm
SWL	180 kg

NDIS



Regency Electric Care Chair

The R3900 Series Regency Care Chair follows the success of the R2900 Series model and is designed to give the large client unprecedented comfort by combining the tilt system with an independent multi-positional back rest and pressure area management seating.

Its many features take into account the needs of both carer and user of the chair. The Regency R3900 derives its status from a proven and successful track record across a broad spectrum of the industry.

- Reclining backrest – full recline to supine position
- Gas strut assisted Tilting seat
- Pressure management seating – gel cushion in seat, self-air adjusting back, profiled adjustable head cushion

- 2-position armrests which drop down for sling insertion
- Forward seat tilt for ease to stand, legrest angles under seat
- Elevating contoured leg rest complete with fold out, extendible footplate
- Low skin shear 2-way stretch fabric on seating surfaces
- Drop down armrests and swing-away head wings – side transfer friendly
- Covers are replaceable
- Patient lifter friendly castor base
- 150mm Twin wheel castors

Colours available:



Tangerine



Purple



Green



Mink



Lagoon



Mustard Pattern



Koala Pattern



Arcadia Marina



Bouquet Navy



Specifications:

Product Code	5050GE
Seat Width	620 mm
Seat Depth	570 mm
Seat Height	550 mm
Overall Width	820 mm
Overall Depth	1050 mm
Overall Height	1320 mm
SWL	200 kg

NDIS



Sertain Care Chair Electric Hi-Lo

The Sertain® Care Chair Series is there when you need it, to reposition the occupant, ease the transfer from bed to chair, increased access for wound dressing and procedures, manoeuvrable in small spaces.

- Electric hi-lo ; low enough to stand from with ease, but high enough for a comfortable working height for carers
- Electric reclining backrest – full recline for supine position
- Electric tilting seat and elevating legrest
- Handset operated with preset auto chair and auto layflat functions
- Pressure management seating – gel cushion in seat, air in backrestin seat, air in backrest
- Fold out, extendible footplate

- Low skin shear 2-way stretch fabric on seating surfaces
- Removable two-position armrests
- Swingaway three position headwings
- Adjustable head rest cushion
- 150mm castors, 3 x braking, 1 x directional lock
- Battery back up system
- Quality Australian made product
- Available in seat widths 520mm, 620mm and 720mm

Colours available:



Tangerine



Purple



Green



Mink



Lagoon



Mustard Pattern



Koala Pattern



Arcadia Marina



Bouquet Navy



Specifications:

Product Code	5142
Seat Width	490 mm
Seat Depth	520 mm
Seat Height (Tilt/Normal)	540 / 940 mm
Backrest Height	700 mm
Armrest Height	150 mm
Legrest Length	440 mm
Overall Width	740 mm
Overall Depth	1080 mm
Overall Height	1200 - 1520 mm
SWL	160 kg

NDIS

NDIS



Regency Result3 Electric Recliner

Packed full of features and with a 200kg weight capacity, this sleek recliner boasts a modern design appealing to wards, aged care facilities and clinics. The fixed lounge-style armrest has been uniquely designed to fit right up to most dining tables to increase social interaction.

- Electric 'tilt-in-space' seat and leg-rest action
- Electric Euro-style reclining backrest
- Elevating standard leg-rest
- Incorporating pressure management Multicush gel seating and air cushion backrest system
- Fixed lounge-style armrest
- Auto lay-flat function, auto chair function, M3 customer set memory position
- 4 x 100 swivel braking castors
- 5 year warranty
- Australian made
- Available in a huge range of Healthcare grade upholstery colours

Regency Bahama Comfort Chair

Our Complex Care Chair is a popular choice for many healthcare environments. It's available in gas or electric formats, easy to maneuver and simple to clean for infection control. It's been carefully designed to keep your patients or residents social at meal times, fitting in closely to many dining tables.

- 'Tilt in space' seat and leg-rest action
- Euro style reclining backrest with built-in Bahama style head wings
- Elevating standard leg-rest
- Incorporates pressure management Multicush gel seating with air cushion backrest system
- Fixed lounge style armrest
- Includes push handle
- 180kg safe working load
- 5 year limited warranty
- Available in a huge range of Healthcare grade upholstery colours
- Australian made



Specifications:

Product Code	5025GH
Seat Width	550 mm
Seat Depth	520 mm
Seat Height	550 mm
SWL	200 kg



Specifications:

Product Code	5022GH
Seat Width	550 mm
Seat Depth	570 mm
Seat Height	550 mm
SWL	180 kg

Also Available:

Regency Bahama Comfort Chair with Supports	
Product Code	5023GH



Ashley Mobile Recliner Luxor Optima

The Ashley Luxor range of chairs have been specially developed for those with pressure care needs or those needing increased comfort when seated. This is achieved with layered pressure reducing foam and stretch PU upholstered on the seating surfaces. All the favourite features of the standard Ashley are still included – a waterfall style backrest, sculptured armrests, side pockets and padded leg chaise. Cascading roll back with isolated head rest and lumbar roll to support both these areas

- Upholstered leg rest with calf support
- Manual operated recliner action with lever
- Castors and push handle for mobility
- No sag springs in seat and lumbar for long lasting support
- Quality Australian made
- 12 month warranty on fabric, 24 month warranty on frame
- Adjustable head rest cushion
- 150 mm castors, 3 x braking, 1 x directional lock
- Allow +/- 20mm tolerances

Ashley Mobile Recliner - Standard

Economically priced, the Ashley features a waterfall style backrest, sculptured armrests, side pockets and padded leg chaise. Australian made. Cascading roll back with isolated head rest and lumbar roll to support both these areas.

Available in either a fabric or vinyl upholstery in a pattern or plain. Contact our team to discuss! Vinyl available in Austex Studio Encore colour range.

- Allow +/- 20mm tolerances



Specifications:

Product Code	5519MI	5519	5519MX
Size	Mini	Medium	Maxi
Seat Width	480 mm	480 mm	530 mm
Seat Depth	470 mm	500 mm	520 mm
Seat Height	460 mm	480 mm	500 mm
Backrest Height	600 mm	640 mm	680 mm
Overall Width	780 mm	780 mm	830 mm
Overall Height	1020 mm	1050 mm	1150 mm
SWL	120 kg	120 kg	120 kg

Specifications:

Product Code	5521
Size	Medium
Seat Width	480 mm
Seat Depth	500 mm
Seat Height	480 mm
Backrest Height	640 mm
Overall Width	780 mm
Overall Height	1050 mm
SWL	120 kg

NDIS



Ashley Luxor Optima Lift & Recline Chair - One Motor

The Ashley Luxor Optima range of chairs have been specially developed for those with pressure care needs or those needing increased comfort when seated. This is achieved with layered pressure reducing foam and a specialised 2-way stretch upholstery on the seating surfaces. All the favourite features of the standard Ashley are still included – a waterfall style backrest, sculptured armrests, side pockets and padded leg chaise.

- Upholstered leg rest with calf support
- One motor wall saver action – chair can be positioned close to within 100mm from wall and still recline
- Backrest reclines to half lying down position
- Chair raises to assist standing transfer
- No sag springs in seat and lumbar for long lasting support
- Tasmanian Oak hardwood frame for durability
- Battery back up capability
- Small castors enable chair to be wheeled (not patient transport)
- Quality Australian made product
- 2 year electronic warranty, 5 year frame warranty
- Allow +/- 20mm tolerances



Specifications:

Product Code	5522	5522MI	5522MX
Size	Medium	Mini	Maxi
Seat Width	480 mm	480 mm	530 mm
Seat Depth	500 mm	490 mm	520 mm
Seat Height	500 mm	470 mm	530 mm
Backrest Height	640 mm	600 mm	680 mm
Overall Width	780 mm	780 mm	830 mm
Overall Height	1070 mm	1020 mm	1170 mm
SWL	120 kg	120 kg	120 kg

NDIS



Ashley Luxor Optima Lift & Recline Chair - Two Motor

The Ashley Luxor Optima range of chairs have been specially developed for those with pressure care needs or those needing increased comfort when seated. This is achieved with layered pressure reducing foam and stretch PU upholstered on the seating surfaces. All the favourite features of the standard Ashley are still included – a waterfall style backrest, sculptured armrests, side pockets and padded leg chaise. The two motor action enables independent control of backrest and leg rest operation for ultimate flexibility of chair position.

- Upholstered leg rest with calf support
- Backrest reclines to flat position
- Chair raises to assist standing transfer
- No sag springs in seat and lumbar for long lasting support
- Battery back up capability
- Small castors enable chair to be wheeled (not patient transport)
- Quality Australian made product
- Battery back up system
- Allow +/- 20mm tolerances



Specifications:

Product Code	5532	5532MI	5532MX
Size	Medium	Mini	Maxi
Seat Width	480 mm	480 mm	530 mm
Seat Depth	500 mm	490 mm	520 mm
Seat Height	500 mm	470 mm	530 mm
Backrest Height	640 mm	600 mm	680 mm
Overall Width	780 mm	780 mm	830 mm
Overall Height	1070 mm	1020 mm	1170 mm
SWL	120 kg	120 kg	120 kg



Ashley Lift & Recline Chair - One Motor

The Ashley range of chairs have long been the favourite choice with their comfortable seating and range of options. Economically priced and upholstered in standard fabric, the Ashley features a waterfall style backrest, sculptured armrests, side pockets and padded leg chaise. The Ashley chairs can also be ordered in a variety of sizes to suit the individual. This quality powered rise and recline chair is both comfortable and durable, and provides the user assistance with hassle free push button operation.

This model features a single motor allowing the occupant to recline and be raised to the standing position.

Available in either a fabric or vinyl upholstery in a pattern or plain. Contact our team to discuss! Vinyl available in Austex Studio Encore colour range.

- Allow +/- 20mm tolerances



Specifications:			
Product Code	5524	5524MI	5524MX
Size	Medium	Mini	Maxi
Seat Width	480 mm	480 mm	530 mm
Seat Depth	500 mm	490 mm	520 mm
Seat Height	500 mm	470 mm	530 mm
Backrest Height	640 mm	600 mm	680 mm
Overall Width	780 mm	780 mm	830 mm
Overall Height	1070 mm	1020 mm	1170 mm
SWL	120 kg	120 kg	120 kg



Ashley Lift & Recline Chair - Two Motor

The Ashley range of chairs have long been the favourite choice with their comfortable seating and range of options. Economically priced, the Ashley features a waterfall style backrest, sculptured armrests, side pockets and padded leg chaise. Australian made, The Ashley chairs can also be ordered in a variety of sizes to suit the individual. This quality powered rise and recline chair is both comfortable and durable, and provides the user assistance with hassle free push button operation. The two motor action enables independent control of backrest and leg rest operation for ultimate flexibility of chair position.

Available in either a fabric or vinyl upholstery in a pattern or plain. Contact our team to discuss! Vinyl available in Austex Studio Encore colour range.

- Allow +/- 20mm tolerances




Specifications:			
Product Code	5525	5525MI	5525MX
Size	Medium	Mini	Maxi
Seat Width	480 mm	480 mm	530 mm
Seat Depth	500 mm	490 mm	520 mm
Seat Height	500 mm	470 mm	530 mm
Backrest Height	640 mm	600 mm	680 mm
Overall Width	780 mm	780 mm	830 mm
Overall Height	1070 mm	1020 mm	1170 mm
SWL	120 kg	120 kg	120 kg

Ashley Recliner & Embark Chair Colours

Available in either a fabric or vinyl upholstery in a pattern or plain. Vinyl available in Austex Studio Encore colour range.



Magma Wine

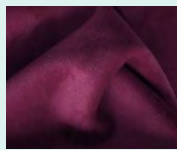


Glamour Scarlet



Raisin




Mazza Charcoal

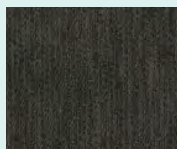


Jazz Walnut



Deluxe Chambray

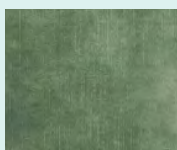

Lovely Jade

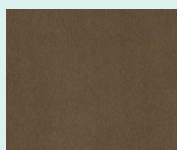

Regis Ginger



Topsuede Denim




Pine




Onyx



Sesame

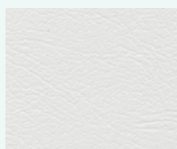


Snow

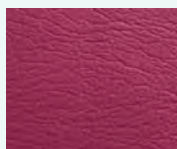

Primrose



Atlantic

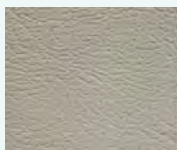


Sand




Grasshopper




Regal

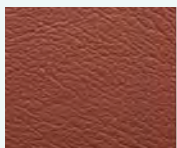


Cameo

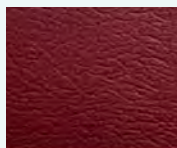

Rhubarb

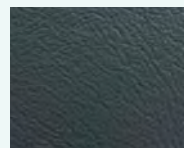

Charcoal

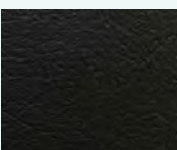

Coffee Bean



Carnival



Smoke

Contact our team to discuss the whole range of fabric or vinyl colours available!

NDIS



Embark Superior Day Chair

The Embark Superior Day Chair was designed and developed by a team of clinicians. It has a strong frame and a high weight capacity of 200 kg. The chair has low resilient foam and 2-way stretch upholstery on seating surfaces. For increased longevity in demanding environments it has 6cm tilt and move wheels and other small improvements.

There have been instances in the past, whereby patients, staff and visitors have sustained cuts and injuries from square tube frames. As a result, the Embark was not only designed for comfort, but to avoid these injuries.

- Low back option also available.
- Comfortable upholstered seat and back
- Padded armrests
- Adjustable seat height and depth
- Powder-coated steel frame
- High backrest incorporating lumbar and lateral support
- Available in different colours! Contact the team for more information



Specifications:

Product Code	5933
Seat Width	500 mm
Seat Depth	450 - 530 mm
Seat Height	450 - 600 mm
Backrest Height	790 mm
Arm Height	230 mm
Overall Width	700 mm
Overall Depth	730 - 800 mm
Overall Height	1200 - 1350 mm
SWL	200 kg

NDIS



Utilite Day Chair

The Utilite is a multi purpose adjustable seat height and angle chair. Ideal for use in an area where adjustable seating is required, for example, lounge, dining room, bedroom, activity room, outdoor area etc. This low back chair provides comfort, stability and practicality at an affordable price.

- Shallow curved low backrest
- Padded armrests
- Comfortable, durable upholstered seat, back and armrests
- Powdercoated steel frame
- Available in different colours! Contact the team for more information



Specifications:

Product Code	5910
Seat Width	460 mm
Seat Depth	460 mm
Seat Height	485 - 640 mm
Backrest Height	455 mm
Arm Height	230 mm
Overall Width	610 mm
Overall Depth	550 mm
Overall Height	900 - 1075 mm
SWL	120 kg

Also Available:

Embark Day Chair Legrest

Product Code	5958
---------------------	------

NDIS



Embark Deluxe Day Chair

An exceptional rehabilitation chair, the Embark Deluxe is available in sturdy quality vinyls for ease of cleaning. The adjustable seat height, seat depth and angle ensures Embark Deluxe can adapt to the widest range of people, making this chair ideal for hospital, home and nursing home. The padded seat and back rest have been designed to enhance the posture, whilst providing maximum comfort.

- High backrest incorporating lumbar and lateral support
- Padded armrests
- Comfortable vinyl upholstered seat and back
- Adjustable seat height and depth
- Powdercoated steel frame
- Standard with tilt and move wheels for ease of manouvrability
- Available in different colours! Contact the team for more information



Specifications:

Product Code	5930
Seat Width	500 mm
Seat Depth	450 - 530 mm
Seat Height	450 - 600 mm
Backrest Height	790 mm
Arm Height	230 mm
Overall Width	700 mm
Overall Depth	730 - 800 mm
Overall Height	1200 - 1350 mm
SWL	200 kg

Also Available:

Embark Day Chair Legrest

Product Code	5958
---------------------	------

WATERLOW PRESSURE ULCER PREVENTION/TREATMENT POLICY
RING SCORES IN TABLE, ADD TOTAL. MORE THAN 1 SCORE/CATEGORY CAN BE USED

BUILD/WEIGHT FOR HEIGHT	◆	SKIN TYPE VISUAL RISK AREAS	◆	SEX AGE	◆	MALNUTRITION SCREENING TOOL (MST) (Nutrition Vol.15, No.6 1999 - Australia)		
AVERAGE BMI = 20-24.9	0	HEALTHY	0	MALE	1	A - HAS PATIENT LOST WEIGHT RECENTLY	B - WEIGHT LOSS SCORE	
ABOVE AVERAGE BMI = 25-29.9	1	TISSUE PAPER	1	FEMALE	2	YES - GO TO B	0.5 - 5kg = 1	
OBESE BMI > 30	2	DRY	1	14 - 49	1	NO - GO TO C	5 - 10kg = 2	
BELOW AVERAGE BMI < 20	3	OEDEMATOUS	1	50 - 64	2	UNSURE - GO TO C AND SCORE 2	10 - 15kg = 3	
BMI=Wt(Kg)/Ht (m) ²		CLAMMY, PYREXIA	1	65 - 74	3		> 15kg = 4	
		DISCOLOURED	2	75 - 80	4		unsure = 2	
		GRADE 1	2	81 +	5	C - PATIENT EATING POORLY OR LACK OF APPETITE	NUTRITION SCORE	
		BROKEN/SPOTS	3			'NO' = 0; 'YES' SCORE = 1	If > 2 refer for nutrition assessment / intervention	
		GRADE 2-4	3					
CONTINENCE	◆	MOBILITY	◆	SPECIAL RISKS				
COMPLETE/ CATHETERISED	0	FULLY	0	TISSUE MALNUTRITION	◆	NEUROLOGICAL DEFICIT		◆
URINE INCONT.	1	RESTLESS/FIDGETY	1	TERMINAL CACHEXIA	8	DIABETES, MS, CVA	4-6	
FAECAL INCONT.	2	APATHETIC	2	MULTIPLE ORGAN FAILURE	8	MOTOR/SENSORY	4-6	
URINARY + FAECAL INCONTINENCE	3	RESTRICTED	3	SINGLE ORGAN FAILURE (RESP, RENAL, CARDIAC,)	5	PARAPLEGIA (MAX OF 6)	4-6	
		BEDBOUND e.g. TRACTION	4	PERIPHERAL VASCULAR DISEASE	5	MAJOR SURGERY or TRAUMA		
		CHAIRBOUND e.g. WHEELCHAIR	5	ANAEMIA (Hb < 8)	2	ORTHOPAEDIC/SPINAL	5	
				SMOKING	1	ON TABLE > 2 HR#	5	
						ON TABLE > 6 HR#	8	
				MEDICATION - CYTOTOXICS, LONG TERM/HIGH DOSE STEROIDS, ANTI-INFLAMMATORY MAX OF 4				

SCORE
10+ AT RISK
15+ HIGH RISK
20+ VERY HIGH RISK

Scores can be discounted after 48 hours provided patient is recovering normally

© J Waterlow 1985 Revised 2005*

Obtainable from the Nook, Stoke Road, Henlade TAUNTON TA3 5LX

* The 2005 revision incorporates the research undertaken by Queensland Health.

www.judy-waterlow.co.uk

REMEMBER TISSUE DAMAGE MAY START PRIOR TO ADMISSION, IN CASUALTY. A SEATED PATIENT IS AT RISK ASSESSMENT (See Over) IF THE PATIENT FALLS INTO ANY OF THE RISK CATEGORIES, THEN PREVENTATIVE NURSING IS REQUIRED A COMBINATION OF GOOD NURSING TECHNIQUES AND PREVENTATIVE AIDS WILL BE NECESSARY ALL ACTIONS MUST BE DOCUMENTED

PREVENTION PRESSURE REDUCING AIDS

- Special**
- Mattress/beds:** 10+ Overlays or specialist foam mattresses.
15+ Alternating pressure overlays, mattresses and bed systems
20+ Bed systems: Fluidised bead, low air loss and alternating pressure mattresses
Note: Preventative aids cover a wide spectrum of specialist features. Efficacy should be judged, if possible, on the basis of independent evidence. No person should sit in a wheelchair without some form of cushioning. If nothing else is available - use the person's own pillow. (Consider infection risk)
10+ 100mm foam cushion
15+ Specialist Gell and/or foam cushion
20+ Specialised cushion, adjustable to individual person.
- Cushions:**
- Bed clothing:** Avoid plastic draw sheets, inco pads and tightly tucked in sheet/sheet covers, especially when using specialist bed and mattress overlay systems
Use duvet - plus vapour permeable membrane.

NURSING CARE

- General** HAND WASHING, frequent changes of position, lying, sitting. Use of pillows
- Pain** Appropriate pain control
- Nutrition** High protein, vitamins and minerals
- Patient Handling** Correct lifting technique - hoists - monkey poles
Transfer devices
- Patient Comfort Aids** Real Sheepskin - bed cradle
- Operating Table** 100mm(4ins) cover plus adequate protection
- Theatre/A&E Trolley**

Skin Care General hygiene, NO rubbing, cover with an appropriate dressing

WOUND GUIDELINES

- Assessment** odour, exudate, measure/photograph position
- WOUND CLASSIFICATION - EPUAP**
- GRADE 1** Discolouration of intact skin not affected by light finger pressure (non-blanching erythema)
This may be difficult to identify in darkly pigmented skin
- GRADE 2** Partial thickness skin loss or damage involving epidermis and/or dermis
The pressure ulcer is superficial and presents clinically as an abrasion, blister or shallow crater
- GRADE 3** Full thickness skin loss involving damage of subcutaneous tissue but not extending to the underlying fascia
The pressure ulcer presents clinically as a deep crater with or without undermining of adjacent tissue
- GRADE 4** Full thickness skin loss with extensive destruction and necrosis extending to underlying tissue.

Dressing Guide Use Local dressings formulary and/or www.worldwidewounds

IF TREATMENT IS REQUIRED, FIRST REMOVE PRESSURE

We are
Here to Help



Let's have a chat!

If you see something you like, you can book in your consultancy with one of our friendly team members. See below information for ways you can get in contact...

1800 267 267 | www.activemedicalsupplies.com.au | hello@activemedicalsupplies.com.au

activeMEDICAL
SUPPLIES > SOLUTIONS > SUPPORT

P. 1800 267 267 | **F.** 07 3899 8825
E. hello@activemedicalsupplies.com.au

www.activemedicalsupplies.com.au