

**active**MEDICAL

# Home Care & Community Portfolio

---

Issue 3







# Contents

<b>01.</b> Bath, Shower & Toilet	<b>6</b>
<b>02.</b> Patient Handling	<b>29</b>
<b>03.</b> Pressure Care	<b>60</b>
<b>04.</b> Wheelchairs	<b>84</b>
<b>05.</b> Seating	<b>94</b>
<b>06.</b> Mobility Solutions	<b>108</b>
<b>07.</b> Bedroom Solutions	<b>118</b>
<b>08.</b> Patient Safety	<b>134</b>

## We help Health Care Providers to Provide Better Patient Outcomes

For over 30 years, Active Medical Supplies has delivered safe and innovative solutions for the care of clients in the Aged Care, Hospital, School and Home Care sectors.

We partner with experts in the field of design and manufacturing to provide the very best resources to ensure that your loved ones are provided with the correct equipment to receive optimum care.

We have a commitment to lower our carbon footprint and protect our environment. We utilise our own delivery vehicles and use sustainable materials for packaging to reduce environmental damage.

Our Home Care team work hard to meet the needs of the patients who require assistive products, but also the carers who support them.

**NDIS Organisation ID# 4050014772**



## Procurement & Distribution

Active Medical has developed a network of manufacturers in Australia and around the world – from Germany and Sweden to the United States of America and Canada. We consult directly with industry professionals and ensure that all our products are ethically sourced and supplied and meet environmental standards.



### \*Disclaimer

Active Medical Supplies does not take responsibility for incorrect application of products. This publication is a guide only to assist in selecting products suitable to your clients. Please use in conjunction with professional advice. Products and information listed in this catalogue are subject to change. Website and digital information supercedes all print. For the most up to date information, please contact us directly.





## Our Commitment to Quality

**All our products at Active Medical Supplies are TGA listed and comply with the latest international standard's where applicable.**

Active Medical Supplies is helping the resellers provide better advice to clients and professionals. This is achieved through product innovation, cost reduction activities and compliance with ISO 9001 requirements.

This Quality Policy seeks to actively encourage stakeholders and members of the company to continually communicate about effective means of improvement for product quality and customer requirements. We will sustain this commitment by the involvement and further training and development of all our employees.

To assist us to do this, Active Medical Supplies undertakes to operate and maintain documented quality management systems which satisfy the appropriate controlling requirements linked to the inspection or certification service being provided.

This quality management system is designed to meet the requirements of AS ISO 9001:2015 Quality Management Systems.

# 01.

## Bath, Shower & Toilet

Shower Chairs	7
Transfer Benches	8
Raised Toilet Seat	10
Toilet Surround	10
Over Toilet Frames	11
Shower Commodes	12
Accessories	20
Change Tables	26
Shower Trolleys	27



DVA

NDIS

MASS



NDIS

MASS



NDIS

MASS



## Staebel Aluminium Shower Chair with Arms

This Shower Chair provides a comfortable seat for those unable to stand whilst showering. The armrests provide security when entering or exiting the shower, and assist the user on and off the seat.

- Strengthened frame for users needing extra room and weight capacity
- Seat height and angle adjusts to suit the individuals requirements
- One piece moulded supportive seat and backrest

## Staebel Extracare Shower Chair with Arms

The Extracare Shower Chair provides a comfortable seat for those unable to stand whilst showering. The armrests provide security when entering or exiting the shower, and assist the user on and off the seat.

- Strengthened frame for users needing extra room and weight capacity
- Seat height and angle adjusts to suit the individuals requirements
- One piece moulded supportive seat and backrest
- Small drainage slots allow water to escape and remove risk of entrapment

## Staebel Shower Stool with Plastic Seat

The aluminium shower stool with arms, provides security when entering or exiting the shower and assist the user on and off the seat.

- Aluminium construction for long life
- Plastic Seat
- Compact size for the smaller shower recess
- Adjustable seat height for the individual's requirements



### Specifications:

<b>Product Code</b>	1315J
<b>Seat Width</b>	440 mm
<b>Seat Depth</b>	430 mm
<b>Seat Height</b>	400 - 550 mm
<b>Bath Rim Clear</b>	460 mm
<b>Backrest Height</b>	350 mm
<b>Overall Width</b>	570 mm
<b>Overall Depth</b>	540 mm
<b>SWL</b>	200 kg



### Specifications:

<b>Product Code</b>	1320J
<b>Seat Width</b>	440 mm
<b>Seat Depth</b>	430 mm
<b>Seat Height</b>	400 - 550 mm
<b>Width Between Arms</b>	500 mm
<b>Backrest Height</b>	350 mm
<b>Overall Width</b>	620 mm
<b>Overall Depth</b>	570 mm
<b>SWL</b>	200 kg



### Specifications:

<b>Product Code</b>	1565J
<b>Seat Width</b>	370 mm
<b>Seat Depth</b>	290 mm
<b>Seat Height</b>	520 - 670 mm
<b>Width Between Arms</b>	460 mm
<b>Armrest Height (From Floor)</b>	640 - 840 mm
<b>Overall Width</b>	550 mm
<b>Overall Depth</b>	370 mm
<b>SWL</b>	175 kg

### Also Available:

<b>Shower Stool with Padded Seat</b>	
<b>Product Code</b>	1567

NDIS



## Staelbel Shower Stool with Arms & Back

This shower stool with arms and backrest, provides security when entering or exiting the shower and assist the user on and off the seat.

- Compact size for the smaller shower recess
- Padded seat and backrest for user comfort
- Adjustable seat height for the individual's requirements



### Specifications:

<b>Product Code</b>	1569
<b>Seat Width</b>	370 mm
<b>Seat Depth</b>	290 mm
<b>Seat Height</b>	520 - 670 mm
<b>Width Between Arms</b>	640 - 840 mm
<b>Armrest Height (From Floor)</b>	640 - 840 mm
<b>Overall Width</b>	570 mm
<b>Overall Depth</b>	450 mm
<b>SWL</b>	125 kg

DVA

NDIS

MASS



## Standard Transfer Bench

Designed for individuals with restricted mobility, the Standard Transfer Bench provides safety and independence while transferring in and out of the bath tub. This durable aluminium transfer bench is adjustable height.

- 3 piece plastic non-slip seat with drain holes
- Hand rail for added support while transferring
- Reversible backrest to suit every situation
- Welded frame



### Specifications:

<b>Product Code</b>	1018J
<b>Seat Width</b>	670 mm
<b>Seat Depth</b>	410 mm
<b>Seat Height</b>	445 - 570 mm
<b>Bath Rim Clear</b>	370 - 495 mm
<b>Seat-Back Height</b>	410 mm
<b>Overall Width</b>	790 mm
<b>Overall Depth</b>	540 mm
<b>SWL</b>	140 kg

### Options Available:

<b>Leg Extension With Suction Cups</b>	
<b>Product Code</b>	1018J - 3
<b>Leg Extension With Stoppers</b>	
<b>Product Code</b>	1018J - 4





NDIS

MASS



NDIS

MASS

## Padded Transfer Bench

Designed for individuals with restricted mobility, the padded Transfer Bench provides safety and independence while transferring in and out of the bath tub. This durable aluminium transfer bench is adjustable in height.

- 3 piece padded vinyl seat
- Hand rail for added support while transferring
- Reversible padded backrest to suit every situation
- Welded frame

## Juvo Heavy Duty Transfer Bench

The Juvo Heavy Duty Transfer bench is both compact and sturdy providing safe and easy access for bathing or showering.

With a weight limit of 200kg this suits a vast range of clients whilst remaining lightweight for easy manoeuvrability. The Heavy Duty bench features height adjustable legs, non slip feet, the option of suction cups and reversible backrest to suit all bath types. With a textured non-slip seat with drainage holes and side rail the transfer bench provides stable seating for our clients.

- Leg Extensions available.
- Suitable for Bariatric Users



### Specifications:

Product Code	1018 J - 1
Seat Width	670 mm
Seat Depth	410 mm
Seat Height	460 - 585 mm
SWL	140 kg



### Specifications:

Product Code	1785
Seat Width	680 mm
Seat Depth	410 mm
Seat Height	470 - 570 mm
Bath Rim Clear	390 - 490 mm
Back Height	360 mm
Overall Width	780 mm
Overall Depth	550 mm
SWL	200 kg

NDIS

DVA

NDIS



## Staebel Raised Toilet Seat - 50mm

- Increases seat height of existing toilet simply and easily.
- Ideal for people with hip or knee injuries or those who find it difficult sitting down or getting up from the toilet.
- Easily installed and removed using front location lip and two rear adjustable brackets.
- Front cutout for easier personal cleaning.

**Order Only\***



## Toilet Surround

This adjustable height frame fits around the toilet and provides assistance in getting on and off the toilet. The toilet surround is manufactured from mild steel, which is zinc treated and epoxy coated making it both strong and durable.

- Ergonomic shaped side frame makes it easy to stand



### Specifications:

Product Code	1737
Thickness	50 mm
SWL	190 kg



### Specifications:

Product Code	1900
Arm Height	625 - 775mm
Width Between Arms	460 mm
Overall Width	540 mm
Overall Depth	490 mm
SWL	125 kg

### Also Available:

#### Toilet Surround

Product Code	1901
SWL	130 kg

#### Toilet Surround Heavy Duty

Product Code	1902
SWL	150 kg

NDIS

MASS



DVA

NDIS



## Stabel Aluminium Over Toilet Frame

This aluminium over toilet frame is adjustable in height and fits over the toilet. This frame has been designed to assist the user who has difficulty sitting down or standing up by raising the toilet seat height.

- Aluminium construction for long life
- No risk of rust
- Removable splash guard included
- Lift up seat that removes for cleaning
- Ideal for those with hip or knee injuries or who require an increased seat height
- Unique stacking design saves space when storing more than one frame
- Multi D leg adjustment ensures stability on uneven surfaces
- Push button height adjustment
- Non slip guides on base of legs



### Specifications:

Product Code	1755J
Width Between Arms	445 mm
Seat Height	460 - 610 mm
Overall Width	550 mm
Armrest Height	700 - 840 mm
Footprint (W x D)	550 x 630 mm
SWL	150 kg

## Folding Aluminium Over Toilet Frame

This Over Toilet frame is ideal for those who need to remove or transport it often.

- Powder coated aluminium frame - light to carry
- Seat is able to be removed for replacement and/or cleaning
- Adjustable seat height to suit the individual's requirements




### Specifications:

Product Code	1760
Width Between Arms	460 mm
Seat Height	480 - 630 mm
Overall Height	680 mm
Armrest Height	680 - 830 mm
Footprint (W x D)	530 x 510 mm
SWL	90 kg

Shower Commodes Innovation by Design

# Leading designer, manufacturer & supplier of Stainless Steel Shower Commodes

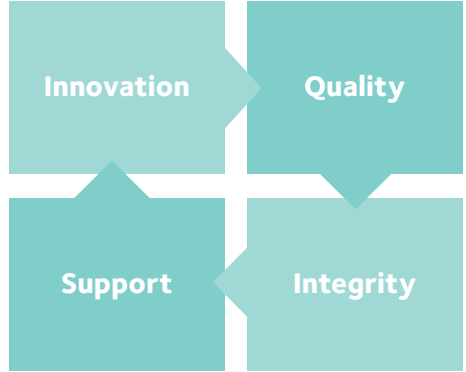

## Build your Standard Shower Commode!

Choose your Commode	Mobile	Static	
Fixed or Tilted Commode	Fixed	Tilted	
Attendant or Self Propelled	Attendant Propelled	Self Propelled	
Footplate Options	Sliding Footplate	Swing-Away Footplate	
Choose your seat width	360 mm	410 mm	460 mm
	510 mm	560 mm	610 mm
	660 mm	710 mm	910 mm (Trix)



Why not enhance your clients situation with high quality accessories! Options include customisable armrests, backrests, footplates & so much more!

**See page 20 for all accessories!**



NDIS

MASS



## Trix 3 In 1 Shower Chair

- Assists patient with showering & toileting
- Clear rear open space enabling positioning over the toilet
- Easy access for hygiene with open armrest design
- Deeper access front trough for easy cleaning
- Integral pan carrier glides. Back saver design, assists with emptying the bowl, can be removed half way, self-supporting while the bowl/pan is removed
- Grade 304 stainless steel heavy wall thickness tubing forming robust frame to withstand extensive wet area use
- Reinforced seat/back frame and armrests – weight bearing
- Easy removable seat, using the Juvo Registered Design "Klipps" System
- Transportation wheels for effortless manoeuvring. Interchangeable so chair can be pulled by back or armrest
- Adjustable seat height for improved client positioning
- Fully adjustable seat depth to assist with client positioning
- Selection of seamless Polyurethane padded seats & pressure care seats (open or closed fronts)
- Armrests project beyond the seat to assist patient in sitting down and standing up
- Splayed armrest design to eliminate pinch point to patient when sitting down and standing up
- Flexible impervious non-woven back straps conforming to client shape for support and comfort. Antibacterial and nonstiction designed straps to avoid skin tear when client is wet
- Curved Backrest design to allow for natural prominence of the gluteal shelf




### Specifications:

<b>Product Code</b>	1783
<b>Seat Width</b>	460 mm
<b>Seat Height (25mm increments)</b>	410 - 560 mm
<b>Seat Depth (25mm increments)</b>	470 - 570 mm
<b>Width between Arms</b>	510 mm
<b>Overall Width</b>	630 mm
<b>Overall Depth</b>	650 - 750 mm
<b>Toilet Max Clearance H x W</b>	390 x 500 mm
<b>SWL</b>	400 kg

### Also Available:

	1779	1780	1781	1782
<b>Seat Width</b>	510 mm	610 mm	710 mm	910 mm
<b>Seat Height (25mm increments)</b>	410 - 560 mm	410 - 560 mm	410 - 560 mm	410 - 560 mm
<b>Seat Depth (25mm increments)</b>	470 - 570 mm	470 - 570 mm	470 - 570 mm	495 - 595 mm
<b>Width between Arms</b>	560 mm	660 mm	760 mm	960 mm
<b>Overall Width</b>	680 mm	780 mm	880 mm	1080 mm
<b>Overall Depth</b>	650 - 750 mm	650 - 750 mm	650 - 750 mm	710 - 810 mm
<b>Toilet Max Clearance H x W</b>	440 x 500 mm	540 x 500 mm	640 x 500 mm	840 x 500 mm
<b>SWL</b>	400 kg	400 kg	400 kg	400 kg

- DVA
- NDIS
- MASS

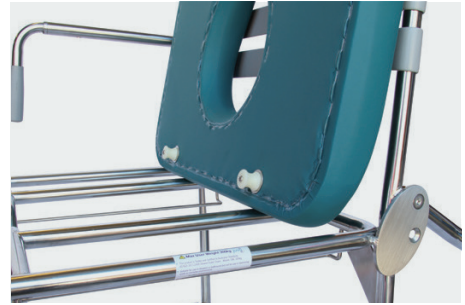

## Juvo Attendant Propelled Shower Commode with Sliding Footplate

The ultimate in strength, quality and stability this Juvo mobile shower chair is ideal for the Petite client! The wheelbase design prevents forward tipping and has proven to be an effective solution in all areas from home to aged care to hospital.

- Attendant propelled
- Quality stainless steel frame for longevity in a demanding environment
- Flip up armrests for ease of side transfers
- Armrests can be flipped so that they are positioned in front of occupant to reduce risk of slipping forward when being wheeled

- Smooth, slatted backrest design to facilitate easy cleaning
- Generous size sliding stainless steel footplate, weight bearing
- To utilise as a commode you must source the slide-in pan rack

**\*Please note that the seat for this product is sold separately**



Seat colours can be varied:



Lemon Twist

Lime

Lipstick

Regal

Emerald

### Specifications:

<b>Product Code</b>	1904J
<b>Seat Width</b>	<b>410 mm</b>
<b>Seat Depth</b>	435 mm
<b>Width Between Arms</b>	440 mm
<b>Overall Width</b>	490 mm
<b>SWL</b>	200 kg

### Sizes Available:

1910J	1914J	1912J	1930J	1936J	1971J
<b>460 mm</b>	<b>510 mm</b>	<b>560 mm</b>	<b>610 mm</b>	<b>660 mm</b>	<b>710 mm</b>
435 mm	435 mm	435 mm	435 mm	435 mm	435 mm
490 mm	540 mm	590 mm	640 mm	690 mm	740 mm
540 mm	590 mm	640 mm	690 mm	740 mm	790 mm
200 kg	300 kg	300 kg	300 kg	300 kg	400 kg



## Juvo Mobile Shower Commode with Swing Away Legrests

The ultimate in strength, quality and stability the Juvo mobile shower commode is built for heavy use having a 200kg SWL! The long wheelbase stops forward tipping and has proven to be an effective solution in all areas from home to hospital.

- Attendant propelled
- Assists transportation of the patient for showering and toileting
- Can be positioned over the toilet
- Durable total stainless steel frame
- Selection of seamless PU padded seats & pressure care seats (open or closed front)
- Removable seat
- Reinforced seat/back frame

- Swing away, removable and height adjustable legrests with flip up footplates
- Patented design, no protrusions
- Stainless steel 125mm castors (4 braking)
- Swing back independent armrests to assist in side transfer and patient access
- Two position height adjustable armrests
- Safety arms provide user security and safety
- Stainless steel pan rack guides
- This product is tested and certified to Australian Standards

**\*Please note that the seat for this product is sold separately**

Seat colours can be varied:



Lemon Twist

Lime

Lipstick

Regal

Emerald



### Specifications:

Product Code	1905J
Seat Width	410 mm
Anterior Seat Height	595 mm
Posterior Seat Height	505 mm
Width Between Arms	440 mm
Overall Width	490 mm
SWL	200 kg

### Sizes Available:

1915J	1916J	1917J	1935J
460 mm	510 mm	560 mm	610 mm
590 mm	590 mm	590 mm	590 mm
565 mm	565 mm	565 mm	565 mm
465 mm	520 mm	570 mm	660 mm
540 mm	610 mm	640 mm	710 mm
200 kg	300 kg	300 kg	300 kg

NDIS

MASS



## Juvo Self-Propelled Shower Commode with Swing Away Legrests

- Assists transportation of the patient for showering and toileting
- Can be positioned over the toilet
- Durable total stainless steel frame
- Selection of seamless PU padded seats & pressure care seats (open or closed fronts)
- Removable seat
- Reinforced seat/back frame
- Swing away, removable & height adjustable legrests with flip-up footplates.
- Patented TEARSAFE design, no protrusions
- Now with braking stainless steel 125mm castors
- Swing back independent armrests to assist in
- Side transfer and patient access
- Two position height adjustable armrests
- Safety arms provide user security & safety
- Stainless steel pan rack guides
- This product is tested and certified to Australian Standards

**\*Please note that the seat for this product is sold separately**



### Seat colours can be varied:



Lemon Twist

Lime

Lipstick

Regal

Emerald

### Specifications:

<b>Product Code</b>	1925J
<b>Seat Width</b>	<b>460 mm</b>
<b>Anterior Seat Height</b>	590 mm
<b>Posterior Seat Height</b>	565 mm
<b>Width Between Arms</b>	465 mm
<b>SWL</b>	200 kg

### Sizes Available:

1926J	1927J	1945J	1939J
<b>510 mm</b>	<b>560 mm</b>	<b>610 mm</b>	<b>660 mm</b>
595 mm	595 mm	595 mm	595 mm
565 mm	565 mm	565 mm	565 mm
540 mm	540 mm	540 mm	540 mm
300 kg	300 kg	300 kg	300 kg





NDIS

MASS

Minimum  
HeightMaximum  
Height

## Juvo 'Tilting in Space' Shower Commode

- Assists transportation of the patient for showering and toileting
- Gas strut tilt mechanism for smooth and easy tilting
- Multi-position 35 deg. tilt
- Can be positioned over the toilet
- Durable total stainless steel frame
- Selection of seamless PU padded seats & pressure care seats (open or closed front)
- Removable seat
- Reinforced seat/back frame
- Swing away, removable and height adjustable legrests with flip up footplates
- Patented TEARSAFE design, no protrusions
- Stainless steel 125mm castors (4 braking)
- Swing back independent armrests to assist in side transfer and patient access
- Adjustable seat height
- Two position height adjustable armrests
- Safety arms provide user security and safety
- Stainless steel pan rack guides

**\*Please note that the seat for this product is sold separately**  
**Electric upgrade available upon request. Call the team to discuss**

Seat colours can be varied:



Lemon Twist

Lime

Lipstick

Regal

Emerald



### Specifications:

<b>Product Code</b>	1950J
<b>Seat Width</b>	<b>410 mm</b>
<b>Overall Width</b>	560 mm
<b>Armrest Height</b>	200 - 240 mm
<b>Back Height from Seat</b>	900 mm
<b>Width Between Arms</b>	440 mm
<b>SWL</b>	200 kg

### Also Available:

1951J	1952J	1956J	1961J
<b>460 mm</b>	<b>510 mm</b>	<b>560 mm</b>	<b>610 mm</b>
570 mm	590 mm	660 mm	590 mm
200 - 240 mm	200 - 240 mm	200 - 240 mm	200 - 240 mm
900 mm	900 mm	800 mm	800 mm
470 mm	520 mm	590 mm	610 mm
200 kg	300 kg	300 kg	300 kg



NDIS

MASS



## Juvo Tilt & Recline Shower Chair

- Manual Reclining Shower Chair which can be positioned over the toilet
- Durable total stainless steel frame
- Selection of seamless PU padded seats & pressure care seats (open or closed front)
- Removable seat
- Reinforced seat/back frame
- Patented TEARSAFE design, no protrusions
- Stainless steel 125mm castors (4 braking)
- Swing away, removable and height adjustable legrests with flip up footplates
- Patented TEARSAFE design, no protrusions
- Swing back independent armrests to assist in side transfer and patient access
- Adjustable seat height
- Two position height adjustable armrests
- Safety arms provide user security and safety
- Stainless steel pan rack guides

**\*Please note that the seat for this product is sold separately**



### Seat colours can be varied:



Lemon Twist



Lime



Lipstick



Regal



Emerald

### Specifications:

Product Code	1976J
Seat Width	410 mm
Overall Width	560 mm
Overall Height	1465 – 1615 mm
Overall Depth	940 – 1600 mm
Width Between Arms	490 mm
Armrest Height	200 – 240 mm
Backrest Height from Seat	800 mm
Seat Tilt Angle	-5° - 38°
SWL	200 kg

### Sizes Available:

1977J	1978J	1980J
460 mm	510 mm	610 mm
560 mm	610 mm	730 mm
1465 – 1615 mm	1465 – 1615 mm	1465 – 1615 mm
940 – 1600 mm	940 – 1600 mm	940 – 1600 mm
490 mm	540 mm	610 mm
200 – 240 mm	200 – 240 mm	200 – 240 mm
800 mm	800 mm	800 mm
-5° - 38°	-5° - 38°	-5° - 38°
200 kg	200 kg	300 kg

DVA

NDIS

MASS

NDIS

MASS



## Juvo Folding Self Propelled Commode with Swingaway Footplates

- Assists transportation of the patient for showering and toileting which also can be positioned over the toilet
- Durable total stainless steel frame
- Selection of seamless PU padded seats & pressure care seats (open or closed fronts)
- Removable seat
- Reinforced seat/back frame
- Swing away, removable & height adjustable legrests with flip-up footplates.
- Patented design, no protrusions
- Now with braking stainless steel 125mm castors
- Swing back independent armrests to assist in side transfer and patient access and provide user security & safety
- Two position height adjustable armrests
- Stainless steel pan rack guides
- Australian Standards Approved
- Patented folding design into standard suitcase size

\*Please note that the seat for this product is sold separately.

## Juvo Paediatric Shower Commode

- Adjustable Backrest design - Adjust the backrest angle from 90 to 120 degrees to improve client comfort and access for the carer.
- Depth adjustable seat - Adjust the seat depth to suit smaller client by up to 75mm.
- Height Adjustable seat up to 90mm
- Self-propelled option available - frame ready to accept self-propelling kit
- Width adjustable



### Specifications:

Product Code	1923J
Seat Width	460 mm
Anterior Seat Height	500 - 590 mm
Posterior Seat Height	475 - 565 mm
Width Between Arms	465 mm
Overall Width	520 mm
SWL	200 kg



### Specifications:

Product Code	1903J
Seat Width	360 mm
Width Between Arms	390 mm
Anterior Seat Height	595 mm
Posterior Seat Height	565 mm
Overall Width	440 mm
SWL	136 kg



### Juvo Angle Adjusted Backrest - 3 Strap

- Shower Commode Accessory
- Self Propelled with 3 Straps
- 340mm high to top strap
- 41cm, 46cm, 51cm, 56cm, 61cm and 66cm widths available

(Contact the team for more information)



### Juvo Angle Adjusted Backrest - 4 Strap

- Shower Commode Accessory
- Self Propelled with 4 Straps
- 415mm high to top strap
- 41cm, 46cm, 51cm, 56cm, 61cm and 66cm widths available

(Contact the team for more information)



### Backrest Extension - 4 Strap

- Suitable for Tilt in Space and Self Propelled with 4 Straps
- 415mm high to top strap
- 41cm, 46cm, 51cm, 56cm, 61cm and 66cm widths available

(Contact the team for more information)

#### Specifications:

**Product Code** 11991 - 46 - 3

#### Specifications:

**Product Code** 11991 - 46 - 4

#### Specifications:

**Product Code** 11945



### Short Backrest Style

- Shower Commode Accessory
- Tilt in Space with 1 Strap
- 41cm, 46cm, 51cm and 61cm widths available

(Contact the team for more information)



### Wheelchair Style Backrest- 3 Strap

- Shower Commode Accessory
- Self Propelled with 3 Straps
- 340mm high to top strap
- 41cm, 46cm, 51cm, 56cm, 61cm and 66cm widths available

(Contact the team for more information)



### Wheelchair Style Backrest - 4 Strap

- Suitable for Tilt in Space and Self Propelled with 4 Straps
- 415mm high to top strap
- 41cm, 46cm, 51cm, 56cm, 61cm and 66cm widths available

(Contact the team for more information)

#### Specifications:

**Product Code** 11944 - 46 - 1

#### Specifications:

**Product Code** 11944-41-3

#### Specifications:

**Product Code** 11944 - 41 - 4



### Backrest Extension - 5 Strap

- Suitable for Tilt in Space and Self Propelled with 5 Straps
- 490mm high to top strap
- 41cm, 46cm, 51cm, 56cm, 61cm and 66cm widths available

(Contact the team for more information)



### Backrest Extension - 6 Strap

- Suitable for Tilt in Space and Self Propelled with 6 Straps
- 570mm high to top strap
- 41cm, 46cm, 51cm, 56cm, 61cm and 66cm widths available

(Contact the team for more information)



### Extended Length Safety Flip Arms

50mm Extended Length Safety Flip Arms for Juvo Commodes.

#### Specifications:

Product Code 11946

#### Specifications:

Product Code 11947

#### Specifications:

Product Code 11989 - EXT



### Transfer Bar / Safety Arm

- Supplied as a pair
- Ergonomic Handgrip
- Spring and button



### Slip-On Lock Down Arm Kit

- Comes in a Pair
- No fixings required



### One Piece Footplate

- One piece footplate for swing-away shower commode models
- Different widths for wider commodes

#### Specifications:

Product Code 11950 - T

#### Specifications:

Product Code 19620

#### Specifications:

Product Code 11963 - 46



### Chest / Pelvic Strap

- Adjustable Chest and Pelvic Strap suitable for Juvo Commodes.



### Arm Trough

The arm trough is of benefit for the shower chair user who has upper limb weakness or paralysis. It attaches to the chair armrest via a sturdy clamp.



### Waterfall Arm Rest

- Juvo waterfall armrest
- Comes in a pair

#### Specifications:

Product Code 11933



### Paediatric Legrest

- For paediatric models
- Left and Right side available (specify when ordering)

#### Specifications:

Product Code	Side
11988 - L	Left
11988 - R	Right

#### Specifications:

Product Code 19613



### Stump Support Leg Extension

- Adjustable stump support for Juvo mobile shower commodes. Installed in replacement for normal design legrest on a needs basis.
- Adjustable angle and length to suit the individual user
- Left and Right side available (specify when ordering)

#### Specifications:

Product Code 11943 - E

#### Specifications:

Product Code 19611 - 2



### Contoured Headrest

- Contoured adjustable headrest to suit Juvo. Comes complete with mounting accessories.
- Fixing Clamp
- Small option available

#### Specifications:

Product Code 19652





## Lateral Support

- Quick release thoracic support
- Swing away
- Left and Right side available (specify when ordering)

### Specifications:

Product Code	Side
11999 - L	Left
11999 - R	Right



## Calf Support Strap

Calf Support now has a colour matched, clip-on calf support for a complete, protective, hygienic solutions for adults and children using shower chairs.

### Specifications:

Product Code	18B216 - 22

NDIS



## Lateral Support - Cover

Cover Only for Lateral Support

### Specifications:

Product Code	11999 - C

NDIS



## Dropdown Pommel

- Height and depth adjustable

### Specifications:

Product Code	11985



## Juvo Commode Bowl & Lid

- Bowl with Lid and Handle
- Can be used with all Shower Commode types

### Specifications:

Product Code	11936



## Pan/Bowl Carrier

Pan and Bowl Carrier that can be used for different commodes. Different sizing available.

### Specifications:

Product Code	Size
19603	410 mm
19607	560 mm

NDIS



## Commode Bowl with Handle

It is suitable for use with the current 7952PA bedside commode, and also with the current style of overtoilet frames and aids. Includes handle.

- Autoclavable
- Top diameter 290 mm
- Base diameter 210 mm
- Height 120 mm

DVA

NDIS



## Standard Bed Pan

These bedpans are made from polypropylene and are autoclavable.

The standard bedpan is heavy duty and ideal for institutional use.

DVA

NDIS



## Male Urinal

A clear polypropylene bottle with integrated handle. Calibrated for output recording, bottle is autoclavable.

Capacity of 1 litre - 100 ml graduation.

### Specifications:

Product Code 1989

DVA

NDIS



## Female Urinal

Clear polypropylene bottle with integral handle, anatomically shaped at opening. Calibrated for output recording and autoclavable.

Capacity of 1 litre - 100 ml graduation.

### Specifications:

Product Code 1990

DVA

NDIS



## Spil - Pruf Urinal

These urinals will not spill or leak in any position. Spil-Pruf urinals help reduce odour problems while increasing confidence contents are not visible.

Autoclavable. With the addition of the 1986 adaptor it is also suitable for women.

### Specifications:

Product Code 1980

NDIS



## Spil - Pruf Urinal Adaptor

This adaptor allows the Spil-Pruf Urinal to also be suitable for women. These urinals will not spill or leak in any position. Spil-Pruf urinals help reduce odour problems while increasing confidence contents are not visible.

### Specifications:

Product Code 1983

### Specifications:

Product Code 1985

### Specifications:

Product Code 1986



## Elevated Legrest

Adjustable elevated legrest for Juvo mobile shower commodes. Installed in replacement for normal design legrest on a needs basis.

- Adjustable angle and length to suit the individual user
- Left and Right side available (specify when ordering)



## Stump Support

Adjustable stump support for Juvo mobile shower commodes. Installed in replacement for normal design legrest on a needs basis.

- Adjustable angle and length to suit the individual user
- Left and Right side available (specify when ordering)



## Juvo Angle Adjusted Footrest

- Juvo Commode Part - Angle adjustable footrest
- Left and Right side available (specify when ordering)

### Specifications:

Product Code	Side
11940 - L	Left
11940 - R	Right

### Specifications:

Product Code	Side
11943 - L	Left
11943 - R	Right

### Specifications:

Product Code	Side
11931 - L	Left
11931 - R	Right

NDIS

NDIS



## Change Table - Electric Hi-Lo

This Change Table includes the following features making the ideal option for multiple environments. With a 450–900mm height adjustment range and a 250kg weight limit, it suits a larger range of clients.

The vinyl padded top is easy to clean and the base has been designed for easy hoist access. Fold down side rails are an optional extra.

- Vinyl upholstered top 700mm wide x 1880mm long
- Electric hi-lo (450mm – 900mm)
- 125mm braking castors
- Designed for hoist access
- Handset operated
- Quality Australian made product



\* Available with side rails

### Specifications:

Product Code	7805
Width	700 mm
Length	1880 mm
Min Height	450 mm
Max Height	900 mm
SWL	250 kg

### Also Available:

Staebel Wall Mounted Fold up Couch	
Product Code	780910

## Change Table - Electric Hi-Lo

Ideal for children, this compact change table features electric hi-lo, padded top, foot control operation. Optional extras include vertical side rails, lap strap.

- Vinyl upholstered top 700mm wide x 1400mm long
- Electric hi-lo (630mm – 930mm)
- 125mm braking castors
- Designed for hoist access
- Foot control operation
- Quality Australian made product



### Specifications:

Product Code	7830
Width	700 mm
Length	1400 mm
Min Height	480 mm
Max Height	880 mm
SWL	90 kg



## Staebel Shower Trolley

- Shower trolley with one-section top
- Central braking system
- Easy steering system (for ease of maneuvering trolley's front left castor can be locked in direction of steering)
- Side rails – foldable, height: 185 mm with safety release mechanism (Material – powder coated stainless steel)
- Metal base cover
- Water-resistant shower mattress (height: 270 mm)
- Mattress integrated drain with 645 – 2000 mm long flexible hose
- Water-resistant head pillow – semi – cylinder

## Staebel Bariatric Shower Trolley

- Hydraulic lowering/raising of the stretcher top via a console lift
- Lowering/raising range: 540 mm to 970 mm
- Castors with a diameter of 150 mm
- Central braking easy steering system (for ease of maneuvering the trolley's front left castor can be locked in the direction of steering directional wheel)
- Side rails – foldable, height: 185mm with safety release mechanism (Material – powder coated stainless steel)
- Metal base cover
- Water resistant shower mattress (height: 270mm)
- Mattress integrated drain with a 645 up to 2000mm long flexible hose
- Water resistant head pillow – semi – cylinder



### Specifications:

Product Code	1209ST
Shower Trolley Top Width	650 mm
Shower Trolley Top Length	1900 mm
Shower Trolley Top Height	540 - 970 mm
Height with Castors	600 - 1030 mm
Trolley Top Inclination	0° - 5°
Length	2070 mm
Width	1070 mm
SWL	170 kg

### Specifications:

Product Code	1216ST
Shower Trolley Top Width	900 mm
Shower Trolley Top Length	2000 mm
Shower Trolley Top Height	540 - 970 mm
Height with Castors	600 - 1030 mm
Trolley Top Inclination	0° - 5°
Length	2170 mm
Width	1070 mm
SWL	250 kg



## Staebel Junior Shower Trolley

- Shower trolley with one-section trolley top
- Trolley top length: 150 cm, trolley top width: 60 cm
- Maximum dimensions: 167 cm (length) and 76 cm (width)
- Adjustable 5° trolley top inclination
- Hydraulic lowering/raising of trolley top via a console from 54 cm to 97 cm
- Castors with diameter of 150 mm
- Central braking system
- Easy steering system (for ease of maneuvering trolley's front left castor can be locked in direction of steering)
- Side rails – foldable, height: 18,5 cm with safety release mechanism (material – powder coated stainless steel)
- Metal base cover
- Water-resistant shower mattress (height: 27 cm)
- Mattress integrated drain from 645 – 2000 mm long flexible hose
- Water-resistant head pillow – semi-cylinder

## Staebel Shower Trolley with Backrest

- Shower trolley with two-section top
- Maximum dimensions: 207 cm (length) and 81 cm (width)
- Trolley top length: 190 cm, trolley top width: 65 cm
- Adjustable 5° trolley top inclination
- Electrically lowering/raising of trolley top via console from 54 cm to 97 cm
- Electrically adjustable head section inclination 45° (head section length: 70 cm)
- Castors with diameter of 150 mm
- Central braking system
- Easy steering system (for ease of maneuvering trolley's front left castor can be locked in direction of steering)
- Side rails – foldable, height: 18,5 cm with safety release mechanism (material – powder coated stainless steel)
- Metal base cover
- Water-resistant shower mattress (height: 27 cm)
- Mattress integrated drain with 645 – 2000 mm long flexible hose
- Water-resistant head pillow



### Specifications:

Product Code	1214ST
Length	1670 mm
Width	760 mm
Castors	150 mm
SWL	170 kg



### Specifications:

Product Code	1211ST
Length	2070 mm
Width	810 mm
Castors	150 mm
SWL	170 kg



# 02.




## Patient Handling

Patient Hoists	30
Hoist Accessories	35
Ceiling Hoists	37
Ceiling Hoist Accessories	41
Slings	42
Transferring Products	45
Positioning Products	50
Supporting Products	56






## What is a Hoist?

Hoists are designed to assist with the manual handling of residents & patients. They're here to help reduce the risk of injury to carers, nurses and doctors. They are also portable and easy to move which allows for the safe movement of patients from one place to another. They are commonly used for shorter transfers such as bed to chair.

## Pick your Hoist!

 Choose Model	Staebel Forte 320	Staebel Salsa 200	
	Staebel Alto 200	Staebel Neos 150	
	Staebel Sonata 150		
 Choose attachment style	Four Point Yoke	Standard Yoke	Tri Hook Yoke
	Powered Pivot Frame	Standard Yoke with Weigh Scale	Bariatric 4 Point Yoke
	Pivot Frame with Weigh Scale	Powered Pivot Frame with Weigh Scales	Bariatric 4 Point Yoke with Weigh Scales
 Accessories	Plastic Leg Protector	Sling Bag	
	Staebel Salsa Handle Kit	Quick Release Adaptor	
	Low Leg Height	Weighing Device	
	Battery Replacement	Handset Replacement	

## Pick your Sling!

 Choose your material		Mesh
		Fabric
		Parasilk
 Choose Sling style	Depending on the model of hoist, choose the appropriate sling style. <b>Contact the team for more information on choosing the right sling!</b>	



## Could not find the information you were after?

Contact our team for a consultation for more hoist options!

NDIS



Featuring Linak Jumbo Care advanced diagnostic control system.



## Staebel NEOS 150 Folding Hoist

The NEOS 150 Folding is a compact ergonomic designed medium sized heavy duty mobile patient lifting hoist with a safe working load of 150kg. This lifter is ideally suited as a compact and lightweight hoist in a home care environment. The NEOS 150 Folding whilst being compact provides superior patient space. The NEOS 150 Folding in fact offers more patient space than most other larger size lifters. The swivel feature of the detachable yoke allows a patient to be rotated through 360° at full lift height, thus making the carer's job easier and safer.

- Folding Aluminium Lifter with a Safe Working Load (SWL) of 150kg
- Quick release folding structure enables lift to be conveniently folded for storage and transfer
- Electric raise and lower powered by Linak
- Foot pedal control for leg adjustment
- Equipped with 75mm front castors and 100mm locking rear castors
- Quick release yoke



### Specifications:

<b>Product Code</b>	2878	<b>Lift Range - Yoke</b>	410 - 1650 mm
<b>Base Width Min</b>	615 mm	<b>Castor Size</b>	75 / 100 mm
<b>Base Width Max</b>	965 mm	<b>Under Bed Clearance</b>	135 mm
<b>Folded Overall Length</b>	1230 mm	<b>Lifter Weight</b>	35 kg
<b>Folded Overall Width</b>	400 mm	<b>SWL</b>	150 kg

DVA

NDIS

MASS



## Staebel Forte 320 Hoist

A hoist designed for the diverse resident portfolio and to enhance future proofing for your facility with mitigating manual handling risk and reducing the risk of skin integrity complications.

Featuring a class-leading lift height with extra low base for better access and maximised patient lift space. The large open cradle design allows the patient to be comfortably transferred without feeling pinched or squeezed. The concealed electric leg spread actuator promotes better hygiene practices and more space to walk behind hoist.

A long lasting lithium battery pack is a standard addition to the Forte Lifter, enabling more lifts from a single charge.

- Smooth mobility and enhanced manoeuvrability.
- Utilises the SPACELINK™ concept, a breakthrough for optimised patient space, that enhances floor lift capability and increases maximum lifting height.
- Emergency stop button & easy to use manual lowering on actuator
- Generous patient space from the Motor to the connection point
- Easy to service and spare parts readily available
- Additional hand control buttons on control box. Can be used when the handheld handset can't be used
- Other attachment models are available



### Specifications:

<b>Product Code</b>	2898	<b>Overall Height when boom is raised</b>	2243 mm
<b>Base Length</b>	1420 mm	<b>Maximum Hoist Reach</b>	810 mm
<b>External Base Width Closed</b>	715 mm	<b>Lifter Weight</b>	63 kg
<b>External Base Width Open</b>	1200 mm	<b>Minimum Storage Height</b>	1445 mm
<b>Under Bed Clearance</b>	111 mm	<b>SWL</b>	320 kg



## Staebel Sonata 150 Lifter

A lightweight, compact mobile patient lifting hoist with innovative and ergonomic design the Allegro Sonata 150 is an easy to use everyday lifter. With a safe working load of 150kg the Sonata offers an outstanding lift capacity in a small patient lifter. It utilises the SPACELINK™ concept, a breakthrough for optimised patient space, that enhances floor lift capability and increases maximum lifting height. With maintenance free bushes in all the major moving parts, this hoist is an Australian Made example of quality and strength that is hard to find in other patient lifters.

The Allegro range of patient lifting equipment, including lifters and accessories, has been designed and tested to comply with: AS ISO 10535-2002 Hoists for the transfer of disabled persons – Requirements and test methods.

- Manual leg adjustment is a standard feature
- Supplied complete with battery, in-built battery charger and instruction manual
- The yoke is the standard attachment for general purpose slings. A unique sling retainer has been designed into the hook of the yoke to prevent the sling loop from coming detached during lifting.
- Other lifter and attachment models are available



### Specifications:

Product Code	2880	Castor Size	100 mm
Base Width Min	605 mm	Under Bed Clearance	126 mm
Base Width Max	1010 mm	Lifter Weight	39 kg
Lift Range - Yoke	414 - 1730 mm	SWL	150 kg

DVA

NDIS

MASS



## Staebel Alto MKII 200 Lifter

Alto Lift 200: A medium sized heavy duty mobile patient lifting hoist with a safe working load of 200kg and a compact ergonomic design. The Alto Lift 200, whilst being compact, provides superior patient lift space through the use of SPACE LINK™ technology.

The Alto Lift 200 is the lifter of choice for institutional use or home care where a high capacity lift is required. Designed and tested to comply with: AS ISO 10535-2002 Hoists for transfer of disabled persons - Requirements and test methods.

- Electric leg adjustment is a standard feature
- Battery life indicator on lifter as well as handset
- Additional up and down buttons on control box
- Emergency stop button
- Supplied complete with battery, battery charger and instruction manual



### Specifications:

<b>Product Code</b>	2891	<b>Maximum Height Yoke</b>	1920 mm
<b>Base Length</b>	1150 mm	<b>Average Lift Range</b>	1620 mm
<b>External Base Width Closed</b>	700 mm	<b>Minimum Storage Height</b>	1350 mm
<b>External Base Width Open</b>	1130 mm	<b>Under Bed Height</b>	55 mm
<b>Internal Base Width Closed</b>	600 mm	<b>Hoist Reach</b>	610 - 1060 mm
<b>Internal Base Width Open</b>	1030 mm	<b>Hoist Reach at Maximum Height</b>	515 mm
<b>Base Width Adjustment</b>	Electric	<b>Mass of Heaviest Part (Base Assembly)</b>	24 kg
<b>Underbed Clearance</b>	130 mm	<b>Lifter Weight</b>	46 kg
<b>Minimum Height Yoke</b>	300 mm	<b>SWL</b>	200 kg



## Staebel Four Point Yoke

The Allegro four point yoke can also be used with the general purpose slings and is particularly useful in the care and comfort of larger patients. A unique sling retainer has been designed into the hook of the yoke to prevent the sling loop from coming detached during lifting.



## Staebel Standard Yoke

The yoke is the standard attachment for general purpose slings.

A unique sling retainer has been designed into the hook of the yoke to prevent the sling loop from coming detached during lifting.



## Staebel Pivot Frame

The pivot frame is used with the pivot sling. The frame spreads the sling attachment points to ensure the patients hips are not flexed excessively, avoiding discomfort some patients experience with standard slings. The ability to position a person in a sitting or lying position, assisting in providing correct seating and positioning.

### Specifications:

Product Code 2886

### Specifications:

Product Code 2882

### Specifications:

Product Code 2884



## Staebel Pivot Frame with Weigh Scale

Complete as one attachment, this allows for standard lifting without losing any critical lift height with separate attachment. Weighing scale Safe Weight Limit 300kg. Suitable for allegro Sonata & Alto Hoists.



## Staebel Standard Yoke with Weigh Scale

Complete as one attachment, this allows for standard lifting without losing any critical lift height with separate attachment. Weighing scale Safe Weight Limit 300kg. Suitable for allegro Sonata & Alto Hoists.



## Staebel Bariatric Four Point Yoke

Complete as one attachment, this allows for standard lifting without losing any critical lift height with separate attachment. Allegro Patient Lifter Option - Bariatric 4-point Yoke spreader bar to suit Forte 320

### Specifications:

Product Code 2895

### Specifications:

Product Code 2894

### Specifications:

Product Code 2873



## Staebel Tri Hook Yoke

These three heavy duty hoists are available with a choice of 6 different sling attachment options. Achieve quality patient transfers by selecting the fitting which best suits your client or healthcare facility.

- Can be used with general purpose (loop attachment) slings



## Staebel Bariatric 4 Point Yoke with Weigh Scales

Allegra Patient Lifter Option - Bariatric 4-point Yoke spreader bar with integrated weighing device and cover to suit Forte 320.



## Staebel Powered Pivot Frame

A powered pivoting frame makes the hoist the perfect choice for lifting larger patients. With the push of a button repositioning and tilting the patient is safe and strain free. A medium sized heavy duty mobile patient lifting hoist with a safe working load of 200kg.

### Specifications:

Product Code 2883



## Staebel Powered Pivot Frame with Scales

A powered pivoting frame makes this hoist the perfect choice for lifting larger patients. With the push of a button repositioning and tilting the patient is safe and strain free. A medium sized heavy duty mobile patient lifting hoist with a safe working load of 200 kg.

### Specifications:

Product Code 2897

### Specifications:

Product Code 2899



## Staebel Lifting Frame with Gliders

The Lifting Frame is an essential piece of equipment for the safe lifting and handling of patients requiring spinal or orthopaedic immobility.

The specialised frame made from a durable metal comes complete with 12 removable nylon gliders and can be dismantled for ease of storage. Strap kit sold separately. \*Pictured being utilised with 4-point yoke. SWL 200kg.

### Specifications:

Product Code 227201

### Specifications:

Product Code 2896



## Staebel Salsa Handle Kit

- Salsa Handle Kit to aid in increased patient security whilst lift is in motion
- Suitable for the Salsa Standing Lifter
- Optional Extra

### Specifications:


Product Code 27THK01






# What is a ceiling hoist?



Ceiling Hoists are designed to assist with the manual handling of residents & patients. They're here to help reduce the risk of injury to carers, nurses and doctors. A ceiling hoist is an overhead hoisting system which can either be hung from the ceiling itself or hung from wall to wall. They are an effective alternative to a floorline hoist and are widely used in home care environments.

## Build your Ceiling Hoist System

 <p><b>Choose Transfer Type</b></p>	Bed to Chair	Bed to Bathroom
	Bed to complete house	Full room access to bathroom
	Bariatric System	Access to other rooms
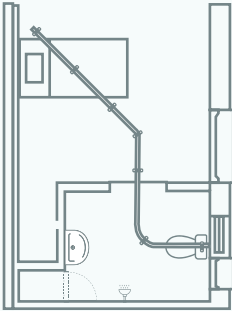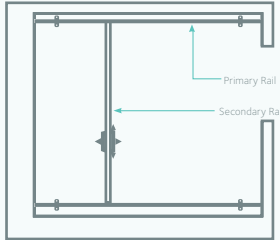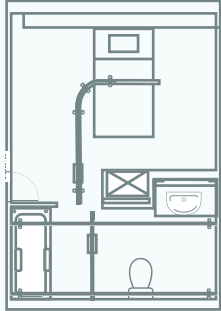
## Choose Tracking System

<p><b>Basic Rail Systems</b></p> 	<p><b>Full Room Coverage System - XY Traverse</b></p> 	<p><b>Fully Integrated Rail System</b></p> 
--	--	---

 <p><b>Weight Limit (SWL)</b></p>	200 kg	270 kg
	350 kg	450 kg +
 <p><b>Choose Motor Style</b></p>	Portable Hoist	Fixed Powered Traverse
	Fixed Manual Traverse	Powered XY Motor
 <p><b>Rail Options</b></p>	Ceiling Mount	Wall Mount

Continue Next Page

# Ceiling Hoist Case Studies

XY system (H System) with Fixed Motor.



XY transitio system nshowing door cutout.



Curved single track transitioned from XY.



Straight track with power traverse motor.



XY system to transition gate showing door cutout.



Straight track with portable motor.



Single Gantry System with Fixed Motor.



Curved single track and straight track with in rail trolley.



**Could not find the information you were after?**

Contact our team for a consultation for more ceiling hoist options!



## SugarGlider Portable Motor

The Monarch Portable Ceiling Hoist is one of the lightest on the market weighing only 5kg, whilst lifting the industry standard of 200kg. Ideal for occasional or multiple room use, the product is robust and charges with an AC Charger.

- Weight: 5KG
- Lifts per Charge: 30 to 40 cycles
- Track Options: H system, Turntable, Gate, Curves, Room-to-Room, Freestanding Gantry: 2m, 3m, 4.5m, 6m, 8m
- Track Compatibility: Compatible with all Monarch tracks



### Specifications:

Product Code	272358000
Length	550 mm
Width	140 mm
Lifts per Charge	30 - 40 Cycles
Total Weight	4.8 kg
SWL	200 kg

## SugarGlider In track Motor

The Monarch Fixed Ceiling Hoist is a complementary product to other products for individuals with mobility challenges. It is designed for use by caregivers and health professionals to lift and transfer a person in homecare, assisted-living or hospital environment.


- Unique Spreader Bar with Secure Sling Mechanism
- Fast - lifts up to 30% more quickly than other models
- Low lift clearance - slim profile means more room to manoeuvre patient
- Li-ion battery - recharges twice as fast as average competitor products
- Ultra lightweight - Weighing 8kg 35% lighter than competitors
- 200kg safe working load
- Ideal for Nursing Homes, Hospitals or Domestic Use



### Specifications:

Product Code	272358001
Length	640 mm
Width	300 mm
Lifts per Charge	30 - 40 Cycles
Total Weight	7.58 kg
SWL	200 kg

NDIS

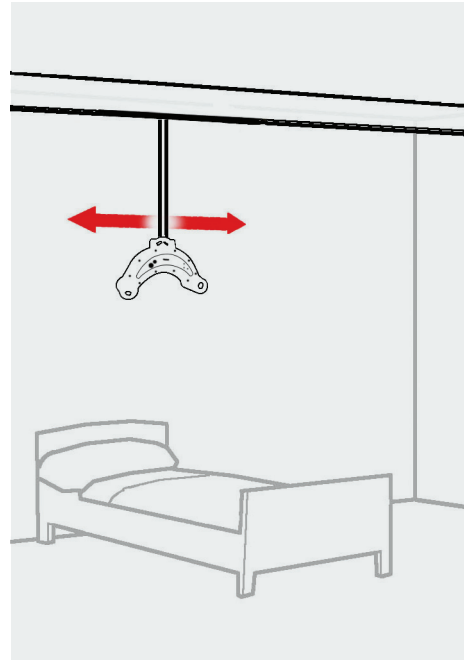

## SugarGlider Nomad Portable Ceiling Hoist

The Nomad is a flexible and user friendly ceiling hoist. It is intuitive and easy to operate. The integrated handles make it truly portable. The curved shape is designed to afford user's head ample space while the integrated 4-point yoke gives ultimate safety, comfort and dignity during the lift and transfer. The user will feel safe and stable without being squeezed.

The combination with a compatible sling ensures that the user will achieve a correct and comfortable sitting position. If rail installations through door openings are not possible, the Nomad can climb from room to room, i.e. from one rail system to another. The Nomad has an emergency lowering and an emergency stop, as well as a soft stop-start function to prevent abrupt movements.

- Lift speed 35 mm/s with soft start/stop function
- High battery performance (70+ lifts per charge)
- Easy overnight charge system
- Built in electronic lift counter
- Splash proof casing
- Lifting Height 2m



### Specifications:

Product Code	7324
Length	520 mm
Width	198 mm
Height	330 mm
Lift Range (Length of Lift Strap)	1980 mm
Lifter Weight	6.8 kg
SWL	230 kg

### Also Available:

Product Code	7325
Length	520 mm
Width	198 mm
Height	330 mm
Lift Range (Length of Lift Strap)	1980 mm
Lifter Weight	6.8 kg
SWL	255 kg



## SugarGlider Transport Trolley

The SugarGlider Transport Trolley makes it a simple task to transport the portable ceiling hoist unit, reach stick and slings from one location to another.

- Two large shelves
- Features a unique column to store reach stick
- Four high quality swivel castors for ease of manoeuvrability
- Robust, heavy duty design favoured by many aged care facilities

## Flexi-Link Arm Reacher

Efficient, safe and simple to operate. Our patented Flexi-Link connection system is today's most advanced device for connecting the portable lift to the trolley. Its removable wand allows carers to move around freely as there are no handles in place and travel of the hoist is not impeded by lights, fans or doors.

### Specifications:


Product Code	7426
Width	520 mm
Depth	198 mm
Height	330 mm

### Specifications:



Product Code	7422FL
Length	600 mm

### Also Available:

Product Code	7423FL
Length	900 mm

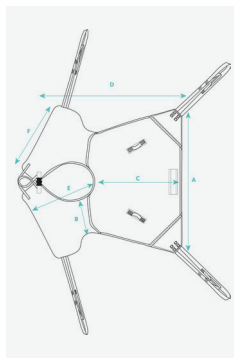

NDIS

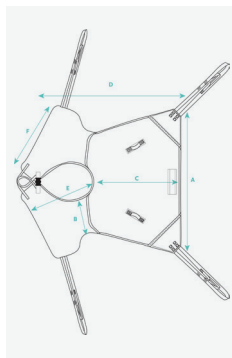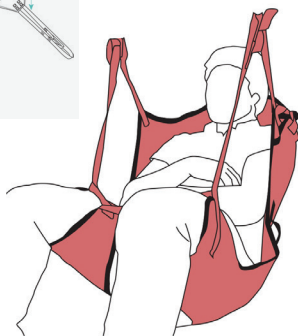

NDIS

MASS



## Staebel Universal Deluxe Pivot Sling

The Universal Deluxe Pivot Sling is designed for use with the pivot/tilt frame and for clients with poor muscle tone or with an amputation. It attaches with the keyhole tags engineered to the highest standards in strength and size tolerances and features auto clip on action (no safety straps required). It also features sling centre marking and positioning handles to assist correct positioning and woven label that will not wash out or fade.

Available in sizes from x-small through to x-large and in mesh or parasilk styles.

## Staebel Universal Deluxe Sling with Head Support

This sling is good for all standard transfers and is exceptionally comfortable for the patient. The leg sections are shaped to provide good support to the hip and thigh areas. It is designed to hold the leg more securely than standard sling leg sections, as the leg is held from both sides instead of just the inside of the leg. The leg sections also spread the load over a larger area on the patient's legs so tend to reduce the pressure on the underside of the legs.



### Mesh

<b>XS</b>	227430P
<b>S</b>	227431P
<b>M</b>	227432P
<b>L</b>	227433P
<b>XL</b>	227434P

### Parasilk

<b>XS</b>	227440P
<b>S</b>	227441P
<b>M</b>	227442P
<b>L</b>	227443P
<b>XL</b>	227444P

	XS	S	M	L	XL
<b>A</b>	715	800	900	1000	1105
<b>B</b>	160	185	210	235	260
<b>C</b>	685	815	945	1075	1205
<b>D</b>	1000	1175	1355	1535	1710
<b>E</b>	295	340	385	430	475
<b>F</b>	375	440	500	560	625

### Mesh

<b>XS</b>	227430
<b>S</b>	227431
<b>M</b>	227432
<b>L</b>	227433
<b>XL</b>	227434
<b>XXL</b>	227450


### Fabric

<b>XS</b>	227435
<b>S</b>	227436
<b>M</b>	227437
<b>L</b>	227438
<b>XL</b>	227439
<b>XXL</b>	227451

### Parasilk



<b>XS</b>	227440
<b>S</b>	227441
<b>M</b>	227442
<b>L</b>	227443
<b>XL</b>	227444
<b>XXL</b>	227212

	XS	S	M	L	XL	XXL
<b>A</b>	715	800	900	1000	1105	1205
<b>B</b>	160	185	210	235	260	285
<b>C</b>	685	815	945	1075	1205	1270
<b>D</b>	1000	1175	1355	1535	1710	1820
<b>E</b>	295	340	385	430	475	520
<b>F</b>	375	440	500	560	625	690

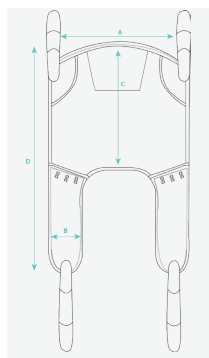

NDIS

MASS

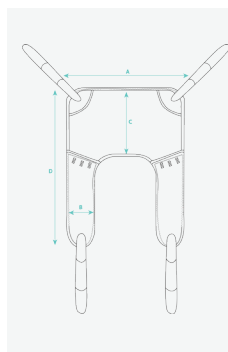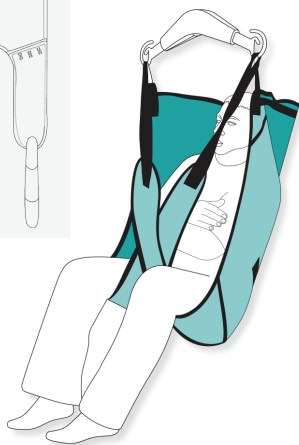

NDIS

MASS



## Stabel General Purpose with Head Support Sling

The General Purpose Head Support sling, as the name suggests, is a full body support style sling and is recommended for general lifting and in situations of lifting from the floor, and for the less secure patient. The loop style attachments have three lengths at each connection point and are used on a patient lifter with a yoke style attachment. Available in sizes from x-small through to XX-large, providing a popular choice for individuals and care facilities. Slings are manufactured from polyester fabric or mesh (suitable for bathing).

### Also available:

Activecare General Purpose Hygiene Slings With Head Support

## Stabel General Purpose Sling

The General Purpose sling is a standard sling providing good support for those with some upper body control. Suitable for general lifting for patient lifters with a yoke style attachment, and available in sizes from X-small to XX-large. Also with the option of fabric or mesh construction, this is a versatile sling that is simple and safe to use.



### Mesh

<b>XS</b>	27GPHS02XS	<b>XS</b>	27GPHS01XS
<b>S</b>	27GPHS02S	<b>S</b>	27GPHS01S
<b>M</b>	27GPHS02M	<b>M</b>	27GPHS01M
<b>L</b>	27GPHS02L	<b>L</b>	27GPHS01L
<b>XL</b>	27GPHS02XL	<b>XL</b>	27GPHS01XL
<b>XXL</b>	27GPHS02XXL	<b>XXL</b>	27GPHS01XXL

	XS	S	M	L	XL	XXL
<b>A</b>	820	790	920	1110	1270	1360
<b>B</b>	210	230	240	260	270	290
<b>C</b>	750	780	920	1000	1020	1050
<b>D</b>	1310	1280	1410	1630	1680	1780
<b>E</b>	585	650	700	740	780	810

### Fabric

<b>XS</b>	27GPHS01XS
<b>S</b>	27GPHS01S
<b>M</b>	27GPHS01M
<b>L</b>	27GPHS01L
<b>XL</b>	27GPHS01XL
<b>XXL</b>	27GPHS01XXL



### Mesh


<b>XS</b>	27GP02XS
<b>S</b>	27GP02S
<b>M</b>	27GP02M
<b>L</b>	27GP02L
<b>XL</b>	27GP02XL
<b>XXL</b>	27GP02XXL

	XS	S	M
<b>A</b>	790	790	890
<b>B</b>	210	230	240
<b>C</b>	470	460	530
<b>D</b>	1120	1130	1190



### Fabric

<b>XS</b>	27GP01XS
<b>S</b>	27GP01S
<b>M</b>	27GP01M
<b>L</b>	27GP01L
<b>XL</b>	27GP01XL
<b>XXL</b>	27GP01XXL

	L	XL	XXL
<b>A</b>	1110	1270	1400
<b>B</b>	260	270	290
<b>C</b>	650	730	750
<b>D</b>	1390	1520	1590

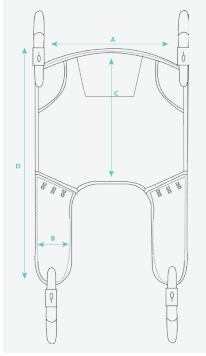

NDIS  
MASS

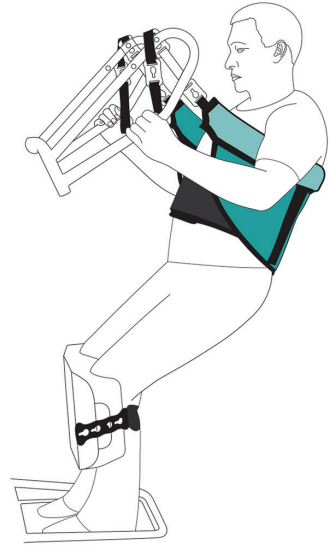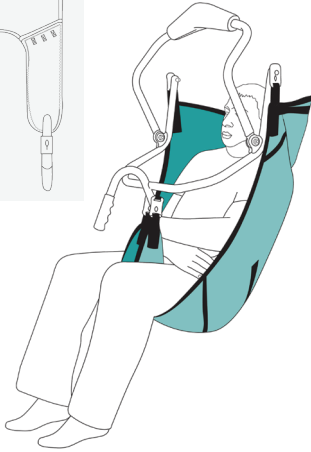

NDIS

### Stael Pivot Frame Sling

The Pivot sling is designed for use with the pivot/tilt frame. It attaches with the keyhole tags engineered to the highest standards in strength and size tolerances and feature auto clip on action (no safety straps required). It also features sling centre marking and positioning handles to assist correct positioning and woven label that will not wash out or fade.

Available in sizes from x-small through to XX-large and in fabric or mesh styles.

### Stael Salsa Sling

The Salsa sling is for use with the Salsa stand up hoist. It features optimum back support and under arm padding for a safe and comfortable lift and also allows for the easy removal of clothing for toilet transfers.

The Allegro Salsa sling features a torso band for security and stability with Buckle (b) or Velcro (v) fastening.




#### Mesh

<b>XS</b>	27P02XS
<b>S</b>	27P02S
<b>M</b>	27P02M
<b>L</b>	27P02L
<b>XL</b>	27P02XL
<b>XXL</b>	27P02XXL

#### Fabric

<b>XS</b>	27P01XS
<b>S</b>	27P01S
<b>M</b>	27P01M
<b>L</b>	27P01L
<b>XL</b>	27P01XL
<b>XXL</b>	27P01XXL

	XS	S	M	L	XL	XXL
<b>A</b>	720	870	1000	1080	1170	1320
<b>B</b>	190	230	270	260	270	290
<b>C</b>	580	830	920	1010	1090	1180
<b>D</b>	870	1220	1310	1410	1540	1710

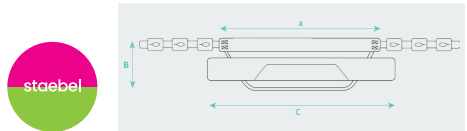

### Specifications:

Buckle Fastening				
Product Code	227701	27702	27703	27704
Size	Small	Medium	Large	Extra Large
<b>A</b>	730 mm	880 mm	1060 mm	1220 mm
<b>B</b>	260 mm	290 mm	370 mm	410 mm
<b>C (b)</b>	900 mm	960 mm	1020 mm	1220 mm
<b>C (v)</b>	1140 mm	1180 mm	1320 mm	1440 mm

Velcro Fastening				
Product Code	27720	27705	27706	27707
Size	Small	Medium	Large	Extra Large
<b>A</b>	730 mm	880 mm	1060 mm	1220 mm
<b>B</b>	260 mm	290 mm	370 mm	410 mm
<b>C (b)</b>	900 mm	960 mm	1020 mm	1220 mm
<b>C (v)</b>	1140 mm	1180 mm	1320 mm	1440 mm





DVA

NDIS

MASS



## Staebel Salsa 200 Stand Up Lifter

The Salsa 200 was designed as an active standing aid which allows a single carer to raise a patient into a comfortable standing position. The Salsa is ideal for transferring between bed and chair, or chair and commode. Allied with the Salsa, our generously padded slings allow easy transfers whilst still providing superior access to client to remove clothing for toiletry purposes. This patient hoist has keyhole sling attachment points. The Salsa200 has a removable narrow foot plate (width 350 mm) which guarantees ease of use when utilised in conjunction with shower commodes etc. With a small footprint the Salsa is highly manoeuvrable and ideal for use in an institutional/aged care environment.

Featuring the Linak Jumbo Care advanced diagnostic control system, with the latest developments in data collection.

This includes service indicator light, and with service data tool a service technician can view the history of the lifters usage including number of lifts, overloads and also program service intervals.

The Allegro Salsa lifter has been designed and tested in house to comply with AS ISO 10535 - 2002 Hoists for the transfer of disabled persons - Requirements and test methods.

The Salsa200 has the following features and benefits :

- Australian made quality product
- Electric leg adjustment
- Patient friendly lift arms enable multiple grip positions during lift
- Removable footplate for ease of transfers
- Height and angle adjustable contoured knee pad with removable calf strap
- Comes complete with battery pack and charger
- T.G.A. 136 065



### Specifications:

<b>Product Code</b>	27700	<b>Lift Arm Height</b>	880 - 1685 mm
<b>Base Length</b>	1010 mm	<b>Under Base Height</b>	55 mm
<b>Base Width Min</b>	655 mm	<b>Under Bed Clearance</b>	130 mm
<b>Base Width Max</b>	1030 mm	<b>Lifter Weight</b>	51 kg
<b>Storage Height</b>	1190 mm	<b>SWL</b>	200 kg



NDIS

MASS

## Staebel Patient Transfer

The Staebel Patient Transfer Device is the perfect solution for safe transfer of a patient or resident within the home or facility. This is your all-in-one solution without the need for additional lifting straps and belts, and promoting independence. The width adjustable base allows good access to most seating, whilst also designed for good access beneath beds etc.



### Specifications:

Product Code	212145
Seat Width	450 mm
Overall Width	640 mm
Overall Length	930 mm
SWL	150 kg

### Also Available:

Product Code	212147
SWL	220 kg



NDIS

MASS

## Staebel Transfer Platform

Staebel Transfer Platform is a manually-operated device designed as a platform for use by healthcare professionals to assist in the physical moving (transfer) of a patient in an upright (standing) position over a short distance.

The device consists of a wheeled platform, a surrounding framework, kneepad and handlebar. It can include safety belt, to hold the person being transferred in the standing position.



### Specifications:

Product Code	212150
Width	620 mm
Height	1180 mm
SWL	150 kg

NDIS

MASS

NDIS



## Stabel Transfer Belt

Stabel Transfer Platform Belt is used with the Stabel Transfer Platform when more support and stability is needed, for instance when the user's weight-bearing ability is questionable or varied. Stabel Transfer Platform Belt can also be used stand alone for sit-to-stand transfers.

- Anti-slip inner surface
- 170kg SWL
- Sizes available, XSmall, Small, Medium, Large, XLarge and XXLLarge
- Each size is colour coded for easy identification



### Specifications:

<b>Product Code</b>	212152M
<b>Size</b>	Medium
<b>Width</b>	210 mm
<b>SWL</b>	170 kg

### Also Available:

<b>Product Code</b>	212152XS	212152S	212152L	212152XL	212152XXL
<b>Size</b>	Extra Small	Small	Large	Extra Large	Extra Extra Large
<b>SWL</b>	170 kg	170 kg	170 kg	170 kg	170 kg

## Stabel Anti Shear Mattress Cover

A simple answer to the problem of constantly putting slide sheets in and out – lightweight and easy to use. You will love the Anti Shear Mattress Cover's simplicity as it makes it so much easier to adjust someone's position without disturbing them and really reduces the risk of injury for staff.

Thinner than a bed sheet, this product fits over the mattress and does not impinge on any pressure reducing qualities the mattress may have. You can sit the client/patient up in bed without the risk of them slipping down, and it only becomes fully slippery when you introduce a flat slide sheet on top to perform the manoeuvre that needs to be performed.

Just place the Anti Shear Mattress Cover over the bed mattress and, if you want to, put the normal bed sheet on top, tucking it in. When you need to reposition someone, untuck the bed sheet and unroll a flat slide sheet in between the Anti Shear Mattress Cover and bed sheet (without disturbing the patient) and use this to achieve a new position. Then slip out the slide sheet easily, retuck the bed sheet and your patient will be comfortable.


- Specialist technical fabric is fully breathable, lightweight and quiet
- Fabric is thinner than a bed sheet and will not impinge on pressure reducing qualities of mattresses
- Light-weight and compact
- Washable at 40°C to 74°C
- When using the base sheet to lie on, independent people can turn more easily in bed
- Breathable fabric means that people will not sweat or get too hot – as can happen with other systems



### Specifications:

<b>Product Code</b>	210290
<b>Width</b>	900 mm
<b>Length</b>	2000 mm

NDIS

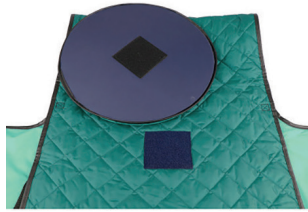

NDIS



NDIS

## Staelbel Transfer Board Slide

The Transfer Board Slide is a simple low-cost device that allows for a sitting transfer across the board without risk of shearing/induced tissue damage.

This is a wonderful little product to use with any sitting transfer board and makes sure that the slide across is what it should be – effortless. Just pop it on top of your board, tuck it under your bottom a little and sliding from one surface to another will be plain sailing.

Simple and highly effective, this Transfer Board Slide is only as wide as the board itself and has a panel of antislip that grips the top of the board.

- Transferring from one sitting surface to another
- Helping people move, both independently or with assistance



### Specifications:

Product Code	210360
Width	300 mm
Length	700 mm

## Staelbel Bed Transfer Assist

The Bed Transfer Assistant has been designed to assist in supporting the lower limbs during functional transfers.

Simple and brilliant, our Bed Transfer Assistant allows one carer to safely raise a person's legs and swivel them around into bed and vice versa!

Cotton upper for comfort.

- Raising legs and swivelling in a sitting position by a single carer



### Specifications:

Product Code	210315
Width	600 mm
Length	640 mm
Thickness	30 mm
SWL	100 kg

## EasyGlide Transfer Board

The flexible three-part sliding board is easy to place in front of the wheels of a wheelchair to protect against soiling during a sliding transfer. The smooth top surface makes it easy to move and provides for easy transfers while the double anti-slip tape on the underside prevents the board from sliding during transfer. It is especially suitable for levelling out height differences or for bridging short distances. EasyGlide is easy to transport and convenient to always have close at hands.

EasyGlide is made of strong polyethylene and is available in different sizes. The longer versions are especially suitable for transfers from wheelchair to car.

- Anti-slip protection on the underside provides for safety and security.
- The long sides can be folded to protect clothing from the wheelchair's wheels.



### Specifications:

Product Code	215020
Width	330 mm
Length	600 mm
Thickness	5 mm
SWL	135 kg

### Also Available:

Product Code	215050
Width	750 mm
Length	330 mm
Thickness	6 mm



## Staebel Car Assist

Staebel Car Assist is simple and easy to use device to allow access into and out of the car seat. There are handles on each side of the Staebel Car Assist to allow the carer to provide assistance if required. Comes standard with a seat belt to secure to the Staebel Car Assist to the car seat, to prevent sliding while in motion.

- Easy to install.
- Robust handles enable carers to provide assistance.
- Secure fitting in the car seat prevents users from sliding forward.
- Soft, padded surface provides comfort



### Specifications:

<b>Product Code</b>	210490
<b>Width</b>	450 mm
<b>Length</b>	450 mm



## Staebel Disposable Rolling Slide

The Disposable Procedure Slide is a simple slide sheet that allows health professionals to easily move patients into the lithotomy position, ready for surgery. It is ideal for head-down procedures where the patient is positioned at an incline of up to 25 degrees. Those who work in Gynaecology will find the Disposable Procedure Slide hugely beneficial as it allows medical practitioners to:

- Slide a patient into a precise position (fine adjustments achievable), usually the Trendelenburg position
- Lock patient position with the sheet's unique and integral locking panels
- Position patients at inclines up to 25 degrees, facing head-down
- Adjust patient position mid-surgery
- Slide patients back up the operating table post-surgery
- Disposable slides help teams uphold excellent infection control
- Slippery material enables easy patient repositioning
- Integral locking pads ensure absolute safety and security
- Can remain under the patient
- 25 per pack

### Specifications:

<b>Product Code</b>	210134
<b>Width</b>	1000 mm
<b>Length</b>	2000 mm

NDIS



## Camel Emergency Lifting Device

The Camel (Complete Air Moving Elevating Lift) is a state of the art air powered cushion designed to sit up and lift a fallen person from the floor.

Operating with low air pressure, the Camel offers an extremely stable and stressfree lift that can be performed either independently or with assistance from a carer.

With a maximum user weight of 320 kg (50 stone), the Camel easily lifts people of any size and age but is particularly useful for lifting bariatric people, minimising the risks associated with moving and handling. Backrest provides extra support.



### Specifications:

Product Code	2635
Height	Flat - 560 mm
Seat Depth	470 mm
Total Length	1120 mm
Width	700 mm
Weight	6.5 kg
SWL	320 kg

NDIS



## Elk Emergency Lifting Device

The ELK (Emergency Lifting Cushion) is an innovative air powered cushion designed to lift a fallen person from the floor. Operating with low air pressure the Elk offers a stable and stress-free lift that can be performed with the assistance of one carer only. The product is extremely compact making it ideal when someone has fallen in a confined space both indoors and out.

With a maximum user weight of 450 kg (70 stone) the Elk can lift people of any size, weight and age and is particularly useful for lifting very large and Bariatric people, minimising the risks associated with moving and handling for both the carer and client.



### Specifications:

Product Code	2650
Height	Flat - 560 mm
Seat Depth	500 mm
Total Length	N/A
Width	570 mm
Weight	3.6 kg
SWL	450 kg

NDIS



## Stabel Slidesheet

The ideal solution when transferring or bed repositioning patients across, up or down a bed. Reduces risk injury by assisting sliding, rather than lifting. Made from strong, non flammable, nylon fabric our slide sheets will not use their slippery quality or colour even after numerous washes. Can be machine washed to 50 degrees Celcius, tumble dried to 60 degrees Celcius and even ironed on a cool setting.

### Available Colours:

- Blue 2434B
- Green 2434G
- Orange 2434O
- Red 2434R



### Specifications:

Product Code	2434
Width	1450 mm
Length	2000 mm

NDIS



NDIS



## Staebel ComfortCare Bed Top Sheet

The ComfortCare is an in-bed sliding system for the care of the long term or 24 hour patient. The amazing ComfortCare offers patients 100% comfort and carers the ability to manoeuvre and reposition them in all directions with ease and then secure that position when needed.

- Moves a patient in all directions
- The locking panel will hold the patient steady
- Offers patients comfort and dignity
- Patients can remain on this indefinitely, either directly or under a cotton bed sheet
- Carers can re-position and manoeuvre patients without the need to disturb or handle directly
- Reduces risk of injury to carers
- Enhances patient care and clinical outcomes
- Suitable for use with hoists
- Easy to use
- Polycotton system, with low-friction inner linings
- Two panel system enables up/down manoeuvres, turning, swivelling, sitting up safely, transferring and turning into prone position
- Integral locking panel means safe position can be secured, even in sitting
- Suitable for single carers
- Washable at 40°C to 80°C and low temperature tumble dry

## Staebel ComfortCare Bed Base Sheet

The ComfortCare is an in-bed sliding system for the care of the long term or 24 hour patient. The amazing ComfortCare offers patients 100% comfort and carers the ability to manoeuvre and reposition them in all directions with ease and then secure that position when needed.

- Moves a patient in all directions
- The locking panel will hold the patient steady
- Offers patients comfort and dignity
- Patients can remain on this indefinitely, either directly or under a cotton bed sheet
- Carers can re-position and manoeuvre patients without the need to disturb or handle directly
- Reduces risk of injury to carers
- Enhances patient care and clinical outcomes
- Suitable for use with hoists
- Easy to use
- Polycotton system, with low-friction inner linings
- Two panel system enables up/down manoeuvres, turning, swivelling, sitting up safely, transferring and turning into prone position
- Integral locking panel means safe position can be secured, even in sitting
- Suitable for single carers
- Washable at 40°C to 80°C and low temperature tumble dry



### Specifications:

Product Code	210295
Width	1400 mm
Length	2000 mm
SWL	200 kg



### Specifications:

Product Code	210296
Width	900 mm
Length	2000 mm





## Staebel In Bed Turner 4-way with Handles

InBed Turner 4-Way with Handles allows safe and easy repositioning for people in bed – enabling them to turn from side to side and move up and down, with or without the aid of a hoist. They are designed for more dependent users who need assistance in bed, often from a single carer, and are designed to be used on top of one of our InBed Turner Base Sheets. They do not work on their own.

- Can remain under the patient at all times.
- Great for those who need multi positioning in bed –up / down and turning.
- Reduces risk of tissue and nerve damage.
- Helps to reduce incidence of pressure sores developing.
- The combination of both silky smooth satin panels working with each other mean easy, quiet and safe movements.
- Cotton and polyester mix for great comfort and repeated washing.
- Washable to 90C.
- Four-way movement in all directions for patient care.
- Allows more dependent people to be moved in comfort and safety.
- Cotton and satin fibres make for comfortable in-bed experience.
- Choice of torso or full length models.
- Main sizes are all 200cm wide, allowing the sheet to be tucked in after repositioning, increasing patient security.



Scan the QR code to watch the **video demonstration!**



### Specifications:

Product Code	210322
Width	2000 mm
Length	1400 mm
SWL	300 kg

### Also Available:

Product Code	210324
Width	2000 mm
Length	2000 mm
SWL	300 kg






## Staebel In Bed Turner Base Sheet

The InBed Turner Base Sheet is part of our renowned satin range for in-bed repositioning and patient care. All InBed Turner Base Sheets allow safe and easy repositioning for people in bed enabling them to turn easily from side to side but not slip down. Extra safety is provided by the striped cotton edges on all of our Base Sheets. They can be used on their own or combined with one of our Satin Top Sheets for those who need more in bed repositioning care.

- Reduces risk of tissue and nerve damage
- Helps reduce pressure sores developing
- No training required
- Easy turning and repositioning with little effort
- Great for those having difficulty turning in bed
- Washable
- Central satin panel is smooth for ease of movement
- Anti-slip cotton polyester on edges reduces risk of slipping
- Cotton polyester for enhanced comfort and repeat washing
- Single and king options are fitted sheets. The double is a flat sheet.

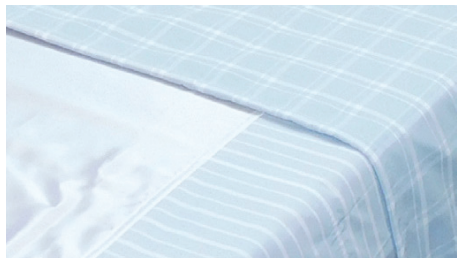

### Specifications:

<b>Product Code</b>	210321
<b>Size</b>	Single
<b>Width</b>	900 mm
<b>Length</b>	2000 mm
<b>SWL</b>	300 kg

### Also Available:

<b>Product Code</b>	210323	210325
<b>Size</b>	King Single	Double
<b>Width</b>	1050 mm	2000 mm
<b>Length</b>	2000 mm	2300 mm
<b>SWL</b>	300 kg	300 kg

NDIS



## Staelbel One Way Slide & Lock

The One Way Slide & Lock are designed to allow carers to help a user slide to the back of a chair or wheelchair and not be able to slide down.

Handles are positioned equally along the length of the product and carers can use the handles directly or attach extension straps to the handles, to prevent them having to bend down.

- To help people to sit back properly in chairs or wheelchairs without sliding back down. People independently, or with assistance by one or two carers depending upon need
- Prevent people sliding down in wheelchairs, chairs and beds
- Improves people's seating posture and comfort
- Allows people to be positioned further back in a chair
- Increases people's independence
- Can be used unassisted or with carer assistance



### Specifications:

Product Code	210190
Width	450 mm
Length	540 mm

NDIS



## Staelbel One Way Slide & Lock Open Ended

The One Way Slide & Lock are designed to allow carers to help a user slide to the back of a chair or wheelchair and not be able to slide down.

Handles are positioned equally along the length of the product and carers can use the handles directly or attach extension straps to the handles, to prevent them having to bend down.

- To help people to sit back properly in chairs or wheelchairs without sliding back down. People independently, or with assistance by one or two carers depending upon need
- Prevent people sliding down in wheelchairs, chairs and beds
- Improves people's seating posture and comfort
- Allows people to be positioned further back in a chair
- Increases people's independence
- Can be used unassisted or with carer assistance



### Specifications:

Product Code	210192
Width	450 mm
Length	1300 mm

### Also Available:

Product Code	210193
Width	550 mm
Length	1300 mm



## Staebel Slideezi with Handles

The Slideezi with Handles are manufactured using an ultra-slippery material developed by pioneering fabric specialists Slideezi®. While one side of the slide sheet is luxuriously soft to offer maximum comfort and quietness, technically advanced fibres give the other side a slippery finish, ideal for repositioning a patient in a bed or trolley.

The Slideezi with Handles allow health professionals to effortlessly manoeuvre patients/clients into a virtually endless list of positions so that they can remain comfortable at all times. Whether patients are using a bed or a trolley, health professionals can use the slide sheets to turn patients into comfortable positions and also transfer them from one surface to another.

- Superb slideability on both sides for ease of use
- Easy to insert and remove
- Superb comfort for patient during repositioning
- Quiet and soft for improved patient experience
- Fabric is manufactured using advanced technical fibres
- Effortless and easy repositioning
- Used in pairs for best performance



### Specifications:

Product Code	210432
Width	900 mm
Length	2000 mm
SWL	300 kg

## Staebel Safemove

The Safemove has been developed with practitioners and allows radiotherapists to position their patients far more precisely than has been possible until now.

Positioning patients for laser precision has always been difficult and achieving that specific position can be uncomfortable and awkward for the patient and extremely hard, tiring and heavy work for staff. It can result in potential skin damage to patients and injury risk to staff. However, by placing the Safemove under the patient (breast or pelvic areas), staff can use the straight and/or curved handles to much more easily reposition their patients to line up precisely with the lasers.

- No need to manually handle patients
- Reduced risk of tissue damage
- Less stress and discomfort for the patient
- Less strain and risk of injury to the staff



### Specifications:

Product Code	210317
Width	370 mm
Length	620 mm

## Staebel Sure Foot Pad

The Sure Foot Pad is a really helpful, inexpensive answer to the problem of feet slipping, on the floor or in bed.

Suitable for use on the floor or in bed, this Sure Foot Pad is designed as a safety aid that helps to prevent feet slipping. It can be used to stand on when sitting or rising from a chair, or it can be used with a slide sheet to aid movement in bed or with a Slide & Lock to aid movement in a chair. A helpful aid for more independent living, it is hard wearing and easy to clean.

Simply place the foot pad on the surface you need to gain more grip on and your feet won't slip. You will love its simplicity because it gives you so much independence and is also a really useful aid for rehabilitation.

Ideal for use with small slide sheets to slide further up the bed or on the floor, on its own or with products such as Slide & Locks. In fact there are so many places that you can use this product, some people even have two; those for the floor and those for use in bed!

- Gaining more floor grip with your feet
- Helping people be more independent
- A rehabilitation aid



### Specifications:

Product Code	210001
Width	1500 mm
Length	2000 mm



NDIS

NDIS



## Staebel Perfect Positioner

Do you hoist people into their wheelchairs and find it difficult to reposition them so that they are sitting properly?

Originally designed to help reposition a tetraplegic gentleman in a moulded wheelchair, this unique handling device has been redesigned and now enables carers to easily and comfortably reposition the user (once hoisted into the wheelchair) accurately at the back of the chair and ensure their pelvic tilt is correct.


- Especially useful in moulded and specialised wheelchairs
- Correct sitting position in wheelchairs
- Accurate pelvic tilt
- One carer or two

## Staebel In Bed Sit Support

This easy-to-attach bed ladder allows patients to maintain control as they sit up and lie down. Patients who would otherwise struggle to manoeuvre themselves gain a whole new level of independence with the bed rope ladder, allowing them to reposition themselves safely and easily.

You will love how simple this multicoloured functional webbing ladder is to attach. No knots to tie or worry about. Just loop it through itself and then use the different coloured rungs to help you as you sit up or lie down. Ideal if you want to maintain your independence and just need a little more help or if you are going through a process of rehabilitation.

- To help people sit up/lie down independently and in a controlled manner

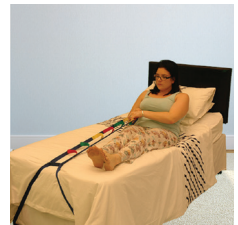

### Specifications:

Product Code	210340
Width	900 mm
Length	2000 mm

### Specifications:

Product Code	210440
Width	160 mm
Length	1880 mm
SWL	300 kg



## Staebel Support Walk Belt

A supportive Belt to help with transferring or walking. The ergonomically designed handles enable the care giver to hold the patient close without the need to grab clothes.

Ranging from Small to Large (adult sizes) this wonderful range of comfortable and supportive belts means that getting close to your patient has never been easier.

Simply place the belt around the person's waist, adjust to fit and then buckle. The ergonomically designed handles are perfectly located for user comfort and functionality. The belt is also washable.

- Helping a person from sitting to standing and vice versa
- Walking
- Supporting
- Helping whilst transferring on a sitting transfer board
- Comfortable to wear
- Give a feeling of added security
- Easily adjusted to fit



### Specifications:

<b>Product Code</b>	210082
<b>Width</b>	250 mm
<b>Length</b>	730 - 910 mm

### Also Available:

<b>Product Code</b>	210083	210084
<b>Width</b>	250 mm	250 mm
<b>Length</b>	880 - 1040 mm	1040 - 1170 mm

## Staebel Knee Cushion

The Knee Cushion will provide a comfortable cushion for those who spend a long time kneeling as part of their job and/or find it painful to kneel on the floor.

Our Knee Cushions can be used on any floor surface, inside or out, and will provide a clean and comfortable surface to kneel on whilst protecting and insulating your knees.

They are completely sealed and have been welded to ensure that they are easy and safe to wipe down.

There is no need to send these Knee Cushions to the laundry or worry about cross-infection anymore.



### Specifications:

<b>Product Code</b>	210250
<b>Width</b>	450 mm
<b>Thickness</b>	40 mm

## Staebel Support Straps - Pair

These Support Straps are used to reduce the risk of over-stretching and reaching in frequent supine transfers. No need to launder reducing the chance of losses.

A new concept in support straps, these can be wiped down and/or washed in soapy water, rinsed and then left to dry. Used to reduce the risk of over-stretching and reaching in frequent supine transfers, they simply clip onto the handles of the flat Slideezy with handles and other patient handling equipment.

- For use in supine transfers
- Reducing the risk of overreaching and stretching



### Specifications:

<b>Product Code</b>	210064
<b>Width</b>	30 mm
<b>Length</b>	750 mm
<b>SWL</b>	100 kg

NDIS



## Staebel Easy Dress Arms

A pair of long slipper sleeves that can slide over a person's arm and ease the putting on of upper body garments.

For best practice: The sleeves work especially well with a flat slide sheet underneath the bed bound person to facilitate the slide.

- Help carers put a shirt, blouse or jumper on a bed bound person/child or someone who requires a lot of assistance

NDIS



## Staebel Easy Dress Legs

A pair of long slipper leg sleeves that can slide over a person's legs and ease the putting on of lower half garments.

For best practice: The sleeves work especially well with a flat slide sheet underneath the bed bound person to facilitate the slide.

- Help carers put trousers/pyjamas on a bed bound person/child or someone who requires a lot of assistance



### Specifications:

Product Code	210159
Size	Medium
Width	280 mm
Length	1070 mm

### Also Available:

Product Code	210161
Size	Large
Width	270 mm
Length	1090 mm

### Specifications:

Product Code	210160
Size	Medium
Width	340 mm
Length	1200 mm

### Also Available:

Product Code	210162
Size	Large
Width	360 mm
Length	1300 mm





## Staebel Easy Dress Kit

Getting dressed in bed has never been easier – the Easy Dress Kit enables dignified dressing for patients who need carer assistance. This Easy Dress Kit is a complete set of in-bed dressing aids that enables one or two carers to dress a patient/client fully in bed. Where dressing is difficult or uncomfortable, they enable easy dressing and promote patient dignity.

- 2 x Easy Dress Legs
- 2 x Easy Dress Arms
- 1 x Flat Sheet (folded in a handy storage bag)
- Easy dressing aid for carers to clothe a person in bed
- Suitable as an aid for an individual who can clothe themselves but require a little help
- Super Slippy fabric for extra ease
- Available in Washable and disposable material
- Available in Standard (one size fits all)
- Washable 40° to 74° and low tumble dry
- Replacement parts available if lost



## Staebel Easy Hands

These silicon coated nylon slipper gloves are designed for carers and can be used to make patients/service users much more comfortable in many ways. For example: skin smoothing, positional adjustments and inserting slings etc.

- Technically designed fabric means the gloves slide easily on the outside but will not slide off the users hands
- Easy to slide under patients for improved patient care
- Can be used for a multitude of patient care tasks:
  - Skin smoothing
  - Tucking in sheets
  - Adjusting patient clothing
  - Moving pillows
  - Inserting slings etc
- Supplied in pairs
- Slippery outer and grippy inner material
- Fits to elbow



### Specifications:

Product Code	210339
Width	200 mm
Length	520 mm

### Specifications:

Product Code 210153 (Complete Kit)

## 03.

# Pressure Care

Action Gel Elbow / Heel Protector	61	 <b>LOW RISK MATTRESS</b>	
Obusforme Back Support – Low Back	61	Viskomattress	69
 <b>LOW RISK CUSHION</b>		Silicone Fibre Mattress Overlay	69
Low Resilient Foam Cushion	62	 <b>MEDIUM RISK MATTRESS</b>	
Relax Trio Cushion	62	Staebel Zonecare Foam Mattress	70
 <b>MEDIUM RISK CUSHION</b>		 <b>HIGH RISK MATTRESS</b>	
Action Commuter Cushion	63	Metzeler Ultratube Firm Feel Mattress	70
Action Centurian Cushion	63	Metzeler Ultratube Soft Feel Mattress	71
Action Gel Pilot Cushion	64	Curocell Area Zone Mattress	71
Action Gel Smart Pad	64	Optimal 5zon Pressure Care Mattress	72
Action Gel Professional Cushion	65	Curocell Sam Pro Air Mattress & Pump	72
Action Gel X-Act Lite	65	CuroCell A4 Cirrus Air Mattress	73
Static Air Cushion	66	CuroCell® A4 CX16 Air Mattress	74
Relax Easy Cushion	66	CuroCell® A4 CX20 Air Mattress	74
 <b>HIGH RISK CUSHION</b>		CuroCell® IQ CX15 Air Mattress	75
Relax Duogel Cushion	67	CuroCell IA Automatic CI17 Mattress	75
Relax GelCell Cushion	67	Curocell UNO Air Mattress	76
Relax Duo Profile Cushion	68	Staebel Action Gel – Segment Pad (3 Per Bed)	76
Relax GelCell G2 Cushion	68	Positioning Products	77





## Action Gel Elbow / Heel Protector

These gel protectors are for use in long term care settings and for assistance in wound management. They allow clients taking medication that affects skin integrity or who are prone to skin tears or abrasions to continue daily activities or to sleep comfortably. They are recommended for people 'at risk' or diagnosed with compromised skin disorders. These protectors require no laundering and need only a wipe over with antiseptic solution or mild detergent to be kept fresh and odour free.

- Unique polymer consistency (not liquid gel) so bottoming out cannot occur
- Helps prevent footdrop
- Quality product of USA



## Obusforme Back Support - Low Back

The ObusForme back support is a rigid framed orthopaedic back support, designed in such a way to relieve back pain by realigning the spine and supporting correct posture. It is supplied with a removable lumbar cushion and an elastic strap to secure to chair or backrest.

- Soft comfortable foam, fabric covered



### Specifications:

Product Code 3170



### Specifications:


Product Code	3030B
Width	450 mm
Height	530 mm

### Also Available:

<b>Obusforme Back Support - High Back</b>	
Product Code	3031B
<b>Obusforme Seat Black</b>	
Product Code	3250B

NDIS

NDIS

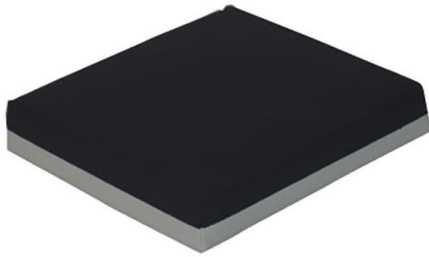

## Low Resilient Foam Cushion

A very popular choice for the low risk user, the low resilient foam cushion flows under pressure to provide a fluid like support. It helps to eliminate and reduce pressure points by transferring support evenly.

- 2-way stretch top cover
- Zipped cover included

## Relax Trio Cushion

The Relax Trio is an ultralight cushion that benefits from four sections of FloGel being positioned into a polyurethane base, with a viscoelastic foam upper layer. It is specially designed to conform to the shape of the user, increasing surface area and relieving pressure. Suitable for medium risk users, the cushion can be rotated, it has no designated front or back.

- 1.5" (38mm) thick
- Clear protective PU film coating is easy to wipe clean and is suitable for incontinent users
- Lightweight and easy to transport
- Outer cover can be machine washed up to 60° C



### Specifications:

Product Code	3210
Width	450 mm
Length	400 mm
Thickness	90 mm
SWL	100 kg

### Specifications:

Product Code	35330
Width	410 mm
Length	410 mm
Thickness	38 mm
SWL	80 kg

### Also Available:

Product Code	35332	35333
Width	450 mm	450 mm
Length	400 mm	450 mm
Thickness	38 mm	38 mm
SWL	80 kg	80 kg

NDIS  
MASS

NDIS  
MASS



## Action Commuter Cushion

The Action Commuter offers basic pressure and shear reduction for low risk users. Combined 1.27cm Akton polymer and 4.45cm foam are sealed in a water proof film for incontinence protection.

- 57mm thick
- Suitable for medium risk clients
- Not restricted by any weight limit
- Easy to care for
- Basic cover included

## Action Centurian Cushion

The Centurian offers pressure and shear reduction for low-to-moderate risk users. Combined layers of 1.6 cm Akton' polymer and 1.6 cm foam are sealed in waterproof film for incontinence protection. Offering significant comfort for active clients who need a light to medium weight cushion with low maintenance, the Centurian is also the first choice of long term care facilities for decubitus prevention and treatment.

- 32 mm thick
- Suitable for medium and low risk clients
- Easy to care for
- Not restricted by any weight limit
- Basic cover included



### Specifications:

Product Code	3560
Width	400 mm
Length	400 mm
Thickness	57 mm
SWL	Unlimited

### Also Available:

Product Code	3561	3562	3563
Width	450 mm	450 mm	500 mm
Length	400 mm	450 mm	450 mm

### Specifications:

Product Code	3570
Width	400 mm
Length	400 mm
Thickness	32 mm
SWL	Unlimited

### Also Available:

Product Code	3571	3572	3573
Width	450 mm	450 mm	500 mm
Length	400 mm	450 mm	450 mm

NDIS  
MASS

NDIS



### Action Gel Pilot Cushion

The Pilot is a low-profile cushion that provides exceptional pressure relief for most users with a pressure injury history. 2.54cm Akton polymer provides immersion for bony prominences and protects clients with fragile skin integrity. This multi-use cushion is recommended for use on manual and power wheelchairs, truck or taxi seats, or for other uses requiring significantly reduced high-pressure areas, vibration, and shear. The low profile allows users to effectively reach the floor for more efficient foot-assisted propulsion. It also conforms to pre-contoured foam or rigid cushions to maximize stability and positioning when required. Custom sizes available upon request.

- 1" (25mm) thick
- Suitable for medium-high risk clients
- Not restricted by any weight limit
- Easy to care for
- Basic cover included



### Specifications:

Product Code	3530
Width	400 mm
Length	400 mm
Thickness	25 mm
SWL	Unlimited

### Also Available:

Product Code	3531	3532	3533
Width	450 mm	450 mm	500 mm
Length	400 mm	450 mm	450 mm

### Action Gel Smart Pad

The SHEAR SMART pad, combined with the breathable double-layer SHEAR SMART cover, significantly reduces the effects of shear while providing superb user comfort. The cube construction of the pad surface allows air circulation, provides cushioning from mild impact, heat dissipation, and the comfort of sitting on 1.4 cm Akton polymer. This cushion will conform to contoured bases for additional support when required. Custom sizes available upon request (maximum of 60.9 x 60.9cm) in 2.54 cm increments.

- 14mm thick
- Suitable for medium risk users
- Not restricted by any weight limit
- Easy to care for
- Shear Smart cover included

### Specifications:

Product Code	3524
Width	400 mm
Length	400 mm
Thickness	14 mm
SWL	Unlimited


### Also Available:

Product Code	3526	3527	3528	35220
Width	400 mm	450 mm	450 mm	500 mm
Length	400 mm	400 mm	450 mm	500 mm

NDIS

MASS

NDIS

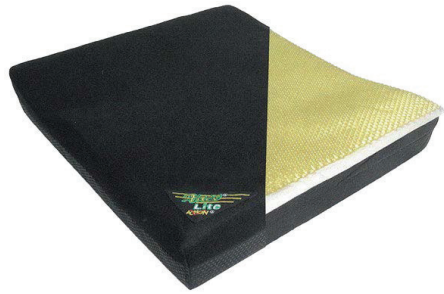

## Action Gel Professional Cushion

The Professional provides exceptional pressure relief for users at risk of skin breakdown or with a history of pressure ulcers. A full 4.1cm thickness of Akton polymer provides effective weight distribution on a contoured or flat surface. The Professional provides an additional 50% cushioning against pressure, shear and vibration. This is significant for those active clients using a power wheelchair who have concerns about bottoming out and potential effects of increased temperatures.

- Suitable for very high risk clients, including medium risk as a preventative measure
- Not restricted by any weight limit
- Easy to care for
- Basic cover included

## Action Gel X-Act Lite

Action Products XACT Cushion range continues to expand with the addition of the XACT LITE, a cushion for the active and mobile patient, including those who require sports-friendly cushions, and can self propel, transfer and shift. The new design of the AKTON polymer reduces shear and pressure and provides maximum surface flexibility with no hammocking effect. Additional features include light leg troughing which improves stability and posture without interfering with transfers, low profile design which helps foot propulsion, pre-ischial ridge to help maintain pelvic neutrality, dual-density foam base with firm rear corners to improve stability, and soft front corners to aid transfers. Cushion weight of 1.45kg.

- 54mm thick
- Easy to care for
- BodySmart cover included



### Specifications:

<b>Product Code</b>	3540
<b>Width</b>	400 mm
<b>Length</b>	400 mm
<b>Thickness</b>	41 mm
<b>SWL</b>	Unlimited

### Also Available:

<b>Product Code</b>	3541	3542	3543
<b>Width</b>	400 mm	450 mm	450 mm
<b>Length</b>	450 mm	450 mm	500 mm

### Specifications:

<b>Product Code</b>	3580
<b>Width</b>	400 mm
<b>Length</b>	400 mm
<b>Thickness</b>	54 mm
<b>SWL</b>	114 kg

### Also Available:


<b>Product Code</b>	3581	3582	3583
<b>Width</b>	450 mm	450 mm	500 mm
<b>Length</b>	400 mm	450 mm	450 mm

NDIS

MASS

NDIS

MASS

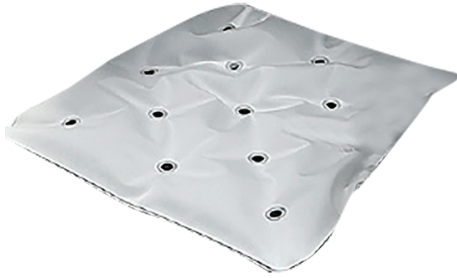

## Static Air Cushion

CuroCell® S.A.C. is a static air seat cushion for comfort, prevention and as an aid in the treatment of pressure injury up to category 2.

The cushion pre-filled and very easy to handle. The durable soft material for long term use provides a superior combination of comfort and support.

- Wipe with cleaning agent and/or disinfectants
- Air holes for venting heat and moisture
- Pre-inflated low profile design
- Fits all chairs
- Durable soft material for long term use
- Patient ready

**\*Pressure Mapping Reports and successful Case Studies Available on request**



### Specifications:

Product Code	3417
Width	450 mm
Length	450 mm
SWL	300 kg

### Also Available:

Product Code	3418
Width	500 mm
Length	500 mm
SWL	300 kg



## Relax Easy Cushion

As the name suggests, this cushion from the Relax range is easy to use and is an ideal pressure relieving cushion for use in the nursing home environment. It has no designated up or down, or front or back, so can be simply placed into position. It features a layer of FloGel sandwiched between an upper and lower layer of viscoelastic foam. This combination of pressure relieving FloGel and viscoelastic foam facilitates effective pressure re-distribution while allowing the patient to immerse themselves comfortably into the cushion. It conforms to the shape of the user, increasing surface area and relieving pressure.

- 2" (50mm) thick
- Suitable for medium - high risk clients
- Clear protective PU film coating is easy to wipe clean and is suitable for incontinent users
- Lightweight and easy to transport
- Includes waterproof wipeable PU stretch outer cover
- Outer cover can be machine washed (60°C)

\* High profile cushion available

### Specifications:

Product Code	35100
Width	400 mm
Length	400 mm
Thickness	50 mm
SWL	130 kg

### Also Available:

Product Code	35312	35313	35315
Width	450 mm	450 mm	500 mm
Length	400 mm	450 mm	450 mm

NDIS  
MASS

DVA  
NDIS



## Relax DuoGel Cushion

The DuoGel is a unique combination of TruGel and FloGel work together effectively to provide excellent pressure relief and weight distribution. This gel and foam combination allows for greater circulation and weight distribution. The DuoGel provides for medial and lateral thigh positioning as well as support for the greater trochanter of the thigh bone. The thigh troughs are lined with TruGel and the ischial area is lined with FloGel which is set on a viscoelastic foam giving greater comfort and seating stability.

- 3" (75mm) thick
- Suitable for high risk clients
- Clear protective PU film coating is easy to wipe clean and is suitable for incontinent users
- Lightweight and easy to transport
- Includes machine washable cover

## Relax GelCell Cushion

A unique combination of AIR and GEL technology effectively provides pressure relief and comfort to the high risk user. The Gelcell consists of 34 interconnected air-filled cells, each covered with a layer of pressure relieving silicone-based TruGel. The action of the air in the cushion enhances capillary function in the skin, maintaining blood flow and oxygen uptake. When pressure is applied, air circulates to each of the cells, distributing weight evenly, thus reducing high pressure concentrations around bony prominences. It ensures enhanced pressure re-distribution and also has a positive effect on skin temperature management whilst minimizing the effects of shear forces.

- 2" (50mm) thick
- Suitable for high risk clients
- Easy to wipe clean and is suitable for incontinent users
- Lightweight and easy to transport
- Easy to inflate and adjust
- Includes outer cover and pump



### Specifications:

Product Code	35300
Width	400 mm
Length	400 mm
Thickness	75 mm
SWL	130 kg

### Also Available:

Product Code	35306	35303	35304
Width	400 mm	450 mm	508 mm
Length	450 mm	450 mm	450 mm

### Specifications:


Product Code	35320
DVA	35300
Width	400 mm
Length	400 mm
Thickness	50 mm
SWL	130 kg

### Also Available:

Product Code	35323	3532320
Width	450 mm	500 mm
Length	450 mm	500 mm

NDIS

NDIS

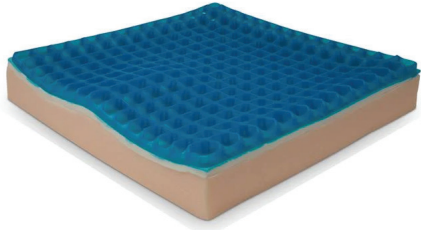

## Relax Duo Profile Cushion

The Duogel is a unique combination of TruGel and FloGel that work together effectively to provide excellent pressure relief. This gel and foam combination allows for greater circulation and weight distribution. The Duogel provides for medial and lateral thigh positioning as well as support for the greater trochanter of the thigh bone. The thigh troughs are lined with TruGel and the ischial area is lined with FloGel which is set on a viscoelastic foam giving greater comfort and seating stability.

- 3" (75mm) thick
- Suitable for high risk clients
- Clear protective PU film coating is easy to wipe clean and is suitable for incontinent users
- Lightweight and easy to transport
- Includes machine washable cover



### Specifications:

Product Code	35380
Width	400 mm
Length	400 mm
Thickness	75 mm
SWL	130 kg



## Relax Gelcell G2 Cushion

The Gelcell G2 consists of 48 interconnected cells, each covered with a layer of pressure relieving silicone-based TruGel

- Enhances pressure re-distribution
- Minimises the effects of shear forces
- Has a positive effect on skin temperature
- Silicone gel is impact dampening

The cells are contoured towards the rear, promoting good postural tilt. When pressure is applied, air circulates to each cell in both chambers distributing weight evenly. This process reduces high pressure concentrations around bony prominences, preventing skin breakdown.

The Gelcell is made up of 3 layers of polyurethane (PU) film making it highly durable and difficult to puncture. It has also easy to wipe clean and suitable for incontinent users.

- 1.5" (38mm) thick
- Suitable for high risk clients
- Easy to wipe clean and is suitable for incontinent users
- Lightweight and easy to transport
- Easy to inflate and adjust
- Includes outer cover and pump




### Specifications:

Product Code	35392
Width	400 mm
Length	400 mm
Thickness	38 mm
SWL	130 kg

### Also Available:

Product Code	35394	35395
Width	450 mm	400 mm
Length	450 mm	450 mm



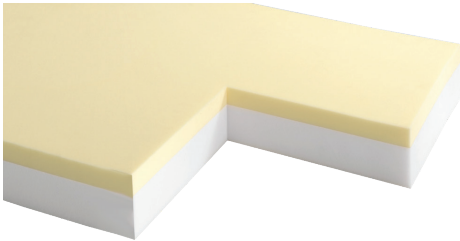

## Viskomattress

The Viskomattress includes a double layer foam mattress with a 40mm top layer of visco elastic foam, suitable for those spending a long time in bed, providing optimum pressure distribution across all body area in contact with the mattress.

- Suitable for medium risk patients
- Two way stretch PU waterproof cover
- Full length flap conceals the zip
- Fully removable zipped cover
- Mattress cover is washable (up to 95°C) and can be tumble dried
- Moisture protection properties mean the cover has a high resistance to bacteria, viruses and fungi
- Mattress does not require turning
- 12 month warranty

## Silicone Fibre Mattress Overlay

This overlay is made of millions of silicone coated cylindrical hollow cored fibres that provide excellent pressure relief while the silicone coating acts as a lubricant to eliminate friction and reduce shear forces. These fibres do not absorb perspiration or urine, making this an excellent product for geriatric and palliative care with elastic corner retention straps.

- Suitable for low risk patients
- Soft breathable cotton cover on top, PVC on bottom
- Machine washable



### Specifications:

Product Code	3755
Width	850 mm
Length	1980 mm
Thickness	140 mm
SWL	120 kg

### Also Available:

Product Code	3755W
Width	1000 mm
Length	2000 mm
SWL	120 kg




### Specifications:

Product Code	3740
Width	920 mm
Length	1900 mm

NDIS

MASS

NDIS

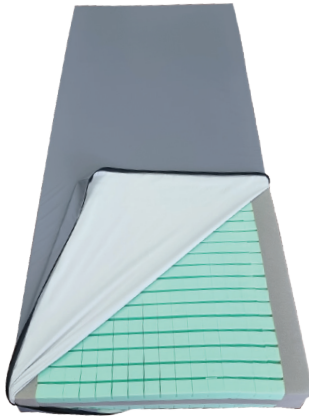

## Stäbel Zonecare Mattress

The Zonecare foam mattress is a high quality foam mattress featuring a profile-cut pressure relieving foam insert. It provides pressure reduction and assists in the treatment of pressure ulcers, whilst enhancing patient comfort levels in homecare and residential care environments. The Eurocare's specialised inner cover aids in shear reduction on a profiling bed.

- Unique reinforced edges designed with longevity in mind
- Mattress turning is not required
- Waterproof zip for infection control
- Mattress Cover is vapour permeable, water resistant and fire retardant
- Mattress cover has welded seams to meet strict infection control standards
- Inner anti-shear cover



### Specifications:

Product Code	3728Z
Width	850 mm
Length	1980 mm
Thickness	150 mm
SWL	150 kg

### Also Available:

Product Code	3729Z
Width	950 mm
Length	1980 mm
SWL	150 kg



## Metzeler Ultratube Firm Feel Mattress


The heart of the Metzeler Ultratube Mattress is the newly developed Rubex Tubes System. The cold formed foam tubes in the mattress form a vertical tube system for perfect air circulation. The tubes, which differ in colour according to their degree of hardness, are integrated into the complete surface. This guarantees the user an ergonomically perfect lying position. Up to 200 litres of perspiration a year leave the body and go into the mattress. The Rubex Tubes System dispels this moisture through the vertical air channels.

- The tubes make up a system of s for perfect air circulation during sleep. The tubes of varying degrees of firmness are arranged to provide an ergonomically ideal sleeping posture reducing the risk of injury caused by shear and friction.
- Reinforced edge makes it easy to get in and out of bed and helps prevent roll-out.
- Standard with a bi-elastic, two-way stretch, waterproof, cover, designed to reduce shear and friction forces.
- Stain resistant, fire-resistant, waterproof, antiallergenic, antimicrobial, hygienically clean. Can be machine washed up to 95°C. Antibacterial treated.
- The Rubex Tubes System is symmetrically structured so that you always lie straight.
- Fully washable in temperatures up to 60°C



### Specifications:

Product Code	3892
Width	1050 mm
Length	2000 mm
Thickness	180 mm
SWL	250 kg

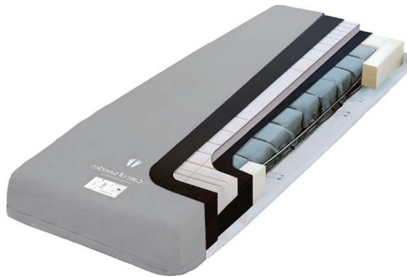

## Metzeler Ultratube Soft Feel Mattress

This mattress can be used on both sides offering the same pressure relief and comfort level; regular turning will assist with long-term mattress performance.

- The Metzeler Mattress provides the perfect night's sleep. Incorporating patented Rubex technology, the mattress provides superior comfort and support offering optimal pressure distribution, ventilation and spinal alignment.
- Reinforced foam edges assist carers and physios with patient transfer, while minimizing the potential for rolling back into or off the mattress.
- The core of this mattress has a unique tubular structure that is zoned to support the spine and provides optimal pressure redistribution.
- The waterproof and breathable visco-elastic cover assists with moisture management, while being easily removable and washable up to 95°C.
- The tubular structure and a choice of soft or medium foam castellated overlay on top and bottom of the inner core provide optimal horizontal and vertical cross-ventilation thermo-regulation needs.
- Unlike most foam mattresses made from petroleum byproducts, Metzeler mattresses use natural, non-toxic materials like sunflower oil and other environmentally-friendly alternatives.



### Specifications:

Product Code	3893
Width	1050 mm
Length	2000 mm
Thickness	180 mm
SWL	250 kg

## Curocell Area Zone Mattress

The AREA Zone mattress replacement system is designed with special valve technology to provide an enhanced mattress surface for the prevention and management of skin breakdown. The AREA provides clinically proven non-powered pressure redistribution and therapeutic support through the use of interconnected pressure redistribution air sectors. The valve system allows for communication between the air chambers which in turn provides a unique, custom low pressure profile for any patient. Upon ingress, a natural chain of events begins in each of the connected air chambers. The slow release valves allow for pressure changes to take place very gently. The mattress top is a layer of visco type memory foam which further enhances the low pressure profile.

- Suitable for high risk patients, and also medium risk as a preventative measure
- Static air replacement mattress system (not an overlay)
- Reinforced sides for improved lateral stability
- Unique integrated heel protection area
- Anti-shear inner cover
- No pump or power cord for hassle free use



### Specifications:

Product Code	3679Z
Width	850 mm
Length	2000 mm
Thickness	180 mm
SWL	230 kg

### Also Available:

Product Code	3679Z	3686Z	3680Z
Width	900 mm	1000 mm	1050 mm
Length	2000 mm	2000 mm	2000 mm
SWL	230 kg	230 kg	230 kg

NDIS

MASS



## Optimal 5zon Pressure Care Mattress

The Optimal 5zon® is used in all types of health care environments for pain relief and as an aid in the prevention and treatment of pressure ulcers up to and including category 2. A soft and flexible mattress replacement system. The Optimal 5zon® offers very good pressure relief properties and is now also equipped with the unique sit-function DuoCore®. This is an insert within the mattress core, that provides support when seated, while offering good comfort in the reclining position. It comes complete with a durable cover that is easy to clean and is equipped with divisible, drip shield zip. The soft and elastic material is vapour permeable, which reduces the risk of skin shear.

- Core: HR-foam 65+50 kg/m3
- Function: as an aid for prevention and treatment up to and including category 2
- Cleaning of cover: wipe with cleaning agent and/or disinfectants.
- Machine wash max 95 °C, tumble drying



### Specifications:

Product Code	3771
Width	850 mm
Length	2000 mm
Thickness	140 mm
SWL	180 kg

### Also Available:

Product Code	3778	3775
Width	1000 mm	1050 mm
Length	2000 mm	2000 mm
SWL	180 kg	180 kg

NDIS

MASS



## Curocell Sam Pro Air Mattress & Pump

The CuroCell SAM PRO is a static air replacement mattress which has been proven to provide exceptional comfort, pressure and pain relief. It is used preventively and in the treatment of pressure ulcers up to and including sore degree 3. The mattress is divided into three zones, which amongst other advantages gives better support when sitting on the mattress. The support pump delivered together with the mattress is easy to use and very quiet. The CuroCell SAM PRO is delivered with a anti-shear inner cover and a soft, elastic, breathable outer cover which has a very well documented hygienic standard.

- Suitable for high risk patients, and also medium risk as a preventative measure
- Dynamic air mattress replacement system (lies directly on bed surface)
- Air cells produced from soft, durable polyurethane Not PVC. Significant noise reduction when cells move compared to other brands
- Comes complete with an anti-shear inner cover which has been proven to reduce pressure by up to 8%
- Easy to use dial pump with table to calculate most appropriate setting
- 2 year manufacturers warranty
- Quality product from Europe



### Specifications:

Product Code	3683
Width	850 mm
Length	2000 mm
Thickness	160 mm
SWL	250 kg

### Also Available:

Curocell Sam Air Mattress Overlay	
Product Code	3670



NDIS



## CuroCell A4 Cirrus Air Mattress

CuroCell® A4 Cirrus is used as an aid in the prevention and treatment of pressure ulcers up to and including category 4 and unclassified pressure ulcers/pressure injuries<sup>(1)</sup>. It is a replacement system suitable for all healthcare environments.

CuroCell® A4 Cirrus combines air with a soft foam surface and provides good user comfort and pressure redistribution<sup>(1)</sup>. The added heel function offers pressure redistribution on the heels.

- Automatic control unit which adapts individual to the user
- User friendly interface
- Audible and visual notifications
- Three program modes: pulsating mode, alternating mode and constant low pressure
- Silent running and vibration free
- Air flow control; Constant temperature, good user comfort, energy reduction
- Soft foam surface
- Smart cable management
- Easy handling and cleaning
- Heel slope, for pressure redistribution of heels
- CPR-function with transport function



### Specifications:

Product Code	3841	
Width	850 mm	
Length	2000 mm	
Thickness	210 mm	
SWL	230 kg	

### Also Available:

Product Code	3843	3844
Width	1000 mm	1050 mm
Length	2000 mm	2000 mm
SWL	230 kg	230 kg

# Case Study

## History

- Female, born 1932
- Operated breast cancer 2003, vertebral compression fractures, hip osteoarthritis with chronic pain
- Arterial fibrillation, hypertension and bilateral leg ulcers.
- Pressure ulcer category 4 in sacrum 7x2 cm with fistulas, R/T ongoing treatment leg ulcer.
- Living at home with help from home health care.
- Mentally clear, no incontinence, normal diet, sitting up during meals.
- Taking regular pain relief, Oxykontin and Oxynorm.
- Points 23, modified Norton scale.
- BMI 21



## Action

- 2010-12-31 prescribed CuroCell Cirrus, alternating pressure mattress with foam core, already have a pressure relieving wheelchair cushion.
- Can manage repositioning by herself, aided move from bed to wheelchair. Wounds redressed daily. Daily supplements of nutrient drinks.
- 2011-01-21 admitted to hospital due to the leg ulcers, mattress goes with the patient.
- 2011-02-18 left leg amputated thigh bone After hospital stay the patient is accommodated short-term at a home and is later placed in a permanent nursing home, mattress goes with the patient. Wounds is decreasing in size.
- 2012-04-04 Wounds completely healed



## CuroCell® A4 CX16 Air Mattress

The CuroCell® A4 CX16 provides simplicity together with user comfort while effectively aiding in prevention and treatment of pressure injuries up to and including category 4, and unclassified pressure injuries. It is designed for excellent user comfort, reduction of shear forces and improvement of the microclimate.

The air mattress is 16cm thick with a foam base for added comfort. The alternation of the cells on the CX16 is every second cell down ensuring 50% of the body is being supported at one time. The pump is intelligent and fully automatic without need for manual settings. It identifies the user and automatically manages all settings to optimise the pressure according to the user's needs.

All mattresses are supplied with a four way stretch hygienic cover which meets standards for hygiene, quality and safety. All covers are manufactured from materials which reduce the risk of skin maceration.

- Unique CPR-function with transport function
- Heel Function
- Silent running and vibration free
- Three program modes: Pulsating, Alternating and CLP (constant low pressure)
- Individually replaceable cells with safety mattress in HR-foam base
- Audible and visual notifications
- Smart cable management



## CuroCell® A4 CX20 Air Mattress

The CuroCell® A4 CX20 provides simplicity together with user comfort while effectively aiding in prevention and treatment of pressure injuries up to and including category 4, and unclassified pressure ulcers. It provides features for excellent user comfort, reduction of shear forces and improvement of the microclimate.

The CX20 mattress is 20 cm thick and the alternation of the cells is every third cell down ensuring 75% of the body is being supported at one time. The pump is intelligent and fully automatic without need for manual settings. It identifies the user and automatically manages all settings to optimise the pressure according to the user's needs.

All mattresses are supplied with a four way stretch hygienic cover which meets standards for hygiene, quality and safety. All covers are manufactured from materials which reduce the risk of skin maceration.

- Unique CPR-function with transport function
- Heel Function
- Audible and visual notifications
- Silent running and vibration free
- Three program modes: Pulsating, Alternating and CLP (constant low pressure)
- Individually replaceable air cells with safety mattress
- Smart cable management



### Specifications:

Product Code	3627
Width	1000 mm
Length	2000 mm
Thickness	160 mm
SWL	200 kg

### Also Available:

Product Code	3626	3628
Width	850 mm	1050 mm
Length	2000 mm	2000 mm

### Specifications:

Product Code	3642
Width	1000 mm
Length	2000 mm
Thickness	150 mm
SWL	250 kg

### Also Available:

Product Code	3641	3643
Width	850 mm	1050 mm
Length	2000 mm	2000 mm




NDIS

MASS



NDIS

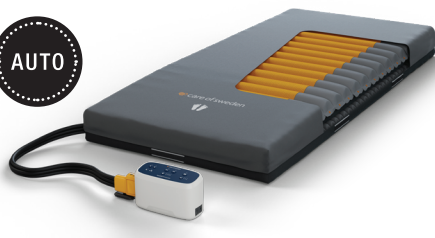

## CuroCell IQ CX15 Air Mattress

The CuroCell® A4 CX15 provides simplicity together with user comfort while aiding in the prevention and treatment of pressure ulcers/pressure injuries up to and including category 4, and unclassified pressure ulcers(1). It is designed for simplicity and provides features for good user comfort, reduction of shear forces and improvement of the microclimate.

CuroCell® A4 CX15 has an automatic control unit with a user-friendly interface. The mattress system regulates independently and without manual adjustment the inner pressure of the mattress to different values based on the user's weight and position. No manual action needs to be done to adjust the internal pressure of the mattress.

All mattresses are supplied with a four way stretch hygienic cover which meets standards for hygiene, quality and safety. All covers are manufactured from materials which reduce the risk of skin maceration.

- Individually replaceable cells with safety mattress
- Static head cells
- Heel function
- Multifunctional inner cover
- CPR-function with transport function
- Smart cable management
- Easy handling and usage
- Well protected zipper

## CuroCell iA Automatic Ci17 Mattress

With Gentle Alternating Low Pressure (GALP), the CuroCell iA® system makes the step from mechanical to digital, from manual to autonomus, and from on-off comfort valves to a fine tuning of total support that we've never been able to deliver before in this segment. GALP has been developed to enable a system to sense, learn and take action whenever the high pressure that can cause ulcers/injuries arises.

CuroCell iA® uses artificial intelligence to constantly monitor pressure across the patient's contact area and calculates its response. The system identifies the patient's weight and height and can detect changes, such as when the patient adopts a new position.

All mattresses are supplied with a four way stretch hygienic cover which meets standards for hygiene, quality and safety. All covers are manufactured from materials which reduce the risk of skin maceration.

- Gentle Alternating Low Pressure (GALP) and Constant Low Pressure (CLP) and Maximum pressure (Caring mode)
- Panel Lock
- Information signals
- CPR connection with transport function
- Air Flow Control - saves energy and creates an even temperature in the support surface
- Low noise-and vibration level
- Individual and replacement cell construction
- Heel Function



### Specifications:

Product Code	3847
Width	850 mm
Length	2000 mm
Thickness	150 mm
SWL	220 kg

### Sizes Available:

Product Code	3848
Width	1000 mm
Length	2000 mm

### Specifications:

Product Code	3928
Width	850 mm
Length	2000 mm
Thickness	170 mm
SWL	210 kg

### Sizes Available:

Product Code	3929	3930	3931
Width	900 mm	1000 mm	1050 mm
Length	2000 mm	2000 mm	2000 mm

NDIS

MASS

NDIS

MASS



## Curocell UNO Air Mattress

The CuroCell UNO® is used as an aid in the prevention and treatment of pressure injuries up to and including category 4. It stands for innovation and modern design and together with quality engineering it lays the foundation for our new generation of fully automated premium mattress systems. It is equipped with a large touch screen panel with a user-friendly interface. The fully automated pump offers automatic adjustment, providing optimal pressure relief and comfort for each individual user while aiding simple use.

- Cell on Cell construction with safety mattress
- Heel Function
- Static head cells
- CPR-function for rapid deflation
- Smart cable management
- Chlorine resistant (10 %) multi stretch cover with welded seams
- Well protected zipper
- Care of Sweden's mattresses are supplied with a four way stretch hygienic outer cover which conforms to relevant standards for hygiene, quality and safety.
- Bariatric capability



### Specifications:

Product Code	3651
Width	850 mm
Length	2000 mm
Thickness	200 mm
SWL	250 kg

### Also Available:

Product Code	3665	3652
Width	1000 mm	1050 mm
Length	2000 mm	2000 mm



## Stabel Action Gel - Segment Pad (3 Per Bed)

The Action bed overlays (pictured) are made of Akton polymer which mimics human fatty tissue in consistency, does not migrate like gel and provides equally distributed pressure management, with varying sizes available, these overlays can be used to cover the entire surface or isolated areas of an existing standard mattress. Patients can use one, two or three segments in combination on their areas of need e.g. shoulders and hips, feet and hips or full body etc.

- Easily cleaned, no maintenance
- Not restricted by any weight limit
- 3 Segments cover length of single bed. Each sold separately



### Specifications:

Product Code	3853
Width	630 mm
Length	690 mm
Thickness	18 mm
SWL	Unlimited






## Stabel Positioning & Support Wedge

This foam positioning support wedge is used to place a user in a side position in order to relieve the sacrum, which is an extra sensitive area. For the side position, two support wedges are used together. The support wedge is easy to use and easy to clean. It is supplied with a zippered, liquid resistant and vapour permeable cover. The cover can be wiped or machine washed.

- Soft comfortable foam, stretch vinyl covered

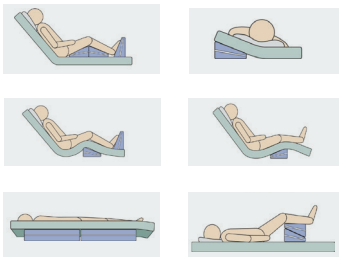

### Specifications:

Product Code	3689
Width	260 mm
Length	700 mm
Height	50 - 150 mm



## Back Pillow

The Back Pillow is a winged-pillow for back relief. The pillow is placed behind the back, and the wing rests under the user to set the position.



### Specifications:

Product Code	32301
Width	500 mm
Length	800 mm

NDIS

NDIS



## Positioning Roll

The positioning roll can give support to the resident in sitting and lying positions. It can be placed in various ways to provide support and relief.

- Supports shoulders, spine, hips, buttocks and heels
- All-round support in seated and reclining positions
- Relieves the knees and supports feet
- Supports in lateral position



## Positioning Pillow

The Positioning Pillow range is an all-round pillow that can be used for different types of patient positioning.

- Support after a hip replacement
- Can be positioned to help breathing problems and general expansion of the lungs
- Positioning after a stroke
- Supports in lateral position

**\* Available in different sizing**



### Specifications:

Product Code	32307
Width	150 mm
Length	600 mm

### Also Available:

Product Code	32308	32309
Width	200 mm	200 mm
Length	1250 mm	2250 mm

### Specifications:

Product Code	32304
Width	250 mm
Length	800 mm



## Multi Pillow

The Multi Pillow is a versatile pillow that can be placed between the knees or under legs and heels for support and relief.

- Supports the feet, counters heel pressure against surfaces
- Locked position in connection with spasticity
- Supports in lateral position
- Relief for heels
- Supports under the knees

## Molly Heel Protection

- Easy to use
- Equalizes pressure
- Faster healing
- Improves circulation
- Makes it easier to freely move the legs
- Prevents crossing the legs after hip surgery
- Soft and comfortable
- Wipe with cleaning agent or disinfect
- Machine wash max. 95°C
- Easy to clean
- Liquid resistant
- Long durability



### Specifications:

Product Code	32302
Width	450 mm
Length	450 mm

### Specifications:

Product Code	32311
Leg Length	780 - 850 mm
Thickness	75 mm

### Also Available:


Product Code	32312	32313
Leg Length	860 - 940 mm	1020 - 1030 mm
Thickness	90 mm	110 mm

# Mattress Matrix

Mattress Name	Size	Comments	Waterlow scale / Risk Score		
			At Risk (10+)	High (15+)	Very High (20+)
<b>Zoncare Pressure Care Foam Mattress</b> 	Single	Standard mattress for all replacements	✓		
	King Single				
<b>Optimal 5zon Pressure Care Mattress</b> 	Single	Highly elastic, high quality material and a unique zone and module system gives effective pressure. Specifically designed to minimise the effects of shear forces. Duocore is designed to give support when sitting and provide good comfort when lying. The cover is vapour permeable, which reduces the risk of skin maceration.  The covers are easy to clean and handle and has a very well documented hygienic standard.	✓		
	King Single				
<b>Curocell S.A.M F16 Pro Static Air Mattress Replacement &amp; Pump</b> 	Single	A protective / preventative strategy for residents with intact skin but identified as at risk of skin integrity issues of pressure injuries.  Beautifully comfortable due to the static nature plus it comes with the comfort cover so you do not create Cat 1-3 pressure injuries.		✓	✓
	King Single				
<b>Curocell A4 CX16</b> 	Single	Cat 1- 4 Pressure Injuries. CX16 is a 3 cell cycle, opposed to the regular 2 cell. 2 cells up, 1 down means the body has 66% support not 50%, however the pulsating is the best program for the greatest comfort.  Great for persons who need the pressure care, but don't like the sensation. Excellent for where high incontinence issues are present.		✓	✓
	King Single				
<b>CuroCell Cirrus</b> 	Single	Cat 1- 4 Pressure Injuries. Uniquely designed to help in pain management, for end of life care. Thick top layer of foam to help create sensory recognition of the user with pulsating mode cells below that create a wave like relief effect opposed to the traditional alternation.  Heel slope, silent pump with sit position function, and option to adjust the timing of the cell change.		✓	✓
	King Single				

Current Pressure				May suit people with particular needs:				
Skin Intact	Skin Broken	Underweight	Bariatric	Client In Pain	Bed Mobility		Bed Transfers	Heels at risk
					Independent	Assisted / Static		
✓			 0-125KG		✓		Firm Transfers	
	 CATEGORY 2	(~30KG)	 0-180KG	✓	✓	✓	YES - the Duo Core adds stability	Zero pressure ability
	 CATEGORY 3	(~30KG)	 0-250KG	✓	May be difficult	✓	Staebel as it is a static mattress	Inbuilt heel slope
	 CATEGORY 4	(~1KG)	 0-200KG	✓	May be difficult	✓	Max function	Inbuilt heel slope
	 CATEGORY 4	(~1KG)	 0-200KG	✓	May be difficult	✓	Foam surround	Inbuilt heel slope

# Cushion Matrix

	Weight Capacity	Cushion Profile Height	Size Range	
			Width - inches	Width - mm
<b>3210 Budget</b> 	100 kg	50 mm		
<b>SAC</b> 	160 kg	25 mm		
<b>Action Pilot</b> 	Unlimited	30 mm		
<b>Action Professional</b> 	Unlimited	45 mm		
<b>Relax Easy Low Profile</b> 	130 kg	50 mm	14 x 16" to 22 x 18"	360 x 400mn to 550 x 450mn
<b>Relax Easy High Profile</b> 	190 kg	85 mm	16 x 16" to 20 x 20"	400 x 400mn to 500 x 500mn
<b>Duogel</b> 	130 kg	85 mm	14 x 14" to 24 x 18"	360 x 360mn to 600 x 450mn
<b>Duo Profile</b> 	130 kg	85 mm	16 x 16" to 18 x 18"	400 x 400mn to 450 x 450mn
<b>Gel Cell</b> 	130 kg	40 mm	15 x 14" to 20 x 18"	380 x 360mn to 500 x 450mn

PI Rating	Thigh Support	Contoured seat	Skin Integrity			Suitable for Lounge chairs
			Skin Intact	Skin Broken	Boney Prominence	
Prevention	✗	✗	✓	✗	✗	✗
Stage 2	✗	✗	✓	✗	✗	✓
Stage 3	✓	✗	✓	✗	✗	✓
Stage 4	✓	✗	✓	✓	✓	✗
Stage 2	✗	✗	✓	✗	✓	✓
Stage 3	✗	✗	✓	✓	✓	✗
Stage 3	✓	✓	✓	✓	✓	✗
Stage 2	✓	✓	✓	✗	✓	✗
Stage 4	✗	✓	✓	✓	✓	✓

## 04.

# Wheelchairs

Staebel Liberation Lightweight Wheelchair	85
Staebel Emperor Wheelchair	85
Staebel Emperor Wheelchair Transit	86
Staebel Thrifty Wheelchair	86
Staebel Sovereign Wheelchair	87
Staebel Deluxe Alloy Wheelchair	88
Staebel Gusto Heavy Duty Wheelchair	89
Accessories	90





## Stabel Liberation Lightweight Wheelchair

The Stabel Liberation Lightweight Wheelchair is a great travel companion. It is a sleek modern design that has been kept compact for easy transport and storage.

- Fully foldable backrest
- Dual brake
- Parking brake push handle
- Fold-up armrests
- Lightweight foldable swing-away footrests



### Specifications:


Product Code	40400
Seat Width	450 mm
Overall Width	540 mm
Overall Height	920 mm
Product Weight	7 kg
SWL	100 kg



## Stabel Emperor Wheelchair

The Stabel Emperor is a comfortable, lightweight aluminium wheelchair, with a really helpful range of high quality features. With its sleek modern frame and lightweight foldable design, the wheelchair is perfect for travel. All features of this chair have been designed to make life easier for the user.

- Foldable backrest
- Lap belt
- Heel strap
- Height adjustable armrests and footplates
- Removable swing-away footrest
- Available in orange (40416) or blue (40411)

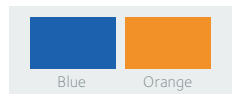

### Specifications:

Product Code	40411
Seat Width	450 mm
Overall Width	670 mm
Overall Height	950 mm
Product Weight	12 kg
SWL	115 kg

### Also Available:

Product Code	40412
Seat Width	500 mm
Overall Width	720 mm
Overall Height	950 mm
Product Weight	15 kg
SWL	115 kg

NDIS


NDIS



## Stael Emperor Wheelchair Transit

The Stael Emperor is a comfortable, lightweight aluminium wheelchair, with a really helpful range of high quality features. With its sleek modern frame and lightweight foldable design, the wheelchair is perfect for travel. All features of this chair have been designed to make life easier for the user.

- Foldable backrest
- Lap belt
- Heel strap
- Height adjustable armrests and footplates
- Removable swing-away footrest
- Available in orange or blue.

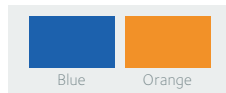

### Specifications:

<b>Product Code</b>	40413
<b>Seat Width</b>	450 mm
<b>Overall Width</b>	670 mm
<b>Overall Height</b>	950 mm
<b>Product Weight</b>	12 kg
<b>SWL</b>	115 kg

### Also Available:

<b>Product Code</b>	40414
<b>Seat Width</b>	500 mm
<b>Overall Width</b>	720 mm
<b>Overall Height</b>	950 mm
<b>Product Weight</b>	15 kg
<b>SWL</b>	135 kg

## Stael Thrifty Wheelchair

The Thrifty is a sturdy and reliable standard weight wheelchair, with removable armrests and leg rests, built for individuals, nursing homes and non-prescription requirements. This wheelchair meets Australian Standards for Manual Wheelchairs in Product Performance. All of these features combine to make the Thrifty of unbeatable value in standard wheelchairs.

- Swing away, nonremovable armrests
- Swing away, detachable legrests
- 24" rear wheels with urethane tyres ( maintenance free )
- Australian Standard AS/NZS 3695.1:2011 approved
- Accessories available ( see page 90 )



### Specifications:

<b>Product Code</b>	4246
<b>Seat Width</b>	450 mm
<b>Seat Depth</b>	400 mm
<b>Overall Width</b>	650 mm
<b>Overall Depth</b>	1080 mm
<b>Overall Height</b>	880 mm
<b>Fold Width</b>	300 mm
<b>Fold Length</b>	790 mm
<b>SWL</b>	115 kg

NDIS

MASS



## Staebel Sovereign Wheelchair

The Sovereign wheelchair can be adjusted to suit almost any user. With much consideration in the design from the perspective of OTs, this chair is suited to the vast range of needs required by our clients.

Rather than numerous script forms for various users, this chair meets the needs of almost any patient. All wheelchair parts suit any adjustments made – no additional parts required.

- Adjustable & removable legrests
- Flip back & removable armrests
- Adjustable armrest height & length
- Seat depth, height & angle adjustment
- Backrest height, tension & angle adjustment
- Accessories available (see page 90)



### Specifications:

Product Code	4340C - 43	4340C - 46	4340C - 49	4352C - 52	4352C - 55	4352C - 58
Seat Width	430 mm	460 mm	490 mm	520 mm	550 mm	580 mm
Anterior Seat Height	420 mm	445 mm	450 mm	470 mm	475 mm	500 mm
Posterior Seat Height	410 mm	430 mm	450 mm	470 mm	490 mm	510 mm
Seat Depth	400 mm	425 mm	450 mm	475 mm	475 mm	500 mm
Backrest Height (2cm increments)	400 - 460 mm	400 - 460 mm	400 - 460 mm	400 - 460 mm	400 - 460 mm	400 - 460 mm
Backrest Angle (5° increments)	69° - 94°	69° - 94°	69° - 94°	69° - 94°	69° - 94°	69° - 94°
Lower Leg Length (2cm increments)	400 - 480 mm	400 - 480 mm	400 - 480 mm	400 - 480 mm	400 - 480 mm	400 - 480 mm
Armrest Height (2cm increments)	180 - 260 mm	180 - 260 mm	180 - 260 mm	180 - 260 mm	180 - 260 mm	180 - 260 mm
Armrest Length (4cm increments)	330 - 450 mm	330 - 450 mm	330 - 450 mm	330 - 450 mm	330 - 450 mm	330 - 450 mm
Footplate Width	105 - 155 mm	105 - 155 mm	105 - 155 mm	105 - 155 mm	105 - 155 mm	105 - 155 mm
Gross Weight	17 kg	17 kg	17 kg	21 kg	21 kg	
SWL	125 kg	125 kg	125 kg	170 kg	170 kg	170 kg

NDIS



## Stael Deluxe Alloy Wheelchair

The Deluxe Ultralight Alloy Wheelchairs have been designed to accurately meet the needs of our clients. They are an ultralight weight wheelchair with the benefit of many additional features. The improvement compared to previous designs (ranging from 25 - 30 kg gross weight) is considerable. Much input from clients and nurses has assisted in the development of these superior quality wheelchairs.

- Flip back armrests
- Swingaway, detachable legrests
- Fold down backrest
- Quick release rear wheels
- 24" rear wheels with urethane tyres
- Aluminium frame
- Padded nylon upholstery
- Complies with AS3695:2008
- Accessories available (see page 90)



### Specifications:

<b>Product Code</b>	4629
<b>Seat Width</b>	460 mm
<b>Seat Depth</b>	400 mm
<b>Seat Height</b>	480 mm
<b>Backrest Height</b>	430 mm
<b>Armrest Height</b>	230 mm
<b>Overall Width</b>	650 mm
<b>Overall Depth</b>	1060 mm
<b>Overall Height</b>	960 mm
<b>Fold Width</b>	320 mm
<b>Net/Gross Weight</b>	11 / 14 kg
<b>SWL</b>	125 kg

### Also Available:

<b>Product Code</b>	45129	45629
<b>Seat Width</b>	510 mm	560 mm
<b>Seat Depth</b>	400 mm	400 mm
<b>Seat Height</b>	490 mm	490 mm
<b>Backrest Height</b>	430 mm	430 mm
<b>Armrest Height</b>	210 - 280 mm	210 - 280 mm
<b>Overall Width</b>	700 mm	750 mm
<b>Overall Depth</b>	1000 mm	1000 mm
<b>Overall Height</b>	960 mm	960 mm
<b>Fold Width</b>	300 mm	300 mm
<b>Net/Gross Weight</b>	14.5 / 17 kg	15 / 17.5 kg
<b>SWL</b>	150 kg	150 kg



## Staebel Gusto Heavy Duty Wheelchair

The Staebel Gusto heavy duty wheelchair caters for the bariatric client. Featuring padded vinyl upholstery, heavy duty folding frame, 24 inch rear wheels and inset foot plates for client comfort, this wheelchair is ideal in the aged care and hospital environment. It comes standard with urethane tyres, eliminating the need to maintain tyre pressure. The swing away leg rests are removable and also adjust in length. The armrests will remove to enable side transfers if necessary. Footplates are alloy for longevity and ease of cleaning.

- Removable, swingaway, adjustable length legrests
- Removable armrests
- Padded vinyl upholstery
- 24" rear mag type wheels, urethane tyres



### Specifications:

<b>Product Code</b>	46108
<b>Seat Width</b>	<b>610 mm</b>
<b>Seat Depth</b>	460 mm
<b>Seat Height</b>	520 mm
<b>Backrest Height</b>	430 mm
<b>Armrest Height</b>	220 mm
<b>Overall Width</b>	850 mm
<b>Overall Depth</b>	1100 mm
<b>Overall Height</b>	930 mm
<b>Fold Width</b>	350 mm
<b>Net/Gross Weight</b>	26 / 30 kg
<b>SWL</b>	240 kg

### Also Available:

<b>Product Code</b>	46608	47108
<b>Seat Width</b>	<b>660 mm</b>	<b>710 mm</b>
<b>Seat Depth</b>	520 mm	550 mm
<b>Seat Height</b>	520 mm	520 mm
<b>Backrest Height</b>	480 mm	520 mm
<b>Armrest Height</b>	235 mm	220 mm
<b>Overall Width</b>	880 mm	930 mm
<b>Overall Depth</b>	1150 mm	1170 mm
<b>Overall Height</b>	1000 mm	990 mm
<b>Fold Width</b>	350 mm	350 mm
<b>Net/Gross Weight</b>	36 / 43 kg	36 / 46 kg
<b>SWL</b>	270 kg	290 kg



## Stabel Anti - Tip Bars

The anti-tip bars are fitted to the wheelchair to prevent the wheelchair tipping backwards. It is particularly beneficial for amputees or active users.

- Attaches to wheelchair tilt bar
- Includes small plastic rollers

Product codes may vary depending on wheelchair model. Please call 1800 267 267 for the correct code. Please note all chairs are supplied with Anti-tip Bars. These are replacements.

### Specifications:

**Product Code** 4982



## Stabel Anti Fold / Anti Theft Device

The anti-fold/ anti-theft device is used on a wheelchair to prevent the frame being folded, and to minimise the risk of theft. It is easily fitted to the frame between the rear wheels.

- Adjustable in width
- Larger sizes available

Product codes may vary depending on wheelchair model. Please call 1800 267 267 for the correct code.

### Specifications:

Product Code	4978	4978L	4978XL
<b>Sizing</b>	400 - 490 mm	510 - 560 mm	610 - 710 mm




## Stabel Calf Strap

This padded calf strap provides comfortable support and keeps the feet safely on the footplates. The Velcro straps simply and easily secure to the legrests of the wheelchair.

- Generous size (suitable for 430-490mm chair width)
- Calf section measures 460mm x 89mm

### Specifications:

**Product Code** 4979

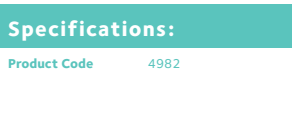

## Stabel Elevating Legrest

The elevating legrest gives the wheelchair user the ability to support their leg in the elevated position. It attaches to the wheelchair in place of the existing legrest.

- Adjustable angle and length to suit the individual user
- Configured left or right - specify when ordering

Product codes may vary depending on wheelchair model. Please call 1800 267 267 for the correct code.

### Specifications:

**Left Leg** 4980 - L  
**Right Leg** 4980 - R



## Stabel Handy Bar

Simple to use and install/remove Transfer Bar wheelchair/car transfer bar.

- Firm secure handle
- Can be used to break glass in emergency
- Inbuilt seat belt cutter
- Fits most vehicles - no modification needed
- Recommended by Doctors and Health Professionals
- Body point belts available. Contact the team for information

### Specifications:

**Product Code** 4998



## Stabel Arm Trough

The arm trough is of benefit for the wheelchair/ shower chair user who has upper limb weakness or paralysis. It attaches to the chair armrest via a sturdy clamp.

- One piece moulded trough
- Supplied with hardware
- Can be mounted left or right side

### Specifications:

**Product Code** 4985



## Stabel Oxygen Bottle Carrier

The oxygen bottle carrier is designed to be attached to the wheelchair to enable the oxygen dependent user to be transported in the wheelchair.

- Suitable for type C oxygen cylinders
- Removes easily for storage and transportation
- Slides onto existing wheelchair tilt bar

Product codes may vary depending on wheelchair model. Please call 1800 267 267 for the correct code.

### Specifications:

**Product Code** 4983



## Stabel Seat Belt

This seat belt is a simple device for safety and positioning of the wheel chair user. It features a push button buckle and is adjustable in length. It has metal eyelets to enable attachment via existing seat upholstery screws on the user's wheelchair.

- Adjustable in length
- Simple and easy fit to most wheelchairs

### Specifications:

**Product Code** 4974



## Stabel IV Pole

The IV pole is attached to the wheelchair and allows the intravenous drip dependent client to be transported in a wheelchair.

- Easy height adjustment with triwheel
- Slides onto existing wheelchair tilt bar

Product codes may vary depending on wheelchair model. Please call 1800 267 267 for the correct code.

### Specifications:

**Product Code** 45129 / 45629



## Stabel Stump Support

The stump support is used on a wheelchair in place of the regular legrest, and is useful for amputees with a below knee amputation.

- Padded stump support
- Adjustable pad height
- Left or right side configured - specify when ordering

Product codes may vary depending on wheelchair model. Please call 1800 267 267 for the correct code.

### Specifications:

**Left Leg** 4981 - L

**Right Leg** 4981 - R



## Stabel Brake Lever Extension

The brake lever extension is used on the wheelchair and is ideal for those who find it difficult to operate standard brakes.

- Lengthens brake handle by approx 10cm
- Increases leverage for easier brake application
- Supplied as each

### Specifications:

**Product Code** 4987



## Stabel Sovereign Power Assist

- Only suitable for Sovereign Wheelchair
- Aluminium frame with PU grey wheels.
- Includes controller, battery and charger.
- Two sizes available: 180W (43428) and 320W (43429).

### Specifications:

**Product Code** 43428



### Staebel Sovereign Arm Trough One Piece

The arm trough is of benefit for the wheelchair/ shower chair user who has upper limb weakness or paralysis. It attaches to the chair armrest via a clamp.

- One piece moulded polyurethane trough
- Supplied with hardware
- Can be mounted left or right side

#### Specifications:

**Product Code** 4985

### Staebel Sovereign Seat Belt

This seat belt is a simple device for safety and positioning. It has a push button buckle and is adjustable in length. It has metal eyelets to enable attachment via existing seat upholstery screws on the wheelchair.

- Adjustable in length
- Simply and easily fitted to wheelchairs

#### Specifications:

**Product Code** 4974

### Staebel Sovereign Drum Brake Kit

- Drum Brake Conversion Kit for the Sovereign Wheelchair
- Includes Brakes, Cables & Wheels

#### Specifications:

**Product Code** 43410



### Staebel Sovereign Stump Support

The stump support is used on a wheelchair in place of the regular legrest, and is useful for amputees with a below knee amputation.

- Adjustable pad height
- Left or right side configured - specify when ordering

#### Specifications:

**Product Code** 43417 - L

**Product Code** 43417 - R

### Staebel Sovereign Elevating Legrest

This seat belt is a simple device for safety and positioning of the wheel chair user. It features a push button buckle and has metal eyelets to enable attachment via existing seat upholstery screws on the wheelchair.

- Adjustable in length

#### Specifications:

**Product Code** 43416 - L

**Product Code** 43416 - R


### Staebel Sovereign Pelvic Strap

- Plastic side release buckle belt, 2-piece design
- Loop tie attachment around wheelchair/ shower chair frame
- Length adjustable 80-155 cm

#### Specifications:

**Product Code** 19636



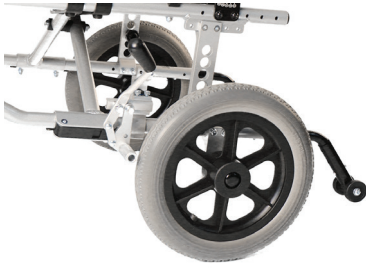

## Staebl Sovereign Transit Wheel Kit

Transit wheel kit for Sovereign Wheelchair  
 Not compatible with 52, 55, 58 cm

### Specifications:

Product Code 43426



## Staebl Sovereign Head Support

Head Support to suit Sovereign Wheelchair  
 • Height Adjustable

### Specifications:

Product Code 49608







## Staebl Sovereign One Arm Drive Kit

The Staebl Sovereign One Arm Drive Kit assists one-handed wheelchair users in controlling the wheelchair. One arm drive systems for manual wheelchairs are designed for people who prefer the advantages offered by lightweight, folding manual wheelchairs.

Double hand rim one arm drives have both two hand rims on the same side of the wheelchair that control the rear wheels on both sides.

One arm drive system with double push rim (Must specify left or right at time of order)

 <p><b>Easy to Use</b></p>	 <p><b>Intelligent</b></p>
 <p><b>Extra Mobility</b></p>	 <p><b>Additional Grip</b></p>

### Outer Hand Rim control

The outer (larger) hand rim controls the wheel on that side of the wheelchair where the kit is located.

### Inner Hand Rim Control

The inner (smaller) hand rim controls the wheel on the opposite side of the wheel chair. The inner hand rim is connected to the wheel by way of a bar between the hand rim and the wheel.

### Specifications:

Product Code 43425

## 05.

## Seating

Ashley Luxor Optima Lift & Recline Chair	95	Regency Result3 Comfort Chair	102
Ashley Lift & Recline Chair	95	Regency Result3 Electric Tilt Chair	103
Ashley Mobile Recliner Luxor Optima	96	Kitchen Perch with Arms	103
Ashley Mobile Recliner - Standard	96	Embark Deluxe Day Chair	105
Cocoon G1 Chair	97	Embark Superior Patient Chair	105
Cocoon G2 Chair	97	Embark Utilite Day Chair	106
Configura Chair	98	Embark Deluxe Extracare Day Chair	106
Regency Care Chair	99	Embark Padded Legrest	107
Regency Care Chair- Electric - Wide	100	Chair / Bed Raiser Biscuits	107
Sertain Care Chair Electric Hi-Lo	101	Chair / Bed Raiser	107
Comflex Bahama Comfort Chair	102		



NDIS



NDIS

## Ashley Luxor Optima Lift & Recline Chair

The Ashley Luxor Optima range of chairs have been specially developed for those with pressure care needs or those needing increased comfort when seated. This is achieved with layered pressure reducing foam and a specialised 2-way stretch upholstery on the seating surfaces. All the favourite features of the standard Ashley are still included – a waterfall style backrest, sculptured armrests, side pockets and padded leg chaise. This model which features a single motor reclines to enable the occupant to rest and will also raise to the standing position.

- Upholstered legrest with calf support
- One motor wallsaver action – chair can be positioned close to within 100mm from wall and still recline
- Backrest reclines to half lying down position
- Chair raises to assist standing transfer
- No sag springs in seat and lumbar for long lasting support
- High quality Dunlop Enduro foams used
- Battery back up capability
- Small castors enable chair to be wheeled
- 2 year electronic warranty, 5 year frame warranty

\* One or two motor available



### Specifications:

Product Code	5522MI	5522	5522MX
Size	Mini	Medium	Maxi
Seat Width	480 mm	480 mm	530 mm
Seat Depth	490 mm	500 mm	520 mm
Seat Height	470 mm	500 mm	530 mm
Backrest Height	600 mm	640 mm	680 mm
Overall Width	780 mm	780 mm	830 mm
Overall Height	1020 mm	1070 mm	1170 mm
SWL	120 kg	120 kg	120 kg

## Ashley Lift & Recline Chair

The Ashley range of chairs have long been the favourite choice with their comfortable seating and range of options. Economically priced and upholstered in standard fabric, the Ashley features a waterfall style backrest, sculptured armrests, side pockets and padded leg chaise. The Ashley chairs can also be ordered in a variety of sizes to suit the individual. This quality powered rise and recline chair is both comfortable and durable, and provides the user assistance with hassle free push button operation.

This model features a single motor allowing the occupant to recline and be raised to the standing position.

Available in either a fabric or vinyl upholstery in a pattern or plain. Contact our team to discuss! Vinyl available in Austex Studio Encore colour range.

\* One or two motor available



### Specifications:

Product Code	5524MI	5524	5524MX
Size	Mini	Medium	Maxi
Seat Width	480 mm	480 mm	530 mm
Seat Depth	490 mm	500 mm	520 mm
Seat Height	470 mm	500 mm	530 mm
Backrest Height	600 mm	640 mm	680 mm
Overall Width	780 mm	780 mm	830 mm
Overall Height	1020 mm	1070 mm	1170 mm
SWL	120 kg	120 kg	120 kg



NDIS



NDIS

## Ashley Mobile Recliner Luxor Optima Ashley Mobile Recliner - Standard

The Ashley Luxor range of chairs have been specially developed for those with pressure care needs or those needing increased comfort when seated. This is achieved with layered pressure reducing foam and stretch PU upholstered on the seating surfaces. All the favourite features of the standard Ashley are still included – a waterfall style backrest, sculptured armrests, side pockets and padded leg chaise. Cascading roll back with isolated head rest and lumbar roll to support both these areas

- Upholstered leg rest with calf support
- Manual operated recliner action with lever
- Castors and push handle for mobility
- No sag springs in seat and lumbar for long lasting support
- Quality Australian made
- 12 month warranty on fabric, 24 month warranty on frame
- Adjustable head rest cushion
- 150 mm castors, 3 x braking, 1 x directional lock
- Allow +/- 20mm tolerances

Economically priced, the Ashley features a waterfall style backrest, sculptured armrests, side pockets and padded leg chaise. Australian made. Cascading roll back with isolated head rest and lumbar roll to support both these areas.

Available in either a fabric or vinyl upholstery in a pattern or plain. Contact our team to discuss! Vinyl available in Austex Studio Encore colour range.

- Allow +/- 20mm tolerances



### Specifications:

Product Code	5519MI	5519	5519MX
Size	Mini	Medium	Maxi
Seat Width	480 mm	480 mm	530 mm
Seat Depth	470 mm	500 mm	520 mm
Seat Height	460 mm	480 mm	500 mm
Backrest Height	600 mm	640 mm	680 mm
Overall Width	780 mm	780 mm	830 mm
Overall Height	1020 mm	1050 mm	1150 mm
SWL	120 kg	120 kg	120 kg

### Specifications:

Product Code	5521
Size	Medium
Seat Width	480 mm
Seat Depth	500 mm
Seat Height	480 mm
Backrest Height	640 mm
Overall Width	780 mm
Overall Height	1050 mm
SWL	120 kg

NDIS



## Cocoon G1 Chair

The lift function of the COCOON G1 armchair is very appreciated by people who have difficulty getting up, but who also find this function a real help when sitting down. Thanks to its coatings that meet fire-retardant standards and are easy to clean, the COCOON G1 can be used both at home and in the community. It consists of removable and washable zipped quilt modules at 30°. This electric lift chair is equipped with removable and / or retractable armrests to facilitate seated – standing transfers (evolving chair in the event of pathological development).

- Optimal hygiene: replaceable, removable and machine washable duvet modules at 30°
- Soft and ultra-comfortable duvets filled with siliconized hollow fiber
- 4 possible positions: Seated, Relaxed (TV), Rest and Lifter for standing transfer
- Different biocompatible coverings available: PVC/PU with M2 fire retardant treatment - mottled fabric - water-repellent velvet
- Different colours variable available depending on the coating
- Guarantees: chair: 2 years - cylinder and mechanics: 5 years
- Adjustable footrests and handle heights
- One or Two Motors Available
- French manufacturing and design.



### Specifications:

Product Code	5572
Seat Width	440 mm
Seat Depth	450 mm
Overall Width	930 mm
Seat Height	520 mm
Overall Height	1100 mm
Chair Weight	59 kg
SWL	160 kg

NDIS



## Cocoon G2 Chair

The COCOON G2 chair is an ultra-comfortable electric lift chair, removable and/or retractable multi-position to facilitate seated – standing transfers (evolving chair in the event of pathological development). The lift function of the COCOON G2 armchair is very appreciated by people who find it difficult to get up, but who also find this function a real help when sitting down. Thanks to its PVC coatings to fire-retardant standards, it can be used both at home and in the community. Using the lifter is very simple. With a simple press of the remote control, the chair returns to a seated position and the chair begins its gentle ascent. You can release the pressure on the remote control at any time, your chair stops.

- Optimal hygiene: replaceable, removable and machine washable duvet modules at 30°
- 4 possible positions: Seated, Relaxed (TV), Rest and Lifter for standing transfer
- Soft and ultra-comfortable duvets filled with siliconized hollow fiber
- Different biocompatible coverings available: PVC/PU with M2 fire retardant treatment - mottled fabric - water-repellent velvet
- 8 possible colors (variable depending on the coating)
- Guarantees: chair: 2 years - cylinder and mechanics: 5 years
- French manufacturing and design.



### Specifications:

Product Code	5575
Seat Width	440 mm
Seat Depth	450 mm
Overall Width	730 mm
Seat Height	520 mm
Overall Height	1100 mm
Chair Weight	59 kg
SWL	160 kg

NDIS



### Configura Chair

The Configura Chair has been designed as a highly adjustable seat that ensures a range of clients are kept comfortable and supported. The chair is adjustable in many ways - change your seat height, seat width and back rest to maximise comfort. Upholstered legrest with calf support.

- Castors for easy manoeuvrability
- Pressure reducing fabric on backrest and armrest
- "Reduced Shear" recline to maximise pelvic stability
- Vapour permeable fabric as standard on seat base, legrest, backrest and armrest
- Adjustable pillow-backrest
- Right hand control with rise setting

#### Optional Features:

- Different levels of pressure care seating cushions
- Different lateral backrest supports for different postural requirements
- Different control options
- Different headrest options



#### Specifications:

<b>Product Code</b>	5550M
<b>Size</b>	Medium
<b>Seat Width</b>	510 mm
<b>Seat Height</b>	460 mm
<b>Armrest Height</b>	180 mm
<b>Overall Height</b>	1140 mm
<b>Overall Width</b>	760 mm
<b>SWL</b>	160 kg

#### Also Available:

	Small	Large
<b>Seat Width</b>	460 mm	560 mm
<b>Overall Width</b>	710 mm	810 mm
<b>SWL</b>	130 kg	160 kg



## Regency Care Chair

The Regency Care Chair is designed to give the patient comfort by combining the original tilt system with an independent multi-positional back rest. Its many features take into account the needs of both carer and user of the chair. The Regency R2900 derives its status from a proven and successful track record across a broad spectrum of the industry.

- Reclining backrest – full recline to supine position
- Pressure management seating – gel cushion in seat, self-air adjusting back, profiled head cushion
- 2-position armrests
- Forward seat tilt for ease to stand, legrest angles under seat

- Elevating contoured legrest complete with fold out, extendible footplate
- Low skin shear 2-way stretch fabric on seating surfaces
- Drop down armrests and swingaway headwings – side transfer friendly
- Covers replaceable
- Patient lifter friendly castor base
- 5 year limited warranty
- 175mm castors, braking on rear, directional lock on front
- Quality Australian made construction



Tangerine



Purple



Green



Pink



Lagoon



Mustard Pattern



Koala Pattern



Arcadia Marina



Bouquet Navy



### Specifications:

<b>Product Code</b>	5050G
<b>Seat Width</b>	490 mm
<b>Seat Depth</b>	520 mm
<b>Seat Height (Tilt/Normal)</b>	570 – 620 mm
<b>Backrest Height</b>	700 mm
<b>Armrest Height</b>	175 mm
<b>Legrest Length</b>	480 – 650 mm
<b>Overall Width</b>	660 mm
<b>Overall Depth</b>	1040 mm
<b>Overall Height</b>	1335 mm
<b>SWL</b>	180 kg

NDIS



### Regency Care Chair - Electric - Wide

The Wide Regency Care Chair follows the success of the R2900 Series model and is designed to give the large client unprecedented comfort by combining the tilt system with an independent multi-positional back rest and pressure area management seating. Its many features take into account the needs of both carer and user of the chair. The Regency derives its status from a proven and successful track record across a broad spectrum of the industry.

- Reclining backrest – full recline to supine position
- Gas strut assisted Tilting seat
- Pressure management seating – gel cushion in seat, self-air adjusting back, profiled adjustable head cushion

- 2-position armrests which drop down for sling insertion
- Forward seat tilt for ease to stand, legrest angles under seat
- Elevating contoured legrest complete with fold out, extendible footplate
- Low skin shear 2-way stretch fabric on seating surfaces
- Drop down armrests and swingaway headwings – side transfer friendly
- Covers are replaceable
- Patient lifter friendly castor base
- 150mm Twin wheel castors



Tangerine



Purple



Green



Pink



Lagoon



Mustard Pattern



Koala Pattern



Arcadia Marina



Bouquet Navy



### Specifications:

Product Code	5050GE - SR
Seat Width	620 mm
Seat Depth	570 mm
Seat Height	550 mm
SWL	300 kg





## Sertain Care Chair Electric Hi-Lo

The Sertain® S4540 Series Pressure Care Seating System is a Team Player! As part of the Care Team the S4540 Series is there when you need it, to reposition the occupant, ease the transfer from bed to chair, increased access for wound dressing and procedures, maneuverable in small spaces - a relatively quiet contributor and achiever but a consistent performer as it adjusts to many of the challenges the Care Team faces - it doesn't take over, it joins in; a real Team needs a real Team Player!

- Electric hi-lo 540-940mm – low enough to stand from with ease, but high enough to achieve a comfortable working height for nursing staff
- Electric reclining backrest – full recline for supine position
- Electric tilting seat and legrest
- Pressure management seating – gel cushion in seat, air in backrest
- Fold out, extendible footplate
- Low skin shear 2-way stretch fabric on seating surfaces
- Removable two-position armrests and head wings
- Battery back up system
- Quality Australian made product
- Grey Frame
- Warranty – 5 year Structural Frame, 2 years Castors, 3 years Electric System



Tangerine



Purple



Green



Pink



Lagoon



Mustard Pattern



Koala Pattern



Arcadia Marina



Bouquet Navy



### Specifications:

<b>Product Code</b>	5141
<b>Seat Width</b>	490 mm
<b>Seat Depth</b>	520 mm
<b>Seat Height (Tilt/Normal)</b>	540 mm / 940 mm
<b>Backrest Height</b>	700 mm
<b>Armrest Height</b>	150 mm
<b>Legrest Length</b>	440 mm
<b>Overall Width</b>	740 mm
<b>Overall Depth</b>	1080 mm
<b>Overall Height</b>	1200 mm - 1520 mm
<b>SWL</b>	160 kg

NDIS



NDIS



### Comflex Bahama Comfort Chair

Our Comflex Care Chair is a popular choice for many healthcare environments. It's available in gas or electric formats, easy to maneuver and simple to clean for infection control. It's been carefully designed to keep your patients or residents social at meal times, fitting in closely to many dining tables.

- 'Tilt in space' seat and leg-rest action
- Euro style reclining backrest with built-in Bahama style head wings
- Elevating standard leg-rest
- Incorporates pressure management Multicush gel seating with air cushion backrest system
- Fixed lounge style armrest
- Includes push handle
- 180kg safe working load
- 5 year limited warranty
- Available in a huge range of Healthcare grade upholstery colours
- Australian made

### Regency Result3 Comfort Chair

The Regency Result3 Comfort Recliner is suited to a huge range of healthcare environments. With a gas-assisted 'tilt in space' seat, reclining backrest and elevating leg rest, it has been designed to keep patients comfortable while being simple to use for healthcare staff. It keeps patients social at mealtimes and in lounge areas.

- Gas-assisted comfort chair
- 'Tilt in space' seat and leg-rest action
- Euro style reclining backrest and elevating standard leg rest
- Incorporating pressure management Muticush gel seating and air cushion backrest system
- Fixed lounge-style armrest
- Includes push handle and swing-away head wings
- 4 x 100mm swivel braking castors
- 5 year limited warranty
- Available in a huge range of Healthcare grade upholstery colours
- Australian made

\* Also available in electric version



Specifications:	
Product Code	5023GH
Seat Width	550 mm
Seat Depth	570 mm
Seat Height	550 mm
SWL	180 kg

Specifications:	
Product Code	5022NV
Seat Width	550 mm
Seat Depth	520 mm
Seat Height	550 mm
SWL	200 kg

NDIS

DVA

NDIS



## Regency Result3 Electric Tilt Chair

Packed full of features and with a 200kg weight capacity, this sleek recliner boasts a modern design appealing to wards, aged care facilities and clinics. The fixed lounge-style armrest has been uniquely designed to fit right up to most dining tables to increase social interaction.

- Electric 'tilt-in-space' seat and leg-rest action
- Electric Euro-style reclining backrest
- Elevating standard leg-rest
- Incorporating pressure management Multicush gel seating and air cushion backrest system
- Fixed lounge-style armrest
- Auto lay-flat function, auto chair function, M3 customer set memory position
- 4 x 100 swivel braking castors
- 5 year warranty
- Australian made

## Kitchen Perch with Arms

An attractive perching stool with an angled padded seat which allow the user to rest and then stand again with minimum effort - particularly helpful with those who tire easily.

Suitable for use in the kitchen, bathroom and any activity where a perched seating position is of assistance. The height of the seat can be adjusted by altering the length of the legs. Zinc treated and powdercoated frame. Seat height 53-69 cm.

- Shallow curved low backrest
- Armrests as part of the frame
- Comfortable vinyl upholstered seat and back
- Adjustable seat height
- Powdercoated steel frame



### Specifications:

Product Code	5025GH
Seat Width	550 mm
Seat Depth	520 mm
Seat Height	550 mm
SWL	200 kg

### Specifications:

Product Code	5941
Seat Width	390 mm
Seat Depth	320 mm
Seat Height	540 - 690 mm
Backrest Height	390 mm
Overall Width	570 mm
SWL	100 kg

# The Embark Range

Specially developed for use in various healthcare environments.



## Choose your Colour!

Emerald

Periwinkle

Smoke

Black

**Custom Colours available upon request.** Contact the team to see custom colours!

NDIS



NDIS



## Embank Deluxe Day Chair

An exceptional rehabilitation chair, the Embark Deluxe is available in sturdy quality vinyls for ease of cleaning. The adjustable seat height, seat depth and angle ensures Embark Deluxe can adapt to the widest range of people, making this chair ideal for hospital, home and nursing home. The padded seat and back rest have been designed to enhance the posture, whilst providing maximum comfort.

- High backrest incorporating lumbar and lateral support
- Padded armrests
- Comfortable vinyl upholstered seat and back
- Adjustable seat height and depth
- Powdercoated steel frame
- Standard with tilt and move wheels for easy manoeuvrability
- Available in a range of different colours. Contact the team for more information



### Specifications:

Product Code	5930
Seat Width	500 mm
Seat Depth	450 - 530 mm
Seat Height	450 - 600 mm
Backrest Height	790 mm
Arm Height	230 mm
Overall Width	700 mm
Overall Depth	730 - 800 mm
Overall Height	1200 - 1350 mm
SWL	200 kg

### Also Available:

Embank Day Chair Legrest	
Product Code	5958

## Embank Superior Day Chair

The Embark Superior Day Chair was designed and developed by a team of clinicians. It has a strong frame and a high weight capacity of 200 kg. The chair has low resilient foam and 2-way stretch upholstery on seating surfaces. For increased longevity in demanding environments it has 6cm tilt and move wheels and other small improvements.

There have been instances in the past, whereby patients, staff and visitors have sustained cuts and injuries from square tube frames. As a result, the Embark was not only designed for comfort, but to avoid these injuries.

- Low back option also available.
- Comfortable upholstered seat and back
- Padded armrests
- Adjustable seat height and depth
- Powder-coated steel frame
- High backrest incorporating lumbar and lateral support



### Specifications:

Product Code	5933
Seat Width	500 mm
Seat Depth	450 - 530 mm
Seat Height	450 - 600 mm
Backrest Height	790 mm
Arm Height	230 mm
Overall Width	700 mm
Overall Depth	730 - 800 mm
Overall Height	1200 - 1350 mm
SWL	200 kg

NDIS

NDIS



### Embark Utilite Day Chair

The Utilite is a multi purpose adjustable seat height and angle chair. Ideal for use in an area where adjustable seating is required, for example, lounge, dining room, bedroom, activity room, outdoor area etc. This low back chair provides comfort, stability and practicality at an affordable price.

- Shallow curved low backrest
- Padded armrests
- Comfortable, durable upholstered seat, back and armrests
- Powdercoated steel frame
- Available in a range of different colours. Contact the team for more information



#### Specifications:

Product Code	5910
Seat Width	460 mm
Seat Depth	460 mm
Seat Height	485 - 640 mm
Backrest Height	455 mm
Arm Height	230 mm
Overall Width	610 mm
Overall Depth	550 mm
Overall Height	900 - 1075 mm
SWL	120 kg

#### Also Available:

Embark Day Chair Legrest	
Product Code	5958

### Embark Deluxe Extracare Day Chair

The Extra Care Embark Deluxe has an impressive 250kg weight capacity and extra wide 60cm width, making it ideal for the heavier client. It still provides all the benefits of the Embark Deluxe, e.g. adjustable seat depth and height, and comes standard with tilt and move wheels.

- High backrest incorporating lumbar and lateral support
- Padded armrests
- Comfortable vinyl upholstered seat and back
- Adjustable seat height and depth
- Powder-coated steel frame with reinforcement for increased load



#### Specifications:


Product Code	5931BN3
Seat Width	600 mm
Seat Depth	450 - 530 mm
Seat Height	470 - 620 mm
Backrest Height	765 mm
Arm Height	220 mm
Overall Width	800 mm
Overall Depth	740 - 820 mm
Overall Height	1180 - 1330 mm
SWL	250 kg

NDIS

NDIS

DVA

NDIS

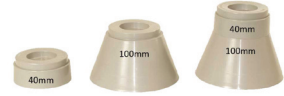

5990

5991

5990 / 5991

## Embark Padded Legrest

This vinyl covered legrest provides comfort for resting legs, with its generously padded top set on a 10 degree angle. Height and angle adjustable using simple push button system. Fitted with non-slip rubber feet.

- Sturdy powdercoated steel frame with vinyl upholstery
- Complements the Embark chair range

## Chair / Bed Raisers Biscuit

This plastic disc is used with the plastic bed raiser, when the bed leg is too large to fit into the raiser. The biscuit is attached with a screw to the base of the bed leg. The biscuit then can then be slotted into the centre top indent of a plastic 40mm or 100mm height raiser.

## Chair / Bed Raiser

Lightweight but extremely strong with a broad base and a slip resistant pad in a recessed top for a chair or bed leg. Can be stacked for greater height increase and storage. Supports up to 250kg.

- Raises height of bed or chair by 40mm
- Can be stacked/combined with the 100mm (5991) raiser to achieve greater height
- Round cavity to support chair leg or bed castor, up to 65mm in diameter



### Specifications:

Product Code	5958
Height	435 - 585 mm
Width	550 mm
Length	430 mm
SWL	200 kg

### Specifications:

Product Code	5992
--------------	------

### Specifications:

Product Code	5990 / 5991
Hole Diameter	65 mm
Height ( 5990 )	40 mm
Height ( 5991 )	100 mm
SWL	250 kg

# 06.

## Mobility Solutions

Walking Stick - Aluminium T-Handle	109	Stael Liger Seat Walker	112
Walking Stick - T-Bar Handle	109	Server HD Seat Walker	112
Walking Stick - Swan Neck	109	Navigator Forearm Walker	113
Easy Reacher	109	Viking Forearm Walker	113
Crutches - Aluminium Height Adjustable	110	Stimulator Rehab Forearm Support Walker	113
Stael Anatomic Grip Elbow Crutches	110	byACRE Carbon Ultralight Walker	114
Lightweight Walking Frame	110	byACRE Carbon Overland Fibre Walker	114
Stael Cruiser Folding Walker	111	byACRE Accessories	115
Grande Seat Walker	111	Rollz Motion Accessories	116
Stael Petite Seat Walker	111	Rollz Motion 2 in 1 Wheelchair Walker	117





NDIS



DVA

NDIS

## Walking Stick - Aluminium T-Handle

One of the best selling walking sticks favoured for its sheer versatility. The shaft is made of lightweight yet durable aluminium which folds into four sections and can be stored neatly in the plastic wallet supplied with it.

## Walking Stick - T-Bar Handle

This adjustable height aluminium walking stick assists the user's mobility. The T-shaped handle is strong and lightweight. It also includes an anti rattle silencer and adjusts in 25mm increments.

- Attractive timber look handle



### Specifications:

Product Code	6241
Height (Med)	790 - 840 mm
SWL	100 kg

### Also Available:

Product Code	6240
Height (Lrg)	840 - 940 mm
SWL	100 kg



### Specifications:

Product Code	6230
Height	760 - 990 mm
SWL	100 kg



NDIS

## Walking Stick - Swan Neck

This adjustable height aluminium walking stick assists the user's mobility. The "swan neck" shaped handle helps to transmit the weight of the user vertically through the body of the stick. It provides exceptional support and comfort, eliminating slipping and encouraging correct posture. It also includes an anti rattle silencer and adjusts in 25mm increments.



### Specifications:

Product Code	62210
Height	730 - 990 mm
SWL	100 kg



NDIS

## Easy Reacher

The Easy Reacher assists in retrieving or reaching a wide variety of items in the home. It is ideal for those who have difficulty bending or reaching, and features a twin gripping, rubber lined jaw for an excellent controlled grip. It is available in different lengths to suit the individual user.

### Specifications:

Product Code	62661
Length	730 mm

### Also Available:

Product Code	62674
Length	10000 mm

DVA  
NDIS



DVA  
NDIS



DVA  
NDIS



### Crutches - Aluminium Height Adjustable

The Aluminium Crutches offer push button height adjustment and are manufactured from a lightweight anodised aluminium. The solid one piece crutch frame gives strength and incorporates a contoured non rotating hand grip. A new anti rattle system makes these crutches silent.

- 3 sizes available (Codes 6011, 6012 and 6013)
- Supplied as a pair

### Staebel Anatomic Grip Elbow Crutches

The Anatomic Elbow Crutches feature contoured moulded handles that evenly distribute pressure across the palm of the hand, contributing to the users ease of operation and comfort. Added to this is the ability to adjust the crutch length from cuff to hand grip and also from hand grip to ground, enabling the user to adjust the crutch to their individual requirements. Manufactured from a larger tubing giving increased stability and weight capacity.

- Supplied as a pair

### Lightweight Walking Frame

This folding lightweight walking frame provides assistance for those who have difficulty walking.

It folds in by pushing a simple release button. When folded this frame will fit into the smallest of cars and is ideal where storage is at a minimum. Other features are soft grip handles, simple push button height adjustment and anti rattle silencers.

- Height adjustable to suit the individuals requirements
- Lightweight aluminium construction
- Folds flat for easy storage and transportation
- Soft handles for comfort and grip



#### Specifications:

Product Code	6011J
Crutch Height	930 - 1130 mm
Patient Height	1372 - 1575 mm
SWL	100 kg

#### Also Available:

Product Code	6012J	6013J
Crutch Height	1140 - 1345 mm	1335 - 1540 mm
Patient Height	1575 - 1778 mm	1778 - 1981 mm
SWL	100 kg	100 kg



#### Specifications:

Product Code	6151J
Cuff to Hand Grip	930 - 1130 mm
Hand Grip Height	1372 - 1575 mm
SWL	180 kg

#### Also Available:

Product Code	6152J
Cuff to Hand Grip	1335 - 1540 mm
Hand Grip Height	1777 - 1981 mm
SWL	180 kg



#### Specifications:

Product Code	6380
Height	690 - 870 mm
SWL	120 kg

#### Also Available:

Product Code	6381
Height	830 - 1010 mm
SWL	120 kg

DVA

NDIS



NDIS

MASS



NDIS

MASS



## Staebel Cruiser Folding Walker

The Cruiser is a popular seat walker manufactured from durable lightweight aluminium. 6" castors make this an ideal outdoor walker. But is also compact enough for use indoors. Features include a padded seat, bag under seat, lockable brakes and padded, folding backrest for comfort while seated. The adjustable height handles are fitted with orthopaedic hand grips, therefore providing the user with excellent support. The walker also folds flat for transport.

- Height adjustable ergonomic shaped handgrips
- Lightweight aluminium construction
- Lockable handbrakes
- 152mm castors

## Grande Seat Walker

- Ideal for those requiring mobility assistance
- Durable lightweight aluminium frame
- Folding for convenient transport and storage
- Height adjustable handles
- Ergonomically designed rubber hand grips
- Fitted with rear braking wheels and front swivel castors
- Locking brake mechanism with reflectors
- Padded seat
- Supplied with reinforced vinyl bag
- Removable backrest

## Staebel Petite Seat Walker

- Ideal for smaller people requiring mobility assistance
- Durable lightweight aluminium frame
- Folding for convenient transport and storage
- Height adjustable handles
- Ergonomically designed rubber hand grips
- Fitted with rear braking and front swivel castors
- Locking brake mechanism with reflectors
- Padded seat
- Supplied with reinforced vinyl bag
- Removable backrest



### Specifications:

<b>Product Code</b>	6730
<b>Seat Width</b>	440 mm
<b>Seat Height</b>	560 mm
<b>Handle Height</b>	820 - 930 mm
<b>Walker Weight</b>	6 kg
<b>SWL</b>	125 kg



### Specifications:

<b>Product Code</b>	6861
<b>Seat Width</b>	440 mm
<b>Seat Depth</b>	300 mm
<b>Seat Height</b>	550 mm
<b>Overall Width</b>	570 mm
<b>Overall Depth</b>	670 mm
<b>Overall Folded WxDxH</b>	570 x 350 x 815 mm
<b>Handle Height</b>	815 - 925 mm
<b>Castor Size</b>	194 x 40 mm
<b>Walker Weight</b>	7.4 kg
<b>SWL</b>	200 kg



### Specifications:

<b>Product Code</b>	6771
<b>Seat Width</b>	380 mm
<b>Seat Depth</b>	300 mm
<b>Seat Height</b>	470 mm
<b>Overall Width</b>	560 mm
<b>Overall Depth</b>	670 mm
<b>Overall Folded WxDxH</b>	570 x 350 x 815 mm
<b>Handle Height</b>	740 - 850 mm
<b>Walker Weight</b>	5.5 kg
<b>SWL</b>	175 kg



DVA  
NDIS



NDIS  
MASS

### Staebel Ligerio Seat Walker

The Ligerio is a popular seat walker manufactured from durable lightweight aluminium. The larger castors make this an ideal outdoor walker, however is also compact enough for use indoors. Features include a padded seat and padded folding backrest for comfort while seated. The adjustable height handles are fitted with Orthopaedic hand grips therefore providing the user with excellent support. The walker also folds for storage and transport. Two seat heights available.

- Seat height of 55cm
- Height adjustable ergonomic shaped hand grips
- Lightweight aluminium construction
- Lockable handbrakes
- 200mm (8") castors



#### Specifications:

Product Code	6865
Seat Width	440 mm
Seat Height	550 mm
Handle Height	790 - 910 mm
Walker Weight	7.5 kg
SWL	120 kg

#### Also Available:

Ligerio Seat Walker - 510 mm	
Product Code	6869

### Server HD Seat Walker

Server HD folds easily and stays that way with a simple locking system that doubles up as an ergonomically shaped carrying handle. Frame and seat are tested and approved for maximum user weight of 200 kg. The brake mechanism is light yet positive. TPEcovered, unbreakable wheels.

- The push handles are height adjustable with 10 positions, from 74 cm to 102 cm.
- User comfort is ensured with an optimal frame design and ergonomic push handles.



#### Specifications:

Product Code	6784
Width	690 mm
Depth	680 mm
Height	740 - 1020 mm
Seat Width	550 mm
Depth (When Folded)	680 mm
Height (When Folded)	860 mm
Width (When Folded)	225 mm
SWL	200 kg

NDIS

MASS



NDIS



NDIS



## Navigator Forearm Walker

The navigator is foldable, staying firmly closed making it easy to transport. The frame and seat are tested and approved on users with a maximum weight of 150 kg. Light and reliable brakes act on the TPE covered unbreakable wheels. The optimal and rigid design of the frame and the fore arm support ensures the best possible user comfort (lightness, turning circle and steadiness).

- The Navigator is available in size large.
- The Navigator is easily height adjustable in the range from 900 – 1150 mm



### Specifications:

Product Code	6782
Overall Width	550 mm
Overall Depth	610 mm
Height	740 - 1020 mm
Seat Height	620 mm
Handle Height	900 - 1150 mm
Walker Weight	10.8 kg
SWL	150 kg

## Viking Forearm Walker

The Viking Forearm Walker presents some unique features that particularly suit a rehab or post surgery environment. Its many adjustments also mean it can be adapted for a wide variety of patients.

- Adjustable angle forearm supports
- Adjustable angle hand grips
- Height adjustable
- Adjustable width base
- Large padded forearm pads for user comfort
- Strong aluminium frame
- Lockable rear castors



### Specifications:

Product Code	6748
Overall Width	700 - 850 mm
Handle Height	880 - 1310 mm
Walker Weight	20 kg
SWL	200 kg

## Stimulator Rehab Forearm Support Walker

The best rehab program and the best rollator will guarantee you the best results – is this what you're looking for? Our Stimulator Rehab Forearm Walker has been specifically designed for the needs of the Aged care and Acute Care/ Hospital environment.

- Very lightweight, easy for the frail, weak patient to use
- High quality castors with super smooth rolling and minimal friction even on carpet, taking away any shuddering or wobbling
- Soft PU multi adjustable arm pads for any patient with any specific position requirement
- Multi adjustable hand grips, for the arthritic, ABI or stroke patients.
- Four swivelling castors for ultimate manoeuvrability
- Individual brakes on each castor, as well as hand brakes
- Modern streamlined design gives an up-to-date image to your rehab and physio area



### Specifications:

Product Code	6755
Overall Width	700 mm
Height of Arm Support	905 - 1310 mm
Width Between Arms	310 - 470 mm
Overall Length	810 mm
SWL	200 kg

NDIS

NDIS



### byACRE Carbon Ultralight Walker

At just 4.8kg this carbon fibre walker / walking frame is the lightest in the world and uses space age technology and the best of Danish design.

- Soft inverted handles to support posture and comfort and reduce strain on fingers, hands and joints when steering
- Unbeatable in terms of the ease to fold, lift and grip when lifting it into the car or up stairs
- Award winning design providing a slick alternative to the traditional looking walking frame
- Brake cables hidden inside the frame for cleaner, less cluttered lines
- Shock absorbing frame and tyres made from a PU/rubber mix to reduce vibration
- Different colours available! Contact the team to see more

### byACRE Carbon Overland Fibre Walker

The byACRE Carbon Overland walker / walking frame is a sturdy, light weight and versatile multi-terrain walker

- Lightest off-road walker at 6.7 kg for easy lifting in y our vehicle
- Easy to fold, easy to grip and easy to lift
- Fully assembled, to make it hassle free to start exploring
- Extremely stable carbon fibre frame – 5 x stronger than steel
- Ultra-durable pneumatic rubber tires for better suspension
- Rigid front forks and sealed bearings
- Aircraft grade aluminium in gunmetal colour, for strength and lightness



Specifications:	
Product Code	6884SR
Width	625 mm
Width when Folded	250 mm
Handle Height	790 - 910 mm
Width between Handles	390 mm
Seat Width	420 mm
Seat Height	610 mm
Length	660 mm
Weight	4.8 kg

Specifications:	
Product Code	6886DG
Handle Height	820 - 950 mm
Length	730 mm
Seat Width	470 mm
Seat Height	610 mm
Width	660 mm
Width when Folded	250 mm
Wheel Diameter	250 mm
Weight	6.7 kg
SWL	120 kg



## byAcre Organiser Bag

The organiser is made to hold your personal belongings, such as keys, wallet, napkins, etc. It can be attached to the back of the walker. With the attachable shoulder strap, the organiser can be carried over your shoulder like a small purse, when you are not using it on the rollator.

- To suit byACRE Walkers.



## byAcre Weekender Bag

Spacious bag and with a small pocket inside to suit byACRE Walkers.

- Wear and weather resistant material with an elegant lining inside
- Hooks on the side so you can attach a strap and use the bag over the shoulder
- Easy to attach and detach and with an extra elastic attachment.

### Specifications:

Product Code

6887

### Specifications:

Product Code

6887 - WB



## byAcre Grocery Bag

Very spacious for a full load of shopping, to suit the byACRE Walkers

- Easy open/close system.
- Room for a tote bag so you can load your groceries easily in to the bag.
- Easy to clean and maintain made from wear and weather resistant polyester that can easily be wiped clean.

### Specifications:

Product Code

6888



## byAcre Cane Holder

- An easy-to-use cane holder to hold your cane in place while walking with your byACRE Walker.
- Easy to attach and detach
- To suit the byACRE Walkers

### Specifications:

Product Code

6889



### Cup Holder for Rollz Motion

The Rollz Motion cup holder attaches conveniently to the side of the frame and allows you to keep hydrated on your journey. Suitable for the Rollz Motion (regular and small) and the Rollz Performance.



### Back Support for Rollz Motion

This back support offers you a comfortable rest while leaning against the soft strap. Suitable for the Rollz Motion (regular and small) and the Rollz Performance.

\*Walker not supplied, just for illustration purposes



### 3-in-1 Accessory for Rollz Motion

This nifty clip provides an a second place on your Rollz Motion frame to slot your wheelchair package. These clips also double as shopping hooks and a walking cane holder. Suitable for the Rollz Motion (regular and small) and the Rollz Performance.

#### Specifications:

Product Code 6880 - 2

#### Specifications:

Product Code 6880 - 3

#### Specifications:

Product Code 6880 - 4



### Basket for Rollz Motion

Need somewhere to place your shopping? This basket sits under the seat of the Rollz Motion and provides another place to store your things possessions.



### Travel Cover for Rollz Motion

The perfect accessory for your Rollz Motion if you like to travel. The lightweight cover slips over the frame and wheelchair package and keeps both safe from any bangs and scratches that sometimes happen on transit. Suitable for the Rollz Motion (regular and small) and the Rollz Performance.



### Umbrella for Rollz Motion

This nifty umbrella can be attach to your Rollz Motion Walker. Suitable for the Rollz Motion or the Rollz Performance.

\*Rollz Walker is not included, shown in image for illustration purposes only.

#### Specifications:

Product Code 6880 - 5

#### Specifications:

Product Code 6880 - 6

#### Specifications:

Product Code 6880 - 7





## Rollz Motion 2 in 1 Wheelchair Walker

The award winning Rollz Motion are one of a kind. Perfect for those who want to walk but are not sure how far they can go. With their large wheels for superior maneuverability and state of the art brakes, these are the ultimate in assisted walking. The Rollz Motion is an award-winning and beautifully designed rollator that converts simply into a wheelchair in seconds.

- Gives you confidence to go further
- Innovative and award-winning design
- Superior brakes suited to wet and dry conditions
- Folds into compact and transportable size
- Strong frame and wheels
- Allows you to stand and walk tall, instead of leaning forward
- Safe with superior manoeuvrability
- Adjustable footrests and handle heights
- Three colours available



### Specifications:





<b>Product Code</b>	6880IB
<b>Seat Height</b>	550 mm
<b>Handle Height Range</b>	870 - 980 mm
<b>Seat Size Rollator</b>	200 x 450 mm
<b>Seat Size Wheelchair</b>	420 x 450 mm
<b>Overall Length</b>	630 mm
<b>Overall Width</b>	310 mm
<b>Width Folded</b>	670 mm
<b>Weight</b>	11.6 kg
<b>SWL</b>	125 kg

# 07.

## Bedroom Solutions

Sentisa SC Bed	120	VT2 Bedside Cabinet	130
Staebel Floorline Bed	121	Staebel Mobile Overbed Table	130
Luxa Homecare Bed	122	Solido 3 Mobile Overbed Tilting Table	130
Nursing Care Bed	123	Staebel Deluxe Bedside Commode	131
Nursing Care Bariatric Bed	123	Staebel Standard Bedside Commode	131
Variana Bed	124	Bed Rail	132
Companion Bed	125	Swivel Hand Control Holder	132
Latera Thema Bed	126	Blanket Support	132
Multicare Bed	127	Bed Pole	132
Nino Adolescent Bed	128	Overbed Frame & Triangle	133
Tom 2 Cot	129	Staebel Bed Rail Protector	133

## Build your Customisable Bed!

 <p>Who is using the bed?</p>	Infant	Child	
	Adult	Elderly	
 <p>Choose your bed width</p>	800 mm	850 mm	1000 mm
	1050 mm	1200 mm	1350 mm
	1500 mm	1800 mm	Companion Bed
 <p>Weight Limit (SWL)</p>	80 kg	150 kg	200 kg
	270 kg	320 kg	450 kg +
 <p>Lowest Bed Height Requirement</p>	100 mm	230 mm	280 mm
	380 mm	450 mm	

## Let's stylise!

 <p>Choose your bed style</p>	Hospital Style		Homecare Style	
				
 <p>Bed Frame Design</p>	Colours	Stains	Laminate	
	Timber Head & Foot Boards	Upholstered Head & Foot Boards	Vinyl Head & Foot Boards	
	Bed Functions	Castors	Controls / Hanset	
 <p>Accessories</p>	Battery	Lighting	Rails / Support Rails	
	Self Help Pole	Clip on Extension	Handset Holders	



**Could not find the information you were after?**

Contact our team for a consultation for more bed options!



## Sentida SC Bed

New concepts for living in old age are in demand. These include the desire to be provided with more homeliness, individuality, safety, comfort and flexibility in the equipment of corresponding living arrangements.


- With the new nursing bed platform sentida sc and customized concepts, we can respond to the different requirements even more flexibly, efficiently and demand-oriented.
- We want to support the maintenance and promotion of individual mobility as an overarching and essential goal in care in a solution-oriented manner.
- SafeFree® side rails – safe mobilization regardless of body size.
- Underbed clearance for mobile patient lifters at any position
- Quality German made construction
- Fully customisable with various models, decor and designs. Contact the team for more information.



NDIS

### Specifications:

<b>Product Code</b>	77754C
<b>External Dimensions</b>	2080 x 1060 mm
<b>Sleeping Surface</b>	870/900 x 2000 mm
<b>Backrest Adjustment</b>	70° ± 5° (retraction by 120 mm)
<b>Upper Backrest Adjustment</b>	30° ± 5° (retraction by 50 mm)
<b>Head/ Foot Swivel</b>	17° / 17°
<b>Height</b>	230 - 800 mm
<b>Bed Weight</b>	115 kg (depending on model)
<b>SWL</b>	270 kg

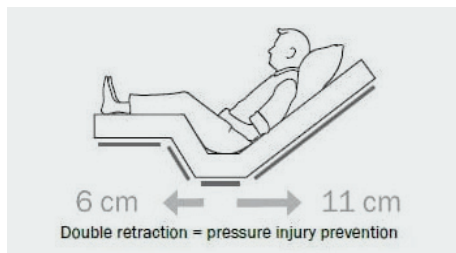




Scan the QR code  
to watch the video  
demonstration!



## Staebel Floorline Bed

The Staebel Floorline Bed is one of the lowest beds available on the market and designed for intensive use. The bed lowers to floor level to prevent injuries and accidents caused by falls from beds. All actuators and mechanisms are concealed beneath the platform, allowing the bed to be constructed with no end pedestals! This allows visibility for users, and nursing staff access for procedures.

- The lifting range of the bed is not compromised by its low level – a platform height of 720 mm is achieved making this bed a superb companion in complex care.
- Electric 4-section profiling
- Auto-regression backrest for minimising shear forces
- Trendelenberg – 15° Reverse – 15°
- 75mm twin wheel castors. Optional central locking
- Various colours and styles available
- 10 year frame / 5 year electric warranty / 1 year Handset

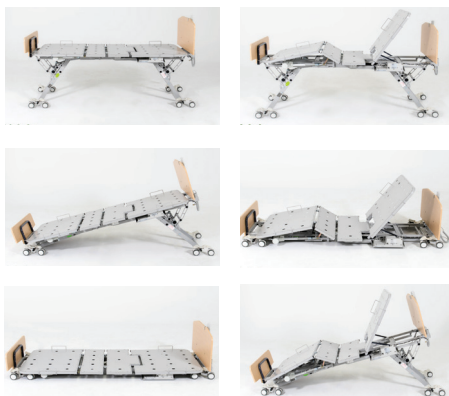

NDIS

### Specifications:

Product Code	7753
Size	Single
Overall Width	905 mm
Overall Length	2080 mm
Height	120 - 720 mm
Sleeping Surface Dimensions	890 x 2000 mm
SWL	200 kg

### Also Available:

Product Code	7755	7758
Size	Wide	King Single
Overall Width	1020 mm	1060 mm
Overall Length	2080 mm	2080 mm
Height	120 - 720 mm	120 - 720 mm
Sleeping Surface Dimensions	1000 x 2000 mm	1060 x 2000 mm
SWL	200 kg	300 kg



## Luxa Homecare Bed

The Luxa is a premium homecare bed that is designed to give the user ultimate comfort and care, while maintaining dignity and style in the bedroom.

This bed is packed with premium features such as an auto-regression backrest, manual vascular support, a high height of 725mm, a low height of 150mm height and optional bluetooth handsets. It also has optional accessories including a sturdy horizontal rail, allowance for self-help poles, underbed lights and more.

- Extremely Low Minimum Position
- Vertical Lift
- Auto-Regression Backrest
- Ash or Ocean Padded Headboard with Full Padded Surround
- Proudly Australian made
- Optional features and styles can be added at time of manufacture (contact the team for more)



NDIS

### Specifications:

Size	Single
External Dimensions	2160 x 1000 mm
Sleeping Surface	2040 x 900 mm
High Position	725 mm
Low Position	150 mm
Backrest / Knee Angle	62° / 35°
TR / ATR Angle	15° / 15°
SWL	250 kg

### Also Available:

Size	King Single	Double	Queen
External Dimensions	2160 x 1170 mm	2160 x 1470 mm	2160 x 1600 mm
Sleeping Surface	2040 x 1070 mm	2040 x 1370 mm	2040 x 1500 mm
High Position	725 mm	725 mm	725 mm
Low Position	150 mm	170 mm	170 mm
Backrest / Knee Angle	62° / 35°	62° / 35°	62° / 35°
TR / ATR Angle	15° / 15°	15° / 15°	15° / 15°
SWL	250 kg	250 kg	250 kg





## Nursing Care Bed

The Nursing Care Bed has been designed for the high care client who is requiring a wider sleeping area than standard single size. This bed offers electric high-low, electric backraise/lower and the use of the electric kneebend assists in preventing sliding down the bed.

- Removable melamine head and foot boards ( shown with optional timber grain boards )
- Fully sealed foot end frame to comply with strict hygiene standards implemented in hospitals
- Hand control has been designed for those who have limited hand function
- Features back lighting for increased visibility at night
- Underbed light to reduce risk of falls in dim conditions
- Powdercoated frame with serial no. record for lifetime tracking
- 5 year electrical warranty, 10 year frame warranty
- Quality Australian made construction with Aust Standard compliance
- 125mm braking castors for ease of manoeuvrability

## Bariatric Nursing Care Bed

The Bariatric Nursing Care bed has been designed for the high care client who is requiring a wider sleeping area than standard single size. This bed offers electric high-low, electric back raise/lower and the use of the electric kneebend assists in preventing sliding down the bed. Other benefits include:

- Removable melamine head and foot boards (shown with optional timber grain boards )
- Fully sealed foot end frame to comply with strict hygiene standards implemented in hospitals
- Powdercoated frame for ease of cleaning
- Hand control has been designed for those who have limited hand function
- 125mm braking castors for ease of manoeuvrability
- 300kg capacity
- Underbed clearance 170mm for ease of patient hoist access
- Quality product, manufactured in Australia
- Angles Back - 65 degrees, Knee 44 degrees




### Specifications:

Product Code	7215
Height	270 - 700 mm
Overall Width	995 mm
Overall Length	2060 mm
Sleeping Surface	1990 x 879 mm
SWL	200 kg

### Specifications:

Product Code	7216
Height	270 - 760 mm
Overall Width	1200 mm
Overall Length	2040 mm
Sleeping Surface	2000 x 1200 mm
SWL	300 kg

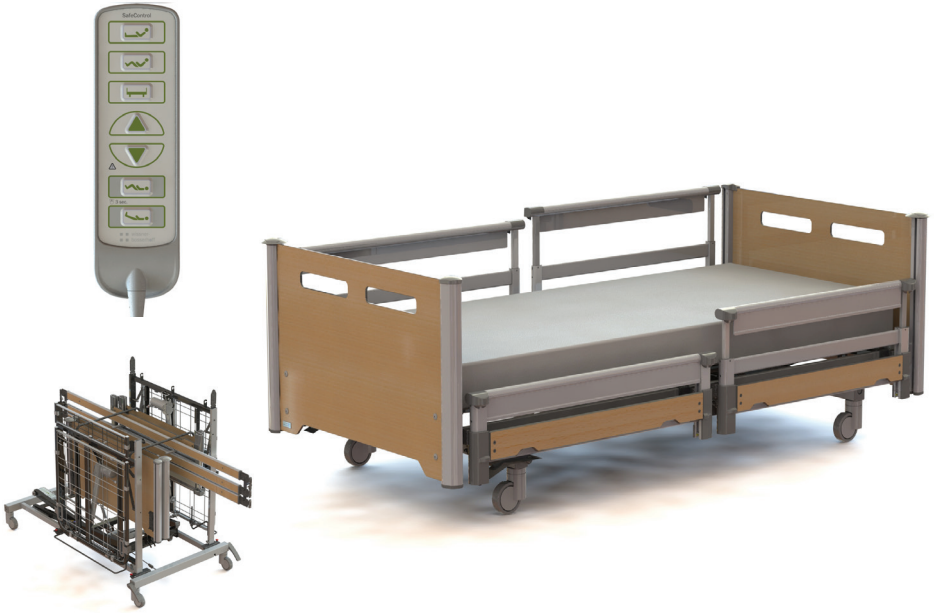

## Variana Bed

The Variana Bed is a demountable and washable nursing bed. The Variana Bed is serially equipped with a Standard washability kit and thus enables flexible and quick cleaning of the bed on a regular basis. However, if everyday care requires more frequent washing cycles, the Variana can also be operated with a Premium kit.

- The bed's handset is an intuitive handset with illuminated buttons making use at night a lot easier. The handset can also be washed to ensure that it is infection control safe.
- The bed also offers what the other beds offer from the WiBo range, with 4-part mobilisation aids or full siderails, dual retraction for pressure injury prevention and a range of colour and décor choices.



### Specifications:

<b>Product Code</b>	77874C
<b>External Dimensions</b>	2180 x 1110 mm
<b>Sleeping Surface</b>	2000 x 900 x 140 mm
<b>Maximum Sleeping Surface Height</b>	340 mm
<b>Minimum Sleeping Surface Height</b>	820 mm
<b>SWL</b>	235 kg





## Companion Bed

The 8000 Series Companion bed allows two patients can be cared for in comfort and dignity. It features electric hi-low, electric knee brake, and a 300kg working load. The actuators are heavy duty, ensuring efficient, reliable function, and are warrantied for five years. The steel bed frame is warrantied for ten years. The 8000 Series bed is available in both Double and Queen size.

- Electric Hi Low
- Dual Electric Backrests running independently.
- Dual Electric Kneebreaks running independently.
- Standard Head & Foot Boards
- Battery Backup
- 125mm Braking & Directional Castors
- 2 x Backlit handsets
- 2 x Underbed lights
- Safe Working Load of 300kg



### Specifications:

<b>Product Code</b>	7756
<b>Overall Length</b>	2060 mm
<b>Overall Width</b>	1390 mm
<b>Height Hi-Lo</b>	450 - 750 mm
<b>Sleeping Surface</b>	1990 x 1370 mm
<b>Underbed Clearance</b>	150 mm
<b>Angles Back</b>	Back 70° / Knee 44°
<b>SWL</b>	300 kg



## Latera Thema Bed

Another product for the Wissner-Bosserhoff range, the Latera Thema was designed after extensive research exposed the need for a system with lateral tilting and reclining that could still be cost-effective. With the Latera Thema this was achieved!

- Fully functional seating and reclining positions.
- Painful friction forces during positioning are reduced
- Significantly reduces stress to care staff when performing clothing changes, wound dressing etc.
- A handset for patient (does not have turning function) and an additional handset for care staff
- Unique 2-part side guards which have multiple positions
- Double retracting sleeping surface to reduce friction & shearing forces
- Vertical raise therefore no holes in walls (not scissor action)
- Underbed clearance for patient lifters at any position
- Hand control has nurse lockout for dementia residents
- Robust hand control with 2 button operation reducing risk of accidental adjustment



NDIS

### Specifications:

<b>Product Code</b>	20 - 0554
<b>External Dimensions</b>	2130 x 1000 mm
<b>Mattress Dimensions</b>	2000 x 860 mm
<b>Lateral Tilt</b>	+ 15° / - 15°
<b>Maximum Height</b>	810 mm
<b>Minimum Height</b>	430 mm
<b>Bed Weight</b>	160 kg
<b>Clearance Height</b>	150 mm
<b>SWL</b>	200 kg

1.



2.



3.



4.





## Multicare Bed

The Multicare is an accomplished member of the critical care team. It simplifies interventions in critical care situations improving patient care through treatment and rehabilitation whilst providing comfortable working conditions for nursing staff. A combination of advanced innovations with functional features ensures that the Multicare exceeds your expectations.

- Automatic Lateral tilting (both sides to 30 degrees)
- Multi Protect 4-part siderails
- Multiboard features an LCD touchscreen (in both siderails)
- Mattress platform tilt and extension adjustment
- Tilt indicators (backrest, lateral tilt, TR/ATR tilt)
- Programmed positions: CPR, orthopaedic chair, 30°HOB, mobilisation, Trendelenburg
- Activation GO button, safety STOP button
- Service data interface, user communication setup, quick reference guide on LCD touchscreen
- Battery back-up indicator (2 batteries)
- The Multicare also includes patient control panel in the side guards
- The Multicare can be customised to suit your specific needs: there are four different kits
- Mobi Lift handles - assists the patient to exit the bed safely.
- I-brake – an automatic brake that is activated after 10 seconds if the bed remains unbraked and is plugged in. It protects against uncontrolled movement, impacts and falls caused by leaning on an unbraked bed.
- Time Saver, Auto tilting bed



NDIS

### Specifications:

<b>Product Code</b>	74200
<b>External Dimensions</b> (Shortest Mattress Platform Position)	1050 x 2150 mm
<b>External Dimensions</b> (Longest Mattress Platform Position)	1050 x 2370 mm
<b>Mattress Dimensions</b>	2080 x 860 mm
<b>Maximum Mattress Height</b>	230 mm
<b>Height Adjustment</b>	440 - 820 mm
<b>Lateral Tilt</b>	±30°
<b>TR / ATR Angle</b>	13° / 16°
<b>Bed Weight</b>	210 kg
<b>SWL</b>	250 kg





## Nino Adolescent Bed

The Nino Series allows children to perceive the world around them while at the same time providing protection and security. The transparent Perspex merges with the room allowing occupants to have an unrestricted view of their immediate surroundings yet still protects against falls and injuries. Also available in different powdercoat colours, or with higher door elements or doors on both sides.

- Electric Hi-low
- Electric Backrest and Kneebreak
- Standard Locking Castors
- Gates one side with Perspex infills
- Fixed Panel one side with Perspex Infills
- 750mm High Sides (from Mattress Base)
- Head & Foot Infills in Natural Oak Matt Finish
- Battery Backup
- Fully customisable Bed



NDIS

### Specifications:

<b>Overall Dimensions</b>	2200 x 1014 mm
<b>Sleeping Surface Dimensions</b>	2000 x 870 mm
<b>NET Bed Weight</b>	140 kg
<b>Underbed Clearance</b>	100 mm (at lowest point)
<b>Low Position</b>	390 mm
<b>High Position</b>	830 mm
<b>SWL</b>	200 kg







## Tom 2 Cot

A state-of-the-art children's cot must be safe and be able to provide excellent access for nursing staff. The electric positioning of TOM 2 children's cot responds to these and further demands in its original concept of telescopic siderails, column construction mechanism and above standard safety parameters.

The Tom 2 bed fully complies with the new standard for hospital children's beds, DIN 32623:2009-11 and its construction guarantees perfect protection against children fall. Also complies with standards: EN 60601-1, EMC 60601-1-2, EMC 60601-1-6

- 3 height positions for side rails, 0cm, 40cm, 80cm
- Backrest retracts when it raises, meaning zero gap
- Integrated collapsible lockable side rails

- 4 x Tente Motion 125mm castors
- Trendelenburg tilt buttons
- Central operated castors (one castor with steer system)
- Undercarriage cover
- Quick removable head/foot ends
- Accessory holders: Pair of accessory rails (located on both sides under side rails)
- Battery Backup
- Horizontal Wall/corner bumpers
- Quality European product



NDIS

### Specifications:

Product Code	76765
Mattress Platform	1370 x 700 mm
Electric Hi-Lo	580 – 830 mm
Electric Tilt/Reverse Tilt	15°
Backrest Adjustment	40°
SWL	75 kg

### Accessory Available:

Tom 2 Cot Mattress

Product Code 76766



NDIS



**VT2 Bedside Cabinet**

The tried-and-tested details of the bedside cabinets are quite impressive. Stable steel frames, high-quality ABS edges on all components of the bedside cabinet, as well as the rounded top board give the bedside cabinet its unmistakable character.

- The soft closing drawers are a standard feature that minimize noise and the danger of entrapment

**Colour options**



Beech Natural

Lindberg Oak

Havana Cherry

**Specifications:**

<b>Width</b>	506 mm
<b>Length</b>	779 mm
<b>Depth</b>	470 mm

NDIS



**Staebel Mobile Overbed Table**

A very popular mobile table with a number of useful features, making it an ideal choice in many situations. With its easy lever operated height adjustment and smooth curved U-design, it fulfils all needs.

- Spring assisted lever height adjustment – lever stays with top eliminating need to bend to adjust table height
- Powdercoated U-shape base, round tube minimises injury risk
- Swivel castors, two are lockable

\* Table top is sold separately. There are a range of colours on offer



**Specifications:**

<b>Product Code</b>	4215
<b>Width</b>	450 mm
<b>Length</b>	750 mm

NDIS



**Solido 3 Mobile Overbed Tilting Table**

The Solido 3 is a solid, stable, robust mobile table, specifically designed to address the demands of a hospital environment. Sourced from Europe, and made from the highest quality materials.

- Moulded table top covers infection control requirements and is easy to clean
- Tilting table top with fixed side piece
- 50mm castors gives low base height of 95mm to clear under beds
- Solid table top with 30kg safe working load
- Tilted positions - 0°, 20°, 40°, -20°, -40°
- Easy to use gas spring adjustment to adjust height (82-117cm)
- Generous table surface 96x45cm with recessed edge



**Specifications:**

<b>Product Code</b>	77821
<b>Table Length</b>	960 mm
<b>Table Width</b>	450 mm
<b>Height</b>	820 - 1170 mm
<b>SWL</b>	30 kg



DVA  
NDIS



DVA  
NDIS  
MASS

### Staebel Deluxe Commode Bedside

This bedside commode offers a full range of features including a slide back seat with durable toilet seat and pan rack. This comes standard with bowl and lid.

A skirt that surrounds the seat conceals the bowl etc. so the chair can be likened to a regular bedroom sitting chair

- Adjustable seat height
- Comfortable upholstered slide back seat
- Comfortable upholstered shallow curve backrest
- Plastic toilet seat, bowl and lid
- Padded armrests
- Powdercoated steel frame
- Vinyl upholstery



#### Specifications:

Product Code	7955
Overall Width	600 mm
Overall Height	870 - 1020 mm
Overall Depth	510 mm
Seat Width	460 mm
Seat Height	420 - 570 mm
Seat Depth	440 mm
Back Height	470 mm
SWL	120 kg

#### Also Available:

Product Code	7960	7961	7962	7963
Seat Width	500 mm	550 mm	600 mm	650 mm
SWL	135 kg	175 kg	195 kg	210 kg

### Staebel Standard Bedside Commode

The standard bedside commode offers padded armrests for comfort and features adjustable height, a removable padded seat cover and a plastic toilet seat with bowl and lid.

- Value for money
- Stackable design



#### Specifications:

Product Code	7952PA
Overall Width	575 mm
Overall Height	940 - 1090 mm
Overall Depth	600 mm
Seat Width	440 mm
Seat Height	500 - 650 mm
Seat Depth	440 mm
Back Height	475 mm
Toilet Seat	475 - 625 mm
SWL	125 kg

NDIS



### Bed Rail

The bed rail is designed to slide under the mattress and act as a gentle barrier to stop the patient falling out the side of the bed. Rail 500 mm high x 600 mm wide

NDIS



### Swivel Hand Control Holder

An accessory suitable for use on any bed from the Wissner-Bosserhoff range this handset holder can be mounted on either side of the bed for easy access to the handset by care staff or the client. Significantly reduces the risk of handset breakage due to handset falling on the floor. Quality German-made and engineered. Must be purchased along with 77611 or 77636 adaptor

### Specifications:

Product Code 7905



DVA

NDIS

### Blanket Support

The bed cradle provides relief from the weight of bed clothes on delicate feet and legs.

- Provides clearance between bed linen and the mattress
- Ideal for keeping bed linen of sensitive areas of the patient
- Slides between mattress and bed base
- Powder coated mild steel construction
- DVA: 2010

### Specifications:

Product Code 7920

### Also Available:

Blanket Support - Adjustable

Product Code 7921

### Specifications:

Product Code 77612S

NDIS



### Bed Pole

This single bed pole is used to help with bed mobility. It simply slides under the mattress. Where a double bed stick is required and the bed is an unusual width two of these single poles can be interlocked under the mattress. Pole height 510 mm.

### Specifications:

Product Code 7900

Height 510 mm

### Also Available:

Double Bed Stick


Product Code 7902



DVA

NDIS

NDIS

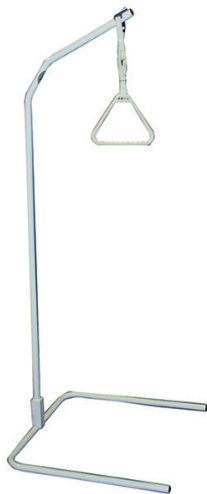

## Overbed Frame & Triangle

The strong steel frame is free standing and does not require any wall or floor fittings. Simply slides under the bed. Provided with an adjustable length strap making the triangle easy to reach.

- Ideal for enhancing bed mobility and self repositioning
- Suitable for users up to 130kg



### Specifications:

Product Code	7935
Overall Width x Length	1150 x 2230 mm
Height	500 - 750 mm
SWL	130 kg



## Staebel Bed Rail Protector

This Bed Rail is made from a polyurethane material, coupled with a bioelastic incontinence cover, making it the ultimate answer for longevity and effectiveness.

No zips, clips, or flimsy foam means a zero risk of skin tears or wounds caused by the patient/resident hitting the side guards. Simply push it down onto the top rail and it is in place.

No time consuming, technical attachment points that foul up with the side guards.

- Skin friendly
- Anti Microbial safety
- Easy to clean and to disinfect
- Easy to use and install
- 2 x protectors can be used together to cover the complete length of a 4-part side guard or full length side rails.
- Sold singularly



### Use with Wi-Bo Bed Mobilisation Aids




### Specifications:

Product Code	3983
Width	130 mm
Length	450 mm
Height	800 mm

# 08.

# Patient Safety

Staebel Socks	135	Cura1 Cordless Bed Pad with Transmitter	137
Posey Non Skid Slippers	135	Cura1 Cordless Chair Pad with Transmitter	137
Cura1 Cordless Floor Mat with Transmitter	135	Cura1 Wristband Alarm	138
iBeam Patient Monitor	136	Cura1 Personal Emergency Transmitter (PET)	138
Ezy Stand	136	Cura1 Magnetic Reed Switch	138
Cura1 6 Channel Cordless Monitor	136	Cura1 Wireless Handset	138
Cura1 LED Pager	137		

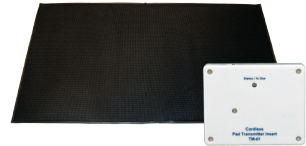

## Staebel Socks

Staebel socks provide safety and prevention of falls.

- Soft terry cloth upper. Latex-free.
- Colour helps care staff identify high fall risk patients
- Grip on top and bottom
- Mesh band on large socks comfortably fits swollen ankles and calves
- Standard size fits up to adult size 13

Available in Green, Orange, Purple, Red or Yellow in sizes Standard and Large

## Posey Non Skid Slippers

Posey fall management slippers provide safety and reduce the risk of falls.

- Soft terry cloth upper
- Latex-free rubber sole
- Colour helps care staff identify high fall risk patients
- Machine washable
- All colours available in five sizes
- One-pair per package

Available in Green, Orange, Purple, Red or Yellow

Please call 1800 267 267 for full range of colour and size product codes.

## Cura1 Cordless Floor Mat with Transmitter

This cordless floor mat is positioned beside a resident's bed to alert staff when the resident has put their feet on the floor, attempting to leave the bed.

- Simple operation
- Easy to clean
- Virtually no false alarms over the lifetime of the pad
- Made from heavy duty non-slip material
- Generous size compared to others on the market
- Quality Australian made product
- 12 month life expectancy
- 50 meter range (\*Subject to environmental conditions)
- Includes transmitter
- Must be used with 92700 Monitor



### Specifications:

<b>Small</b>	9900
<b>Medium</b>	9902
<b>Large</b>	9904

### Specifications:

<b>Product Code</b>	9PSFP	
<b>Size</b>	Ladies (Shoe Size)	Mens (Shoe Size)
<b>Small</b>	5 - 6	N/A
<b>Medium</b>	6 - 7	5 - 6
<b>Medium / Large</b>	7 - 8	6 - 7
<b>Large</b>	8 - 9	7 - 9
<b>Extra Large</b>	10 +	10 - 12

### Specifications:

<b>Product Code</b>	92663
<b>Width</b>	610 mm
<b>Length</b>	1220 mm

NDIS



### iBeam Patient Monitor

The narrow beam creates a precise detection zone. This is especially useful for preventing unwanted alarms if other people frequent the same room, or where only a narrow area requires monitoring.


Whether across a door, stairs or beside a bed, there are countless ways that our iBeam can be used to protect those in your care and to make the job of being a carer easier.

When used with our sturdy, adjustable-height EzyStand, precise positioning of the iBeam could not be simpler.

**Specifications:**

Product Code	927413
--------------	--------

NDIS

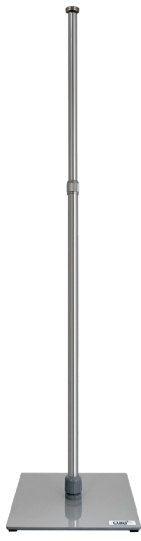

### Ezy Stand

The EzyStand / iBeam combination allows a caregiver to quickly deploy targeted motion detection alarms in a variety of positions such as next to a bed, in front of a chair, near a door or stairs.

Our innovative 2723 EzyStand magnetically attaches to the 2721 iBeam and significantly adds to its ease of use.

The stand is sturdy, attractive and perfect for use in both professional and home environments.

**Specifications:**

Product Code	92723
--------------	-------

NDIS




### Cura1 6 Channel Cordless Monitor

No more cords between monitor and pads. Can be used with up to six Cordless pads, mats or other triggering devices and one corded (Premium) device. Modern ergonomic design

- Supplied with batteries, durable ABS plastic enclosure and Caregiver Reset Key
- Alarm memory and auto re-arm
- Amber LED visual alarm indicator
- Caregiver secure reset option
- Battery operated or optional external plug pack
- Up to 6 months battery life (nominal)
- Easy to access battery compartment
- Can be connected to a nurse call system
- 3 selectable alarm tones
- Easy thumbwheel adjustable alarm volume

**Specifications:**

Product Code	92700
Width	300 mm
Height	770 mm
Length	1220 mm

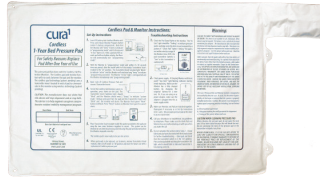

## Cura1 LED Pager

- Can be programmed to monitor up to 3 wireless devices
- Vibration and tone selection for hard of hearing
- Auto-reset when paired with a 92719 transmitter
- Direct signalling from wireless devices for use in home care environments
- Easily cleanable buttons for a hospital setting
- Power requirements: 2 x AA batteries (included)

## Cura1 Cordless Bed Pad with Transmitter

- Pad includes pouch that contains a transmitter
- Guaranteed 1 year life expectancy
- Easy to clean
- Non-allergenic soft PVC
- For use with cura1 cordless bed monitors
- 50 meter range (\*Subject to environmental conditions)
- Vent holes in the neck reduce the risk of moisture -ingress
- Must be used with 92700 Monitor

## Cura1 Cordless Chair Pad with Transmitter

- Simple operation and maintenance
- Includes transmitter
- For use with Cura1 cordless monitor
- Easy to clean
- Virtually no false alarms over the lifetime of the pad
- All bed/chair pads made from non-allergenic soft PVC
- 12 month life expectancy
- Must be used with 92700 Monitor

### Specifications:

Product Code 93508

Width 550 mm

Height 830 mm

Length 260 mm

### Specifications:

Product Code 92662

Width 760 mm

Length 510 mm

### Specifications:

Product Code 92661

Width 380 mm

Length 255 mm

NDIS



## Cura1 Wristband Alarm

- For use with cura1 StaffAlert or LRCS Duress System
- Lightweight and small in size
- Low cost and low maintenance
- Replaceable battery
- Includes buckle strap. Optional lanyard available
- Medical Device Immunity: EN 60601-1-2

NDIS



## Cura1 Personal Emergency Transmitter (PET)

- Large rubber button for easy use by elderly patients
- Standard lanyard for patient to wear around neck
- Works with single use hospital wristbands
- Changeable battery with projected 1 year life
- Splash Resistant

### Specifications:

Product Code 92595



NDIS

## Cura1 Magnetic Reed Switch

- Sends wireless signal when door or window is opened
- Built-in call button extends usage options - eg. use it as a portable panic button
- Double sided tape included for easy mounting
- Battery included
- LED alarm active indication
- Weight: 28 g

### Specifications:

Product Code 93503

### Specifications:

Product Code 92588



NDIS

## Cura1 Wireless Handset

- Green cancel button resets alarm\*
- Large rubber button for easy use by elderly patients
- Reassurance LED when transmission is sent
- Red CALL button to trigger alarm for caregiver
- Lanyard clip allows patient to connect with clothes
- Wall mount bracket for when handset is not in use
- Power requirements: 2 x AAA alkaline batteries

### Specifications:

Product Code 927001











We are  
*Here to Help*



# Let's have a chat!

If you see something you like, you can book in your consultancy with one of our friendly team members. See below information for ways you can get in contact...

1800 267 267 | [www.activemedicalsupplies.com.au](http://www.activemedicalsupplies.com.au) | [hello@activemedicalsupplies.com.au](mailto:hello@activemedicalsupplies.com.au)

**activeMEDICAL**  
SUPPLIES > SOLUTIONS > SUPPORT

**P.** 1800 267 267 | **F.** 07 3899 8825

**E.** [hello@activemedicalsupplies.com.au](mailto:hello@activemedicalsupplies.com.au)

---

**[www.activemedicalsupplies.com.au](http://www.activemedicalsupplies.com.au)**