

Relax DuoGel Cushion



**The Relax Duogel effectively
provides pressure relief &
weight distribution.**



Relax DuoGel Cushion

Product Features:


The Relax DuoGel is a unique combination of TruGel and FloGel working together effectively to provide excellent pressure relief and weight distribution. The DuoGel's anatomical pre-contoured base provides for medial and lateral thigh positioning as well as support for the greater trochanter of the thigh bone. The thigh troughs are lined with TruGel and the ischial area is lined with FloGel which is set on a viscoelastic foam giving greater comfort and seating stability.

Product Benefits:

- Thigh troughs assist with seating stability
- Fast Recovery Foam (Anatomically Shaped)
- Includes machine washable outer cover (not pictured)
- Clear protective PU film coating - easy to clean and suitable for incontinent users
- Anatomical Pre-Contoured Foam Base

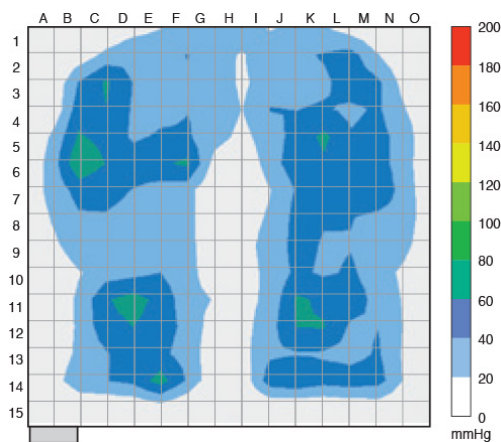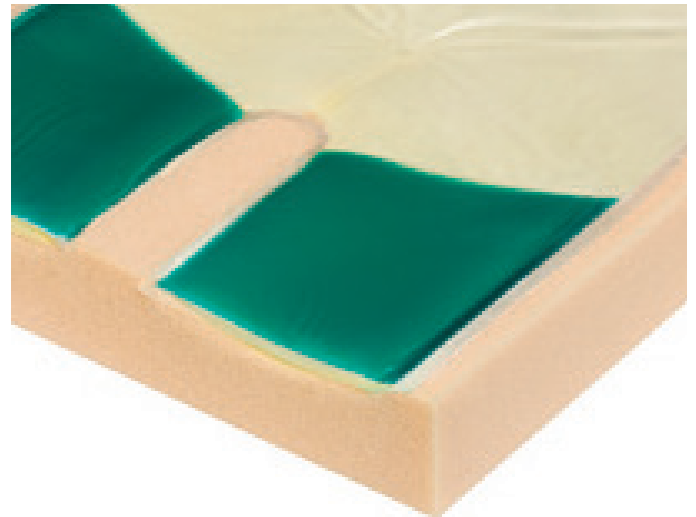

Pressure diagram with Relax DuoGel cushion

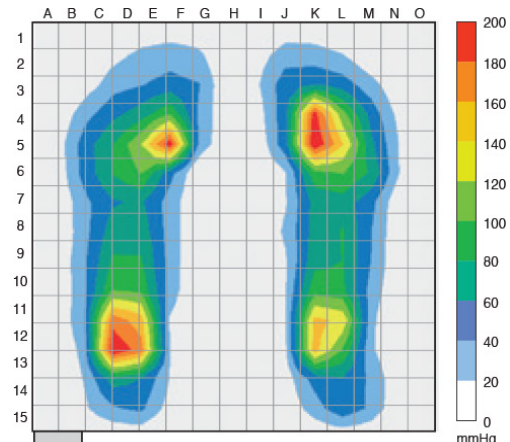

Pressure diagram without wheelchair cushion

Specifications:		Sizes Available:					
Product Code	35302	Product Code	35300	35306	35303	35304	35307
Width	450 mm	Width	400 mm	400 mm	450 mm	500 mm	500 mm
Length	400 mm	Length	400 mm	450 mm	450 mm	450 mm	500 mm
Thickness	80 mm	Thickness	80 mm	80 mm	80 mm	80 mm	80 mm
SWL	130 kg	SWL	130 kg	130 kg	130 kg	130 kg	130 kg